

manuals.plus /

› [SKINACT](#) /

› [SKINACT Bellage Electric Treatment Table \(Chair\) 4 Electric Adjustable Motors \(Gray\) - User Manual](#)

## SKINACT SKU164247

# SKINACT Bellage Electric Treatment Table (Chair)

Model: SKU164247

## INTRODUCTION

Thank you for choosing the SKINACT Bellage Electric Treatment Table. This versatile and durable electric chair is designed for professional use in spa and medical environments, offering exceptional comfort and adjustability. Equipped with four independent electric motors, it allows for precise positioning of the backrest, leg rest, seat tilt, and overall height, ensuring optimal client comfort and practitioner convenience. This manual provides essential information for the safe and effective setup, operation, and maintenance of your new treatment table.



Figure 1: Overview of the SKINACT Bellage Electric Treatment Table with key features labeled.

## SAFETY INFORMATION

Please read all safety warnings and instructions carefully before assembly and operation. Failure to do so may result in injury or damage to the product. Ensure the table is placed on a stable, level surface. Keep hands and loose clothing clear of moving parts during operation. Do not exceed the maximum weight capacity of 350 lbs. Disconnect power before performing any maintenance or cleaning.

# SETUP

---

## 1. Unpacking and Placement

- Carefully remove all components from the packaging.
- Place the base of the treatment table on a firm, level floor. Ensure there is sufficient clear space around the table for full range of motion and safe operation.
- Verify all leveling feet on the base are properly adjusted to ensure stability.

## 2. Component Assembly

- **Attach Headrests:** The Bellage table comes with two different headrests. Select the appropriate headrest for your procedure and securely attach it to the head section of the table. One headrest includes a breather hole for clients lying face down.
- **Attach Armrests:** Slide the removable armrests into their designated slots on the sides of the main body of the chair.
- **Foot Rest Extension:** The foot rest section can be extended as needed.



*Figure 2: The two interchangeable headrests provided with the table.*


### **3. Power Connection**

- Locate the power cord and plug it securely into the power input on the table's base.
- Plug the other end of the power cord into a grounded electrical outlet (100-120 volts, 50-60 hertz).

### **4. Remote Control Connection**

- Connect the remote control cable to the designated port on the underside of the table.

**Adjustable height:  
26"-34"**

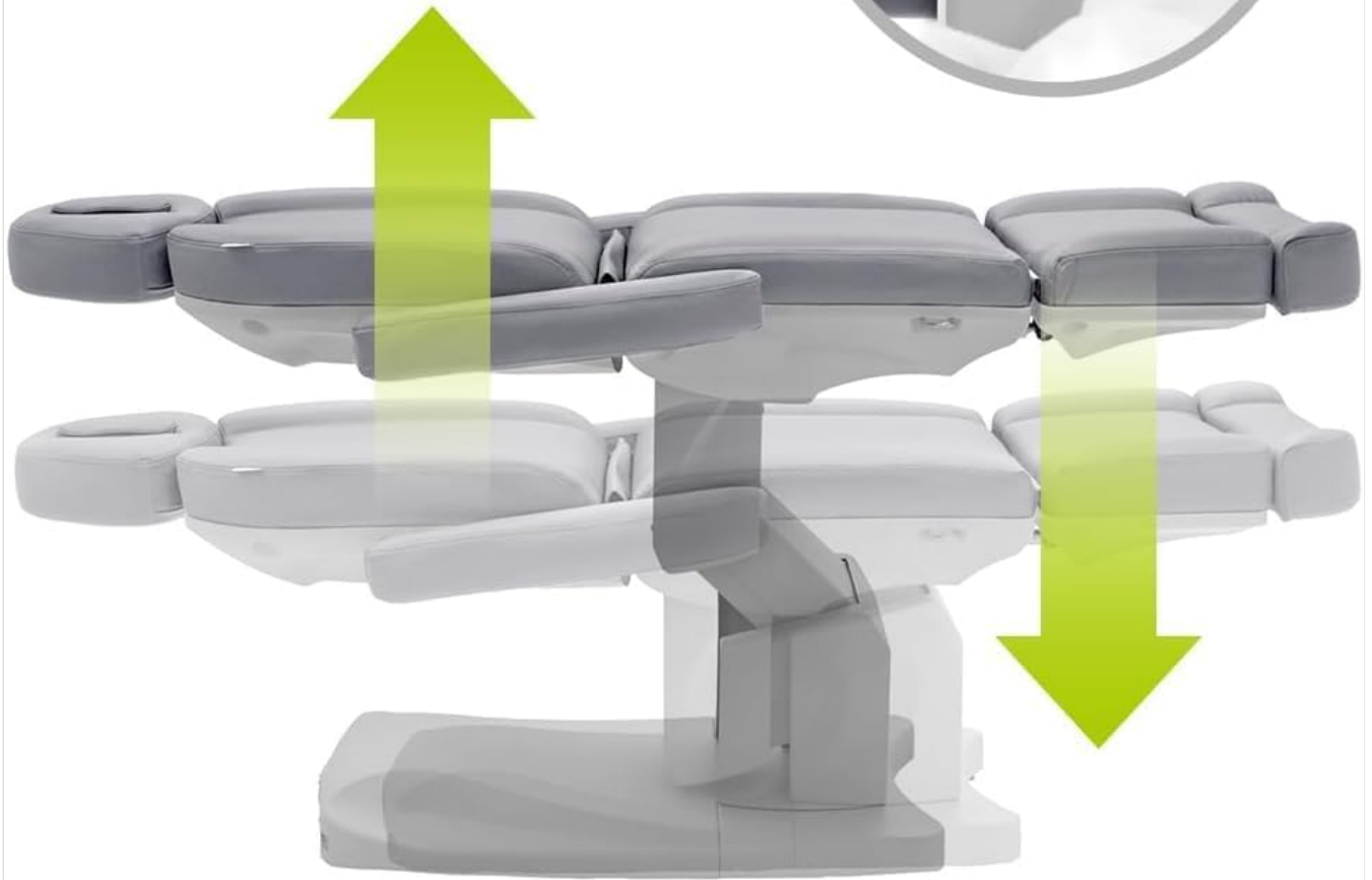

*Figure 3: The wired remote control for the Bellage Electric Treatment Table.*

## Setup Demonstration Video

*Video 1: A demonstration of the Bellage Electric Treatment Table's features and adjustments.*

## OPERATING INSTRUCTIONS

The SKINACT Bellage Electric Treatment Table is operated using the wired remote control. Each button on the remote corresponds to a specific motor, allowing for independent adjustment of different sections of the table.

### 1. Power On/Off

- Press the **Power** button (usually indicated by a circle with a vertical line) on the remote to turn the table on or off.

### 2. Adjusting Height

- Use the **Up/Down arrows** associated with the overall table icon to raise or lower the table's height. The height can be

adjusted between 25.5 and 31.5 inches.

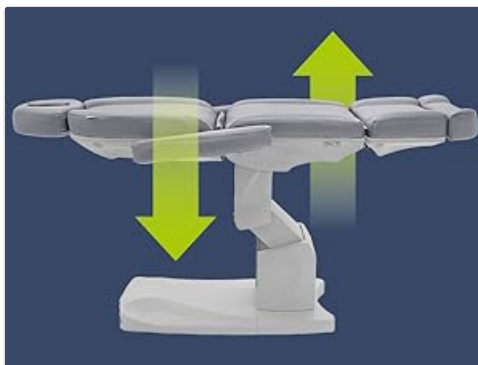


Figure 4: Electric height adjustment.

### 3. Adjusting Backrest Angle

- Use the **Up/Down arrows** associated with the backrest icon to adjust the backrest angle from 0° to 80° incline.



Figure 5: Electric backrest adjustment.

### 4. Adjusting Leg Rest Angle

- Use the **Up/Down arrows** associated with the leg rest icon to adjust the leg rest angle from 0° to 85° decline.



Figure 6: Electric leg rest adjustment.

### 5. Adjusting Bed Tilt Angle

- Use the **Up/Down arrows** associated with the full bed tilt icon to adjust the entire bed's tilt angle from 0° to 10°.



Figure 7: Electric bed tilt adjustment.

## 6. Memory Presets (P1, P2, P3)

- The remote control may feature preset buttons (P1, P2, P3) to save and recall preferred positions. Consult the remote's specific instructions for programming these presets.

## MAINTENANCE

Regular maintenance ensures the longevity and optimal performance of your SKINACT Bellage Electric Treatment Table.

- **Cleaning:** Wipe the surfaces of the table with a dry cloth. For deeper cleaning, use a mild, non-abrasive cleaner and a soft cloth. Avoid harsh chemicals or abrasive materials that could damage the upholstery or frame.
- **Inspection:** Periodically inspect all electrical cords, connections, and moving parts for any signs of wear, damage, or loose fittings. If any issues are found, discontinue use and contact customer support.
- **Lubrication:** The electric motors and moving mechanisms are generally maintenance-free. Do not attempt to lubricate internal components unless specifically instructed by a qualified technician.

## TROUBLESHOOTING

If you encounter any issues with your treatment table, refer to the following common troubleshooting steps:

Problem	Possible Cause	Solution
Table does not respond to remote control.	No power; Loose connection; Remote control malfunction.	Check if the power cord is securely plugged into the outlet and the table. Ensure the remote control cable is firmly connected. Try a different outlet. If the issue persists, contact customer support.
Motor makes unusual noise or movement is jerky.	Obstruction; Motor issue.	Check for any obstructions preventing smooth movement. Ensure the table is not overloaded. Discontinue use and contact customer support if the noise or jerky movement continues.
Table is unstable or wobbles.	Uneven floor; Unadjusted leveling feet.	Ensure the table is on a flat, level surface. Adjust the leveling feet on the base until the table is stable.

## SPECIFICATIONS

Feature	Detail
---------	--------

Feature	Detail
Brand	SKINACT
Model Number	SKU164247
Material	Metal frame, plastic cover, Foam, PU upholstery
Color	Dark Gray
Item Weight	200 pounds
Product Dimensions (L x W x H)	68 x 28 x 32 inches
Length (with headrest)	76.8 inches
Height Range	25.5 / 31.5 inches
Width (with armrests)	35.4 inches (24.8 inches without armrests)
Backrest Angle	0° to 80° incline
Bed Tilt Angle	0° to 10°
Leg Rest Angle	0° to 85° decline
Power Requirements	100-120 volts, 50-60 hertz, 200 watts, 2 amp
Maximum Weight Capacity	350 lbs.

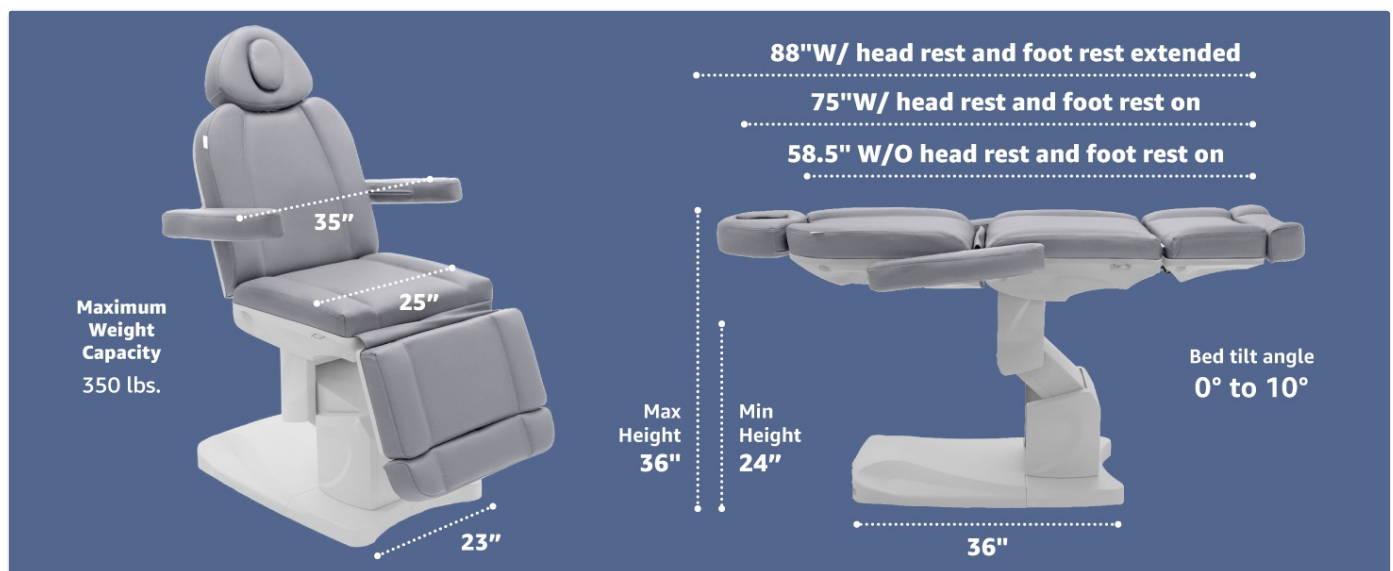


Figure 8: Detailed dimensions of the treatment table.

## WARRANTY & SUPPORT

The SKINACT Bellage Electric Treatment Table comes with a **ONE YEAR WARRANTY** from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use.

For technical support, warranty claims, or any questions regarding your product, please contact SKINACT customer service. Refer to your purchase documentation for specific contact details or visit the official SKINACT website.

## Related Documents - SKU164247

### Slate Metal Cart



SKIN•ACT

Toll-free: 1.800.416.0830	<a href="http://www.800CT.COM">www.800CT.COM</a> <a href="http://www.greenequipment.com">www.greenequipment.com</a>	International: 818.834.6640
---------------------------	--	-----------------------------

## SkinAct Slate Metal Cart - Product Overview and Specifications

Detailed overview and specifications for the SkinAct Slate Metal Cart, a durable all-metal utility trolley designed for medical, spa, and salon environments. Features include ample storage, smooth-rolling casters, and a sleek design.