

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- > [CorDx](#) /
- > [CorDx Ovulation + Early Detection Pregnancy Test Combo \(Model ACT21013\) Instruction Manual](#)

CorDx ACT21013

CorDx Ovulation + Early Detection Pregnancy Test Combo Instruction Manual

Model: ACT21013

1. PRODUCT OVERVIEW

The CorDx Ovulation + Early Detection Pregnancy Test Combo provides a comprehensive solution for family planning. This kit includes both Ovulation Test Strips and Early Detection Pregnancy Test Strips, designed for at-home urine testing. The ovulation tests help identify your most fertile days by detecting Luteinizing Hormone (LH) surges, while the pregnancy tests detect Human Chorionic Gonadotropin (hCG) to confirm pregnancy as early as the first day of a missed period. Both tests offer over 99% accuracy when used correctly.



Image: The product packaging for the CorDx Ovulation + Pregnancy Test Combo, showing the brand name, product type, and quantity of tests included.

2. WHAT'S INCLUDED

Your CorDx Ovulation + Early Detection Pregnancy Test Combo kit contains the following items:

- 50 individually sealed Ovulation Test Strips
- 20 individually sealed Early Detection Pregnancy Test Strips
- Instruction Manuals for both test types
- Sample Cups (quantity may vary)



Image: A visual representation of the kit contents, including the box, individual ovulation and pregnancy test strips, instruction manuals, and sample cups.

3. IMPORTANT SAFETY INFORMATION

- For **in vitro diagnostic use only**. Do not ingest.
- Keep out of reach of children.
- Do not use the test strips after the expiration date printed on the foil pouch.
- Store tests in their original sealed pouches at room temperature (4-30°C or 39-86°F). Do not freeze.
- Do not use if the foil pouch is damaged or opened.
- Each test strip is for single use only. Discard after use.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is for diagnostic purposes.

4. SETUP AND PREPARATION

Before performing either test, ensure you have:

- The CorDx Ovulation or Pregnancy Test Strip.
- A clean, dry urine collection cup.
- A timer or clock.

Allow the test strip to reach room temperature (18-30°C or 64-86°F) before opening the pouch. Collect a urine sample in a clean, dry container. For ovulation tests, it is recommended to test around the same time each day,

but avoid using first morning urine as LH levels may not be concentrated enough. For pregnancy tests, first morning urine is recommended as it typically contains the highest concentration of hCG.

5. OPERATING INSTRUCTIONS

5.1. Ovulation Test Instructions

1. Open the sealed foil pouch and remove the test strip.
2. Immerse the strip vertically into the urine sample with the arrow pointing towards the urine. **Do not immerse past the MAX line.**
3. Hold the strip in the urine for **5 seconds**.
4. Remove the strip from the urine and lay it flat on a clean, dry, non-absorbent surface.
5. Read the results between **3 to 8 minutes**. Do not interpret results after 10 minutes.

INSTRUCTIONS FOR OVULATION TEST

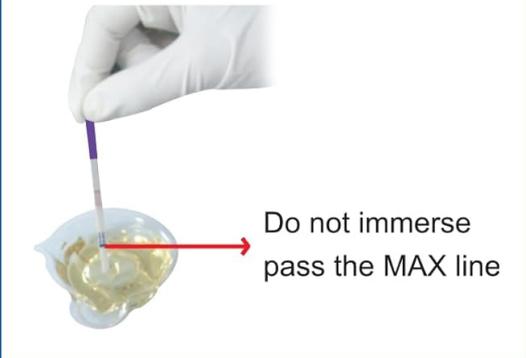
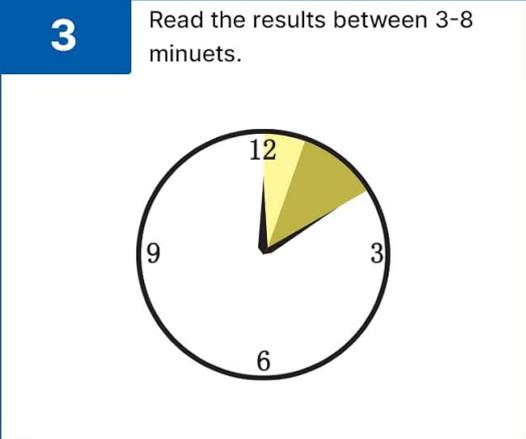
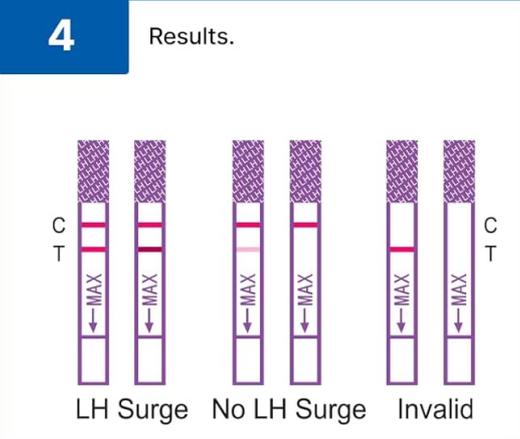
1	Immerse the strip in the sample for 5 seconds.	2	Lay the strip flat.
			
3	Read the results between 3-8 minutes.	4	Results.
			

Image: Step-by-step visual instructions for performing the ovulation test, including immersion time, maximum immersion line, and result reading time.

5.2. Pregnancy Test Instructions

1. Open the sealed foil pouch and remove the test strip.
2. Immerse the strip vertically into the urine sample with the arrow pointing towards the urine. **Do not immerse past the MAX line.**

3. Hold the strip in the urine for **5 seconds**.
4. Remove the strip from the urine and lay it flat on a clean, dry, non-absorbent surface.
5. Read the results at **5 minutes**. Do not interpret results after 10 minutes.

INSTRUCTIONS FOR PREGNANCY TEST

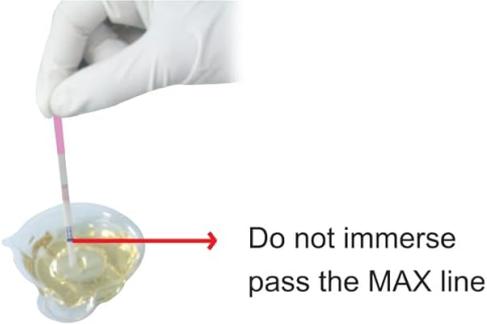
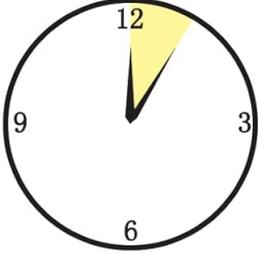
<p>1 Immerse the strip in the sample for 5 seconds.</p> 	<p>2 Lay the strip flat.</p> 
<p>3 Read the results at 5 minutes.</p> 	<p>4 Results.</p> 

Image: Step-by-step visual instructions for performing the pregnancy test, including immersion time, maximum immersion line, and result reading time.

6. INTERPRETING RESULTS

6.1. Ovulation Test Results

- **LH Surge (Positive Result):** If two color bands are visible, and the test band (T) is equal to or darker than the control band (C), an LH surge is detected. Ovulation is likely to occur within 24-36 hours.
- **No LH Surge (Negative Result):** If only the control band (C) appears, or if the test band (T) is lighter than the control band (C), no LH surge is detected.
- **Invalid Result:** If no control band (C) appears, the test is invalid. Repeat with a new test strip.

6.2. Pregnancy Test Results

- **Pregnant (Positive Result):** If two distinct color bands appear, one in the control region (C) and one in the test region (T), the result is positive. The color intensity of the test band (T) may vary depending on the hCG concentration.

- **Not Pregnant (Negative Result):** If only one color band appears in the control region (C), and no band appears in the test region (T), the result is negative.
- **Invalid Result:** If no control band (C) appears, the test is invalid. Repeat with a new test strip.

CLEAR RESULTS, FAST

Get results in just minutes, helping you take the next steps confidently in your conception journey.



Image: Visual representation of the time required to read results for both pregnancy (5 minutes) and ovulation (3-8 minutes) tests.

7. MAINTENANCE AND STORAGE

To ensure the longevity and accuracy of your CorDx test strips:

- Store the sealed pouches at room temperature (4-30°C or 39-86°F).
- Keep away from direct sunlight, moisture, and heat.
- Do not freeze the test strips.
- Use the test strip immediately after opening the foil pouch.

8. TROUBLESHOOTING

8.1. Invalid Result

If no control band (C) appears, the test is invalid. This may be due to insufficient urine sample, incorrect testing procedure, or a damaged test strip. Discard the invalid test and repeat the procedure with a new test strip,

carefully following all instructions.

8.2. Faint Test Line

For pregnancy tests, a faint test line (T) still indicates a positive result. The intensity of the line can vary with the concentration of hCG in the urine. For ovulation tests, a faint test line (T) that is lighter than the control line (C) indicates no LH surge. If you are unsure, consider retesting in 24-48 hours or consulting a healthcare professional.

9. SPECIFICATIONS

Product Name	CorDx Ovulation + Early Detection Pregnancy Test Combo
Model Number	ACT21013
Product Dimensions	7 x 2 x 3 inches
Product Weight	3.2 ounces
Manufacturer	CorDx
Test Type	Ovulation (LH) and Pregnancy (hCG) Urine Test Strips
Accuracy	Over 99%
Results Time	Ovulation: 3-8 minutes; Pregnancy: 5 minutes

10. WARRANTY AND SUPPORT

CorDx products are manufactured to high-quality standards. For any concerns regarding product performance or results, please refer to the contact information provided on the product packaging or visit the official CorDx website.

Return Policy

This product is typically eligible for a 30-day return policy (refund/replacement) from the date of purchase. Please check with your retailer for specific return conditions.

Manufacturer Information

Manufacturer: CorDx

For additional support or inquiries, please visit the [CorDx Store on Amazon](#).