

VERMEIREN VrmJazz20

VERMEIREN JAZZ S20 Economy Foldable Wheelchair User Manual

Model: VrmJazz20

1. INTRODUCTION

This manual provides essential information for the safe and effective use, maintenance, and care of your VERMEIREN JAZZ S20 Economy Foldable Wheelchair. The Jazz S20 is a lightweight, sturdy, and foldable wheelchair designed for mobility assistance, featuring ergonomic handgrips, adjustable and removable armrests, and adaptable leg rests. It is engineered in Belgium and manufactured with high-quality materials to ensure durability and user comfort.

Please read this manual thoroughly before operating the wheelchair to ensure proper handling and to maximize its lifespan.

2. SETUP AND ASSEMBLY

The VERMEIREN JAZZ S20 wheelchair is designed for easy setup. Follow these steps to prepare your wheelchair for use:

- Unfolding the Wheelchair:** Place the folded wheelchair on a flat, stable surface. Gently pull the armrests outwards until the frame fully expands and locks into place. Ensure the seat canvas is flat and taut.
- Attaching Footplates:** Insert the leg rest assemblies into the receivers on the front of the frame. Swing the footplates into position and ensure they are securely locked. Adjust the height of the footplates to ensure comfortable and proper leg positioning for the user.
- Adjusting Armrests:** The armrests are designed for comfort. Ensure they are securely in place. If removable, verify they are properly seated before use.
- Brake Check:** Before first use, test both parking brakes to ensure they engage and disengage smoothly and hold the wheels firmly.



Image: The VERMEIREN JAZZ S20 wheelchair fully unfolded, showcasing its sturdy frame and ready-to-use configuration.



Foldable for **easy**
Transport or **Storage**

Image: Demonstrates the wheelchair's foldable design, making it convenient for transport and storage.

3. OPERATING INSTRUCTIONS

Proper operation ensures the safety and comfort of the user.

3.1. Braking System

- **Engaging Brakes:** To engage the parking brakes, push the brake levers forward until they lock into place against the wheels. Always engage both brakes when the wheelchair is stationary, especially during transfers.
- **Disengaging Brakes:** Pull the brake levers backward to release the brakes before moving the wheelchair.



Advanced features



Image: Detailed view of the high-quality brakes and premium footplates, highlighting key operational components.

3.2. User Transfer

- Always engage both parking brakes before attempting to transfer a user into or out of the wheelchair.
- Ensure the footplates are swung away or removed to clear the path for the user's feet.
- If possible, position the wheelchair as close as possible to the transfer surface (bed, chair, etc.).

3.3. Attendant Operation

- Use the ergonomic handgrips at the rear of the wheelchair to push and maneuver.
- Always be aware of the user's comfort and safety, especially when navigating uneven terrain or inclines.



Additional features



Image: Focus on the ergonomic handgrips for comfortable attendant operation and the fire-retardant nylon material of the seat and backrest.

4. MAINTENANCE AND CARE

Regular maintenance will prolong the life of your wheelchair and ensure its safe operation.

- **Cleaning:** Clean the frame with a damp cloth and mild detergent. The fire-retardant nylon material can be wiped clean. Avoid abrasive cleaners.
- **Wheel and Tire Inspection:** Regularly check the polyurethane tubeless tires for wear and tear. Ensure the stainless steel spokes are intact and the wheels spin freely.
- **Brake Adjustment:** Periodically check the brakes for proper function. If they feel loose or do not engage effectively, consult a qualified technician for adjustment.
- **Fasteners:** Check all nuts, bolts, and screws regularly to ensure they are tight. Do not overtighten.
- **Storage:** Store the wheelchair in a dry place, away from direct sunlight and extreme temperatures. The

foldable design allows for easy storage in compact spaces.



Image: Highlights durable components such as stainless steel spokes, tubeless tires, and the overall powder-coated glossy finish for easy maintenance.

5. TROUBLESHOOTING

This section addresses common issues you might encounter with your wheelchair.

Problem	Possible Cause	Solution
Wheelchair does not roll smoothly.	Debris in wheels, loose axle, worn tires.	Clean wheels, check axle tightness, inspect and replace worn tires if necessary.
Brakes do not hold the wheels.	Brake pads worn, brake mechanism loose.	Inspect brake pads for wear. Tighten or adjust brake mechanism. If unsure, contact support.

Problem	Possible Cause	Solution
Wheelchair is difficult to fold/unfold.	Frame components are stiff or obstructed.	Ensure no obstructions. Apply a small amount of silicone lubricant to moving parts if necessary.
Footplates are loose or wobbly.	Locking mechanism not fully engaged or worn.	Ensure footplates are fully inserted and locked. Check for damage to the locking pins.

6. SPECIFICATIONS

Detailed technical specifications for the VERMEIREN JAZZ S20 Wheelchair:

Feature	Specification
Model	VrmJazz20
Weight Capacity	Up to 120 kg
Item Weight	16.6 kg
Product Dimensions (LxWxH)	112 x 93.5 x 67.5 cm
Rear Wheel Size	24 inches
Front Caster Size	8 inches
Frame Material	Sturdy Frame (Powder coated glossy look)
Seat/Backrest Material	Fire retardant nylon material
Tires	Polyurethane tubeless tires
Spokes	Stainless steel
Country of Origin	Belgium



Product Dimensions



Weight capacity – 110kgs

units in mm

Image: Technical diagram illustrating the precise dimensions of the wheelchair in millimeters, including height, width, and depth.

7. WARRANTY AND SUPPORT

VERMEIREN is committed to providing high-quality products and excellent customer service.

7.1. Warranty Information

A 1-year extended warranty is available for the VERMEIREN JAZZ S20 Wheelchair. Please refer to your purchase documentation for specific terms and conditions regarding warranty coverage.

7.2. Customer Support

For feedback, service assistance, or warranty inquiries, please contact us during business hours (10 AM to 5 PM, Monday to Friday):

- **Phone:** [7499852548](tel:7499852548)

- Email: bussiness@vermeirenindia.com



JAZZ S20

- FEEDBACK
- SERVICE ASSISTANCE
- WARRANTY



10am to 5pm Mon to Friday

contact



7499652548


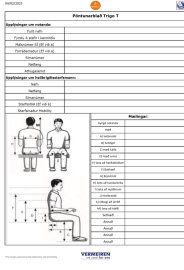

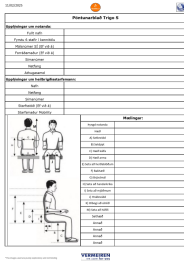


bussiness@vermeirenindia.com

Image: Contact details for customer support, including phone number and email address, for any product-related queries.

© 2024 VERMEIREN. All rights reserved.

Related Documents

 <p>Jazz S50 Standard wheelchair</p>	<p>Vermeiren Jazz S50 Standard Wheelchair: Features, Options, and Specifications</p> <p>Discover the Vermeiren Jazz S50 standard wheelchair, a blend of technicality, comfort, and quality. Learn about its features, customization options like transit models and anti-tipping wheels, and detailed technical specifications.</p>
	<p>Vermeiren Trigo T Wheelchair: Product Catalog and Specifications</p> <p>Detailed product catalog for the Vermeiren Trigo T, a trendy, active, and lightweight wheelchair. Includes specifications, dimensions, options for seat, backrest, armrests, footrests, wheels, brakes, and accessories, along with an order form section.</p>
 <p>User manual</p> <p>Jazz S50 Kids</p>	<p>Vermeiren Jazz S50 Kids User Manual: Operation, Safety, Maintenance, and Specifications</p> <p>This user manual provides comprehensive guidance for the Vermeiren Jazz S50 Kids wheelchair, covering its features, safe operation, maintenance procedures, troubleshooting, and detailed technical specifications.</p>
	<p>VERMEIREN Trigo S: Pöntunarleiðbeiningar fyrir virkan léttan hjólastól</p> <p>Uppgötvaðu VERMEIREN Trigo S, stílhreinan, virkan og léttan hjólastól hannaðan fyrir hámarks hreyfanleika og þægindi. Þessi ítarlegu leiðbeiningar útskýra allar tiltækar stillingar, mál og sérsniðnar valkosti fyrir Trigo S líkanið.</p>