Manuals+

Q & A | Deep Search | Upload

manuals.plus /

- arboleaf /
- arboleaf Wireless Heated Eye Mask ME10A User Manual

arboleaf ME10A

arboleaf Wireless Heated Eye Mask User Manual

Model: ME1OA

1. Introduction

This manual provides essential instructions for the safe and effective use of your arboleaf Wireless Heated Eye Mask. Designed for comfort and relief, this device offers adjustable warmth to soothe tired or dry eyes. Please read this manual thoroughly before first use and retain it for future reference.



Figure 1: arboleaf Wireless Heated Eye Mask and accessories.

2. SAFETY INFORMATION

To ensure your safety and comfort, adhere to the following guidelines:

- Do not use the heated eye mask for more than 20 minutes per session to prevent discomfort or potential skin irritation.
- During charging, observe the power button indicator light. It will flash red while charging and turn blue when fully charged.
- Always read the product label, warnings, and instructions provided on the product packaging before
 use.
- Keep out of reach of children.
- Do not use if the mask appears damaged or if the heating element is exposed.

3. PRODUCT COMPONENTS

The arboleaf Wireless Heated Eye Mask package includes:

• 1 x arboleaf Wireless Heated Eye Mask

- 1 x USB Charging Cable
- 1 x User Manual



Figure 2: Package contents.

4. SETUP

4.1. Charging the Device

Before initial use, fully charge the eye mask. Connect the provided USB charging cable to the mask's charging port and to a compatible USB power source (e.g., computer, wall adapter). The power button indicator will flash red during charging and turn blue once fully charged.

• Battery Capacity: 1000 mAh

- Charging Time: Approximately 1.5 hours
- Usage per Charge: Up to 4 sessions (based on 15-minute usage)

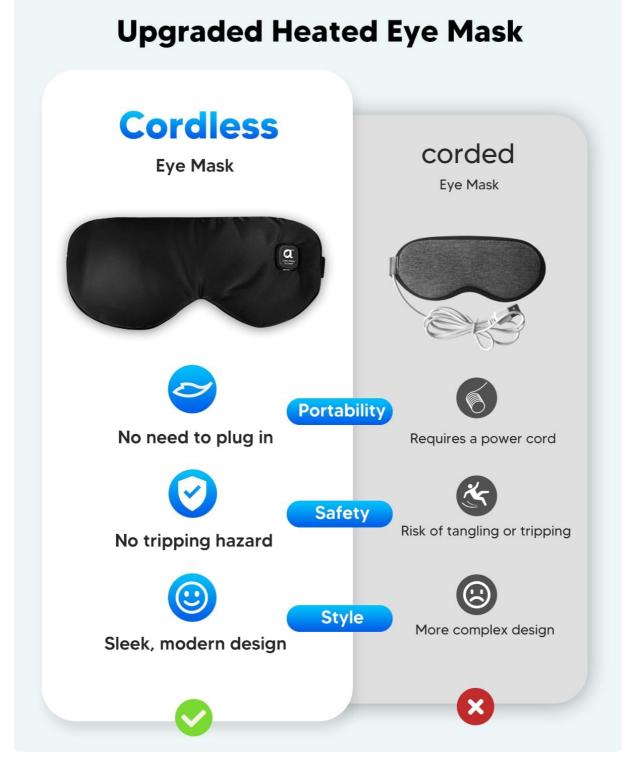


Figure 3: Charging the eye mask and indicator lights.

5. OPERATING INSTRUCTIONS

5.1. Power On/Off

To power on the device, press and hold the power button located on the mask for approximately 3 seconds. To power off, press and hold the power button again for 3 seconds.

5.2. Temperature Modes

The eye mask features three distinct temperature modes for personalized comfort. Once powered on, short press the power button to cycle through the modes:

- Mode 1: Approximately 100°F (38°C) Gentle warmth
- Mode 2: Approximately 108°F (42°C) Moderate warmth
- Mode 3: Approximately 115°F (46°C) Deep warmth



Figure 4: Temperature mode selection.

5.3. Automatic Shut-off

For safety and convenience, the eye mask is equipped with an automatic shut-off feature. It will automatically turn off after approximately 15 minutes of continuous operation.

6. MAINTENANCE

6.1. Cleaning

The eye mask features a removable and washable jacket. To clean, carefully remove the outer jacket from the heating element and hand wash it with mild soap and water. Allow it to air dry completely before reassembling or using. Do not wash the internal heating component.

6.2. Storage

Store the eye mask in a cool, dry place away from direct sunlight and extreme temperatures. Ensure the mask is fully charged before storing for extended periods to maintain battery health.

7. TROUBLESHOOTING

If you encounter any issues with your arboleaf Wireless Heated Eye Mask, consider the following:

- Device not turning on: Ensure the battery is sufficiently charged. Connect to a power source and check the charging indicator light.
- **No heat or insufficient heat:** Verify that the device is powered on and you have selected a temperature mode. Ensure the mask is properly positioned on your eyes.
- **Charging issues:** Check the USB cable and power source for proper connection. Try a different USB cable or power adapter.

If problems persist, please contact arboleaf customer support for assistance.

8. Specifications

Feature	Specification
Model Number	ME1OA
Dimensions	24.64 x 10.4 x 2.3 cm
Weight	80 grams
Battery Type	Lithium-ion (included)
Battery Capacity	1000 mAh
Material	Polyester, Artificial Silk
Color	Black
Special Features	Adjustable strap, Lightweight, Removable washable jacket
Age Range	Adult
Country of Origin	China

9. WARRANTY AND SUPPORT

For warranty information, product support, or any inquiries not covered in this manual, please refer to the official arboleaf website or contact arboleaf customer service directly. Contact details can typically be found on the product packaging or the brand's official online presence.

Related Documents - ME10A



Arboleaf Scalp Massager CM11B Quick Start Guide

Quick start guide for the Arboleaf Scalp Massager, Model CM11B. Learn about product overview, specifications, warranty, troubleshooting, and cautions.



Arboleaf J2 Percussion Massager User Manual

Comprehensive user manual for the Arboleaf J2 percussion massager. Includes detailed specifications, product overview, quick start guide, important safety precautions, and warranty information.



Arboleaf CM40G Massage Gun User Manual

Comprehensive user manual for the Arboleaf CM40G Massage Gun, detailing specifications, safety guidelines, operation instructions, and warranty information.



Arboleaf CM40G Massage Gun User Manual - Operation and Safety Guide

Comprehensive user manual for the Arboleaf CM40G massage gun, covering important safeguards, specifications, product overview, quick start guide, massage areas, and warranty information.



Arboleaf CM20C Massager User Manual

User manual for the Arboleaf CM20C massager, including important safeguards, specifications, product overview, quick start guide, and warranty information.

Arboleaf CM40G Massage Gun Quick Start Guide & Specifications
Comprehensive quick start guide and specifications for the Arboleaf CM40G Massage Gun. Learn about its features, operation, accessories, safety precautions, and warranty.
about its realures, operation, accessories, safety precautions, and warranty.