

[Manuals.plus](#) /

> [SSYA](#) /

> SSYA Dementia Clock User Manual - Model Z071

SSYA Z071

SSYA Dementia Clock User Manual

Model: Z071

Brand: SSYA

1. INTRODUCTION AND OVERVIEW

The SSYA Dementia Clock is designed to provide clear and easy-to-read time and date information, specifically beneficial for seniors, individuals with memory impairment, or those with vision challenges. Its large, non-abbreviated display helps users maintain a sense of time and day, reducing confusion and promoting independence.

Key features include:

- Ultra-large 11.5-inch (10-inch screen) 1280x800 HD display.
- Clearly spells out full DAY of the WEEK, MONTH, and DATE without abbreviations.
- Displays time of day (MORNING, AFTERNOON, EVENING, NIGHT, PREDAWN) with Sun & Moon icons.
- Supports 7 different display modes to suit various preferences.
- Features 12 customizable alarms for reminders like medication or appointments.
- Automatic Daylight Saving Time (DST) adjustment.
- Automatic dimming function for comfortable viewing at night.
- Supports 8 languages.



Figure 1: Front view of the SSYA Dementia Clock displaying the current time and date.

2. PACKAGE CONTENTS

Please check the package contents upon unboxing to ensure all items are present:

- SSYA Dementia Clock Unit (1)
- Power Adapter (1)
- User Manual (1)

3. SETUP

1. **Unpacking:** Carefully remove the clock unit and all accessories from the packaging.
2. **Power Connection:** Connect the power adapter to the DC IN port on the back of the clock. Plug the adapter into a standard electrical outlet. The clock is designed for plug-and-play operation and should power on automatically.
3. **Initial Display:** Upon first power-on, the clock will display the default time and date.
4. **Automatic Time Synchronization:** The clock is equipped with an automatic Daylight Saving Time (DST) adjustment feature. Ensure this setting is enabled if your region observes DST. The clock will automatically adjust for Summer/Winter time changes in March and November.



Figure 2: Overview of clock dimensions and available display modes.

4. OPERATING INSTRUCTIONS

4.1. Setting Time and Date

The clock typically sets its time automatically. If manual adjustment is needed, refer to the buttons on the back of the unit. Use the navigation buttons to select and adjust the hour, minute, year, month, and day. Confirm selections to save changes.

4.2. Display Modes

The clock offers 7 distinct display modes to cater to different visual preferences. These modes can be cycled through using the dedicated button on the back of the clock. Options include various color schemes and layouts, such as analog and digital combinations.



Figure 4: The clock's alarm interface, showing a "Medication time" reminder.

4.4. Language and Time Format

The clock supports 8 different languages. You can select your preferred language through the settings menu. Additionally, you can switch between 12-hour (AM/PM) and 24-hour time formats, and choose between DAY-MONTH-YEAR or MONTH-DAY-YEAR date formats.

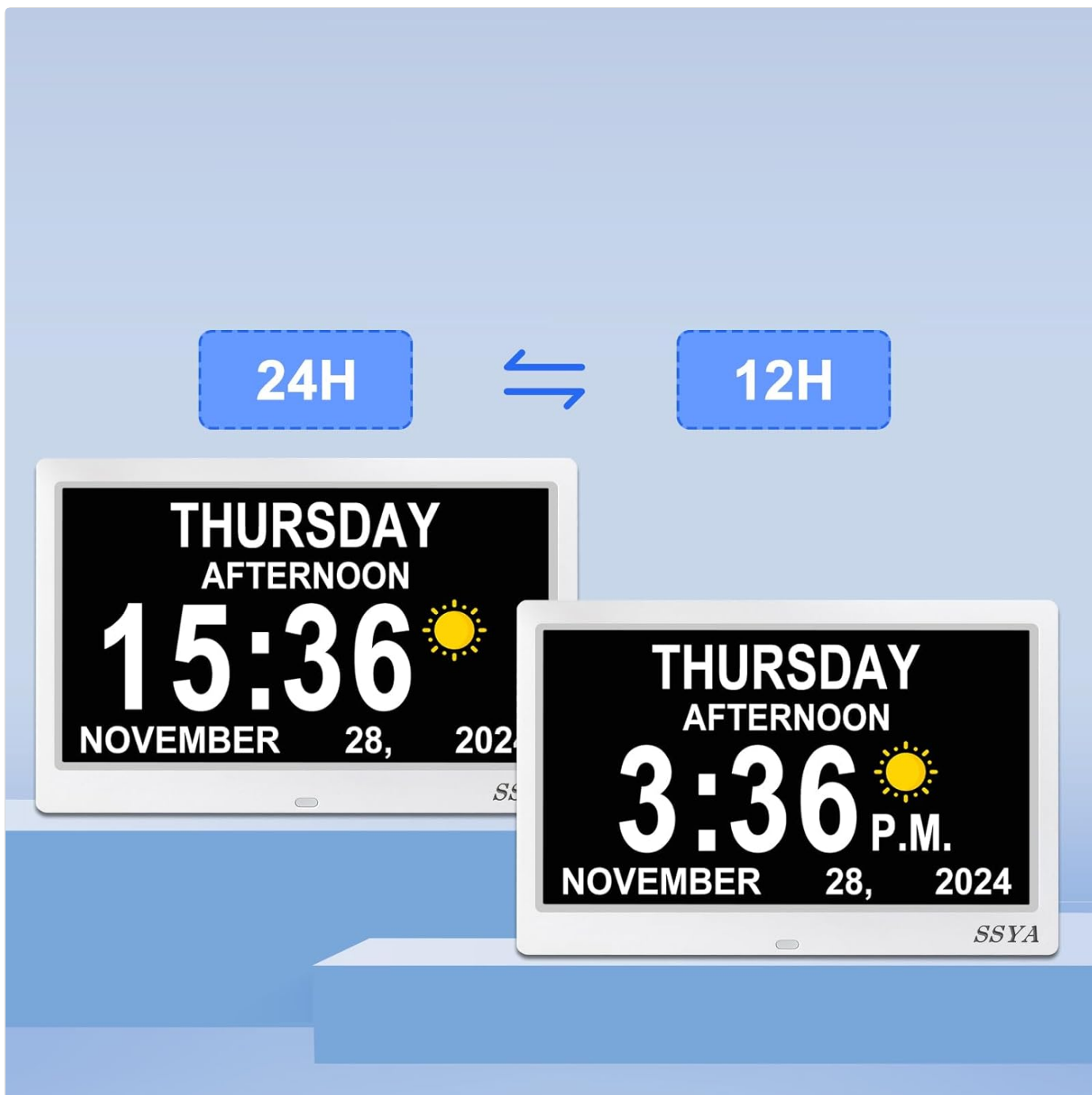


Figure 5: Illustration of 24-hour versus 12-hour time display options.

4.5. Automatic Dimming and Time Cycles

The clock features an automatic dimming function that adjusts screen brightness for optimal viewing during different times of the day, ensuring it's not too bright at night. It also categorizes the day into 5 distinct cycles:

- **MORNING:** 5:00 am - 11:59 am
- **AFTERNOON:** 12:00 pm - 5:59 pm
- **EVENING:** 6:00 pm - 8:59 pm
- **NIGHT:** 9:00 pm - 11:59 pm
- **PREDAWN:** 12:00 am - 4:59 am

These cycles are indicated by corresponding Sun & Moon icons on the display.

5. MAINTENANCE

To ensure the longevity and optimal performance of your SSYA Dementia Clock, please follow these maintenance guidelines:

- **Cleaning:** Use a soft, dry, lint-free cloth to gently wipe the screen and frame. Avoid using abrasive cleaners, solvents, or spray directly onto the screen.
- **Placement:** Place the clock on a stable, flat surface or mount it securely on a wall. Avoid areas with extreme temperatures, direct sunlight, or high humidity.
- **Power:** Always use the provided power adapter. Disconnect the power adapter from the outlet if the clock will not be used for an extended period.

6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Clock does not power on.	Power adapter not connected or outlet is not functioning.	Ensure the power adapter is securely plugged into the clock and a working electrical outlet. Try a different outlet.
Incorrect time or date displayed.	Manual settings are incorrect or DST setting is off.	Access the settings menu to manually adjust the time and date. Verify the Automatic DST setting is enabled if applicable. The clock retains settings after a power outage.
Display is too bright or too dim.	Automatic dimming feature is not set correctly or preference.	Check the automatic dimming settings in the menu. Adjust brightness levels if manual control is preferred.
Alarms are not sounding.	Alarm is not activated or volume is too low.	Ensure the alarm is set to ON and the volume is audible in the alarm settings.

7. SPECIFICATIONS

Feature	Detail
Brand	SSYA
Model Number	Z071
Color	11.5" White
Display Type	Digital
Style	Modern
Special Features	Daylight Saving, Corded, LED Display, Automatic Dimming, Multiple Display Modes, Customizable Alarms
Product Dimensions	9.72"W x 6.36"H
Power Source	Corded Electric
Age Range (Description)	Adult
Room Type	Bedroom, Living Room, Home Office
Shape	Rectangular
Indoor/Outdoor Usage	Indoor

Feature	Detail
Material	Plastic
Frame Material	Plastic
Batteries Included	Yes (for memory retention, not primary power)
Mounting Type	Tabletop
Alarm Clock	Yes
Warranty Type	Limited
Number of Items	1
Operation Mode	Electrical
Item Weight	1.34 pounds
Country of Origin	China

8. WARRANTY AND SUPPORT

The SSYA Dementia Clock comes with a 360-day warranty covering various situations. For warranty claims, technical support, or any inquiries regarding your product, please contact SSYA customer service through the retailer's platform or the contact information provided in your product packaging.

Please retain your purchase receipt as proof of purchase for warranty purposes.