

BEBIRD I35R-bl


Bebird EarSightPlus HD 1920P Flexible Otoscope User Manual

Model: I35R-bl

1. PRODUCT OVERVIEW

The Bebird EarSightPlus is a flexible otoscope equipped with a 1920P HD camera, designed for visual inspection and cleaning of the ear canal. Its flexible tip adapts to the natural curves of the ear, providing a clear view for effective wax removal and examination. This device is also suitable for inspecting the nose, throat, and teeth.

Flexi-Otoskop Passt Besser

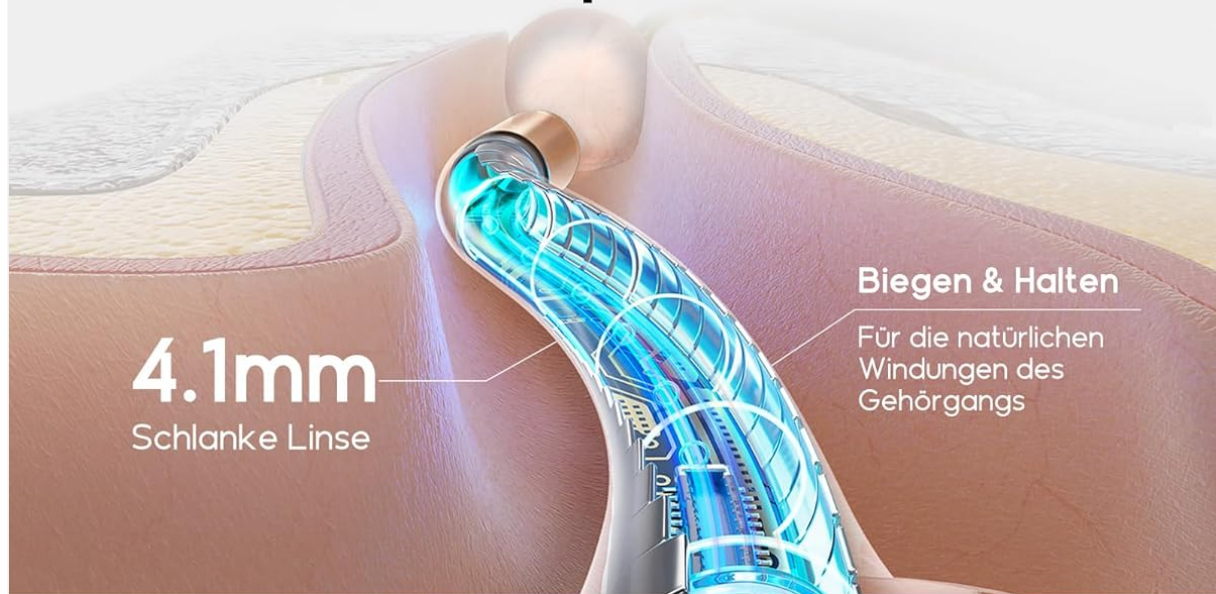

Bietet kristallklare Sicht



Image 1.1: Bebird EarSightPlus Otoscope and included accessories.

2. PACKAGE CONTENTS

- Bebird EarSightPlus Flexible Otoscope
- Flexible Tweezers
- Ear Spoons (various types)
- Scraper Tips
- Sponge Tips
- Protective Cover
- Charging Cable
- Storage Case

- User Manual

3. SAFETY INFORMATION

To ensure safe operation and prevent injury, please adhere to the following guidelines:

- Always use the soft silicone tips to prevent scratches or discomfort to the ear canal.
- The device features a built-in safety indicator and protective cover to prevent excessive insertion, safeguarding the eardrum during cleaning. Do not force the device beyond comfortable limits.
- Ensure all attachments are securely fastened before use.
- Clean and disinfect all tips before and after each use.
- Keep out of reach of children. Adult supervision is recommended when used by children.
- Do not use if you experience pain, discomfort, or bleeding. Consult a medical professional.



Image 3.1: The device is designed for precise and gentle use, suitable for sensitive ears.

4. GETTING STARTED

4.1 Charging the Device

Before first use, fully charge the otoscope using the provided charging cable. Connect the cable to the device and a suitable USB power source. The indicator light will show charging status.

4.2 App Installation and Connection

1. Download the official Bebird application from your smartphone's app store (iOS or Android).
2. Turn on the Bebird otoscope.
3. Open your smartphone's Wi-Fi settings and connect to the Wi-Fi network broadcast by the Bebird device (usually named 'Bebird-XXXX').
4. Open the Bebird application. The device's camera feed should appear on your screen, providing a real-time view.


Ohrhygiene bis zu 100 %

Erhalten Sie immer die richtige Ohrenreinigung!



Pinzette

Entfernt hartnäckigen
Ohrenschmalz

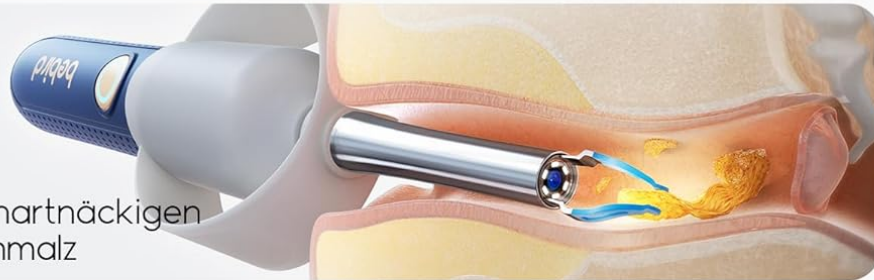

Löffel

Entfernt normalen
Ohrenschmalz



Schaber

Entfernt Ablagerungen



Schwamm

Entfernt feuchten
Ohrenschmalz



Image 4.1: Connecting the otoscope to the Bebird app for real-time viewing.

5. OPERATING INSTRUCTIONS

5.1 Ear Examination and Wax Removal

The flexible otoscope and HD camera allow for a clear and magnified view of the ear canal. Select the appropriate cleaning tip based on the type of earwax:

- **Flexible Tweezers:** For gently removing stubborn earwax without touching the skin.
- **Ear Spoons:** Suitable for general earwax removal.
- **Scraper Tips:** Designed to remove accumulated deposits.
- **Sponge Tips:** Ideal for moist earwax.

Carefully insert the otoscope into the ear canal while observing the real-time feed on your smartphone. Use gentle movements to remove earwax. The flexible design helps navigate the natural curves of the ear canal.



Image 5.1: Various cleaning tips for different types of earwax.

5.2 Multi-functional Use

Beyond ear inspection, the Bebird EarSightPlus can be used for:

- **Nose Examination:** Inspect the nasal cavity.
- **Throat Examination:** View the throat area.
- **Oral Examination:** Check teeth and gums.

The device supports capturing photos and videos, which can be useful for sharing health information with medical professionals.



Image 5.2: The otoscope's versatility for various home examinations.

6. MAINTENANCE AND STORAGE

6.1 Cleaning

After each use, remove the cleaning tip and clean it thoroughly with alcohol or an appropriate disinfectant. Ensure the camera lens is clean and free of debris for optimal viewing.

6.2 Storage

Store the otoscope and all accessories in the provided compact storage case. This keeps the components organized, clean, and protected from damage, making it suitable for both home use and travel.



Image 6.1: Hygienic storage of the otoscope and its attachments.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
No image on screen	Device not powered on; Wi-Fi not connected; Low battery.	Ensure device is on; Reconnect Wi-Fi; Charge the device.
Image is blurry	Camera lens is dirty; Incorrect focus.	Clean the camera lens; Adjust the device's position for better focus.
Device disconnects frequently	Weak Wi-Fi signal; Interference.	Ensure device is close to smartphone; Avoid other Wi-Fi interference.
Tips are loose	Tip not securely attached.	Ensure the tip is firmly screwed or pushed onto the device.

8. SPECIFICATIONS

- **Brand:** BEBIRD
- **Model Number:** I35R-bl
- **Camera Resolution:** 1920P HD
- **Weight:** Approximately 80 grams
- **Package Dimensions:** 19.9 x 13.4 x 3.2 cm
- **Connectivity:** Wi-Fi
- **Compatibility:** iOS and Android devices

9. WARRANTY AND SUPPORT

For warranty information, please refer to the documentation included with your purchase or contact BEBIRD customer support directly. If you encounter any issues or have questions regarding the use of your Bebird EarSightPlus otoscope, please reach out to the manufacturer for assistance.

Manufacturer: BEBIRD