

## FLEXISPOT XL5

# FLEXISPOT XL5 Dual OKIN Motor Power Lift Recliner Chair User Manual

Model: XL5

## 1. IMPORTANT SAFETY INSTRUCTIONS

---

Please read all instructions carefully before using the FLEXISPOT XL5 Power Lift Recliner Chair. Retain this manual for future reference.

- Ensure all components are securely assembled before connecting to power.
- Keep children and pets away from the chair's moving parts.
- Do not allow children to operate the chair without adult supervision.
- Avoid placing hands or feet near the lifting mechanism during operation.
- Do not expose the chair to excessive moisture or liquids.
- Unplug the chair from the power outlet before cleaning or maintenance.
- If the power cord is damaged, it must be replaced by the manufacturer or a qualified service agent to avoid hazards.
- Do not sit on the footrest when it is extended.
- This chair is designed for indoor use only.

## 2. PACKAGE CONTENTS

---

The FLEXISPOT XL5 recliner is shipped in one box. Carefully unpack all components and verify against the list below. If any parts are missing or damaged, contact customer support immediately.

- Recliner Base
- Backrest
- Side Wings (2)
- Power Adapter
- Remote Control (2) - one for lift/recline, one for massage/heat
- User Manual

# ALL IN ONE BOX

Receive & Enjoy it instantly. No more Waiting!



Image: The FLEXISPOT XL5 recliner shown next to a shipping box, illustrating the 'All in One Box' delivery.

## 3. ASSEMBLY INSTRUCTIONS

---

Assembly of the FLEXISPOT XL5 recliner is designed to be straightforward. Please follow these steps:

1. **Unpack Components:** Carefully remove all parts from the shipping box. Place the recliner base on a flat, stable surface.
2. **Attach Backrest:** Slide the backrest onto the metal posts on the recliner base until it clicks securely into place. Ensure both sides are fully engaged.
3. **Connect Side Wings:** Attach the side wings to the designated slots on the backrest.
4. **Connect Power:** Locate the power cable from the chair and connect it to the power adapter. Plug the power adapter into a standard electrical outlet.
5. **Allow for Rebound:** After assembly, allow the cushions approximately 48 hours to fully rebound to their intended plushness.

## 4. OPERATING INSTRUCTIONS

---

The FLEXISPOT XL5 features dual OKIN motors for independent control of the backrest and footrest, along with massage and heating functions.

## 4.1 Lift and Recline Functions

Use the dedicated lift and recline remote control to adjust the chair's position.



Image: Diagram illustrating the buttons on the lift and recline remote control and their corresponding chair movements.

- **Recline Backrest:** Press the 'Recline the Backrest' button (down arrow for backrest) to lean the backrest backward.
- **Raise Backrest:** Press the 'Raise the Backrest' button (up arrow for backrest) to bring the backrest upright.
- **Raise Footrest:** Press the 'Raise the Footrest' button (up arrow for footrest) to elevate the footrest.
- **Lower Footrest & Lift to Standing:** Press the 'Lower the Footrest, then Lift to Standing Position' button (down arrow for footrest) to retract the footrest and then gently lift the chair to assist in standing.
- **Lay Flat Position:** The dual motors allow for a nearly 180° lay-flat recline by independently adjusting the backrest and footrest to your desired comfort.

# POWER LIFT RECLINER CHAIR

Experience the comfort of Laying Flat and  
the ease of stand assistance



Image: The recliner demonstrating its lay-flat capability and the lift-assist function for easy standing.

## 4.2 Massage and Heating Functions

Utilize the separate massage and heat remote control to activate these features.



# MASSAGE & HEATING

A Spa-Like Experience



Lumbar  
Heating



8 Points  
Massage



3 Massage  
Modes



Timing  
Control

Image: The recliner highlighting the lumbar heating area and the 8 massage points across the back, seat, and leg rest.

- **Massage Points:** The chair features 8 vibrating massage points located in the back, lumbar, seat, and leg areas.
- **Massage Modes:** Select from 3 adjustable massage modes to suit your preference.
- **Lumbar Heating:** Activate the lumbar heating function for soothing warmth.
- **Intensity and Timer:** Adjust massage intensity and set a timer for your massage sessions using the remote.

## 4.3 USB Charging Ports and Storage

The recliner is equipped with convenient features for your comfort.

# PRODUCT DETAILS



**Built-in Cup Holders**



**Side Pockets**



**USB Charging Ports**

Image: Detailed view of the built-in cup holders, side storage pockets, and the USB charging port located on the lift/recline remote.

- **USB Charging Port:** A USB port is integrated into the lift/recline remote control, allowing you to charge your electronic devices while relaxing.
- **Cup Holders:** Two built-in cup holders are provided for beverages.
- **Side Pockets:** Dual side pockets offer convenient storage for remotes, magazines, or other small items.

## 5. CARE AND MAINTENANCE

Proper care will extend the life and appearance of your recliner.

- **Cleaning Fabric:** For general cleaning, use a soft cloth and a mild, non-abrasive fabric cleaner. Test on an inconspicuous area first.
- **Spills:** Blot spills immediately with a clean, dry cloth. Do not rub.
- **Vacuuming:** Regularly vacuum the fabric to remove dust and debris.
- **Moving Parts:** Periodically check all bolts and connections to ensure they are tight. Do not lubricate moving parts unless specifically instructed by the manufacturer.
- **Power Cord:** Inspect the power cord for any damage. If damaged, discontinue use and contact customer support.

## 6. TROUBLESHOOTING GUIDE

---

If you encounter issues with your recliner, refer to the following common problems and solutions.

Problem	Possible Cause	Solution
Chair does not operate (no power)	Power cord unplugged, outlet not working, power adapter faulty.	Ensure power cord is securely plugged into the chair and a working wall outlet. Test the outlet with another device. Check for damage to the power adapter.
Motor makes unusual noise	Loose connection, obstruction in mechanism.	Check for any obstructions under the chair. Ensure all connections are secure. If noise persists, contact customer support.
Massage/Heat not working	Massage/heat remote not connected, power issue.	Ensure the massage/heat remote is properly connected to the chair's power system. Verify the main power supply to the chair.
Chair does not lift or recline smoothly	Obstruction, weight exceeding capacity, motor issue.	Check for obstructions. Ensure user weight does not exceed 350 lbs. If issues persist, contact customer support.

If the problem persists after attempting these solutions, please contact FLEXISPOT customer support.

## 7. PRODUCT SPECIFICATIONS

---

Detailed specifications for the FLEXISPOT XL5 Dual OKIN Motor Power Lift Recliner Chair.



# PRODUCT DIMENSIONS



Image: Technical drawing displaying the width, depth, height, seat depth, and seat width of the recliner.

- **Brand:** FLEXISPOT
- **Model:** XL5
- **Color:** Blue
- **Material:** Textile
- **Product Dimensions:** 29"D x 36.8"W x 41.5"H
- **Seat Depth:** 21 inches
- **Special Feature:** Dual OKIN Motor, Lay Flat, Massage, Heat, USB Ports, Cup Holders
- **Weight Capacity:** Up to 350 lbs (as indicated in product images)
- **Power:** Electric (Power Adapter included)

## 8. WARRANTY AND CUSTOMER SUPPORT

FLEXISPOT is committed to providing quality products and excellent customer service.



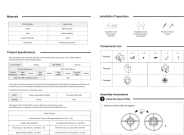

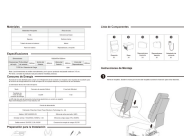

- **Warranty:** For detailed warranty information regarding your FLEXISPOT XL5 recliner, please refer to the warranty card included with your purchase or contact FLEXISPOT customer support directly.
- **Customer Support:** If you have any questions, require assistance with assembly, operation, or



troubleshooting, FLEXISPOT's customer service team is available to provide prompt support.

**Contact Information:** Please visit the official FLEXISPOT website or refer to your purchase documentation for the most current customer support contact details.

## Related Documents - XL5

	<p><a href="#">FlexiSpot DT3B-EU User Manual and Warranty Information</a></p> <p>Comprehensive user manual for the FlexiSpot DT3B-EU height adjustable desk frame, including detailed assembly instructions, parts list, and warranty coverage information.</p>
	<p><a href="#">FlexiSpot Waffle Recliner User Manual</a></p> <p>Comprehensive user manual for the FlexiSpot Waffle Recliner, detailing assembly, operation, safety warnings, and troubleshooting. Learn how to set up and use your new recliner.</p>
	<p><a href="#">FlexiSpot Power Swivel Rocker Recliner Assembly and User Guide</a></p> <p>Detailed assembly instructions, product specifications, usage guidelines, and troubleshooting for the FlexiSpot Power Swivel Rocker Recliner. Learn how to set up and operate your new recliner chair safely and efficiently.</p>
	<p><a href="#">FlexiSpot XP6 Pushback Recliner Chair User Manual</a></p> <p>User manual for the FlexiSpot XP6 Pushback Recliner Chair, providing assembly instructions, usage guidelines, maintenance tips, warnings, and troubleshooting information.</p>
	<p><a href="#">FlexiSpot XR5 Electric Recliner Chair: Assembly, Specifications, and Operation Guide</a></p> <p>Comprehensive guide for the FlexiSpot XR5 electric recliner chair. Learn about its materials, specifications, energy efficiency, step-by-step assembly, switch functions, and troubleshooting. Includes power supply details and safety warnings.</p>
	<p><a href="#">FlexiSpot XP6 Pushback Recliner Chair User Manual</a></p> <p>Comprehensive user manual for the FlexiSpot XP6 Pushback Recliner Chair, covering assembly instructions, product specifications, usage guidelines, maintenance tips, warnings, and troubleshooting.</p>