



[Manuals.plus](#) /

> [Ceretone](#) /

> Ceretone Fusion OTC Rechargeable Hearing Aids User Manual (Model A61)

Ceretone A61

Ceretone Fusion OTC Rechargeable Hearing Aids User Manual

Model: A61

INTRODUCTION

This manual provides essential instructions for the safe and effective use of your Ceretone Fusion OTC Rechargeable Hearing Aids. Please read this manual thoroughly before using the device and retain it for future reference. These hearing aids are designed to assist individuals with mild to moderate hearing loss by amplifying sounds and reducing ambient noise.

WHAT'S IN THE BOX

Upon opening your Ceretone Fusion package, you should find the following items:

- 1x Pair of Fusion Hearing Aids
- 1x Portable Charging Case
- 1x Charging Cable
- 3x Pairs of Eartips (Small, Medium, Large)
- 1x Cleaning Tool



Image: The complete packaging list, showing the hearing aids, charging case, cable, eartips, and cleaning tool.

SETUP

1. Initial Charging

Before first use, fully charge the portable charging case and the hearing aids. Connect the charging cable to the case and a power source. The case indicators will show charging status. A full charge typically takes approximately 2 hours.

First time using a hearing aid?

No need to worry,
we provide a comprehensive user manual,
video tutorials and after-sales support,
Providing you with comprehensive hearing support.



Instruction Manual



Quick Start Guide

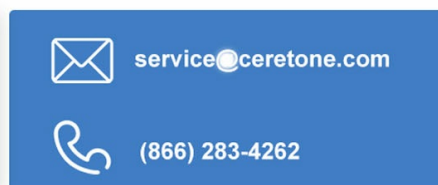


Image: The hearing aids and charging case connected to a power adapter for charging, illustrating the 2-hour rapid charge and extended battery life.

2. Selecting Eartips

Choose the eartips that provide the most comfortable and secure fit in your ear canal. The package includes small, medium, and large sizes. A proper fit is crucial for optimal sound quality and to prevent feedback.



Image: A visual guide to selecting the correct eartip size for all-day comfort and proper fit.

3. Inserting the Hearing Aids

Gently insert the hearing aid into your ear canal. Ensure it sits snugly and comfortably. The ergonomic design is intended to fit the natural curve of the ear.

OPERATING INSTRUCTIONS

1. Power On/Off

The hearing aids automatically power on when removed from the charging case and power off when placed back into the case.

2. Volume and Mode Adjustment

Each hearing aid features an ergonomic one-click button for adjustments:

- **Short Press:** Adjusts the volume through 6 available levels.
- **Long Press:** Switches between 5 different hearing modes.

The device includes a memory function, retaining your last used volume and hearing mode settings for convenience.

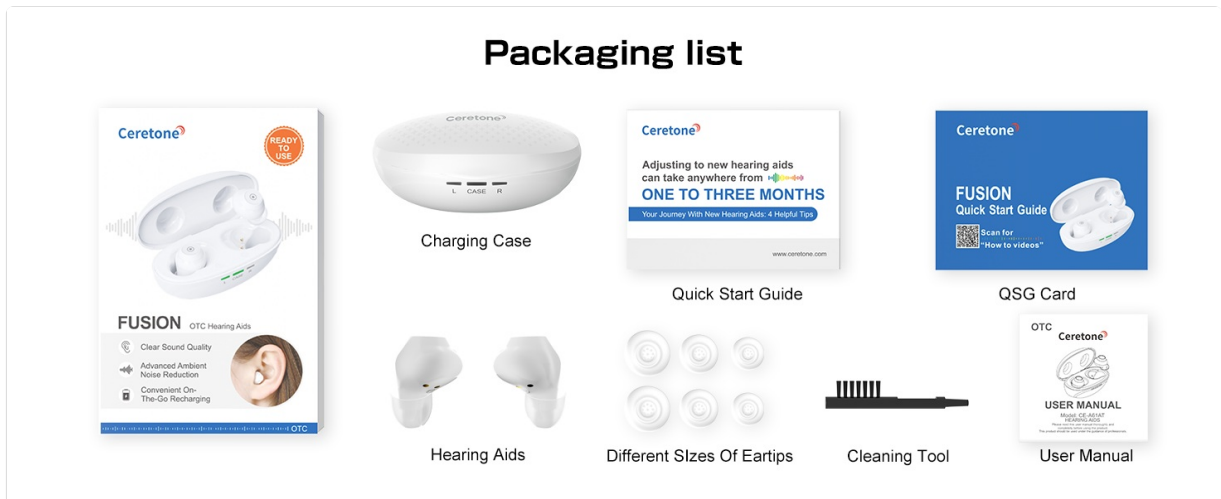


Image: Detailed view of the hearing aid, highlighting the microphone, charging contacts, eartip, and the multi-function button for volume and mode switching.

3. Hearing Modes

The Ceretone Fusion offers versatile hearing modes to optimize performance in various environments:

- **Standard Mode:** Ideal for everyday listening.
- **Restaurant Mode:** Enhances speech clarity in noisy dining environments.
- **Outdoor Mode:** Optimizes hearing in outdoor settings with varying background noise.
- **Music Mode:** Delivers enhanced sound quality for enjoying music.

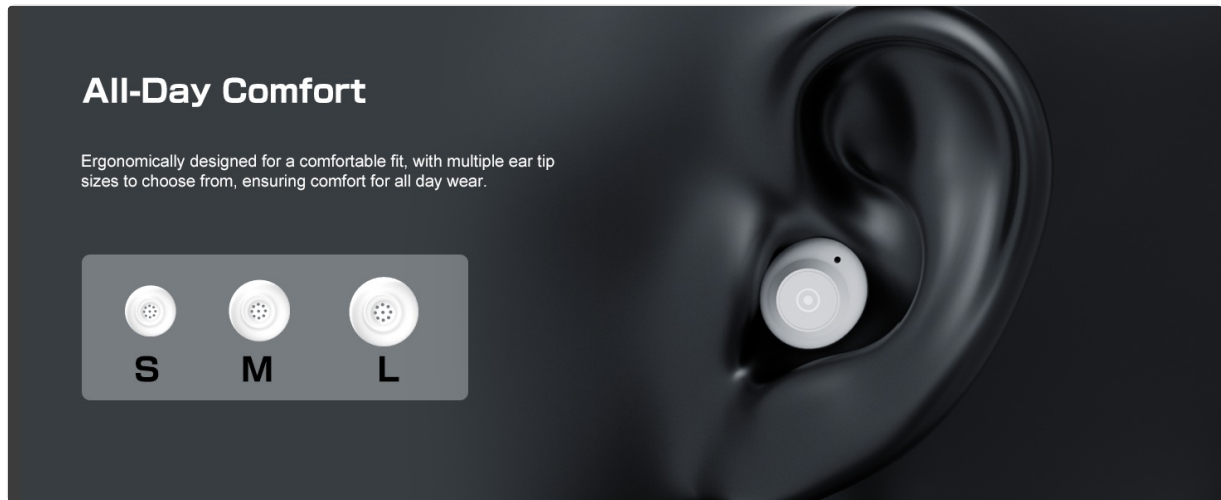


Image: Illustrations depicting the Standard, Restaurant, Outdoor, and Music modes, showing how the hearing aids adapt to different soundscapes.

MAINTENANCE

1. Cleaning

Regular cleaning is essential for maintaining optimal performance and hygiene. Use the provided cleaning tool to gently remove any earwax or debris from the eartips and sound outlets. Wipe the body of the hearing aids with a soft, dry cloth. Do not use liquids or harsh chemicals.

2. Storage

When not in use, store the hearing aids in their charging case. This protects them from dust and damage, and keeps them charged for your next use.

3. Water and Dust Resistance

The hearing aids are rated IP54 for dust and sweat resistance. While they can withstand splashes and light exposure to dust, they are not waterproof. Avoid submerging them in water or exposing them to heavy rain.

TROUBLESHOOTING

If you encounter any issues with your Ceretone Fusion Hearing Aids, please refer to the following common problems and solutions:

- **No Sound:**
 - Ensure the hearing aids are fully charged.
 - Check if the volume is set to an audible level.
 - Verify that the eartips are not blocked by earwax or debris. Clean if necessary.

- **Whistling/Feedback:**

- Ensure the eartips are properly seated in your ear canal and form a tight seal. Try a different size eartip if needed.
- Reduce the volume level.

- **Sound is Distorted or Unclear:**

- Clean the eartips and sound outlets.
- Ensure the hearing aids are fully charged.
- Try switching to a different hearing mode.

- **Hearing Aids Not Charging:**

- Ensure the charging cable is securely connected to both the case and the power source.
- Check the charging contacts on both the hearing aids and the case for any debris.

- **Low Battery Warning:**

- The hearing aids will provide a low battery warning. Place them in the charging case promptly.

If these steps do not resolve the issue, please contact customer support.

SPECIFICATIONS

Item Model Number	A61
Package Dimensions	6.73 x 4.57 x 2.36 inches
Item Weight	10.23 ounces
Battery Life (Single Charge)	Up to 25 hours
Total Battery Life (with Charging Case)	Up to 100 hours
Charging Time	Approximately 2 hours
Dust and Water Resistance	IP54

WARRANTY AND SUPPORT

Ceretone offers comprehensive customer support for your hearing aids. We are committed to ensuring your satisfaction and providing assistance whenever needed.

Lifetime After-Sale Service

We provide lifetime after-sale service. Our U.S.-based customer support team is available to assist you with any questions or concerns.

Contact Information

- **Phone:** (866) 283-4262
- **Email:** service@ceretone.com

Easier Operation Than Ever, One Press Solves Everything



Image: Information on Ceretone's customer service, including phone number and email for support.

© 2025 Ceretone. All rights reserved.