

[Manuals.plus](#) /

› [Snailax](#) /

› Snailax SL-54C6 Shiatsu Foot Massager User Manual

## Snailax SL-54C6

# Snailax SL-54C6 Shiatsu Foot Massager User Manual

Model: SL-54C6

## PRODUCT OVERVIEW

The Snailax SL-54C6 Shiatsu Foot Massager is designed to provide a comprehensive massage experience for your feet, calves, and even arms. It combines deep kneading shiatsu massage, air compression, heat therapy, and vibration to help relax muscles and promote circulation. This unit features an intuitive control panel and a convenient remote control for ease of use.



Figure 1: Snailax SL-54C6 Shiatsu Foot Massager. The main image displays the massager unit with its remote control. Inset images demonstrate its versatility for foot, calf, and arm massage.

## IMPORTANT SAFETY INSTRUCTIONS

When using electrical appliances, basic precautions should always be followed, including the following:

- Always unplug the appliance from the electrical outlet immediately after using and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for an appliance that has fallen into water. Unplug it immediately.
- An appliance should never be left unattended when plugged in.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not

recommended by the manufacturer.

- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on or near eyes or other highly sensitive areas.
- To disconnect, turn all controls to the "OFF" position, then remove plug from outlet.
- This appliance is designed for household use only.
- Do not attempt to repair the appliance. All servicing of this appliance must be performed by authorized Snailax service personnel only.

*Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is not intended for medical use. Consult your physician before using this product if you have any medical conditions.*

## SETUP AND ASSEMBLY

---

The Snailax SL-54C6 Foot Massager comes pre-assembled and ready for use. Follow these steps to set up your device:

1. **Unpack:** Carefully remove the massager from its packaging. Retain packaging for future storage or transport.
2. **Placement:** Place the massager on a flat, stable surface. Ensure there is enough space around the unit for comfortable use and proper ventilation.
3. **Power Connection:** Plug the power adapter into the massager's power input port, then plug the other end into a standard electrical outlet (110-240V).
4. **Adjustable Stand:** The massager features a 360° adjustable handle/stand. Adjust it to your desired angle for foot, calf, or arm massage.



## 360° Adjustable Stand



Figure 2: The 360° adjustable stand allows for versatile positioning of the massager for different body parts, such as feet, calves, or arms.

## OPERATING INSTRUCTIONS

The Snailax SL-54C6 can be operated using either the integrated LCD smart touch panel or the included remote control.

### Control Panel and Remote Control Overview

## Panel Description



Power



Vibration



Knead



Compression



Heat



Figure 3: Control Panel Description. Buttons include Power, Vibration, Knead, Compression, and Heat, each with multiple intensity levels.

## Remote Control Description



Figure 4: Remote Control Description. The remote provides convenient access to all massage functions, including Mute, Knead, Compression, Heat, Vibration, Pause, and Foot & Calf mode selection.

### Massage Modes and Settings

The massager offers a variety of customizable settings:

- **Power On/Off:** Press the power button on the control panel or remote to turn the unit on or off.
- **Shiatsu Kneading Massage:** Features 8 compression units for deep kneading. Select from 3 intensity levels (Low, Medium, High).
- **Air Compression:** Provides soothing compression massage. Choose from 3 intensity levels (Low, Medium, High). **Note:** For first-time users, it is recommended to start with compression level 1.
- **Heat Therapy:** Offers gentle warmth to enhance relaxation. Select from 3 heat settings.

- **Vibration Massage:** Two vibration massage modes are available.
- **Foot/Calf Mode:** An optional mode to optimize massage for either feet or calves.
- **Automatic Shut-off:** The unit is equipped with a 15-minute automatic shut-off timer for safety and convenience.



Figure 5: Customizable Comfort. This image highlights the various massage settings available, including Shiatsu, Compression, Vibration, and two optional Auto Modes, allowing users to tailor their experience.

# Multifunctional Foot & Calf Massager

Relax Your Tiredness of The Day



Vibrating



Heating



Kneading



Compressing



Figure 6: Multifunctional Massage. This visual emphasizes the four primary massage functions: Vibrating, Heating, Kneading, and Compressing, all integrated into the device.

# Professional Air Compression



Figure 7: Professional Air Compression Benefits. This image illustrates the positive effects of the massager's air compression feature, including improved blood circulation, enhanced relaxation, muscle relief, and promotion of better sleep.

## CARE AND MAINTENANCE

Proper care and maintenance will ensure the longevity of your Snailax SL-54C6 Foot Massager.

- **Cleaning:** The removable and washable cover ensures easy cleaning. Unzip and remove the fabric cover, then hand wash with mild detergent and air dry. Do not machine wash or tumble dry.

- **Unit Cleaning:** Unplug the unit before cleaning. Wipe the surface of the massager with a soft, damp cloth. Do not use abrasive cleaners, brushes, gasoline, kerosene, glass furniture polish, or paint thinner to clean the massager.
- **Storage:** Store the massager in a cool, dry place when not in use. Avoid direct sunlight and extreme temperatures. Do not store with the power cord wrapped tightly around the unit.
- **Avoid Damage:** Do not stand on the unit. Do not use sharp objects near the massager.

## TROUBLESHOOTING GUIDE

If you encounter any issues with your Snailax SL-54C6 Foot Massager, please refer to the following common problems and solutions:

| Problem                                  | Possible Cause                                                  | Solution                                                                                                                                                                                |
|------------------------------------------|-----------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Massager does not turn on.               | Not plugged in correctly; Power outlet issue; Unit malfunction. | Ensure the power cord is securely plugged into both the massager and a working electrical outlet. Try a different outlet. If the issue persists, contact customer support.              |
| Massage intensity feels too strong/weak. | Incorrect intensity setting; User preference.                   | Adjust the massage intensity levels (Shiatsu, Compression, Vibration) using the control panel or remote control. Remember to start with lower compression levels if you are a beginner. |
| Heat function is not working.            | Heat function not activated; Unit malfunction.                  | Ensure the heat function is turned on via the control panel or remote. Allow a few minutes for the heat to become noticeable. If no heat is felt, contact customer support.             |
| Remote control is not responding.        | Dead batteries; Obstruction between remote and unit.            | Replace the batteries in the remote control. Ensure there is no obstruction between the remote and the massager's sensor.                                                               |

If the problem persists after attempting these solutions, please contact Snailax Customer Support for further assistance.


## PRODUCT SPECIFICATIONS

| Feature      | Detail  |
|--------------|---------|
| Model Number | SL-54C6 |
| Brand        | Snailax |

| Feature            | Detail                                                                                                                                                 |
|--------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------|
| Product Dimensions | 19.29"L x 13.58"W x 9.76"H                                                                                                                             |
| Item Weight        | 15.25 Pounds                                                                                                                                           |
| Material           | Acrylonitrile Butadiene Styrene (ABS)                                                                                                                  |
| Color              | Grey                                                                                                                                                   |
| Special Features   | Remote Control, 3-in-1 Massage (Shiatsu, Compression, Vibration), Heat Function, 360° Adjustable Stand, Removable/Washable Cover, 15-min Auto Shut-off |



## Space Saving

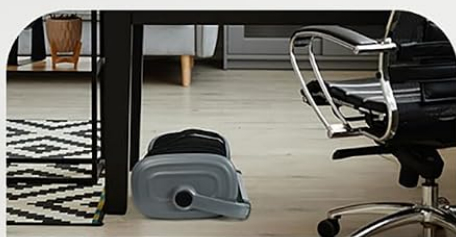

Office



Home

Figure 8: Product Dimensions and Space-Saving Design. This image provides a visual representation of the massager's dimensions and suggests its suitability for various environments like office and home due to its compact design.

## WARRANTY AND SUPPORT

Snailax products are designed to offer lasting comfort and peace of mind. For warranty information and customer support, please refer to the contact details provided with your product packaging or visit the official Snailax website.

For any questions or concerns regarding your Snailax SL-54C6 Foot Massager, please contact Snailax

Customer Service. Please have your model number (SL-54C6) and purchase information ready when contacting support.

*Note: Warranty terms may vary by region and retailer. Please consult your purchase documentation for specific details.*