

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Momcozy](#) /

› [Momcozy M5 Hands-Free Wearable Breast Pump Instruction Manual](#)

## Momcozy M5

# Momcozy M5 Hands-Free Wearable Breast Pump Instruction Manual

Brand: Momcozy | Model: M5

## 1. INTRODUCTION

The Momcozy M5 Hands-Free Wearable Breast Pump is designed for convenience and comfort, allowing mothers to express milk discreetly and efficiently. Its compact and lightweight design, combined with advanced pumping technology, aims to provide a comfortable and effective pumping experience.

Key features include:

- **Lightweight and Discreet:** Weighing only 230g, it is designed to be worn inside a bra without being overly noticeable.
- **Comfortable Lactation:** Features a 'Baby Mouth' concept with a double-sealed flange and ergonomic breast arc design for a painless and comfortable experience.
- **Efficient Milk Expression:** Incorporates micro-vibration technology to promote milk flow. Offers 3 modes (Stimulation, Expression, Mixed) and 9 suction levels for personalized pumping.
- **Tailored Fit:** Comes with various flange sizes (17/19/21/24mm) to ensure a proper fit for different breast shapes. Made from BPA-free, food-grade silicone.
- **Easy to Clean:** Designed with minimal parts for straightforward assembly and cleaning.

## 2. IMPORTANT SAFETY INFORMATION

**WARNING:** Do NOT use this device if you have a cardiac pacemaker, implanted defibrillator, or any other metallic or electronic implant, as it may cause electromagnetic interference or serious injury.

Always consult with a healthcare professional regarding any concerns about breastfeeding or using a breast pump.

## 3. PRODUCT COMPONENTS

The Momcozy M5 Wearable Breast Pump package includes the following components:

# What's Your Get - 1 Pack

## Parts needed to pump



1 Motor



1 Flange (24mm)



1 Milk Collector



1 Diaphragm



1 Valve

## Also included



1 Valve



1 Insert (17mm)



1 Insert (19mm)



1 Insert (21mm)



1 Flange Cover



1 Type-C Cable



1 Size Ruler



1 Pumps Case

Tips: Suggest replacing every 1-3 months for performance

Image: Components included in a single pack of the Momcozy M5 Wearable Breast Pump.

- 1 x Motor Unit
- 1 x Flange (24mm)
- 1 x Milk Collector Cup
- 1 x Diaphragm
- 1 x Valve
- 1 x Insert (17mm)

- 1 x Insert (19mm)
- 1 x Insert (21mm)
- 1 x Flange Cover
- 1 x Type-C Charging Cable
- 1 x Nipple Size Ruler
- 1 x Pump Case

*Note: It is suggested to replace parts like diaphragms and valves every 1-3 months for optimal performance.*

## 4. ASSEMBLY AND SETUP

---


Proper assembly is crucial for the efficient operation of your Momcozy M5 breast pump. Follow these general steps:

1. Ensure all parts are clean and dry before assembly.
2. Insert the valve securely into the bottom of the milk collector.
3. Place the diaphragm onto the milk collector, ensuring a tight seal.
4. Attach the flange (or appropriate insert) to the milk collector.
5. Connect the assembled milk collector and flange to the motor unit, ensuring it clicks into place.
6. If using, attach the flange cover for storage or transport.

# Like a baby sucking naturally

Painless Lactation

2nd Gen Motor  
More Power

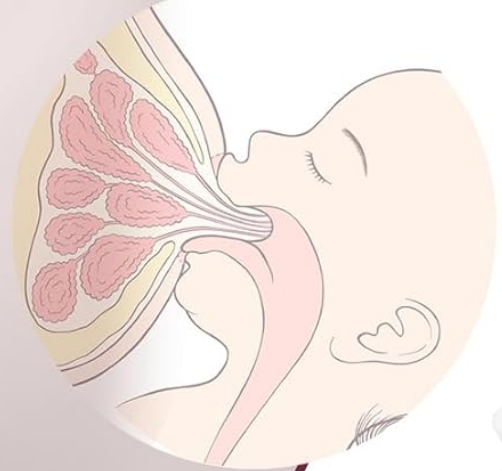

Increased Milk Output

More Comfort

No Pull Pain

Image: An assembled Momcozy M5 Wearable Breast Pump, ready for use.

## 5. FLANGE SIZING GUIDE

Selecting the correct flange size is essential for comfort and effective milk expression. Use the provided nipple ruler to determine your size:



## "I LOVE my M5!

Great pump for a working mom on the go  
As far as the convenience and ease of use go- 10/10. I am very comfortable using these in public.

As a full-time working professional who drives a lot for work and with a child who is primarily fed breast milk this wearable pump is a lifesaver. I would not say that it is not a quality pump, but this pump allows for more quantity pumps, which will eventually be quality. I get the amount of milk my baby needs, and it's more discreet, I can wear it around the office and I only have to clean it once a day. So all in all a good purchase for me, making my life better! 🙌



-- Momcozy Pump User

Image: Visual guide demonstrating how to measure nipple size and select the appropriate flange.

1. Gently stimulate your nipple until it becomes erect.
2. Use the nipple ruler to measure the diameter of your nipple at its base (not including the areola).
3. Refer to the sizing chart on the ruler to find the corresponding flange size. A correct fit ensures your nipple moves freely within the tunnel without excessive friction or compression.

## 6. OPERATING INSTRUCTIONS

The Momcozy M5 offers three distinct pumping modes and nine suction levels to customize your experience. To begin, ensure the pump is fully charged.

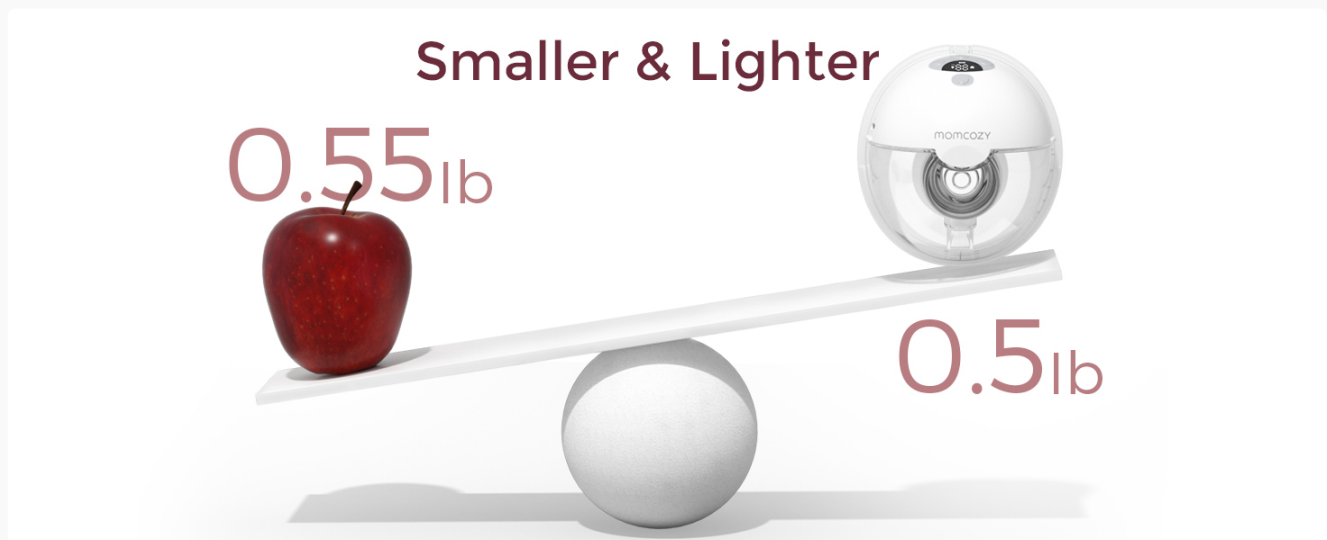


Image: Display of the Momcozy M5 pump illustrating the three modes and nine suction levels.

### Pumping Modes:

- **Stimulation Mode (Heart Icon):** This mode uses fast and shallow suction to stimulate milk letdown, mimicking a baby's initial suckling.
- **Expression Mode (Water Droplet Icon):** Once milk flow begins, switch to this mode for slower, deeper suction to efficiently express milk.
- **Mixed Mode (Heart + Water Droplet Icons):** This mode alternates between stimulation and expression, providing a balanced approach for comfort and efficiency.

### Adjusting Suction Levels:

Use the '+' and '-' buttons on the motor unit to adjust the suction intensity within each mode. There are 9 levels available. Start with a lower level and gradually increase until you find a comfortable and effective setting.

### Hands-Free Operation:

Place the assembled pump securely inside your nursing bra. Ensure the flange is centered on your nipple. The pump is designed to be discreet and allow you to perform other tasks while pumping.

Your browser does not support the video tag.

Video: Demonstrates the Momcozy M5 Wearable Breast Pump in various hands-free scenarios, highlighting its convenience and discreet design.

The pump features an auto shut-off after 30 minutes to protect breast health and conserve battery life.

## 7. CLEANING AND MAINTENANCE

Regular cleaning of your breast pump parts is essential for hygiene and optimal performance. The Momcozy M5 is designed for easy disassembly and cleaning.



Image: Illustrates the minimal number of parts for the Momcozy M5, making cleaning simple.

1. **Disassembly:** Carefully detach the motor unit from the milk collector. Separate all individual parts: flange, diaphragm, valve, and milk collector.
2. **Washing:** Wash all parts that come into contact with breast milk (flange, diaphragm, valve, milk collector) immediately after each use. Use warm, soapy water and a dedicated brush. Rinse thoroughly with clean water.
3. **Sterilization:** Sterilize parts regularly according to manufacturer guidelines (e.g., boiling, steam sterilizer, or dishwasher top rack if compatible).
4. **Drying:** Allow all parts to air dry completely on a clean surface or drying rack. Avoid using cloth towels, which can transfer germs.
5. **Motor Unit:** The motor unit should not be submerged in water. Wipe it clean with a damp cloth.

## 8. TROUBLESHOOTING

If you encounter issues with your Momcozy M5 breast pump, consider the following common troubleshooting steps:

- **No Suction or Weak Suction:**

- Check that all parts are correctly assembled and securely sealed. A loose diaphragm or valve can cause suction loss.
- Ensure the flange is properly fitted to your breast, creating a good seal.

- Verify the motor unit is fully charged.
  - Inspect the valve and diaphragm for any tears, damage, or wear. Replace if necessary.
- **Pump Not Turning On:**
    - Ensure the battery is fully charged. Connect the pump to the charging cable and allow sufficient time to charge.
    - Press and hold the power button for a few seconds to ensure it registers.
  - **Pain or Discomfort During Pumping:**
    - Re-evaluate your flange size using the sizing guide. An incorrect flange size is a common cause of discomfort.
    - Reduce the suction level. Pumping should not be painful.
    - Ensure your nipple is centered in the flange tunnel.

If problems persist, please refer to the warranty and support section for further assistance.

## 9. SPECIFICATIONS

Feature	Specification
Product Dimensions	2 x 3 x 4 inches
Item Model Number	M5-W1
Batteries	1 Lithium Polymer battery required (included)
Item Weight	1.46 pounds
Brand	Momcozy
Color	Cozy Red
Size	1 Count
UPC	810184052739
Manufacturer	Momcozy

## 10. WARRANTY AND SUPPORT

For warranty information, product support, or any questions regarding your Momcozy M5 Wearable Breast Pump, please contact Momcozy customer service directly. You can typically find contact details on the official Momcozy website or through the retailer where you purchased the product.

Keep your proof of purchase for any warranty claims.



