

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

> [Momcozy](#) /

> [Momcozy M5 Smart APP Controlled Wearable Breast Pump Instruction Manual](#)

## Momcozy M5

# Momcozy M5 Smart APP Controlled Wearable Breast Pump Instruction Manual

Model: M5 | Brand: Momcozy

## INTRODUCTION

---

Thank you for choosing the Momcozy M5 Smart APP Controlled Wearable Breast Pump. This manual provides essential information for the safe and effective use of your breast pump. Please read it thoroughly before first use and keep it for future reference.

## SETUP

---

### 1. Unboxing and Components

Carefully unpack all items from the carrying case. The Momcozy M5 Smart comes with the following components:

- 2 Pump Motors
  - 2 Milk Collectors
  - 2 24mm Flanges
  - 2 Flange Covers
  - 2 Diaphragms
  - 2 Valves
  - 2 17mm Flange Inserts
  - 2 19mm Flange Inserts
  - 2 21mm Flange Inserts
  - 2 USB Type-C Cables
  - 1 Flange Size Ruler
  - 1 User Manual
  - 1 Quick Start Guide
  - 1 Storage Case
-



Image: All components of the Momcozy M5 Smart Breast Pump, including two pump units, various flange inserts, charging cables, and documentation, neatly arranged within its storage case.

## 2. Assembly

Assemble the breast pump by carefully connecting the components. Ensure all parts are securely fitted to maintain proper suction.

1. Insert the valve into the milk collector.
2. Place the diaphragm onto the milk collector, ensuring a tight seal.
3. Attach the flange (or appropriate insert) to the milk collector.
4. Secure the pump motor onto the assembled milk collector.



Image: A visual guide demonstrating the simple assembly of the Momcozy M5 Smart Breast Pump, highlighting the four main parts.

## 3. Flange Sizing

Correct flange sizing is crucial for comfort and efficiency. Use the provided nipple measure tool to determine your ideal flange size. The M5 Smart comes with four flange sizes (17/19/21/24mm) to ensure a personalized fit.

Only **4 Parts** to Clean  
Sterilizes in **5 Minutes**

\*Heat safe up to 212°F (100°C)

Image: A visual representation of different flange sizes (15mm, 17mm, 19mm, 21mm, 24mm, 27mm) and a nipple measure tool to help users select the correct size for optimal comfort and milk expression.

#### 4. Charging

Charge the pump motors using the provided USB Type-C cables. A full charge takes approximately 2 hours and provides enough power for up to 6 pumping sessions.

Made With High-Quality,  
Baby-Safe Materials



Image: The Momcozy M5 Smart Breast Pump and its charging cables, illustrating the charging process and battery life.

#### 5. App Connection

Download the Momcozy APP from your smartphone's app store. Follow the in-app instructions to pair your M5 Smart breast pump for remote control and tracking features.

Control Your  
Pump Effortlessly  
From Your Phone

Adjust modes and suction  
levels with ease

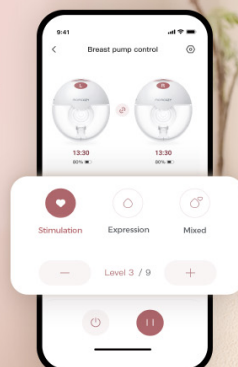


Image: A smartphone displaying the Momcozy app, which allows users to control the M5 Smart breast pump and monitor their lactation data.

## OPERATING INSTRUCTIONS

### 1. Modes and Suction Levels

The M5 Smart offers 3 customizable modes and 9 suction levels to personalize your pumping experience:

- **Stimulation Mode:** Quick, gentle suction to initiate milk flow.
- **Expression Mode:** Deeper, slower suction for efficient milk expression.
- **Mixed Mode:** Combines stimulation and expression for optimal comfort and milk output.

Adjust modes and suction levels directly on the pump unit or via the Momcozy APP.

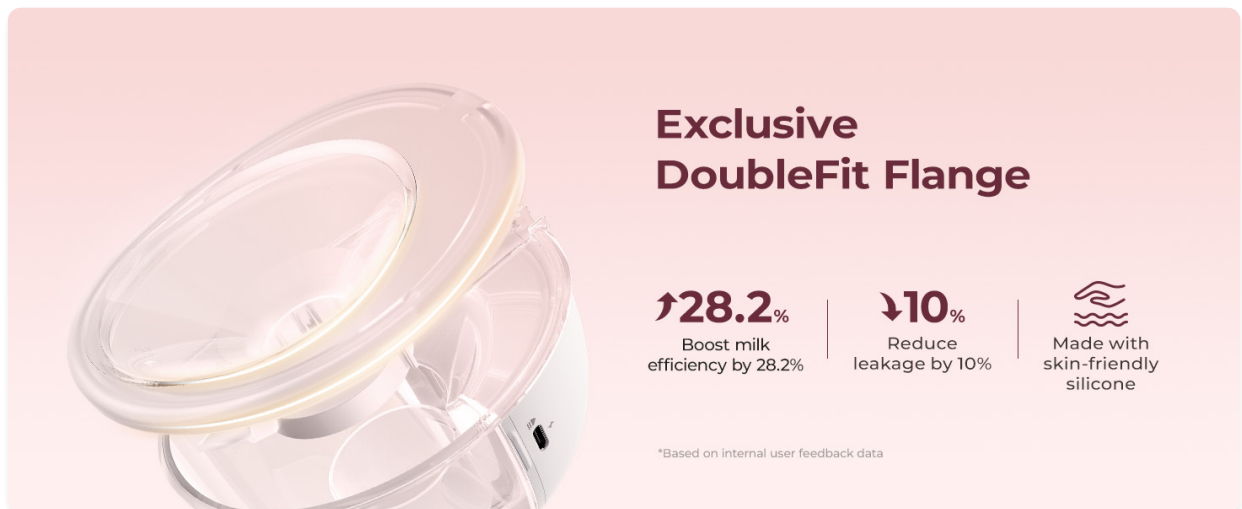


Image: A graphic depicting the three distinct pumping modes (Stimulation, Expression, Mixed) available on the Momcozy M5 Smart, explaining their purpose in the pumping cycle.

### 2. Micro-vibration Technology

The M5 Smart incorporates micro-vibration technology to gently encourage milk flow and clear milk ducts, potentially leading to more efficient pumping sessions.




Image: An animated diagram illustrating the micro-vibration technology within the breast pump, designed to stimulate milk flow and clear ducts.

### 3. Hands-Free Use

The compact and lightweight design of the M5 Smart allows for discreet, hands-free pumping. Simply place the

assembled pump inside your nursing bra and operate it remotely via the app or directly on the unit.



Supports **6 Sessions** Per Charge

Full Charge in **2 Hours**

⚠ Use ONLY 5V/1A adapter. Incompatible with fast/high-output chargers.

Image: A woman is shown comfortably wearing the Momcozy M5 Smart Breast Pump while performing daily activities, demonstrating its hands-free and discreet nature.


Your browser does not support the video tag.

Video: This video demonstrates the Momcozy M5 Smart Breast Pump's hands-free functionality, showing how it can be worn discreetly and used while performing various tasks.

## MAINTENANCE

### 1. Cleaning

The M5 Smart is designed with minimal parts for easy cleaning. Disassemble all milk-contact parts immediately after each use. Wash with warm soapy water and rinse thoroughly.

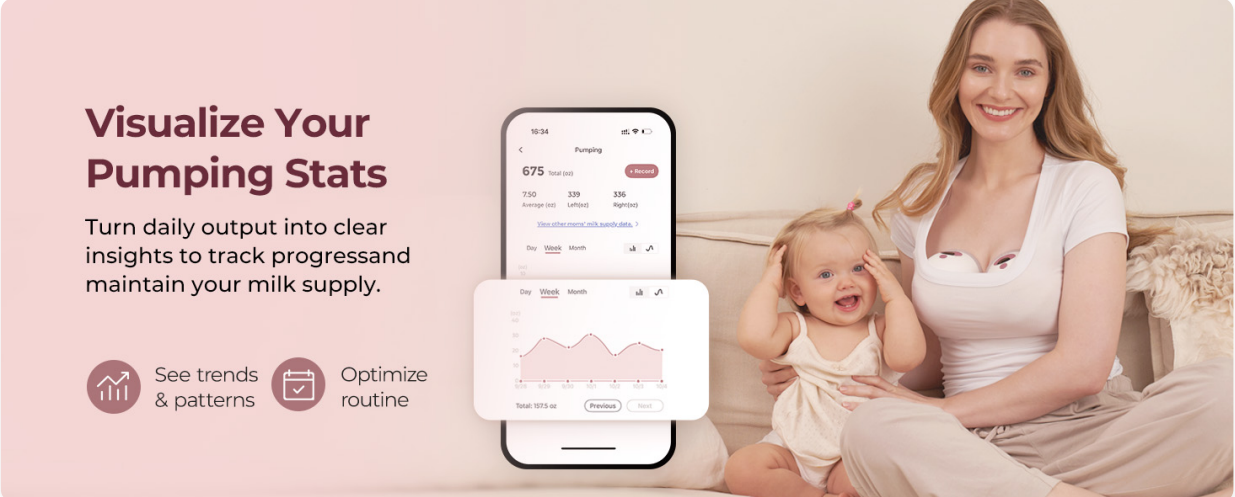

**Visualize Your Pumping Stats**

Turn daily output into clear insights to track progress and maintain your milk supply.

- See trends & patterns
- Optimize routine

Pumping	
675	Total (oz)
750	Average (oz)
339	Left (oz)
336	Right (oz)

Total: 1275 oz

Image: A person is shown cleaning the individual components of the Momcozy M5 Smart Breast Pump under a kitchen faucet, emphasizing its easy-to-clean design.

### 2. Sterilization

All milk-contact parts are made from BPA-free, food-grade silicone and are heat safe up to 212°F (100°C). Sterilize parts by boiling for 3-5 minutes or using a steam sterilizer.



Image: The Momcozy M5 Smart Breast Pump components are displayed alongside icons indicating they are BPA-free, made of food-grade silicone, and safe for sterilization.

### 3. Part Replacement

For optimal performance, it is suggested to replace certain parts like valves and diaphragms every 1-3 months, depending on usage.

## TROUBLESHOOTING

- **No Suction or Weak Suction:** Ensure all parts are correctly assembled and form a tight seal. Check that the valve and diaphragm are not damaged or worn.
- **Pump Not Turning On:** Verify the battery is charged. Connect the pump to the charger and try again.
- **App Connectivity Issues:** Ensure Bluetooth is enabled on your phone and the pump is within range. Restart both the app and the pump if issues persist.
- **Milk Leakage:** Check your flange size using the nipple measure tool. An incorrect flange size can lead to leakage. Ensure the milk collector is not overfilled.

## SPECIFICATIONS

<b>Brand</b>	Momcozy
<b>Model Number</b>	BP380
<b>Package Dimensions</b>	11.14 x 7.2 x 4.49 inches
<b>Item Weight</b>	2.65 Pounds
<b>Modes</b>	Stimulation, Expression, Mixed
<b>Suction Levels</b>	9
<b>Flange Sizes Included</b>	17mm, 19mm, 21mm, 24mm
<b>Battery Life</b>	Up to 6 sessions per full charge
<b>Charging Time</b>	Approx. 2 hours
<b>Material</b>	BPA-Free, Food-Grade Silicone
<b>Additional Features</b>	App Control, LCD Display, Lightweight, Portable, Micro-vibration Technology

## WARRANTY

---

For detailed warranty information, please refer to the official Momcozy website or contact their customer support directly. Warranty terms may vary based on region and purchase date.

## SUPPORT

---

If you have any questions, concerns, or require assistance with your Momcozy M5 Smart Breast Pump, please visit the [Momcozy Official Store](#) or contact their customer service team.

**Your browser does not support the video tag.**

Video: An official Momcozy video showcasing the M5 Smart Breast Pump's features and benefits, emphasizing its role in providing comfortable and convenient care for mothers.