

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Pride Mobility](#) /
- › [Pride Mobility Jazzy Ultra Light Power Wheelchair Instruction Manual \(10 AH Battery\)](#)

### Pride Mobility Jazzy Ultra Light

# Pride Mobility Jazzy Ultra Light Power Wheelchair Instruction Manual

Model: Jazzy Ultra Light (10 AH Battery)

## 1. INTRODUCTION

---

The Pride Mobility Jazzy Ultra Light Power Wheelchair is designed to provide enhanced mobility and convenience. Featuring an all-carbon fiber frame and a lithium battery, this wheelchair is built for durability and portability. It weighs approximately 33 pounds without the battery, making it suitable for transport and storage. The design incorporates a mesh back and seat with a standard cushion for user comfort, along with front suspension for a smooth ride across various terrains.

This manual provides essential information for the safe and effective setup, operation, and maintenance of your Jazzy Ultra Light Power Wheelchair. Please read it thoroughly before using the device.



Figure 1: Front view of the Jazzy Ultra Light Power Wheelchair.

## 2. SETUP

### 2.1 Unfolding the Wheelchair

1. Ensure the wheelchair is on a flat, stable surface.
2. Grasp the backrest and gently pull upwards until the frame locks into its open position.
3. Lower the seat and ensure it is securely in place.
4. Adjust the footrest to the desired position. The footrest can be folded up or down for user comfort and transfer.

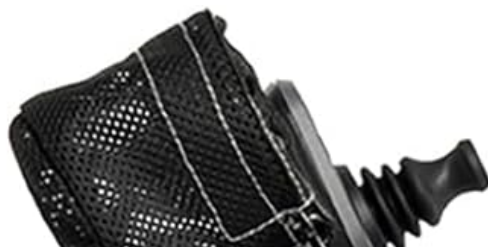




Figure 2: The Jazzy Ultra Light Power Wheelchair in its folded state.

## **2.2 Battery Installation and Charging**

The Jazzy Ultra Light is powered by a lithium-ion battery. The 10 AH battery provides an approximate range of 8 miles at 200 lbs. (A 20 AH battery option is available for up to 16.4 miles range).

1. Locate the battery compartment, typically beneath the seat.
2. Insert the battery pack firmly until it clicks into place.
3. To charge, connect the off-board charger to the battery pack and then to a standard electrical outlet. Ensure the battery is fully charged before first use.

## **2.3 Joystick Controller Attachment**

The programmable joystick controller can be attached to either the left or right armrest, depending on user preference and mobility needs.

1. Slide the joystick mounting bracket onto the desired armrest.
2. Secure the bracket using the provided fasteners.
3. Connect the joystick cable to the corresponding port on the wheelchair.

# **3. OPERATING INSTRUCTIONS**

---

## **3.1 Powering On/Off**

Press the power button on the joystick controller to turn the wheelchair on or off. The battery indicator lights will illuminate upon powering on.

## **3.2 Driving the Wheelchair**

Use the joystick to control the direction and speed of the wheelchair. Push the joystick forward to move forward, pull back to move backward, and push left or right to turn. The maximum speed is approximately 3.69 miles per hour.



# LIGHTWEIGHT TRAVEL REDEFINED

Figure 3: Key features including joystick control and mesh backrest.

### 3.3 Speed Adjustment

The joystick controller features buttons to adjust the speed settings. Use the '+' button to increase speed and the '-' button to decrease speed. Start at a low speed and gradually increase as you become comfortable with the controls.

### 3.4 Safety Features

- **Rear Lights:** The wheelchair is equipped with rear lights that illuminate when the unit is in motion, enhancing visibility.
- **Braking System:** Features an intelligent braking electronic regenerative system and a disc park brake for reliable stopping.
- **Cup Holder:** A mesh cup holder is standard and can be attached to either armrest.

Your browser does not support the video tag.

Video 1: This video demonstrates the key features of the Jazzy Ultra Light Power Wheelchair, including its lightweight design, easy folding mechanism, and joystick controls.

## 4. MAINTENANCE

### 4.1 Battery Care

To maximize the lifespan of your lithium-ion battery:

- Charge the battery regularly, even if the wheelchair is not in frequent use.
- Avoid completely draining the battery before recharging.
- Store the battery in a cool, dry place.

### 4.2 Cleaning

Clean the carbon fiber frame and mesh components with a damp cloth and mild soap. Avoid using abrasive cleaners or solvents that could damage the finish or materials.

### 4.3 Regular Inspections

Periodically check all connections, fasteners, and moving parts for wear or damage. Ensure tires are in good condition. Contact customer support for any necessary repairs or replacement parts.

## 5. TROUBLESHOOTING

If you encounter issues with your Jazzy Ultra Light Power Wheelchair, consider the following common solutions:

- **Wheelchair not powering on:** Check if the battery is properly installed and fully charged. Ensure the power button on the joystick is pressed.
- **Loss of power during operation:** The battery may be low. Recharge the battery immediately.
- **Unresponsive joystick:** Ensure all cables are securely connected. If the issue persists, contact customer support.
- **Unusual noises:** Inspect for any loose parts or obstructions around the wheels or motor. If the noise continues, discontinue use and contact support.

For complex issues or repairs, always consult with authorized service providers.

## 6. SPECIFICATIONS



Figure 4: Infographic summarizing key features and specifications.

Feature	Specification
Model	Jazzy Ultra Light

Feature	Specification
Weight Capacity	300 lbs.
Item Weight (without battery)	33 lbs.
Turn Radius	33.19 Inches
Maximum Speed	Up to 3.69 mph
Ground Clearance	2.17" at motor bracket
Front Wheels	8" solid casters
Drive Wheels	10.71" flat free
Drivetrain	Two motor, rear-wheel drive
Braking System	Intelligent Braking electronic regenerative, disc park brake
Seat Size	18" x 18"
Battery Requirements	10AH, 25.2V Lithium Ion (or 20Ah, 25.2V Lithium Ion)
Standard Battery Charger	Off-board
Frame Material	Carbon Fiber

## 7. WARRANTY AND SUPPORT

The Jazzy Ultra Light Power Wheelchair comes with an extended warranty provided exclusively by Living Well Stores, Inc. This warranty includes batteries, battery charger, wiring, joystick, motor/gearbox, brakes, controller, PC board, throttle, metal connection pins, hardware, nuts and bolts, fuses, lightbulbs, bearings and bushings, seat swivel mechanism, and armrests.

Living Well Stores, Inc. offers a nationwide service network with over 400 local, factory-authorized service providers across the country to support your product. For any questions, service needs, or warranty claims, please contact Living Well Stores, Inc. directly.

The Jazzy Ultra Light is an FDA Class II Medical Device, designed to aid individuals with mobility impairments.