

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Beurer](#) /

› [Beurer SI31 Steam Inhaler Instruction Manual](#)

Beurer SI31

Beurer SI31 Steam Inhaler Instruction Manual

Model: SI31

IMPORTANT SAFETY INFORMATION

Always read the instruction manual completely before first use. Failure to follow these instructions may result in injury or damage to the device.

- Keep the steamer out of reach of children.
- Do not leave the device unattended while in operation.
- Allow the steamer to cool down completely before handling or refilling the water tank.
- Avoid direct contact with hot steam to prevent burns or scalding.
- Use only clean tap water. Do not add essential oils, medications, or other substances to the water tank unless specifically instructed by a healthcare professional and compatible with the device.

PRODUCT OVERVIEW

The Beurer SI31 Steam Inhaler is designed to provide fast, drug-free relief from nasal and sinus congestion, allergies, colds, coughs, and flu symptoms. It features adjustable steam control and a universal mask suitable for all ages. Its filter-free design allows for easy cleaning and maintenance.

Fast Symptom Relief

Ready to use in 3-5 minutes



Breathe easy in 5 minutes



Image: The Beurer SI31 Steam Inhaler, fully assembled and ready for use.

Key Features:

- **Quick & Effective Relief:** Delivers soothing steam in 3-5 minutes for 5-15 minute treatment times.
- **Personalized Steam Comfort:** Adjustable steam output for customized relief.
- **Universal Mask:** Flexible mask fits both mouth and nose, suitable for children and adults.
- **Easy Maintenance:** Removable parts and filter-free design for simple cleaning.
- **Natural Congestion Relief:** Utilizes steam technology to alleviate congestion.
- **Menthol Pad Compatible:** Can be used with menthol inhalant pads (sold separately).

PACKAGE CONTENTS

Verify that all components are present in the package:

- Beurer SI31 Steam Inhaler Base Unit
- Universal Mask

- Water Tank/Reservoir
- Measuring Cup
- Instruction Manual

SETUP

Follow these steps to prepare your Beurer SI31 Steam Inhaler for first use:

1. **Place the Device:** Position the inhaler on a flat, stable surface. Ensure the power cord is fully extended and not obstructed.
2. **Remove Mask and Top:** Gently twist and lift the universal mask and the top cover from the base unit.
3. **Fill with Water:** Using the provided measuring cup, fill the water tank with clean tap water up to the MAX line. Do not overfill. Do not add any other liquids or substances.
4. **Reassemble:** Securely place the top cover back onto the base unit, ensuring it clicks into place. Then, attach the universal mask firmly.



Image: Visual guide demonstrating the 5 steps for quick setup and use of the steam inhaler.



Easy to use

Works with tap water

Filter free & comes apart for easy cleaning

Image: A person filling the water reservoir of the Beurer SI31 Steam Inhaler with tap water using a measuring cup.

OPERATING INSTRUCTIONS

Once the device is set up, follow these steps for operation:

1. **Power On:** Plug the power cord into a suitable electrical outlet. Press the power button on the front of the unit. The indicator light will illuminate.
2. **Heating Time:** Allow 3-5 minutes for the device to heat up and begin producing steam.
3. **Adjust Steam Flow:** Use the adjustable steam control dial located on the side of the unit to regulate the intensity of the steam output to your comfort level.
4. **Inhale:** Position your face comfortably over the universal mask, ensuring a good seal around your nose and mouth. Breathe normally through the mask.
5. **Treatment Duration:** Inhale the steam for 5-15 minutes, or as recommended by your healthcare professional.
6. **Power Off:** After your treatment, press the power button to turn off the device. Unplug the power cord from the outlet.

Natural & Effective Relief

Adjustable steam levels

5-15 minute treatment time

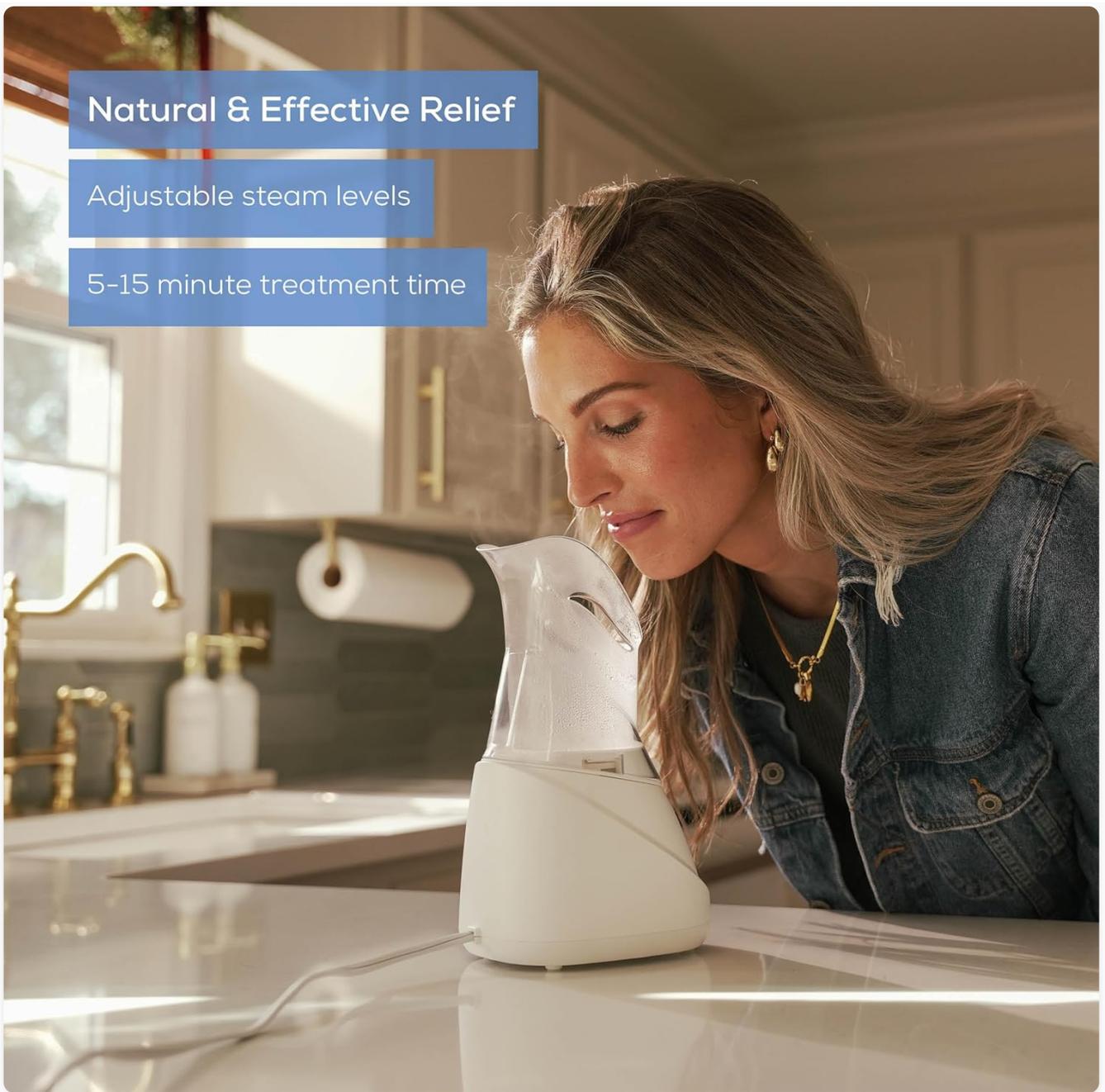


Image: A woman using the Beurer SI31 Steam Inhaler, demonstrating proper inhalation technique.

Soothing Steam Therapy

Customizable steam control

Ultra-quiet operation



Image: A close-up view of the Beurer SI31 Steam Inhaler producing steam, highlighting its customizable steam control.

Using with Menthol Inhalant Pads (Optional):

If desired, menthol inhalant pads (sold separately) can be used with the SI31. Refer to the instructions provided with the menthol pads for proper placement and usage. Do not place menthol pads directly into the water tank.



Compatible With Menthol Inhalant Pads (sold separately)

Image: The Beurer SI31 Steam Inhaler shown alongside menthol inhalant pads, indicating compatibility.

Instructional Video:

Your browser does not support the video tag.

Video: A comprehensive guide demonstrating the setup and operation of the Beurer SI31 Steam Inhaler.

MAINTENANCE AND CLEANING

Regular cleaning ensures optimal performance and hygiene of your steam inhaler. Always unplug the device and allow it to cool completely before cleaning.

1. **Disassemble:** Remove the universal mask and the top cover from the base unit.
2. **Empty Water Tank:** Carefully pour out any remaining water from the tank.
3. **Clean Components:** Wash the universal mask, top cover, and water tank with warm, soapy water. Rinse thoroughly with clean water.
4. **Wipe Base Unit:** Wipe the exterior of the base unit with a damp cloth. Do not immerse the base unit in water.

5. **Dry Thoroughly:** Ensure all parts are completely dry before reassembling or storing the device.
6. **Descaling (if necessary):** If mineral deposits build up, mix a solution of equal parts white vinegar and water. Fill the water tank with this solution and let it sit for 20-30 minutes. Empty and rinse thoroughly. Do not turn on the device with the descaling solution.



Image: The Beurer SI31 Steam Inhaler disassembled, illustrating its filter-free design and ease of cleaning.

TROUBLESHOOTING

If you encounter issues with your Beurer SI31 Steam Inhaler, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Not plugged in; Power button not pressed; Power outage.	Ensure the power cord is securely plugged into a working outlet. Press the power button firmly. Check your household's power supply.

Problem	Possible Cause	Solution
No steam or very little steam produced.	Insufficient water; Device not heated up; Mineral buildup.	Check water level and refill if necessary. Allow 3-5 minutes for heating. Perform descaling procedure if mineral buildup is suspected.
Steam is too hot or too cold.	Steam control setting.	Adjust the steam control dial to increase or decrease steam intensity to your preferred comfort level.
Unusual odor during operation.	First use; Residue from previous use; Mineral buildup.	A slight odor may occur during first use. Ensure all parts are clean. Perform descaling if necessary. If odor persists, contact customer support.

If the problem persists after attempting these solutions, please contact Beurer customer support for further assistance.

SPECIFICATIONS

- **Model:** SI31
- **Brand:** Beurer
- **Product Dimensions:** 7.6 x 6 x 9.7 inches
- **Item Weight:** 1.86 Pounds
- **Material:** Plastic
- **Color:** White
- **Power Source:** Corded Electric
- **UPC:** 810110200265

WARRANTY AND SUPPORT

Beurer products are manufactured with quality and reliability in mind. For specific warranty information, please refer to the warranty card included with your product or visit the official Beurer website.

For technical support, troubleshooting assistance, or to inquire about replacement parts, please contact Beurer customer service. Contact details can typically be found on the Beurer website or in the product packaging.

Online Resources: [Visit the Beurer Store on Amazon](#)

