

## BEBIRD Bebird Earsight Plus D39R

# Bebird Earsight Plus D39R Instruction Manual

Your comprehensive guide to safe and effective ear care.

### 1. PRODUCT OVERVIEW

The Bebird Earsight Plus D39R is an innovative ear care device designed to provide a clear, visual method for ear wax removal and general ear canal inspection. It combines a high-definition camera with various cleaning tools, ensuring precision and safety for users of all ages.



Figure 1: Bebird Earsight Plus D39R with all included accessories.

## 2. SETUP AND INITIAL USE

### 2.1 Unboxing and Component Identification

Carefully remove all components from the packaging. The main device, accessory case with various tips, charging cable, and protective plugs should be present.

Video 1: Unboxing the Bebird Earsight Plus D39R and identifying its components.

### 2.2 Charging the Device

Connect the provided USB-C charging cable to the device's charging port and a compatible USB power source. The indicator light will show charging status.

### 2.3 App Download and Connection

Download the official Bebird app from your smartphone's app store. Power on the device by pressing and holding the power button until a blue light illuminates. Open the Bebird app, connect to the device's Wi-Fi network (named 'bebird-xx-xxxxxx'), and then select 'Begin to Clean Ears' to start using the device.

## 2.4 Attaching Tips and Protective Plugs

Select the appropriate tip (ear scoop, tweezers, scraper, sponge) for your needs. Attach it securely to the device. For safety, especially with children or first-time users, utilize the ear canal protective plug to prevent over-insertion. The flexible hose can be bent to fit your ear canal, ensuring the scale side faces forward.

# Flexi-otoscope Fits Better



# Delivers Crystal Clear Vision

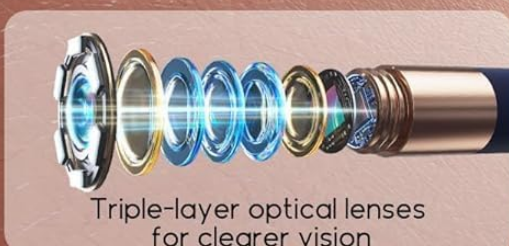


Figure 2: The flexible otoscope can be bent to fit various ear canal shapes.

Video 2: Demonstration of bending the flexible otoscope for optimal use.

### 3. OPERATING INSTRUCTIONS

#### 3.1 Ear Wax Removal

Once connected to the app, gently insert the device into the ear canal while observing the real-time view on your smartphone. Use the appropriate tip for the type of earwax:

- **Scooper Tips:** Ideal for removing regular earwax.
- **Tweezers:** Effective for stubborn earwax or debris stuck to the ear canal walls. These are designed to be gentle on sensitive ear canals.
- **Scraper:** Used for extracting debris.
- **Sponge:** Suitable for removing wet earwax.



# Home Check-Up Must-Have

Check ears, nose, tooth and throat for your family



Figure 3: Illustration of various ear cleaning methods using different tips.

# Ear Hygiene up to 100%

Always get ear cleaning right!

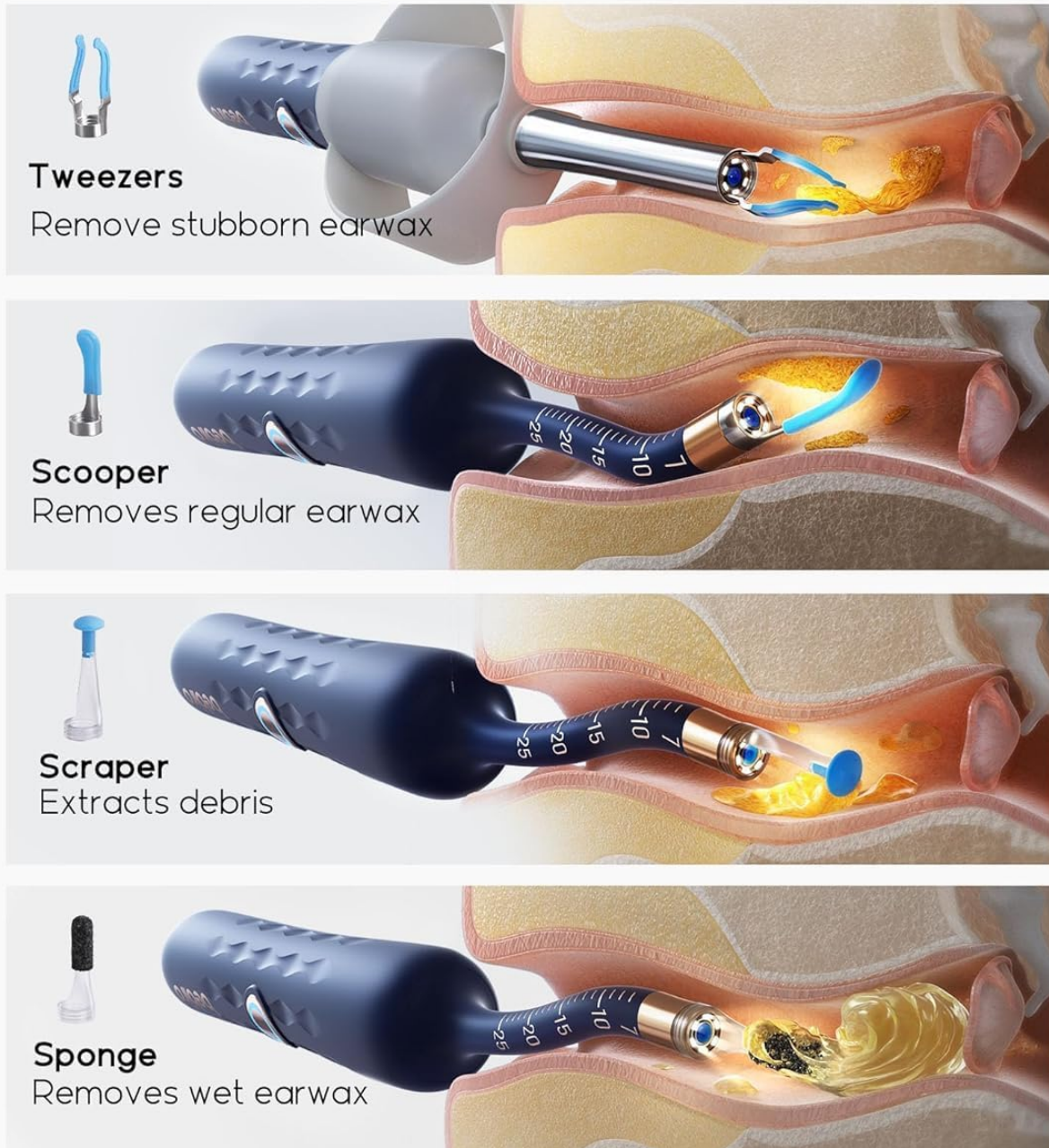


Figure 4: The painless ear tweezers are designed for gentle and precise removal of earwax.

## 3.2 Multi-purpose Inspection

The Bebird Earsight Plus can also be used for inspecting other body parts. Attach the appropriate speculum (nasal or oral) for visual checks of teeth, nasal passages, throat, or scalp. This allows for confident and safe ENT checks at home.

# Clean and Hygienic Storage

Organize all ear cleaning tips and get them ready!



Figure 5: The device can be used for ear, nasal, and oral examinations.

## 4. MAINTENANCE

### 4.1 Cleaning the Device and Tips

After each use, carefully remove the tip and clean it thoroughly with the provided lens brush or a soft, damp cloth. Ensure no debris remains on the camera lens. The tips can be washed with mild soap and water.



## 4.2 Storage

Store the Bebird Earsight Plus D39R and all its accessories in the provided accessory case. This ensures hygienic storage and keeps all components organized and ready for next use.



Figure 6: The compact case provides clean and hygienic storage for the device and all tips.

## 5. TROUBLESHOOTING



### 5.1 Connection Issues

If the device fails to connect to your smartphone, ensure the device is fully charged and powered on. Verify that your phone's Wi-Fi is enabled and connected to the Bebird device's network. Restarting both the device and the app may resolve the issue.

### 5.2 Image Clarity Problems

If the image appears blurry or unclear, ensure the camera lens is clean and free of any earwax or debris. Use the provided lens brush for cleaning. Also, check your phone's display settings for optimal brightness and resolution.

## 6. SPECIFICATIONS

Item Model Number	Bebird Earsight Plus D39R
Package Dimensions	7.91 x 5.31 x 1.34 inches
Item Weight	9.59 ounces
Batteries	2 AA batteries required (not included)
Manufacturer	BEBIRD
First Available Date	November 27, 2024

## 7. WARRANTY AND SUPPORT

### 7.1 Warranty Information


The Bebird Earsight Plus D39R comes with a 1-year warranty from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use conditions.

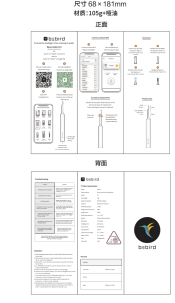
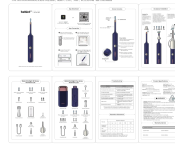


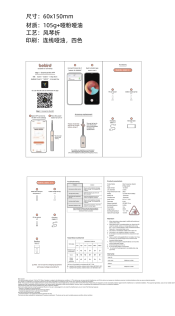
### 7.2 Customer Support

For any questions, technical assistance, or warranty claims, please contact our customer service team. Refer to the product packaging or the official Bebird website for the most up-to-date contact information.

© 2024 BEBIRD. All rights reserved.

### Related Documents - Bebird Earsight Plus D39R

	<a href="#">bebird Home 30S Smart Ear Cleaning Device User Manual</a> Comprehensive user manual for the bebird Home 30S Smart Ear Cleaning Device, covering setup, app connection, accessory usage, troubleshooting, and safety guidelines.
---	--

	<p><a href="#">Bebird C3 Smart Visual Ear Cleaning Rod - User Manual and Specifications</a></p> <p>Comprehensive guide for the Bebird C3 smart visual ear cleaning rod, covering setup, usage, troubleshooting, product specifications, and warranty information. Learn how to connect, use accessories, and maintain your device.</p>
	<p><a href="#">Bebird Smart Visual Ear Cleaner SE &amp; Plus Series User Manual</a></p> <p>Comprehensive guide for Bebird Smart Visual Ear Cleaner SE and Plus series, covering device overview, app connection, accessory installation, troubleshooting, product specifications, and safety warnings.</p>
	<p><a href="#">bebird Ear-care: Your Guide to Visual Ear Inspection &amp; Cleaning</a></p> <p>Explore the bebird Ear-care digital endoscope with this user manual. Learn app connection, device operation, accessory installation, troubleshooting, and safety guidelines for effective ear care.</p>
	<p><a href="#">Bebird C3 Smart Visual Ear Cleaning Rod User Guide</a></p> <p>This document provides instructions and specifications for the Bebird C3 Smart Visual Ear Cleaning Rod, including setup, usage, troubleshooting, and product details.</p>
	<p><a href="#">Bebird K3 Visual Ear Cleaner User Manual and Specifications</a></p> <p>Detailed user guide for the Bebird K3 and K3 Pro visual ear cleaners, covering app setup, Wi-Fi connection, accessories, troubleshooting, technical specifications, safety warnings, and FCC compliance.</p>