

BEBIRD NOTE5

BEBIRD Note 5 Ear Camera Wax Removal Tool User Manual

Model: NOTE5

1. INTRODUCTION

Thank you for choosing the BEBIRD Note 5 Ear Camera Wax Removal Tool. This innovative device combines a high-definition ear camera with precision tweezers and ear spoons, offering a visual and effective solution for earwax removal and ear canal inspection. Designed for safety and ease of use, the Note 5 provides a clear view of the ear canal, allowing for precise and gentle cleaning. Please read this manual thoroughly before use to ensure proper operation and to maximize the benefits of your device.



Image 1.1: The BEBIRD Note 5 device, showing its camera, tweezers, and included ear spoons, alongside a smartphone displaying the live ear canal view.

2. IMPORTANT SAFETY INFORMATION

Your safety is paramount. Please adhere to the following guidelines:

- **Do not insert the device too deeply into the ear canal.** The device is equipped with a safety cap to prevent excessive penetration.
- Always use the device gently. Avoid applying excessive force, especially when using the tweezers or ear spoons.
- Ensure the ear canal is dry before use.
- If you experience any pain, discomfort, or dizziness during use, stop immediately and consult a medical professional.
- Keep the device out of reach of children and pets when not in use. Adult supervision is recommended when used on children.

- Clean the ear spoons and tweezers thoroughly before and after each use to prevent cross-contamination.
- Do not use the device if you have an ear infection, perforated eardrum, or any other ear injury.
- The device is designed for personal use. For medical diagnoses, consult a healthcare professional.

Safe Ear Cleaning

The safety cap design, which limits the depth of entry into the ear, ensures safer and more reliable ear canal cleaning.



Image 2.1: Illustration of the safety cap feature, designed to prevent over-insertion into the ear canal.

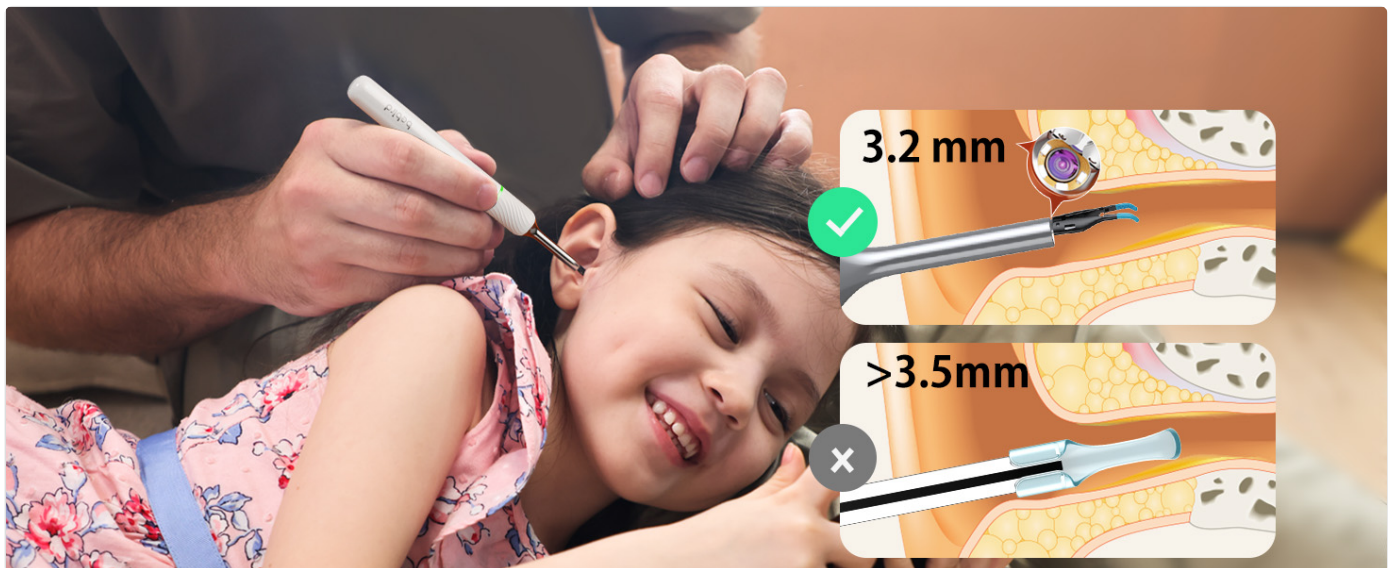


Image 2.2: The BEBIRD Note 5 features a 3.2mm ultra-fine lens, suitable for narrower ear canals, including those of children and pets.

3. PACKAGE CONTENTS

Please check the package to ensure all items are present:

- BEBIRD Note 5 Ear Camera Device
- Assorted Ear Spoons (various colors and types)
- Charging Cable (Type-C)
- Storage Case
- User Manual (this document)



Image 3.1: All components included in the BEBIRD Note 5 package.

4. SETUP GUIDE

1. **Download the BEBIRD App:** Search for "Bebird" in the Apple App Store (for iOS devices) or Google Play Store (for Android devices) and download the official application.
2. **Power On the Device:** Remove the magnetic cap from the BEBIRD Note 5. The device will automatically power on.
3. **Connect via Wi-Fi:** On your smartphone or tablet, go to your Wi-Fi settings. Locate and connect to the Wi-Fi network named "bebird-Note5-XXXXXX" (where XXXXXX is a unique device identifier). No password is required.
4. **Launch the App:** Open the BEBIRD app. The live view from the ear camera should now be displayed on your screen.

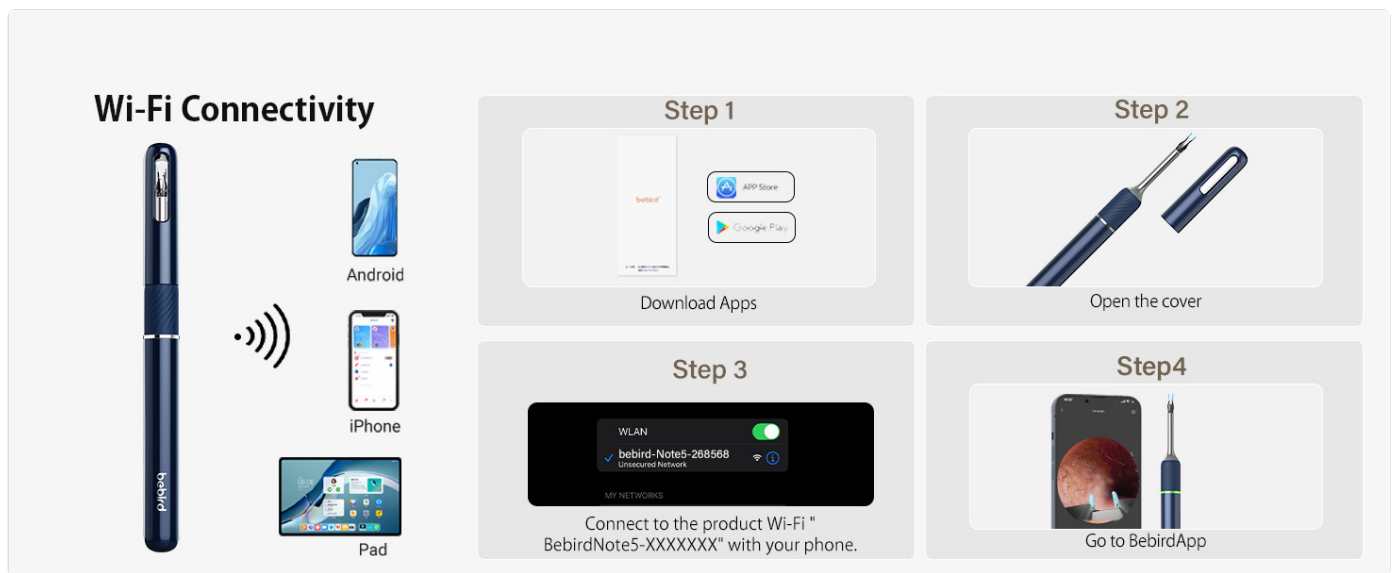


Image 4.1: Step-by-step instructions for connecting the BEBIRD Note 5 to your smartphone via Wi-Fi.

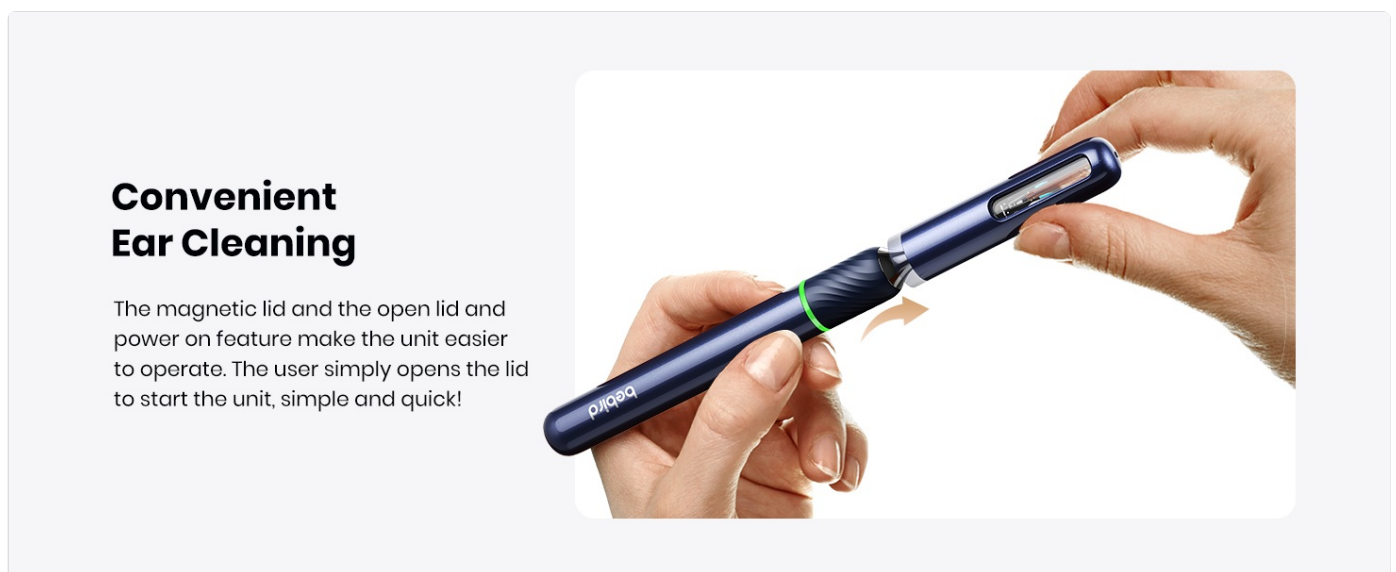


Image 4.2: The magnetic lid provides convenient automatic power-on when opened.

5. OPERATING INSTRUCTIONS

5.1. Understanding the Modes

The BEBIRD Note 5 offers three primary modes for ear care:

- **Tweezer Mode:** For gripping and removing larger earwax or foreign objects. The high-precision tweezers have a 0.05mm error margin and a maximum clamping force of 0.98 newtons.
- **Earscoop Mode:** For general ear cleaning and gentle removal of earwax using the provided silicone ear spoons.
- **Observe Mode:** For visual inspection of the ear canal without active cleaning, useful for monitoring ear health or showing a doctor.

3-in-1 Function



Image 5.1: The 3-in-1 functionality of the BEBIRD Note 5, showing its tweezers, earscoop, and observation capabilities.

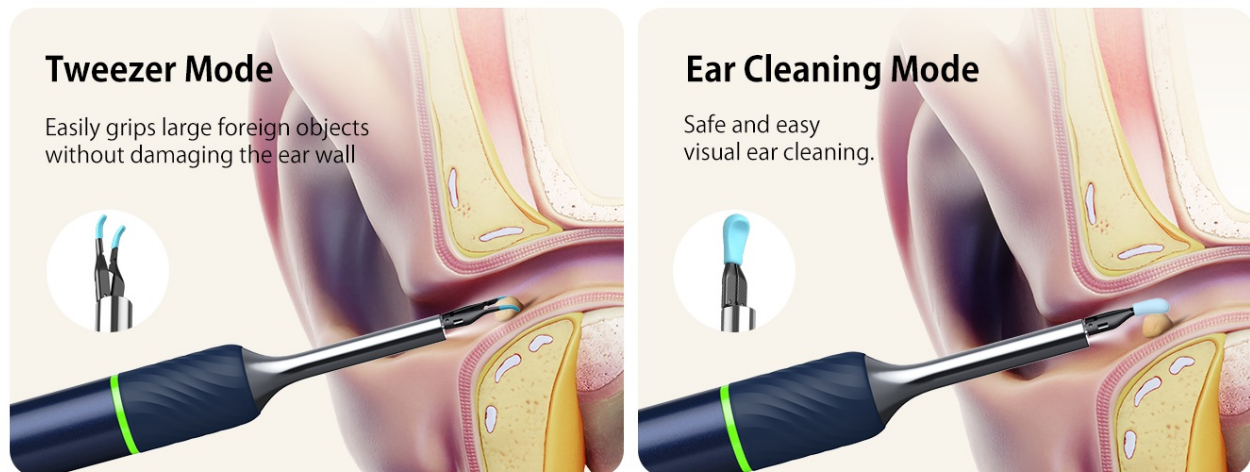


Image 5.2: Visual representation of Tweezer Mode for foreign object removal and Ear Cleaning Mode for wax removal.

5.2. Using the Tweezers and Ear Spoons

- To use the tweezers, gently press the 360-degree ring button on the device handle to control the opening and closing of the tweezers.
- For earscoop functionality, attach one of the food-grade silicone ear spoons onto the tip of the tweezers. Ensure it is securely fastened.
- Always operate with a steady hand and refer to the live view on your connected device for precise movements.

Technology and Painless Ear Cleaning

The 360-degree ring switch is easy to grasp and adopts an ergonomic design. The tweezers can be easily controlled by pressing manually to achieve painless cleaning.



Image 5.3: The ergonomic 360-degree ring button allows for precise control of the tweezers.

5.3. App Features

The BEBIRD app enhances your ear cleaning experience:

- **Live Visualization:** View your ear canal in real-time with the 10-megapixel camera.
- **Photo & Video Recording:** Capture images and videos of your ear canal for personal records or to share with a medical professional.
- **Remote Consultation:** The app facilitates sharing images and videos, allowing for online consultations with doctors, potentially saving time and medical expenses.

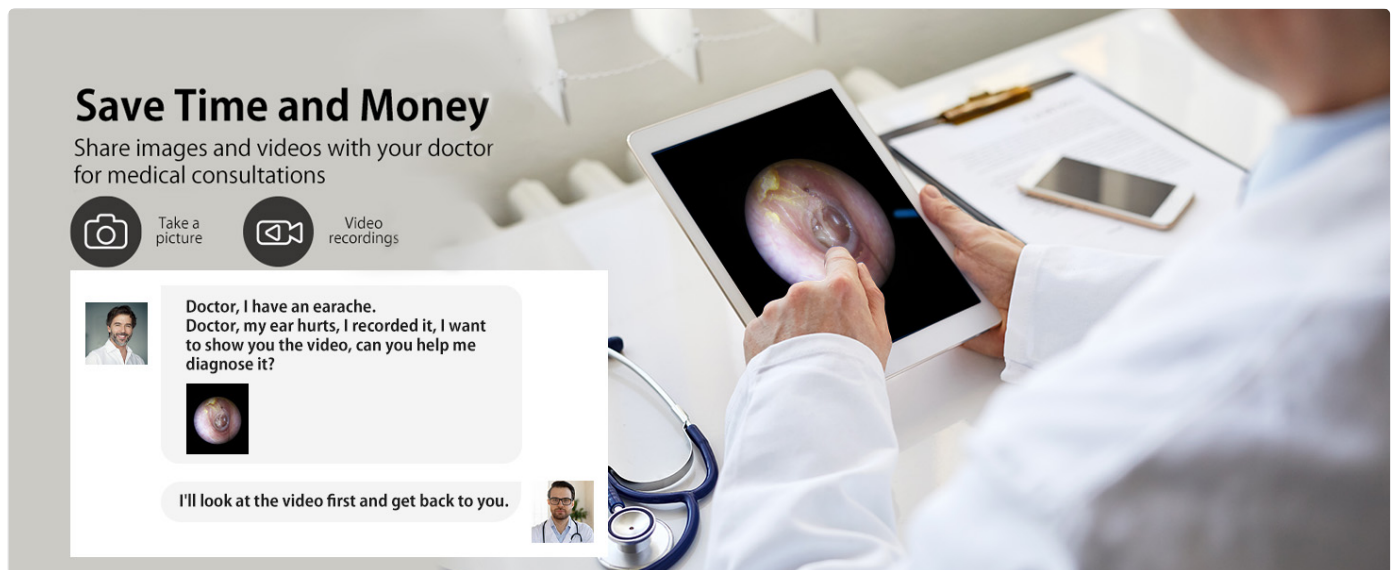


Image 5.4: The BEBIRD app allows users to capture images and videos for medical consultation, saving time and money.

Use App to Consult with Your Doctor Anytime, Anywhere to Make Sure Your Ears are Healthy!

With APP video transmission and remote doctor consultation, this time and money saving app also allows you to assess your ear health online for your family and ear health.



Image 5.5: The BEBIRD app facilitates remote consultations with medical professionals for ear health assessment.

5.4. Temperature Control

The BEBIRD Note 5 features dual-layer temperature control to ensure a comfortable experience. The device maintains a temperature close to body temperature (lens at 25°C, handle at 32°C) to prevent discomfort during use.

Comfortable Ear Cleaning

Comfortable ear cleaning with lens temperature control ensures user comfort and safety and improves the overall performance of the visual ear spoons

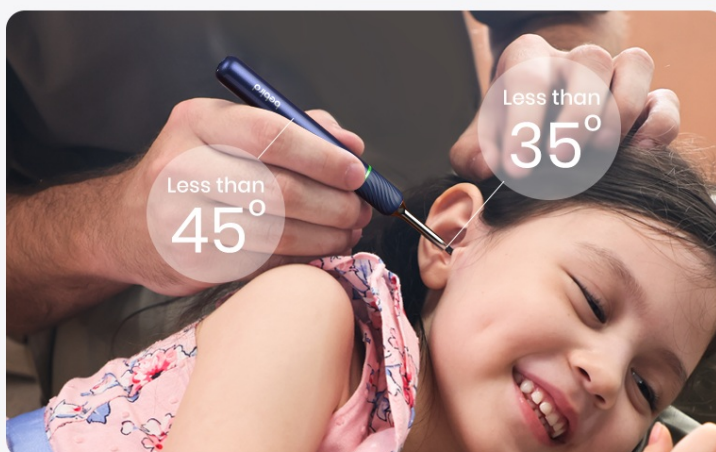


Image 5.6: The device's temperature control ensures a comfortable and safe ear cleaning experience.

6. MAINTENANCE

- **Cleaning:** After each use, carefully remove the ear spoon or clean the tweezers. Wipe the lens and tips with a soft, dry cloth or an alcohol wipe. Ensure no debris obstructs the camera lens.
- **Charging:** When the battery is low, connect the device to a power source using the provided Type-C charging cable. The device has a 230 mAh battery.
- **Storage:** Store the BEBIRD Note 5 in its protective case in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Ear Spoon Replacement:** The multi-color replacement ear spoons are designed for individual family members to prevent cross-contamination. Replace them as needed.

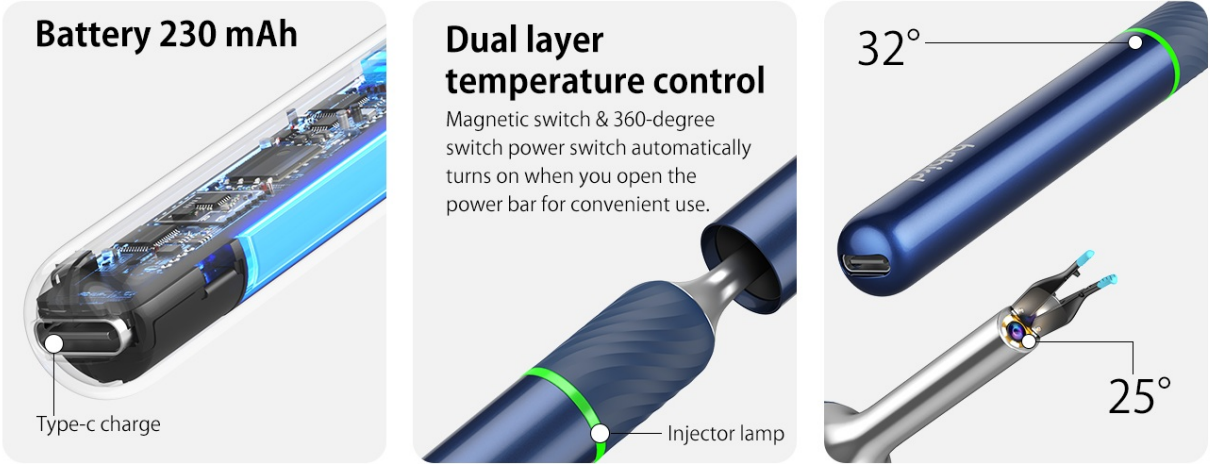


Image 6.1: Overview of the device's battery, charging port, and temperature control features.

Real Family Time



Three-colored 6 tips replacement and an observation cap for each family member


Image 6.2: The multi-colored ear spoons help prevent cross-contamination when used by multiple family members.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not power on.	Low battery or magnetic cap not removed.	Charge the device using the Type-C cable. Ensure the magnetic cap is fully removed.
Cannot connect to Wi-Fi.	Device not powered on, incorrect Wi-Fi selected, or app not launched.	Ensure the device is on. Select the correct "bebird-Note5-XXXXXX" Wi-Fi network. Relaunch the BEBIRD app.
Image is blurry or unclear.	Camera lens is dirty or fogged.	Clean the camera lens with a soft, dry cloth. The device is dustproof and fog-resistant, but cleaning is still recommended.
Tweezers not gripping effectively.	Misalignment or damage to the tweezer jaws.	Ensure the tweezers are properly aligned. If damaged, contact customer support for assistance.

8. SPECIFICATIONS

- Model:** NOTE5
- Camera Resolution:** 10 Megapixels
- Lens Diameter:** 3.2 mm
- Battery Capacity:** 230 mAh
- Charging Port:** Type-C
- Connectivity:** Wi-Fi
- Tweezer Precision:** 0.05mm error
- Max Clamping Force:** 0.98 Newton
- Temperature Control:** Dual-layer (Lens 25°C, Handle 32°C)
- Dimensions:** 7.83 x 3.39 x 0.79 inches (Package)
- Weight:** 4.16 ounces (Package)
- Manufacturer:** BEBIRD

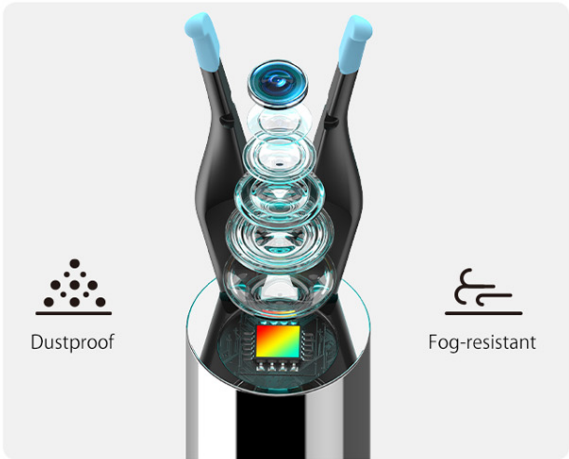

Dustproof

Fog-resistant

Power Mela with 10 Megapixels & Led Light

High-definition endoscope allows you to clearly distinguish the intra-aural situation at a glance

10 million pixels ✓

BEBIRD


3 million pixels ✗

Other

Image 8.1: The 10-megapixel camera provides superior image quality compared to lower resolution models.

9. WARRANTY AND SUPPORT

BEBIRD provides a one-year warranty for the Note 5 Ear Camera Wax Removal Tool. If you encounter any issues or have questions regarding the use of your device, please do not hesitate to contact our customer service team.

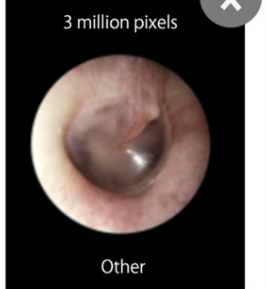
Customer Support: For assistance, please contact us directly via email. Refer to the product packaging or the BEBIRD official website for contact details.


10. PRODUCT VIDEOS

No specific video URLs with creator type "Seller" were provided in the product data. If available, relevant instructional videos would be embedded here to demonstrate product features and usage.



Related Documents - NOTE5

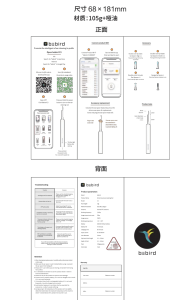
	Bebird Note5 Smart Visual Ear Cleaner User Manual and Specifications Comprehensive guide to the Bebird Note5 smart visual ear cleaner, covering app connection, product usage, troubleshooting, specifications, and safety information. Learn how to use the device for effective and safe ear cleaning.
	Bebird C3 Smart Visual Ear Cleaning Rod - User Manual and Specifications Comprehensive guide for the Bebird C3 smart visual ear cleaning rod, covering setup, usage, troubleshooting, product specifications, and warranty information. Learn how to connect, use accessories, and maintain your device.
	Bebird C3 Smart Visual Ear Cleaning Rod User Guide This document provides instructions and specifications for the Bebird C3 Smart Visual Ear Cleaning Rod, including setup, usage, troubleshooting, and product details.
	bebird Home 30S Smart Ear Cleaning Device User Manual Comprehensive user manual for the bebird Home 30S Smart Ear Cleaning Device, covering setup, app connection, accessory usage, troubleshooting, and safety guidelines.
	bebird Ear-care: Your Guide to Visual Ear Inspection & Cleaning Explore the bebird Ear-care digital endoscope with this user manual. Learn app connection, device operation, accessory installation, troubleshooting, and safety guidelines for effective ear care.

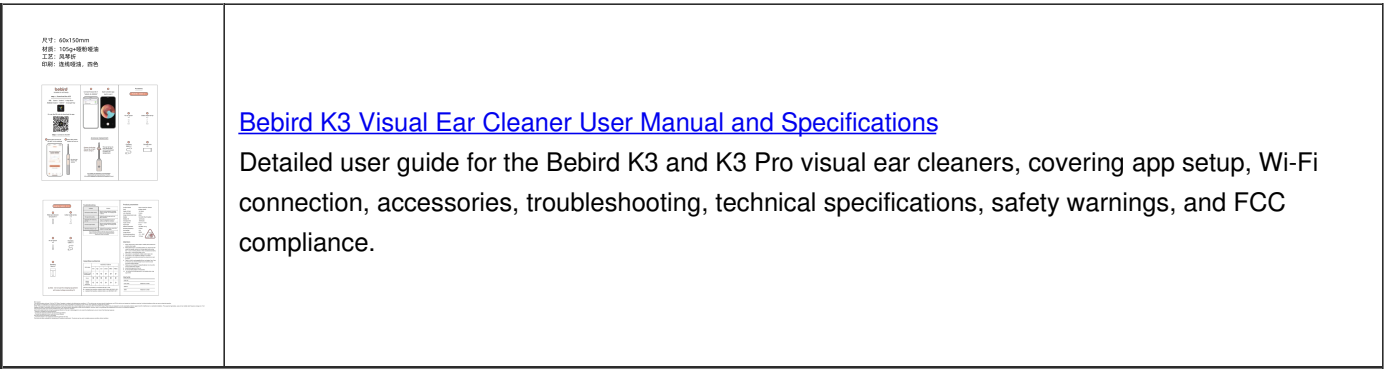

[Bebird K3 Visual Ear Cleaner User Manual and Specifications](#)

Detailed user guide for the Bebird K3 and K3 Pro visual ear cleaners, covering app setup, Wi-Fi connection, accessories, troubleshooting, technical specifications, safety warnings, and FCC compliance.

Documents - BEBIRD – NOTE5

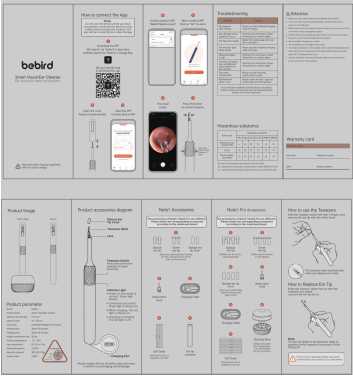

[\[pdf\] User Manual Installation Guide Troubleshooting Guide](#)

Note5 Pro 20230727 User manual eng Bebird Note 5 otoszkóp fültisztításhoz kamerával kép 39660 36729 79b5ef59308f4272e70869800f3bb333 puputeve hu media wysiwyg attachments ||| Smart Visual Ear Cleaner Ear cleaning for better ear protection Recommended charging equipment with ... age sensor IEEE 802.11b/g Name CMOS Model Inspectors Production Date Smart Visual Ear Cleaner **NOTE5** Series Product accessories diagram Silicone Ear Tip Cover Tweezers Head Lens **Note5** Accessori... lang:en score:38 filesize: 8.09 M page_count: 1 document date: 2023-12-18



[\[pdf\] Declaration of Conformity](#)

Declaration Letter of Model Difference Heifeng Zhizao Shenzhen Technology Co Ltd NOTE5 Smart visual ear clean Rod 2AU63 2AU63NOTE5 note5 Heifeng Zhizao Shenzhen Technology Co.,Ltd Add: Floor 1, Building 6, Gangzai Street, Furong Industrial Zone, Xiangshan Community, Xinqiao Street, Bao'an District, Shenzhen China Date: 2022-8-3 FCC ID: 2AU63-**NOTE5** Dear Sir/Madam, We, Heifeng Zhizao Shenzhen Technology Co.,Ltd Hereby declares the f... lang:tl score:29 filesize: 120.92 K page_count: 1 document date: 2022-08-10



[\[pdf\] Declaration of Conformity SAR Rating](#)


tcIsevers RF Exposure Evaluation Heifeng Zhizao Shenzhen Technology Co Ltd NOTE5 Smart visual ear clean Rod 2AU63 2AU63NOTE5 note5 Page 1 of 4 FCC ID: 2AU63-**NOTE5** FCC RF Exposure Evaluation 1. Product Information FCC ID Product name Test Model Additional Model No. Model Declaration Power supply Bluetooth Channel Number Modulation Type Antenna Type Antenna Gain Hardware version Software version Exposure category EUT Type Dev... lang:en score:28 filesize: 340.6 K page_count: 4 document date: 2022-08-11



[pdf] Test Report


TELEC LEO RADIO TEST REPORT For Heifeng Zhizao Shenzhen Technology Jul 26 2022 · Note5
bebird Pro Max Plus 5 for 2 4G WIFI RF Test Data Note SA set to from 4965GHz 219 229169 02 005
tele soumu go jp giteki SearchServlet2 PageID jt01 ATF 96499005 |||

Page 1 of 24 Report No.: LCSA071222049EA RADIO TEST REPORT For Heifeng Zhizao Shenzhen Technology Co.,Ltd Smart visual ear-clean Rod Test Model: **Note5**
Additional Model No.: Please Refer to Page 6 Prepared for Address : Heifeng Zhizao Shenzhen Technology Co.,Ltd : Floor 1, Building 6, No.3, Gangz...
lang:en score:28 filesize: 941.47 K page_count: 24 document date: 2022-08-12

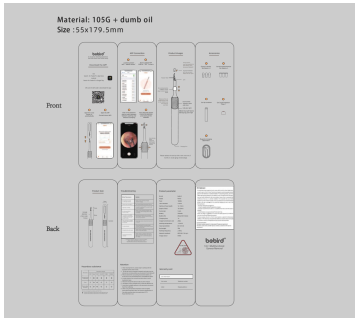

Bebird Note5 Smart Visual Ear Cleaner User Manual and Specifications


Comprehensive guide to the Bebird Note5 smart visual ear cleaner, covering app connection, product usage, troubleshooting, specifications, and safety information. Learn how to use the device for effective and safe ear cleaning.
lang:en score:28 filesize: 1.4 M page_count: 1 document date: 2022-08-03

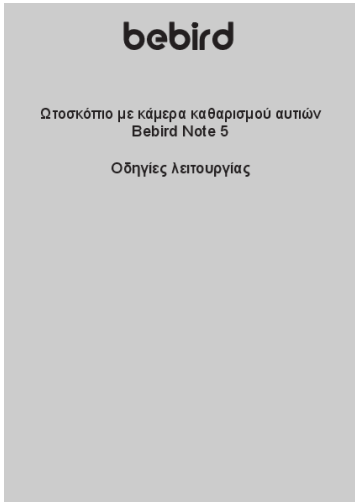

[pdf]

Bebird Note 5 PL docId 487F3BC8DA8E529EA31E2871CB747B0C Οδηγίες λειτουργίας gre otoszkóp
fültisztításhoz kamerával kém 49920 36729 f469aec7f9193693e96ace2ecad91161 puputeve hu media
wysiwyg attachments |||

Bebird Note 5 1: iOS: Android: bebird Google Play QR WiFi. WiFi. bebird**Note5**-
xxxxxxx bebird GO . bebird-**Note5**-323369 bebird-**Note5** 99 . . . LED . . . , . . .
lang:tl score:26 filesize: 1.02 M page_count: 12 document date: 2024-04-08

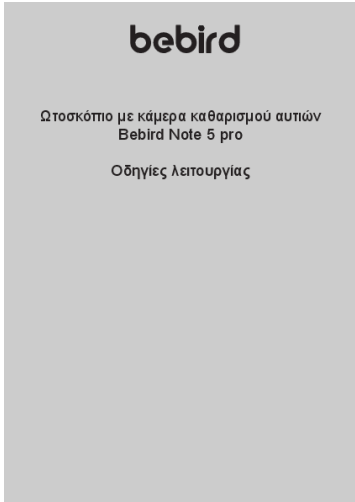

[pdf]

Bebird Note 5 Pro PL docId B5695A743B0A9FAE3B3DD757B17D64CA Οδηγίες λειτουργίας gre pro
otoszkóp fültisztításhoz kamerával fehér 49908 36730 10434ecd0e7cf46fead4dd3bc4b84a62 puputeve
hu media wysiwyg attachments |||

Bebird Note 5 pro 1: iOS: Android: bebird Google Play QR WiFi. WiFi.
bebird**Note5**-xxxxxxx bebird GO . bebird-**Note5**-323369 bebird-**Note5** 99 . . . , . 12
bits Glowic...
lang:tl score:26 filesize: 1.02 M page_count: 12 document date: 2024-04-08

bebird

Otoskop z kamerą do czyszczenia uszu
Bebird Note 5

Instrukcja obsługi

[\[pdf\]](#)

Bebird Note 5 PL Instrukcja obsługi pol otoszkóp fültisztításhoz kamerával kék 39659 36729

4129c2a65f10ab6d8d4d6cd493155060 puputeve hu media wysiwyg attachments |||

Otoskop z kamer do czyszczenia uszu Bebird Note 5 Instrukcja obsługi Krok 1:

Pobranie aplikacji iO ... ie wlcza si automatycznie. Otwrz aplikacj. Polcz urzdzenie z WiFi. Polcz produkt z WiFi. bebird-**Note5**-xxxxxx Powrt do aplikacji bebird Kliknij GO , aby jej uy. bebird-**Note5**-323369 bebird-Note...

lang:pl score:26 filesize: 8.18 M page_count: 12 document date: 2023-12-21

bebird

Otoskop z kamerą do czyszczenia uszu
Bebird Note 5 pro

Instrukcja obsługi

[\[pdf\]](#)

Bebird Note 5 Pro PL Instrukcja obsługi pol pro otoszkóp fültisztításhoz kamerával fehér 39657 36730

e5e7cb2a36429233b21cb918cec530db puputeve hu media wysiwyg attachments |||

Otoskop z kamer do czyszczenia uszu Bebird Note 5 pro Instrukcja obsługi Krok 1:

Pobranie aplikacj ... ie wlcza si automatycznie. Otwrz aplikacj. Polcz urzdzenie z WiFi. Polcz produkt z WiFi. bebird-**Note5**-xxxxxx Powrt do aplikacji bebird Kliknij GO , aby jej uy. bebird-**Note5**-323369 bebird-Note...

lang:pl score:26 filesize: 8.37 M page_count: 12 document date: 2023-12-21

bebird

Otoskop s kamerou na čistenie uší Bebird
Note 5 pro

Návod na obsluhu


[\[pdf\]](#)

Bebird Note 5 Pro PL docId 55A1EB880521B95F268D43602794D9A6 Návod na obsluhu slo pro otoszkóp fültisztításhoz kamerával fehér 49912 36730 c16f032c9696826f123aee3eb2373443 puputeve hu media wysiwyg attachments |||

Otoskop s kamerou na cistenie us Bebird Note 5 pro Nvod na obsluhu Krok 1:

Stiahnutie aplikcie iOS ... zapne. Otvorte aplikciu. Pripojte zariadenie k sieti Wi-Fi. Pripojte vrobok k sieti WiFi. bebird**Note5**-xxxxxx Sp na aplikciu Bebird Kliknutm na tlačidlo GO ju mzete pouzva. Bebird-**Note5**-323369...

lang:hr score:25 filesize: 978.73 K page_count: 12 document date: 2024-04-08

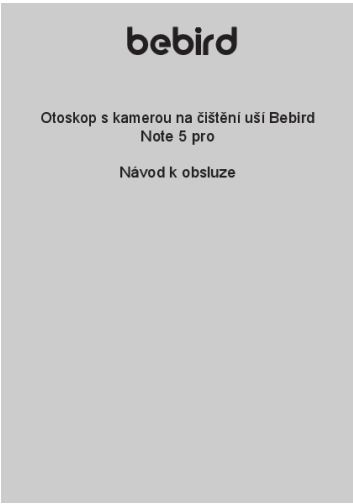

[pdf]

Bebird Note 5 Pro PL docId E05BF0CC612ABE648E8B23CF0655DB33 Návod k obsluze cze pro otoszkóp fültisztításhoz kamerával fehér 49907 36730 c54d032e9d069ef1670e218d3dc53c5a puputeve hu media wysiwyg attachments |||

Otoskop s kamerou na cistn us Bebird Note 5 pro Nvod k obsluze Krok 1: Stazen aplikace iOS: Androi ... tomaticky zapne. Otevete aplikaci. Pipojte zazen k sti Wi-Fi. Pipojte vrobek k sti WiFi. bebird-**Note5**xxxxxxx Zpt do aplikace Bebird Kliknutm na tlactko GO ji mzete pouzt. bebird-**Note5**-323369 ...

lang:hr score:25 filesize: 1019.72 K page_count: 13 document date: 2024-04-08




[pdf] Test Report

FCC PART 15 news Test Report 2 4GWLAN Heifeng Zhizao Shenzhen Technology Co Ltd NOTE5 Smart visual ear clean Rod 2AU63 2AU63NOTE5 note5

Page 1 of 35 FCC ID: 2AU63-**NOTE5** Report No.: LCSA071222058EA FCC TEST REPORT For Heifeng Zhizao Shenzhen Technology Co.,Ltd Smart visual ear-clean Rod Test Model: **Note5** Additional Model No.: Please Refer to Page 6 Prepared for Address : Heifeng Zhizao Shenzhen Technology Co.,Ltd : Floor 1, Buil...

lang:en score:25 filesize: 1.56 M page_count: 35 document date: 2022-08-10

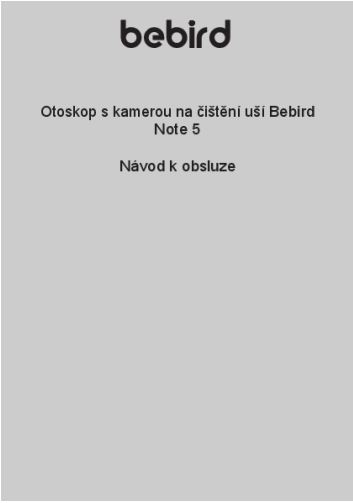

[pdf]

Bebird Note 5 PL docId 4CBDA8817F30C5ED45D520B1E0C2C6EA Návod k obsluze cze otoszkóp fültisztításhoz kamerával kék 49918 36729 da6f3c81cb415963cb97160924fa3d9a puputeve hu media wysiwyg attachments |||

Otoskop s kamerou na cistn us Bebird Note 5 Nvod k obsluze Krok 1: Stazen aplikace iOS: Android: V ... tomaticky zapne. Otevete aplikaci. Pipojte zazen k sti Wi-Fi. Pipojte vrobek k sti WiFi. bebird-**Note5**xxxxxxx Zpt do aplikace Bebird Kliknutm na tlactko GO ji mzete pouzt. bebird-**Note5**-323369 ...

lang:hr score:25 filesize: 1010.76 K page_count: 13 document date: 2024-04-08

bebird

Otoskop s kamerou na čistenie uší Bebird
Note 5

Návod na obsluhu

[\[pdf\]](#)

Bebird Note 5 PL docId F55DD7E42A788251E2754C9838B1AB6B Návod na obsluhu slo otoszkóp
fültisztításhoz kamerával kék 49922 36729 af9dc054217bcef6bb859c8e08995a6f puputeve hu media
wysiwyg attachments |||

Otoskop s kamerou na cistenie us Bebird Note 5 Nvod na obsluhu Krok 1: Stiahnutie
aplikcie iOS: Vy ... zapne. Otvorte aplikciu. Pripojte zariadenie k sieti Wi-Fi. Pripojte
vrobok k sieti WiFi. bebird**Note5**-xxxxxx Sp na aplikciu Bebird Kliknutm na tlačidlo GO
ju mzete pouzva. Bebird-**Note5**-323369...

lang:hr **score:25** filesize: 969.76 K page_count: 12 document date: 2024-04-08

bebird

Otoskop cu cameră de curățare a urechii
Bebird Note 5 pro

Instrucțiuni de utilizare

[\[pdf\]](#)

Bebird Note 5 Pro PL docId 373FF422F6185405404F84D8186511CC Instrucțiuni de utilizare rum pro
otoszkóp fültisztításhoz kamerával fehér 49911 36730 62a0dbddcf9949988e7e4e18bbeb8b25 puputeve
hu media wysiwyg attachments |||

Otoscop cu camer de curare a urechii Bebird Note 5 pro Instrucuni de utilizare Pasul
1: Descrcare ... automat. Deschidei aplicaia. Conectai-v dispozitivul la WiFi. Conectai
produsul la WiFi. bebird-**Note5**xxxxxx napoi la aplicaia bebird Facei clic pe GO pentru
a o utiliza. bebird-**Note5**-323369 ...

lang:es **score:25** filesize: 1.04 M page_count: 12 document date: 2024-04-08

bebird

Otoskop cu cameră de curățare a urechii
Bebird Note 5

Instrucțiuni de utilizare

[\[pdf\]](#)

Bebird Note 5 PL docId 878D79797D09CE5414B3B3B64A49520E Instrucțiuni de utilizare rum otoszkóp
fültisztításhoz kamerával kék 49921 36729 d7be4f4b7702a71e769064cdc22efe4f puputeve hu media
wysiwyg attachments |||

Otoscop cu camer de curare a urechii Bebird Note 5 Instrucuni de utilizare Pasul 1:
Descrcarea ap ... automat. Deschidei aplicaia. Conectai-v dispozitivul la WiFi.
Conectai produsul la WiFi. bebird-**Note5**xxxxxx napoi la aplicaia Bebird Facei clic pe
GO pentru a o utiliza. bebird-**Note5**-323369 ...

lang:it **score:25** filesize: 1.01 M page_count: 13 document date: 2024-04-08

bebird

Otoszkóp fültisztító kamerával Bebird Note 5
pro

Használati utasítás

[\[pdf\]](#)

Bebird Note 5 Pro PL docId 726E70AB3898AE24FE303C5C6E02A039 Használati utasítás hun pro
otoszkóp fültisztításhoz kamerával fehér 49910 36730 7f28deb252988c081067cce064971427 puputeve
hu media wysiwyg attachments |||

Otoszkoftisztító kamerával Bebird Note 5 pro Használati utasítás 1. lps: Az alkalmazás letöltése
iOS: Android ... mázst. Csatlakoztassa a készüléket a WiFi-hálzathoz. Csatlakoztassa a
terméket a WiFi-hálzathoz. bebird**Note5**-xxxxxxx Vissza a bebird alkalmazáshoz
Kattintson a GO gombra a használatához. bebird-**Note5**-...

lang: id **score:24** filesize: 1.13 M page_count: 16 document date: 2024-04-08

bebird

Otoszkóp fültisztító kamerával Bebird Note 5

Használati utasítás

[\[pdf\]](#)

Bebird Note 5 PL docId 3DA78CC09C9E6675FAA41D9C3913EAC9 Használati utasítás hun otoszkóp
fültisztításhoz kamerával kék 49919 36729 5a96927f13a5e3d84b7fe32f68edebda puputeve hu media
wysiwyg attachments |||

Otoszkóftisztító kamerával Bebird Note 5 Használati utasítás 1. lps: Az alkalmazás letöltése
iOS: Android ... mázst. Csatlakoztassa a készüléket a WiFi-hálzathoz. Csatlakoztassa a
terméket a WiFi-hálzathoz. bebird**Note5**-xxxxxxx Vissza a bebird alkalmazáshoz
Kattintson a GO gombra a használatához. bebird-**Note5**-...

lang: no **score:24** filesize: 1.12 M page_count: 16 document date: 2024-04-08



[\[pdf\]](#) User Manual

Note5 Pro 20221212q User manual eng Bebird Note 5 pro otoszkóp fültisztításhoz kamerával fehér
39658 36730 5ab370e3275ccc5c2430d7c607143a93 puputeve hu media wysiwyg attachments |||

...

lang: no **score:19** filesize: 1.55 M page_count: 1 document date: 2023-12-18