

[Manuals.plus](#) /

› [Vive Mobility](#) /

› Vive Mobility Ultra-Lightweight Transport Wheelchair Instruction Manual

Vive Mobility B0DPRBTP18

Vive Mobility Ultra-Lightweight Transport Wheelchair Instruction Manual

Model: B0DPRBTP18

1. PRODUCT OVERVIEW

The Vive Mobility Ultra-Lightweight Transport Wheelchair is designed for easy mobility and transport. It features a compact, foldable design and is suitable for individuals between 4'10" and 5'7" in height, supporting up to 300 pounds.

Key features include:

- **Ultra-Lightweight Design:** Weighs approximately 23 pounds for effortless maneuverability and transport.
- **Foldable for Compact Storage:** Quickly folds for easy storage in vehicles or small spaces.
- **Pre-Assembled:** Ready for immediate use with no assembly required.
- **Sturdy Aluminum Frame:** Provides stability and supports users up to 300 pounds.
- **Narrow 22" Frame:** Allows for easy navigation through tight spaces.
- **Adjustable Armrests:** Facilitate easier transfers for the user.

Lightweight Transport Chair for Small Adults

Weighs only 23 lbs



Image 1.1: The Vive Mobility Ultra-Lightweight Transport Wheelchair folded for compact storage in a car trunk.

2. SETUP AND UNFOLDING

The Vive Mobility Ultra-Lightweight Transport Wheelchair arrives pre-assembled and ready for use. No tools or complicated assembly steps are required.

2.1 Unfolding the Wheelchair

1. Remove the wheelchair from its packaging.
2. Place the folded wheelchair on a flat, stable surface.
3. Gently pull the armrests outwards to begin unfolding the frame.
4. Push down firmly on the seat rails until the seat is fully flat and the frame is locked into its open position. Ensure all locking mechanisms are engaged.
5. Lower the footrests into position.

Folds Quickly & Easily

For easy travel & compact storage



Image 2.1: Demonstrating the quick and easy folding and unfolding process of the wheelchair for travel and storage.

Pre-Assembled Frame

Ready to go out of the



Image 2.2: The wheelchair is pre-assembled and ready for immediate use out of the box.

3. OPERATING INSTRUCTIONS

This transport wheelchair is designed to be pushed by an attendant. It is not intended for self-propulsion.

3.1 Using the Hand Brakes

The wheelchair is equipped with hand brakes located on the push handles. These brakes allow the attendant to control the speed and secure the wheelchair in a stationary position.

- **To slow down or stop:** Squeeze the brake levers gently.
- **To lock the brakes:** Push the brake levers down until they click into the locked position. This will prevent the wheels from rolling. Always engage the locking brakes when the user is transferring in or out of the chair, or when the chair is stationary.
- **To release the brakes:** Pull the brake levers up until they disengage from the locked position.

Sturdy & Long Lasting

Supports up to 300 lbs



Image 3.1: Detail of the locking hand brake mechanism for secure stopping and parking.

3.2 Adjustable Armrests

The armrests can be adjusted or flipped back to facilitate easier transfers for the user into and out of the wheelchair.

- Locate the release mechanism near the base of the armrest.
- Engage the mechanism to allow the armrest to pivot upwards and backwards.
- Return the armrest to its original position and ensure it is securely locked before use.

Adjustable Armrests

Let you get into & exit your wheelchair more easily



Image 3.2: The adjustable armrests allow for easier transfers for the user.

3.3 Smooth-Rolling Wheels

The wheelchair features solid, swivel wheels designed for smooth movement on level surfaces. The front wheels are 6 inches in diameter, and the rear wheels are 8 inches in diameter.

Smooth-Rolling Wheels



Image 3.3: The solid, swivel wheels provide smooth movement. Front wheels are 6 inches, rear wheels are 8 inches.

3.4 Product Video Demonstration

Your browser does not support the video tag.

Video 3.1: An official product video demonstrating the features and use of the Vive Ultra-Lightweight Transport Wheelchair, including its compact design, wide seat, flip-back armrests, and overall ease of use for confident mobility.

4. MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your transport wheelchair.

- **Cleaning:** Wipe down the frame and seat with a damp cloth and mild soap as needed. Avoid abrasive cleaners.
- **Brakes:** Periodically check the hand brakes for proper function. Ensure they engage and disengage smoothly and hold the wheels securely when locked. Adjust if necessary (refer to a qualified technician if unsure).
- **Wheels:** Inspect wheels for wear and tear. Ensure they spin freely and are free of debris.
- **Frame:** Regularly check the aluminum frame for any signs of damage, cracks, or loose fasteners. Tighten any loose screws or bolts.
- **Fabric Components:** Inspect the seat and backrest fabric for tears or excessive wear.

5. TROUBLESHOOTING

If you encounter any issues with your Vive Mobility Ultra-Lightweight Transport Wheelchair, refer to the following common troubleshooting tips:

- **Wheelchair does not roll smoothly:**
 - Check for debris (e.g., hair, lint) wrapped around the wheel axles.
 - Ensure the brakes are fully disengaged.
 - Inspect wheels for damage or excessive wear.
- **Brakes are not holding or engaging properly:**
 - Verify that the brake levers are pushed down completely into the locked position.
 - Check for any obstructions preventing the brake mechanism from fully engaging the wheels.
 - If the brakes feel loose or ineffective, they may require adjustment by a qualified technician. Do not use the wheelchair if brakes are not functioning correctly.
- **Frame feels unstable or twists:**
 - Ensure the wheelchair is fully unfolded and the seat is pushed down firmly to lock the frame.
 - Check all visible bolts and screws for tightness.
 - If the frame appears bent or damaged, discontinue use immediately and contact customer support.
- **Difficulty maneuvering in tight spaces:**
 - This model has a narrow 22" frame designed for tight spaces. However, sharp turns require careful handling by the attendant.
 - Ensure the swivel wheels are functioning correctly.

For issues not resolved by these steps, please contact Vive Mobility customer support.

6. SPECIFICATIONS

Feature	Specification
Model Name	Ultra-Lightweight Transport
ASIN	B0DPRBTP18
Product Dimensions (D x W x H)	30" x 22" x 40"
Product Weight	Approximately 23 pounds
Maximum Weight Recommendation	300 pounds
Ideal User Height	4'10" - 5'7"
Seat Width	18 inches
Frame Material	Aluminum Alloy Steel
Wheel Type	Solid, Swivel (Front: 6", Rear: 8")
Brake Style	Hand brake
Handle Type	Push Handles

Feature	Specification
UPC	810178881253

Comfortable Seat & Backrest

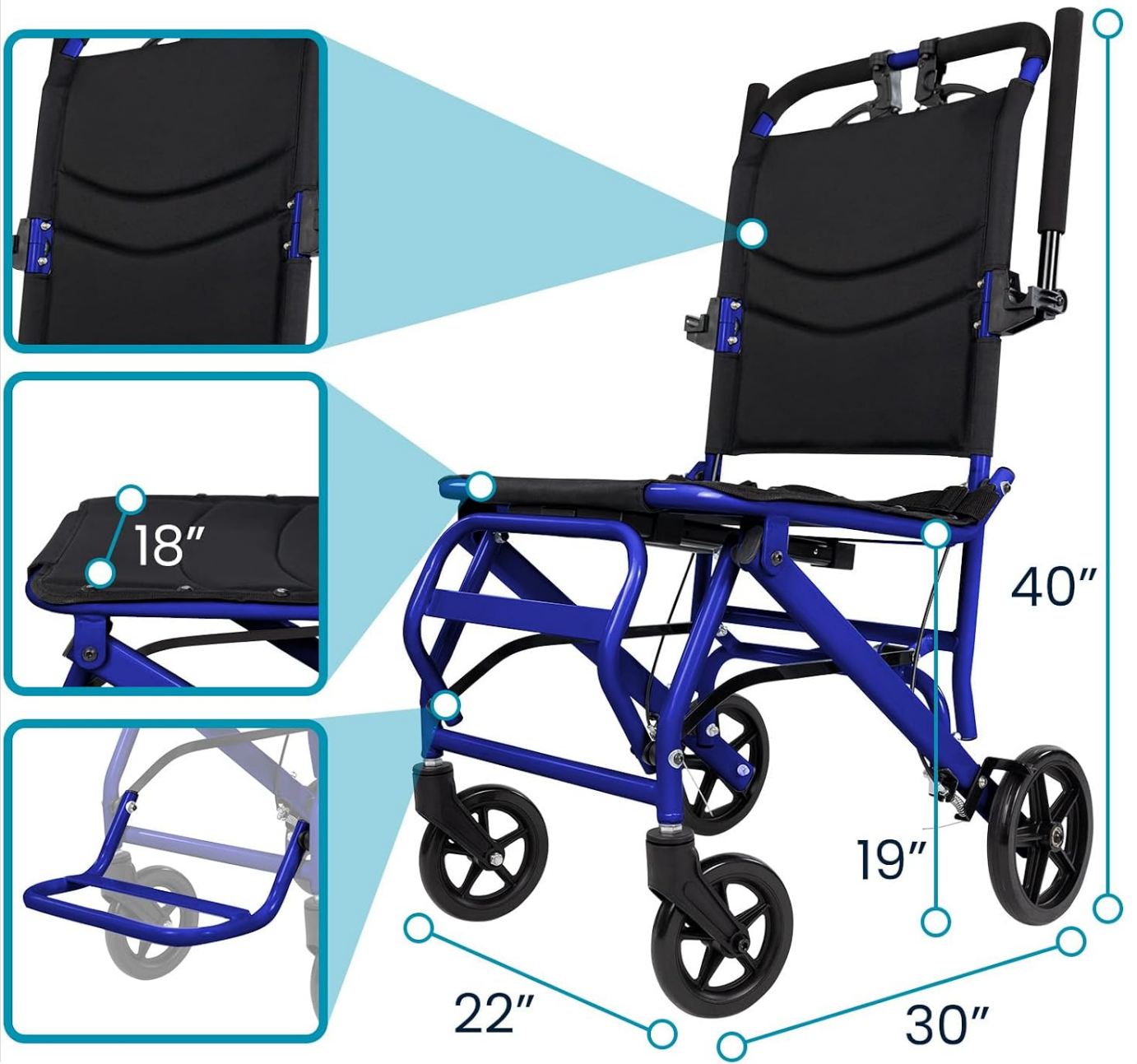


Image 6.1: Detailed dimensions of the wheelchair, including height, width, depth, and seat width.

7. WARRANTY AND SUPPORT

The Vive Mobility Ultra-Lightweight Transport Wheelchair comes with a **1-Year Guarantee**.

For warranty claims, technical assistance, or any questions regarding your product, please contact Vive Mobility customer support directly.

Manufacturer: Vive Health

1-Year Guarantee

Purchase with confidence



Image 7.1: The Vive Health 1-Year Guarantee ensures product quality and customer satisfaction.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question