

UNIWA WPEC-bai

UNIWA Electric Ear Cleaner Instruction Manual

Model: WPEC-bai

Your guide to safe and effective ear wax removal.

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the UNIWA Electric Ear Cleaner. Failure to follow these guidelines may result in injury or damage to the device.

- **Water Temperature:** Always use warm water close to body temperature. Cold water can cause dizziness or discomfort.
- **Ear Tips:** Use only the provided soft, flared ear tips. Ensure they are securely attached before use. Replace tips regularly, especially when used by multiple individuals, to prevent infection. The flared design helps protect the eardrum.
- **Medical Conditions:** Do not use if you have ear infections, a perforated eardrum, ear tubes, recent ear surgery, or any other ear condition without consulting a medical professional.
- **Children:** This device is not intended for use by children without adult supervision.
- **Discomfort:** If you experience pain, severe dizziness, or any unusual discomfort during use, stop immediately and consult a doctor.
- **Device Integrity:** Do not attempt to modify or repair the device. Contact customer support for assistance.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This device is for ear wax removal only.

2. PRODUCT OVERVIEW AND CONTENTS

The UNIWA Electric Ear Cleaner is designed for safe and effective ear wax removal at home. It features an intelligent control system for stable water pressure and a dual-tank design for clean and waste water.

Package Contents:

- Electric Ear Cleaner System (Product Body) x 1

- Replaceable Ear Tips x 2
- Type-C Charging Cable x 1
- Cotton Swabs x 1 pack
- Temperature Sticker x 1
- User Manual x 1



Image: Contents of the UNIWA Electric Ear Cleaner kit, including the main unit, ear tips, charging cable, cotton swabs, temperature sticker, and user manual.

Key Features:

- **4 Water Pressure Modes:** Normal, Gentle, Pulse, Strong.
- **High Capacity Battery:** 2000mAh for extended use.
- **IPX7 Waterproof Rating:** Ensures durability and safety during use.
- **Dual Water Tanks:** 600ml clean water tank and 600ml dirty water tank for continuous cleaning without interruption.
- **One-Handed Operation:** Designed for ease of use.
- **Automatic Waste Water Discharge:** Ensures a mess-free experience.

High Quality

Effortlessly removes
earwax build-up at home



4 Mode
Pressure Settings



2000mAh
Battery



IPX6
Waterproof



600ML
Water Tank



Image: The UNIWA Electric Ear Cleaner unit, illustrating its key features including 4 pressure settings, 2000mAh battery, IPX6 waterproof rating, and 600ml water tank capacity.

3. SETUP INSTRUCTIONS

Follow these steps to prepare your UNIWA Electric Ear Cleaner for use:

1. **Charge the Device:** Connect the provided Type-C charging cable to the device and a suitable power source. Ensure the device is fully charged before first use.
2. **Prepare Water:** Fill the clean water tank (marked 'FT') with warm water. The water should be close to body temperature to prevent dizziness. You can use the included temperature sticker to verify.
3. **Soften Earwax (Optional but Recommended):** If you have hardened earwax, consider using earwax softening drops a few minutes before irrigation, as recommended by a healthcare professional.
4. **Install Cleaning Tip:** Securely attach one of the replaceable ear tips to the device's nozzle. Ensure it clicks into place.

How To Use?



Image: Visual guide showing steps for softening earwax, checking water temperature, installing the cleaning tip, and preparing the device for use.

4. OPERATING INSTRUCTIONS

Once the device is set up, follow these steps for ear cleaning:

1. **Switch On:** Press the power button to turn on the device.
2. **Select Pressure Mode:** Choose one of the four pressure modes (Soft, Normal, Pulse, Strong) using the mode selection button. For first-time users or sensitive ears, start with the 'Soft' mode and gradually increase if needed.
3. **Position the Device:** Gently place the flared ear tip against the entrance of your ear canal. Do not force it deep into the ear. The device is designed for one-handed operation.
4. **Begin Cleaning:** Press the spray button to start the water flow. The device will spray water into your ear and automatically collect the waste water into the dirty water tank (marked 'DT').
5. **Continue as Needed:** Move the device slightly to ensure thorough cleaning. Continue until you feel your ear is clean or the clean water tank is empty.
6. **Dry Your Ear:** After cleaning, gently dry your ear with a clean cotton swab or soft towel.

4 Cleaning Pressure Mode Settings



Image: Diagram illustrating the four cleaning pressure modes: Soft, Normal, Pulse, and Strong, with a note to start with 'Soft' mode.



Image: A user demonstrating the one-handed operation of the UNIWA Electric Ear Cleaner, highlighting its ease of use compared to traditional methods.

5. MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your UNIWA Electric Ear Cleaner:

- **Empty and Clean Tanks:** After each use, empty both the clean and dirty water tanks. Rinse them thoroughly with clean water and allow them to air dry completely.
- **Clean Ear Tips:** Remove the ear tip and rinse it under running water. For deeper cleaning, tips can be soaked in a mild disinfectant solution. Ensure they are completely dry before storage or next use. Replace ear tips as needed.
- **Wipe Down Device:** Wipe the main body of the device with a soft, damp cloth. Do not submerge the main unit in water.
- **Storage:** Store the device in a clean, dry place away from direct sunlight and extreme temperatures. The device can be mounted on a wall or placed on a flat surface.

Two Storage Options Perfectly sized so you can mount it on the wall or place it anywhere you like



Image: Two storage options for the UNIWA Electric Ear Cleaner, showing it mounted on a wall and placed on a countertop.

6. TROUBLESHOOTING

If you encounter issues with your UNIWA Electric Ear Cleaner, refer to the following common solutions:

- **Device Not Turning On:**

- Ensure the device is fully charged. Connect it to the Type-C charging cable and allow sufficient time to charge.

- **No Water Flow or Weak Flow:**

- Check if the clean water tank is filled correctly.
- Ensure the ear tip is securely attached and not clogged.
- Verify that the water intake tube inside the clean water tank is properly positioned.
- Select a higher pressure mode if the current mode is too gentle.

- **Water Leakage:**

- Ensure the ear tip is firmly attached.
- Check that the water tanks are properly seated and sealed.

- **Discomfort or Loud Noise During Use:**

- Ensure the water temperature is warm, not cold.
- Start with the 'Soft' pressure mode.
- Ensure the ear tip is positioned correctly at the ear canal entrance and not pressed too hard.
- If discomfort persists, discontinue use and consult a medical professional.

- **Waste Water Not Collecting:**

- Ensure the dirty water tank is empty and properly installed.
- Check for any blockages in the waste water collection pathway.

If the problem persists after attempting these solutions, please contact UNIWA customer support.

7. SPECIFICATIONS

Model Number	WPEC-bai
Manufacturer	UNIWA
Battery Capacity	2000mAh
Water Tank Capacity	600ml (Clean), 600ml (Dirty)
Water Pressure Modes	4 (Soft, Normal, Pulse, Strong)
Waterproof Rating	IPX7
Charging Interface	USB-C

8. WARRANTY AND SUPPORT

UNIWA is committed to providing quality products and customer satisfaction.

Warranty Information:

For detailed warranty terms and conditions, please refer to the warranty card included in your product packaging or visit the official UNIWA website. Keep your proof of purchase for warranty claims.

Customer Support:


If you have any questions, concerns, or require assistance with your UNIWA Electric Ear Cleaner, please do not hesitate to contact our customer support team. We are dedicated to providing 100% quality product and service.

Contact information can typically be found on the product packaging, the official UNIWA website, or through

your retailer.

Related Documents - WPEC-bai

	<p>UNIWA V909T 4G Flip Phone User Manual</p> <p>User manual for the UNIWA V909T 4G flip phone, covering installation, operation, settings, and specifications. This feature phone supports single Nano SIM and offers features like FM radio, SOS, and speed dial.</p>
	<p>UNIWA F41 Rugged Android Smart Public Network Intercom User Manual</p> <p>Comprehensive user manual for the UNIWA F41, a rugged Android smart public network intercom designed for industries requiring instant communication, featuring PoC technology, GPS, and durable construction.</p>
	<p>UNIWA V909T Manual de Usuario: Guía Completa del Teléfono Plegable</p> <p>Descubre el manual de usuario del UNIWA V909T, un teléfono plegable 2G/3G/4G fácil de usar. Aprende sobre sus características, instalación, llamadas, mensajes, configuraciones y especificaciones técnicas.</p>
	<p>UNIWA F400 User Manual: Smart Handheld Intercom PDA</p> <p>Comprehensive user manual for the UNIWA F400, an IP65-rated Android smart handheld intercom PDA. Features include PTT communication, GPS, durable design, 4350mAh battery, and Android 15 OS. Details on setup, functions, maintenance, and specifications. Document metadata: Paper: 100G coated paper, Printing: Double-sided printing, four-color printing, Folding method: Spiral fold, Size: 87*150mm.</p>
	

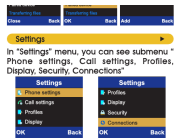

- 07

► Phone settings
Time& date, language settings, shortcut settings, auto power on/off, restore factory settings.
In "Phone settings" menu, you can restore the phone to factory settings. **Please backup your data before restore factory setting, the data will be lost after restore factory. The default restore factory password is "0000"**

► Call settings
Call divert, Call waiting, Call barred, Hide ID, others. See figure on the right:



► Profiles
Normal, Silent, Meeting, Indoor, Outdoor. If you want to change "Profiles" setting, press **[M]** key to enter "Options" menu, and then press **[M]** key to select "Customize". You will see "Rename, Ring settings, Adjust volume, Ring type, Other alert rings". See figure below:

- 08



You can also press the volume +/- keys on the left side of the phone to adjust volume.

► Display
Animation effect, Wallpaper settings, Screen saver settings, Contrast, Keypad backlight time, Auto keypad lock.

► Security
PIN Modify PIN2, Phone locked, Modify the cellphone passwords, Privacy, Fixed Dialing number, blacklist.
In order to avoid forgetting password and lock the phone forever, we do not recommend user to set password. Please set and keep password carefully if you use this function.

► Connections
Network account, Data service, Data Roaming, Network selection, Network type, VoLTE.


- 09

► Speed Dial
For some frequent contacts, "Speed dial" can dial them quickly. Key 2-9 can be set as speed dial keys. See figure below:



After the frequent contact is set, press the corresponding key for 5 seconds, the phone will call the corresponding person automatically.

► SOS Settings
When user is in an emergency, helpline can press **[M]** button on the back of the phone to call and send SMS to pre-set emergency contact by pressing **[M]** button for 5 seconds. You can pre-set 5 contacts and 1 SOS message. See figure below:

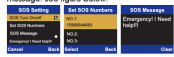

- 10

► Torch light
You can slide up/down key on left side of the phone to turn on/off torch light.



Product specifications

Model No.	UNWA V1000
Dimensions	122.5*56.3*15.1mm
Weight	114g (include battery)
Material	Plastic
Style	Bar
Service	Unlocked
SIM	Single Nano SIM card single standby
Band	2G: GSM 850/900/1800/1900MHz 3G: WCDMA 850/900/1700/1900/2100MHz 4G: TD-LTE B3/B4/B5/B7/B12/B13/B17/B44

- 11

Phonebook Number	1000
LED Indicator	RGB
Message	SMS/MMS
USB Port	Micro USB
Earphone Port	3.5mm audio jack
Speaker	2030 BOX loudspeaker
Screen	2.31 Inch QVGA TFT display screen, 320*240px
Camera	Single camera, 0.3MP back camera
Audio formats	MP3
Video formats	MP4 etc.
Radio	Wireless radio
SOS Function	Support preset 5 SOS emergency contacts and 1 emergency SMS
Torch	One key to turn on/off torch

- 12

BT	V3.0
Battery	Removable 1700mAh Li-Ion battery
Language	English, Arabic, French, Hindi, Indonesian, Malay, Portuguese, Russian, Spanish, Thai, Vietnamese, Italian, Turkish, Deutsch, Greek, Hebrew, Czech, Polish, Swedish, Ukrainian, Dutch, Simplified/Traditional Chinese.
Accessories	1* Battery 1* Charger 1* Charging dock 1* User manual



- 13



[UNIWA V202T 4G Flip Feature Phone User Manual](#)

User manual for the UNIWA V202T 4G flip feature phone, covering unboxing, basic functions, contact management, settings, SOS emergency calls, product specifications, and technical support.