



Manuals.plus /

- › VEVOR /
- › VEVOR Power Lift Recliner Chair Instruction Manual (Model SX-81541-2)

## VEVOR SX-81541-2

# VEVOR Power Lift Recliner Chair Instruction Manual

Model: SX-81541-2

## INTRODUCTION

Thank you for choosing the VEVOR Power Lift Recliner Chair. This manual provides essential information for the safe and efficient operation, assembly, and maintenance of your recliner. Please read these instructions thoroughly before use and retain them for future reference.

## IMPORTANT SAFETY INFORMATION

To reduce the risk of injury, electric shock, fire, or damage to the product, always follow basic safety precautions.

- **Weight Capacity:** Do not exceed the maximum weight capacity of **350 lbs (159 kg)**.
- **Stable Surface:** Ensure the chair is placed on a level and stable surface before operation.
- **Clearance:** Maintain adequate clearance around the chair to prevent obstructions during reclining or lifting.
- **Children and Pets:** Keep children and pets away from the chair's moving parts.
- **Power Cord:** Do not place the power cord under the chair or in a position where it could be pinched or damaged.
- **Anti-Fall Ropes:** Ensure the anti-fall ropes are securely attached and within reach for emergency use.
- **Liquid Contact:** Avoid spilling liquids on the chair, especially near electrical components.
- **Disassembly:** Do not attempt to disassemble or repair the chair yourself. Contact qualified service personnel.

## PACKAGE CONTENTS

Upon unpacking, verify that all components are present and undamaged:

- Recliner Base
- Backrest
- Armrests (Left and Right)
- Power Adapter and Cord

- Remote Control
- Anti-Fall Ropes (2)
- User Manual

If any parts are missing or damaged, do not proceed with assembly. Contact VEVOR customer support.

## ASSEMBLY AND SETUP

Assembly typically involves connecting the main sections of the chair. Follow these general steps:

1. **Unpack Components:** Carefully remove all parts from the packaging. Place them on a soft, clean surface to prevent scratches.
2. **Attach Backrest:** Slide the backrest onto the metal posts located on the recliner base until it clicks securely into place.
3. **Connect Armrests:** Attach the left and right armrests to the sides of the chair, ensuring they are firmly seated.
4. **Connect Wiring:** Locate the electrical connectors from the backrest and armrests and connect them to the corresponding ports on the recliner base. Ensure all connections are secure.
5. **Connect Power:** Plug the power adapter into the chair's power input port, then plug the power cord into a standard electrical outlet.
6. **Attach Anti-Fall Ropes:** Securely attach the two anti-fall ropes to the designated points on the chair, ensuring they are easily accessible.
7. **Initial Test:** Before sitting, test the chair's functions using the remote control to ensure proper operation.

# Massage & Heat Function

Easy control with the remote



A full front view of the VEVOR Power Lift Recliner Chair in dark gray, showcasing its comfortable design.

## OPERATING INSTRUCTIONS

Your VEVOR Power Lift Recliner Chair is equipped with a remote control to manage its various functions.

### Power Lift and Recline Functions

The chair features dual OKIN motors, allowing independent adjustment of the backrest and footrest. Use the remote control buttons to achieve desired positions:

- **Up/Down Arrows (Backrest):** Adjust the angle of the backrest.
- **Up/Down Arrows (Footrest):** Adjust the height and extension of the footrest.
- **Lift Function:** Press the designated lift button to gently tilt the chair forward, assisting you in standing up.
- **Recline Function:** Press the recline button to move the chair into a reclined position.
- **Home Button:** Press the 'Home' button to return the chair to its upright, default position with one click.

The chair supports multiple angles for various activities:

- **Assisted Standing:** Approximately -45° tilt.
- **Reading:** Approximately 90° upright position.
- **Watching TV/Relaxation:** Approximately 165° recline.
- **Full Recline/Sleep:** Up to 180° zero-gravity position.

# Comfortable Experience



## Chenille Fabric

Breathable and  
Comfortable



## 25D High-Density Sponge

Soft Yet Supportive



## Footrest

Long and Comfortable

A diagram demonstrating the independent backrest and footrest adjustment powered by dual OKIN motors, showing various reclining positions from upright to fully reclined.


# Convenient And Thoughtful Design

**2 Remote Control Pockets**

**2 Anti-Fall Ropes**

**2 Side Pockets**

**Anti-Tip Base**

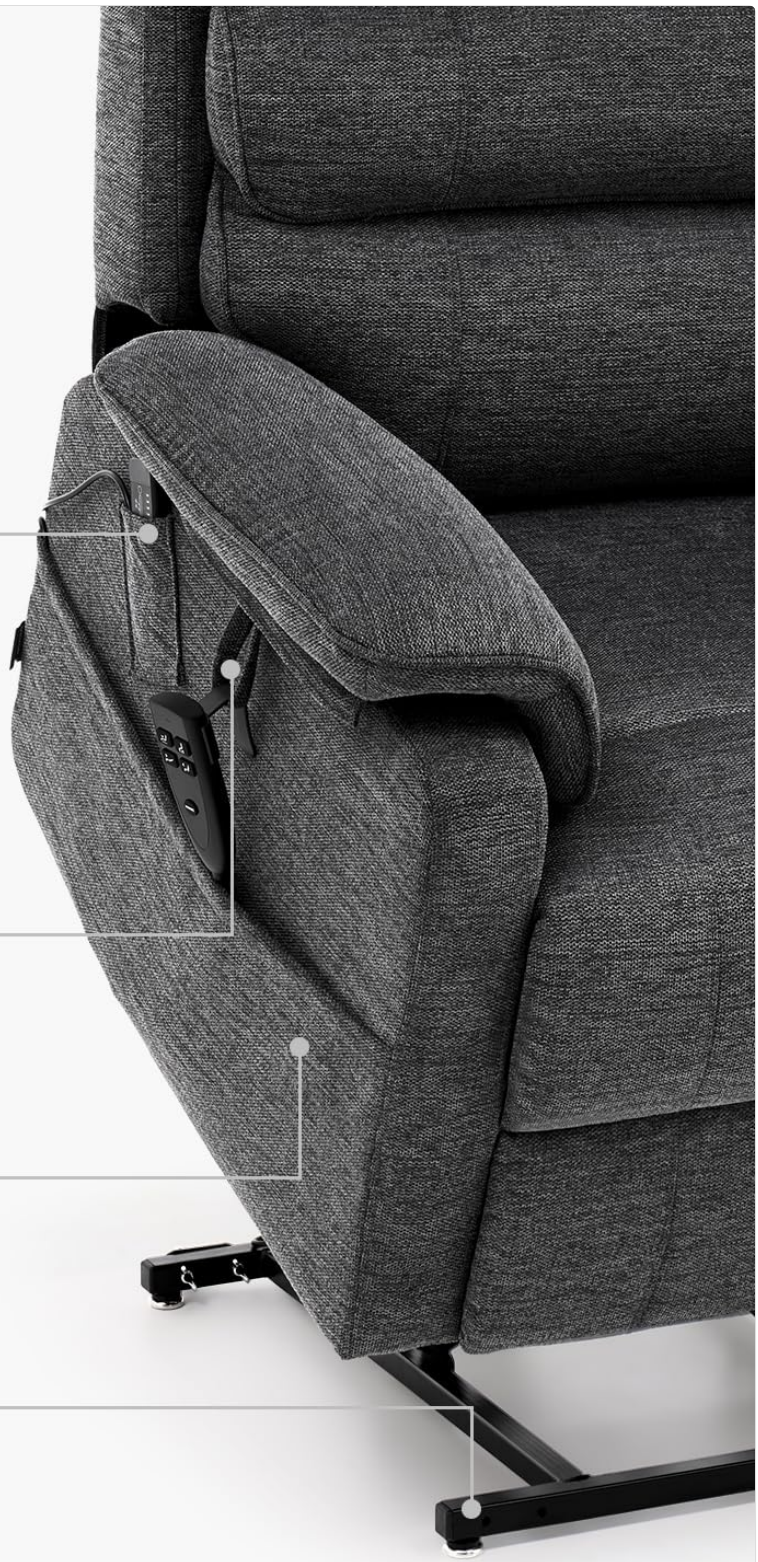

An image depicting the power lift function of the recliner, which gently tilts forward to assist users, particularly seniors, in standing up easily.

## Massage and Heat Functions

The recliner features multiple massage points and an independent lumbar heating function.

- **Massage Points:** Back, lumbar, thighs, and calves.
- **Massage Modes:** Select from 8 different massage modes using the remote.
- **Intensity Levels:** Adjust massage intensity with 5 available levels.

- **Lumbar Heating:** Activate the independent lumbar heating function for added comfort.
- **Auto Shut-off:** Both massage and heating functions automatically turn off after 15 minutes for safety and energy efficiency.

## Power Lift Assistance

Assisting senior or pregnant women to stand up



An illustration highlighting the massage points (back, lumbar, thighs, calves) and independent lumbar heating function of the recliner, along with its remote control for easy operation.

### Convenient Features

- **Side Pockets:** Use the large side pockets for storing remotes, books, or other small items.

- **Remote Control Pouch:** A dedicated pouch keeps your remote control secure and easily accessible.
- **Anti-Fall Ropes:** These ropes provide additional support and stability when using the lift function.



An image highlighting the practical design elements of the recliner, including integrated remote control pockets, anti-fall ropes for safety, spacious side pockets, and a stable anti-tip base.

## MAINTENANCE

Proper care and maintenance will extend the life of your recliner.

- **Cleaning Fabric:** For general cleaning, use a soft, damp cloth with mild soap. Avoid harsh chemicals or abrasive cleaners. Test any cleaning solution on an inconspicuous area first.
- **Vacuuming:** Regularly vacuum the chair to remove dust and debris.
- **Power Cord:** Periodically inspect the power cord for any signs of damage. If damaged, discontinue use and contact customer support.
- **Moving the Chair:** When moving the chair, lift it rather than dragging it to prevent damage to the base or flooring.

Ensure all electrical connections are disconnected before moving.



5'1"-5'9" / 1.55-1.75 m

Suitable Height Range



350 lbs / 159 kg

Weight Capacity



Detailed views showcasing the chair's comfortable chenille fabric, the supportive high-density foam, and the extended footrest for optimal relaxation.

## TROUBLESHOOTING


If you encounter issues with your recliner, refer to the following common problems and solutions:

<b>Problem</b>	<b>Possible Cause</b>	<b>Solution</b>
Chair does not respond to remote.	No power; loose connection; remote issue.	Ensure power cord is securely plugged into outlet and chair. Check all internal wiring connections. Replace remote batteries if applicable (this model uses a wired remote, so check for damage).
Motor makes noise but chair does not move.	Obstruction in mechanism; motor malfunction.	Check for any objects obstructing the chair's moving parts. If no obstruction, contact customer support.
Massage or heat function not working.	Loose connection; function not activated; auto shut-off.	Ensure all electrical connections are secure. Activate the function using the remote. Note that these functions have a 15-minute auto shut-off.
Chair is unstable or wobbles.	Uneven floor; loose assembly.	Ensure the chair is on a level surface. Check all assembly connections for tightness.

If the problem persists after attempting these solutions, please contact VEVOR customer support for assistance.

## **SPECIFICATIONS**

Feature	Detail
Brand	VEVOR
Model Number	SX-81541-2
Color	Medium- Dark Gray
Product Dimensions (D x W x H)	38.4" x 32.7" x 42.3" (975mm x 830mm x 1075mm)
Suitable Height Range	5'1" - 5'9" (1.55 - 1.75 m)
Maximum Weight Capacity	350 lbs (159 kg)
Net Weight	99.2 lbs (45 kg)
Material	Flannel, high-density sponge, metal frame, engineering board, foam
Motor Type	Dual OKIN Motors
Motor Input	AC100V-240V 50/60Hz
Motor Output	DC29V 1.8A
Massage & Heating Input	AC110V-240V 50/60Hz
Massage & Heating Output	DC12V 1A
Massage Modes	8
Massage Intensity Levels	5
Heating Function	Independent Lumbar Heating
Auto Shut-off	15 minutes (Massage & Heat)

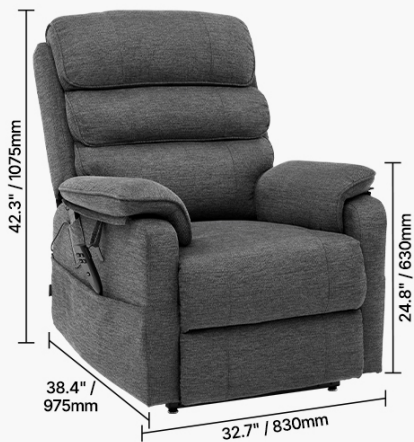

Item Model Number: **SX-81541-1(US), SX-81541-2(EU)**

Motor Input: **AC100V-240V 50/60Hz**

Motor Output: **DC29V 1.8A**

Massage & Heating Input: **AC110V-240V 50/60Hz**

Massage & Heating Output: **DC12V 1A**

Weight Capacity: **350 lbs / 159 kg**

Net Weight: **99.2 lbs / 45 kg**

Product Dimensions: **D38.4\"**

A graphic displaying the item model number (SX-81541-2), motor input/output, massage and heating input/output, weight capacity, net weight, and product dimensions for the VEVOR Power Lift Recliner Chair.



A diagram illustrating the dimensions of the VEVOR Power Lift Recliner Chair, including height, width, depth, and suitable user height range (5'1"-5'9") and weight capacity (350 lbs).

## WARRANTY AND SUPPORT

VEVOR products are designed for durability and performance. For specific warranty details regarding your Power Lift Recliner Chair, please refer to the warranty card included with your purchase or visit the official VEVOR website. If you require technical assistance, have questions about assembly, operation, or need to report a defect, please contact VEVOR customer support. Have your model number (SX-81541-2) and purchase information ready when contacting support.

**VEVOR Customer Service:** Please refer to your purchase documentation or the VEVOR official website for contact information.