

manuals.plus /

› [Frida Baby](#) /

› [Frida Baby Electric NoseFrida Nasal Aspirator User Manual](#)

## Frida Baby 100001971

# Frida Baby Electric NoseFrida Nasal Aspirator

Model: 100001971

### INTRODUCTION

The Frida Baby Electric NoseFrida Nasal Aspirator is designed to provide gentle yet effective relief from nasal congestion for babies and toddlers. Its user-friendly design and multiple suction levels ensure a comfortable experience for your child, helping them breathe easier.

# frida baby® Electric NoseFrida®

The speedy, self-driving SnotSucker®  
that banishes boogers and snot



1 Nasal Aspirator Set + 1 Case



*The Frida Baby Electric NoseFrida Nasal Aspirator, shown with its packaging and included silicone tips. This device is designed for efficient and gentle nasal congestion relief in infants and toddlers.*

## WHAT'S INCLUDED

Your Frida Baby Electric NoseFrida package includes:

- 1 Electric NoseFrida Nasal Aspirator unit
- 2 Silicone Tips (different sizes for baby and toddler)
- 1 Cleaning Brush
- 1 USB Rechargeable Charging Cable (USB-C)
- 1 Storage Case for convenient travel and organization

# Fully-loaded boogie buster



*All essential components of the Electric NoseFrida are included for immediate use and easy maintenance. This image displays the main unit, two silicone tips, a cleaning brush, and the USB-C charging cable.*

## SETUP

1. **Charge the Device:** Before first use, ensure the device is fully charged using the provided USB-C charging cable. The charging port is located on the side of the unit.
2. **Attach Silicone Tip:** Select the appropriate silicone tip size for your child's age and comfort. Gently push the chosen silicone tip onto the top of the aspirator unit until it is securely in place.

## OPERATING INSTRUCTIONS

1. **Power On:** Press the power button on the device to turn it on. It will typically start on the lowest suction setting.
2. **Select Suction Level:** The aspirator features 3 power levels. Press the '+' button to cycle through the suction levels (low, medium, high) to find the most effective and comfortable setting for your child's congestion.



3. **Optional Distraction Light:** Press the light button to activate the color-changing light. This can serve as a distraction to help calm restless children during use.
4. **Aspirate:** Gently insert the silicone tip into your child's nostril. Ensure a snug fit to create effective suction. Move the tip as needed to clear mucus. Repeat for the other nostril.
5. **Power Off:** Press the power button again to turn off the device when finished.

## Self driving suction for instant congestion relief



*The Electric NoseFrida provides self-driving suction for instant congestion relief. Gently place the silicone tip into the nostril to clear mucus effectively.*

# Power through sick days



*The device features an optional color-changing distraction light and three distinct suction levels, indicated by LED lights, allowing for customized and comfortable use.*

## MAINTENANCE & CLEANING

Regular cleaning is essential to maintain hygiene and performance.

1. **Disassemble:** Remove the silicone tip and the clear collection cup from the main unit.
2. **Clean Components:** Wash the silicone tip and collection cup with warm soapy water. The parts are designed to come apart for thorough cleaning and are dishwasher safe.
3. **Use Cleaning Brush:** Utilize the provided cleaning brush to ensure all mucus is removed from the narrow passages of the tip and cup.
4. **Wipe Main Unit:** Wipe the main electric unit with a damp cloth. **Do not submerge the main unit in water.**
5. **Dry & Reassemble:** Allow all parts to air dry completely before reassembling the aspirator.

# This



VS

# That



Comes apart for cleaning

Gets gross + moldy

Dishwasher safe parts

Hides nasty buildup inside

*Unlike traditional bulb syringes that can harbor hidden buildup, the Frida Baby Electric NoseFrida features parts that easily detach for thorough cleaning, ensuring hygiene and preventing mold.*

## TROUBLESHOOTING

If you encounter any issues with your Frida Baby Electric NoseFrida, please ensure the device is fully charged and all components are correctly assembled and clean. Refer to the cleaning instructions to ensure no blockages are present. If problems persist, contact customer support.

## SPECIFICATIONS

Feature	Detail
Model Name	Electric NoseFrida
Item Model Number	100001971






Brand	Frida Baby
Age Range	Baby
Power Source	Battery Powered (USB Rechargeable)
Material Type	Silicone
Color	White
Item Weight	13.4 ounces
Package Dimensions	8.74 x 4.02 x 2.48 inches
Product Care Instructions	Hand Wash Only (for removable parts)

## WARRANTY AND SUPPORT

For warranty information, product registration, or further assistance, please visit the official Frida Baby website or contact their customer support directly. Keep your purchase receipt for any warranty claims.

© 2025 Frida Baby. All rights reserved.


## Related Documents - 100001971

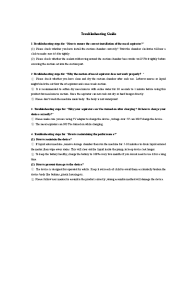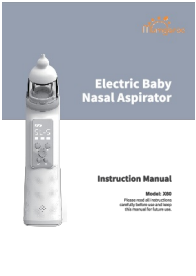	<p><a href="#">Troubleshooting Guide for Baby Nasal Aspirator</a></p> <p>Find solutions to common issues with your baby nasal aspirator, including installation, suction problems, power issues, and maintenance. Learn how to keep your device performing optimally.</p>
	<p><a href="#">Mumgaroo Electric Baby Nasal Aspirator Instruction Manual</a></p> <p>This instruction manual provides comprehensive guidance for the Mumgaroo Electric Baby Nasal Aspirator Model X60, covering its features, assembly, operation, and cleaning for safe and effective infant nasal care.</p>
	<p><a href="#">PLQSX Electric Nasal Aspirator UB-007 User Guide</a></p> <p>Comprehensive guide for the PLQSX Electric Nasal Aspirator (Model UB-007), detailing product information, component descriptions, usage steps, function explanations, important tips, and warranty information for safe and effective infant nasal care.</p>

Product instruction

BC024

© 2017 G.F. 11/16/2017

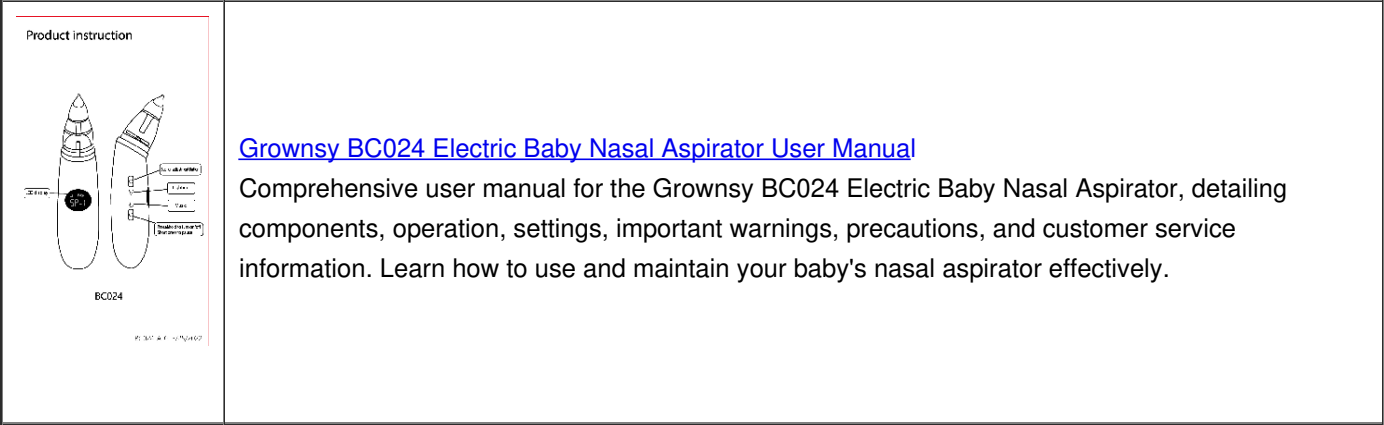

Product instruction


BC024


© 2017 G.F. 11/16/2017

Product instruction

BC024

© 2017 G.F. 11/16/2017

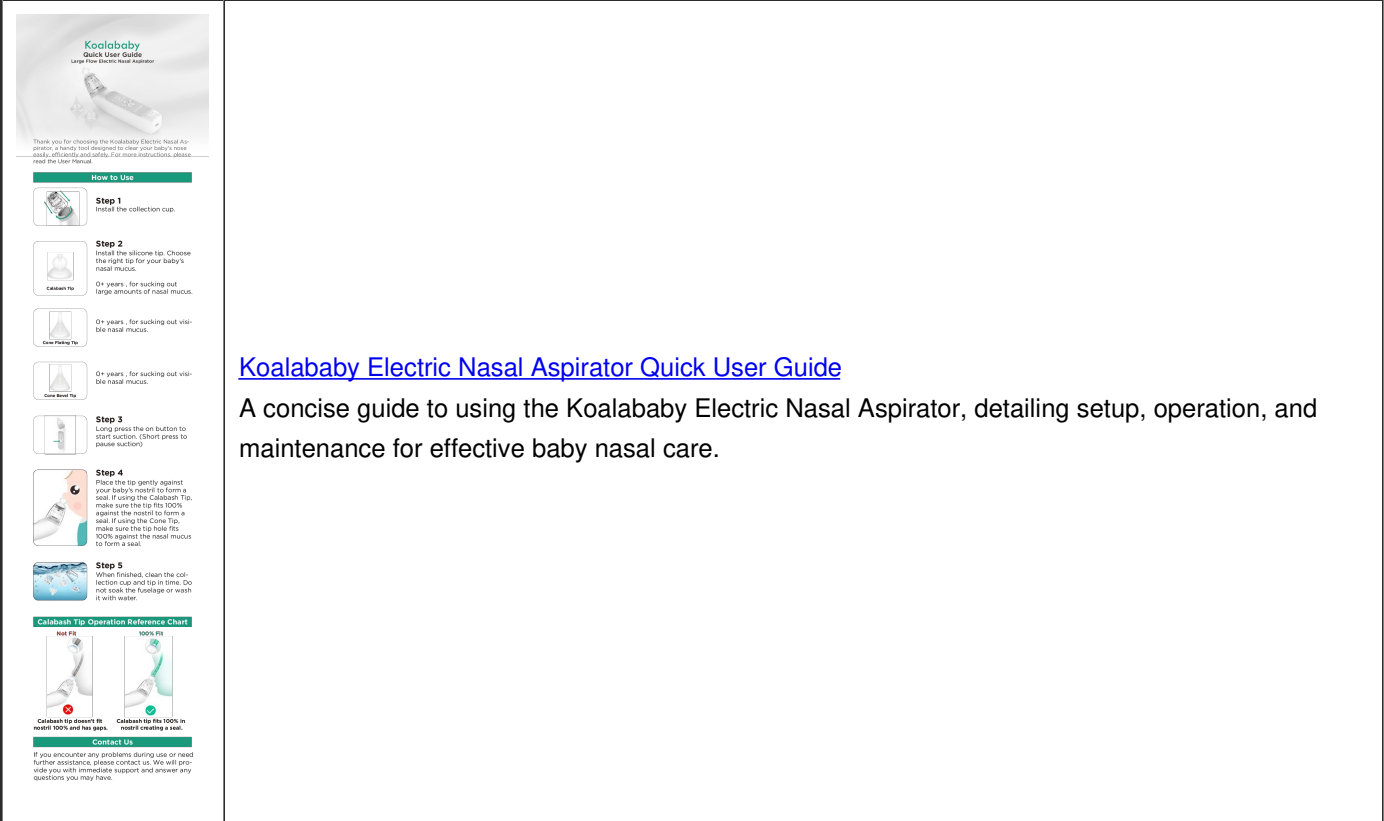





**Koalababy**  
Quick User Guide  
Large Tip Electric Nasal Aspirator

Thank you for choosing the Koalababy Electric Nasal Aspirator, a handy tool designed to clear your baby's nose safely, efficiently and safely. For more instructions, please read the User Manual.

### How to Use



**Step 1**  
Insert the collection cup.




**Step 2**  
Install the silicone tip. Choose the right tip for your baby's nasal mucus.


**Calabash Tip**  
0+ years - for sucking out large amounts of nasal mucus.

**Cone Picking Tip**  
0+ years - for sucking out visible nasal mucus.


**Cone Round Tip**  
0+ years - for sucking out visible nasal mucus.



**Step 3**  
Long press the on button to start suction. (Short press to pause suction)



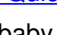

**Step 4**  
Place the tip gently against your baby's nostril to form a seal. If using the Calabash Tip, make sure the tip fits 100% against the nostril to form a seal. If using the Cone Tip, make sure the tip hole fits 100% against the nasal mucus to form a seal.



**Step 5**  
When finished, clean the collection cup and tip in time. Do not cook the Koalababy or wash it with water.

### Calabash Tip Operation Reference Chart


Not Fit	100% Fit
 <p>Calabash tip doesn't fit nostril 100% and has gaps.</p>	 <p>Calabash tip fits 100% in nostril creating a seal.</p>

### Contact Us

If you encounter any problems during use or need further assistance, please contact us. We will provide you with immediate support and answer any questions you may have.

## Koalababy Electric Nasal Aspirator Quick User Guide


A concise guide to using the Koalababy Electric Nasal Aspirator, detailing setup, operation, and maintenance for effective baby nasal care.




**Koalababy**  
Quick User Guide  
Large Tip Electric Nasal Aspirator

Thank you for choosing the Koalababy Electric Nasal Aspirator, a handy tool designed to clear your baby's nose safely, efficiently and safely. For more instructions, please read the User Manual.

### How to Use



**Step 1**  
Insert the collection cup.




**Step 2**  
Install the silicone tip. Choose the right tip for your baby's nasal mucus.


**Calabash Tip**  
0+ years - for sucking out large amounts of nasal mucus.

**Cone Picking Tip**  
0+ years - for sucking out visible nasal mucus.


**Cone Round Tip**  
0+ years - for sucking out visible nasal mucus.



**Step 3**  
Long press the on button to start suction. (Short press to pause suction)



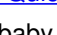

**Step 4**  
Place the tip gently against your baby's nostril to form a seal. If using the Calabash Tip, make sure the tip fits 100% against the nostril to form a seal. If using the Cone Tip, make sure the tip hole fits 100% against the nasal mucus to form a seal.



**Step 5**  
When finished, clean the collection cup and tip in time. Do not cook the Koalababy or wash it with water.

### Calabash Tip Operation Reference Chart


Not Fit	100% Fit
 <p>Calabash tip doesn't fit nostril 100% and has gaps.</p>	 <p>Calabash tip fits 100% in nostril creating a seal.</p>

### Contact Us

If you encounter any problems during use or need further assistance, please contact us. We will provide you with immediate support and answer any questions you may have.

## Koalababy Electric Nasal Aspirator Quick User Guide

A concise guide to using the Koalababy Electric Nasal Aspirator, detailing setup, operation, and maintenance for effective baby nasal care.

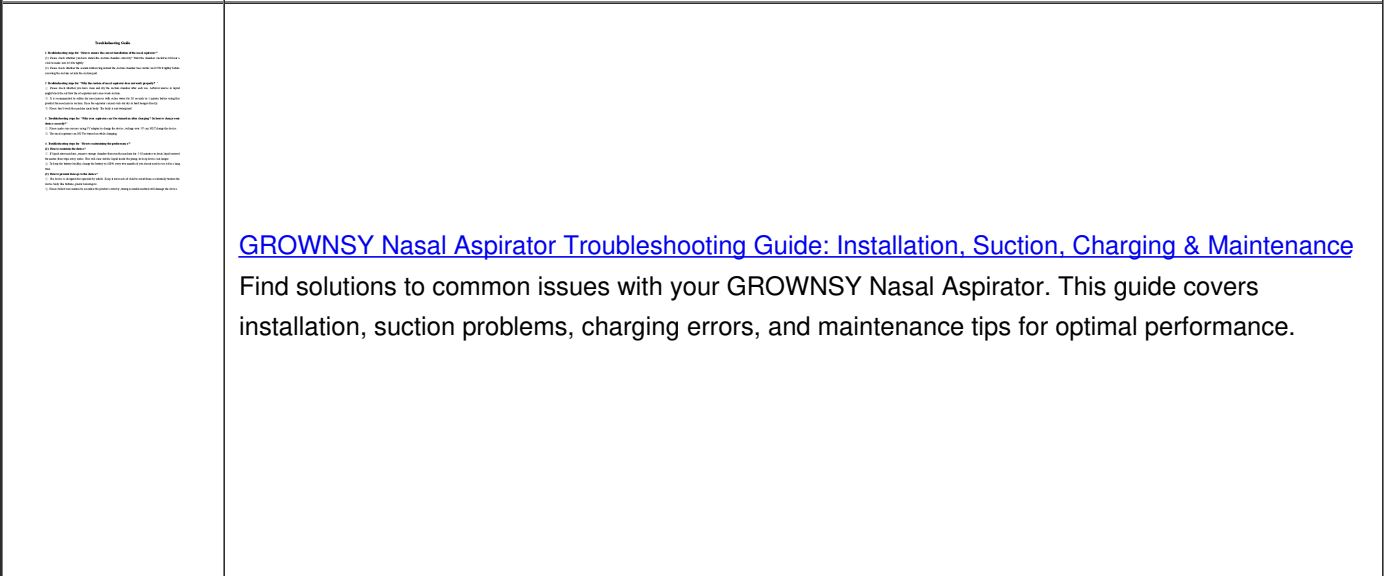

<p><b>Troubleshooting Guide</b></p> <p><b>1. Installation Issues:</b> Ensure the device is correctly assembled. Check the suction tube for leaks and ensure the motor is properly secured.</p> <p><b>2. Suction Problems:</b> Verify the suction tube is not blocked or kinked. Check the patient's airway for obstructions. Ensure the device is properly sealed to the patient's face.</p> <p><b>3. Charging Errors:</b> Check the battery level. Ensure the device is fully charged before use. Verify the charging cable is properly connected.</p> <p><b>4. Maintenance Tips:</b> Clean the device after each use. Replace the suction canister when full. Store the device in a clean, dry place.</p>	<h2><u>GROWNSY Nasal Aspirator Troubleshooting Guide: Installation, Suction, Charging &amp; Maintenance</u></h2> <p>Find solutions to common issues with your GROWNSY Nasal Aspirator. This guide covers installation, suction problems, charging errors, and maintenance tips for optimal performance.</p>
---	---

<p><b>Troubleshooting Guide</b></p> <p><b>1. Installation Issues:</b> Ensure the device is correctly assembled. Check the suction tube for leaks and ensure the motor is properly secured.</p> <p><b>2. Suction Problems:</b> Verify the suction tube is not blocked or kinked. Check the patient's airway for obstructions. Ensure the device is properly sealed to the patient's face.</p> <p><b>3. Charging Errors:</b> Check the battery level. Ensure the device is fully charged before use. Verify the charging cable is properly connected.</p> <p><b>4. Maintenance Tips:</b> Clean the device after each use. Replace the suction canister when full. Store the device in a clean, dry place.</p>	<h2><u>GROWNSY Nasal Aspirator Troubleshooting Guide: Installation, Suction, Charging &amp; Maintenance</u></h2> <p>Find solutions to common issues with your GROWNSY Nasal Aspirator. This guide covers installation, suction problems, charging errors, and maintenance tips for optimal performance.</p>
---	---