

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [PAISEEC](#) /

› [PAISEEC R1 Foldable Rollator Walker with Seat Instruction Manual](#)

PAISEEC R1


PAISEEC R1 Foldable Rollator Walker with Seat Instruction Manual

1. PRODUCT OVERVIEW

The PAISEEC R1 Foldable Rollator Walker with Seat is designed to provide mobility assistance and comfort. It functions as both a rollator walker and a transport chair, featuring adjustable components and a sturdy aluminum frame.

2-IN-1 Rollator

Capacity up to 300lbs

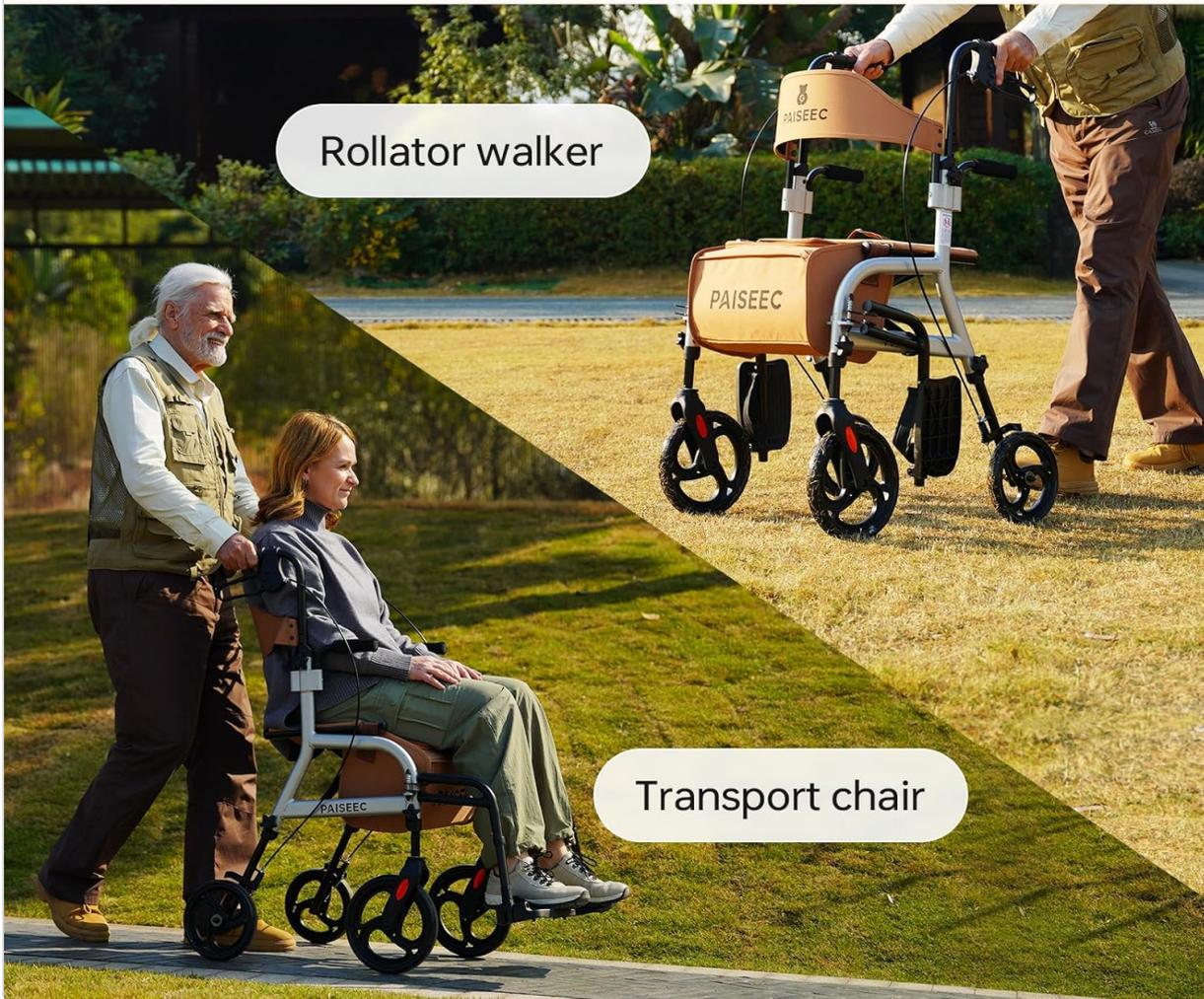

Rollator walker

Transport chair

Image: The PAISEEC R1 Rollator Walker demonstrating its dual functionality as a walker and a transport chair.

2. SAFETY INFORMATION

- Before each use, ensure all screws and knobs are fully tightened.
- If any parts feel loose or make unusual noises after proper installation, contact customer support immediately.
- Always engage the parking brakes when sitting on the seat or when the unit is stationary to prevent unintended movement.
- Do not exceed the maximum weight capacity of 300 lbs.
- When using as a transport chair, ensure the user's feet are properly placed on the footrests.

Great Control and Enhanced Safety



Brake- Provides resistance
Safer for going downhill



Lock- Use when sitting or
to prevent rolling



Image: Detailed view of the dual braking system, illustrating normal brake and parking brake functions for enhanced safety.

3. PACKAGE CONTENTS

The PAISEEC R1 Rollator Walker package includes the following components:

- Main frame with pre-attached seat and backrest
- Front wheels (2)
- Rear wheels (2)
- Adjustable handles with braking system (2)
- Footrests (2)
- Detachable 20L storage bag
- Assembly tools (if applicable)
- Instruction manual

Your browser does not support the video tag.

Video: An official PAISEEC video demonstrating the unboxing and assembly process of the R1 Rollator Walker.

4. SETUP AND ASSEMBLY

Follow these steps to assemble your PAISEEC R1 Rollator Walker. Assembly typically takes approximately 10 minutes.

1. Unfold the main frame of the rollator until it locks into place.
2. Insert the front wheels into the designated slots on the front legs. Secure them by tightening the provided knobs.
3. Insert the rear wheels into the designated slots on the rear legs. Secure them by tightening the provided knobs.
4. Attach the adjustable handles to the main frame, ensuring the braking cables are correctly routed. Secure with the provided fasteners.
5. Attach the backrest to the designated points on the handles.
6. Attach the detachable storage bag under the seat using the straps.
7. For transport chair mode, attach the footrests to the front legs by aligning the pins and securing them.

Your browser does not support the video tag.

Video: An official PAISEEC video providing a step-by-step guide for the installation and assembly of the R1 Rollator Walker.

5. OPERATING INSTRUCTIONS

5.1. Walker Mode

To use the PAISEEC R1 as a walker, stand behind the unit and grasp the adjustable handles. Push the walker forward as you walk, maintaining a comfortable and stable pace. The 10-inch front wheels and 8-inch rear wheels provide smooth movement across various terrains.

Larger All-Terrain Wheels

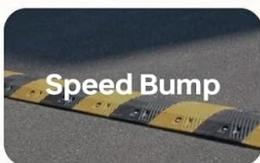
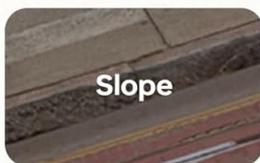


Image: The large all-terrain wheels of the PAISEEC R1, designed for smooth navigation on various indoor and outdoor surfaces.

5.2. Transport Chair Mode

To convert the unit into a transport chair, ensure the footrests are attached and extended. The user can then sit on the padded seat, and a caregiver can push the chair using the handles. The backrest provides support for the user.

Your browser does not support the video tag.

Video: An official PAISEEC video tutorial demonstrating the usage of the R1 Rollator Walker, including conversion to transport chair mode.

5.3. Braking System

The PAISEEC R1 features a dual braking system:

- **Normal Brake:** Squeeze the brake levers on the handles to slow down or stop the rollator.
- **Parking Brake:** Push the brake levers down until they click into place to engage the parking brake. This locks the wheels, providing stability when sitting or when the unit is stationary. Pull the levers up to release.

5.4. Adjustments

The handles and seat height are adjustable to accommodate users from 5'1" to 6'7". To adjust, locate the adjustment pins or knobs on the handle posts and seat frame, press or loosen, set to the desired height, and re-secure.



Image: Illustration of the adjustable handle and seat heights, designed for ergonomic comfort for a wide range of users.

5.5. Folding and Storage


To fold the rollator, lift the seat-mounted pull strap upwards. The unit will collapse for compact storage. Its lightweight design (22 lbs) makes it easy to carry and store in car trunks or tight spaces.

Easy to Fold

Compact Size for Storage



One Hand Folding

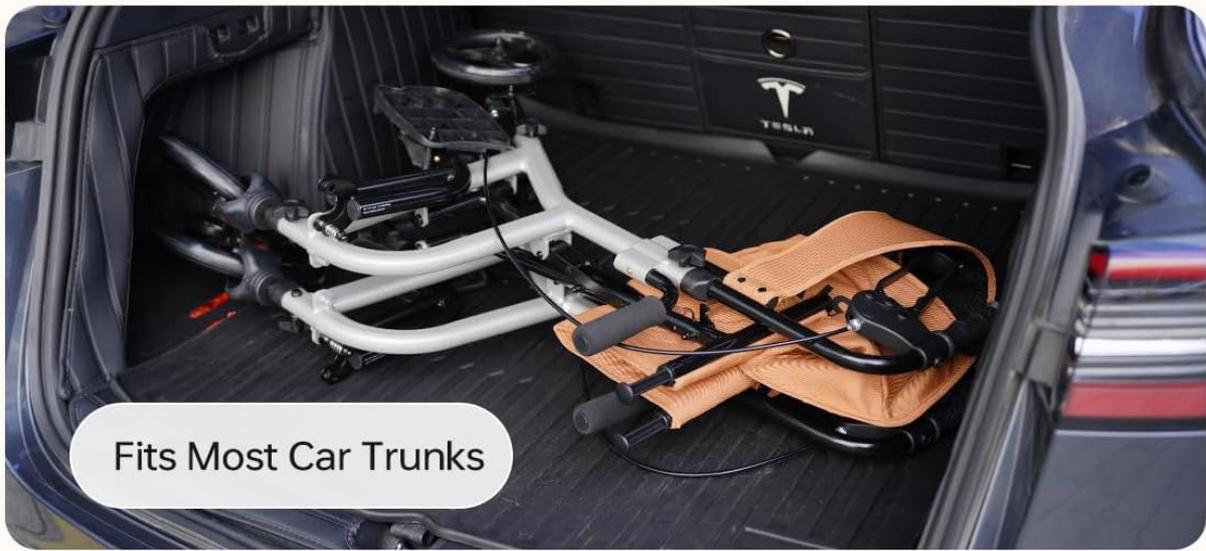

Fits Most Car Trunks

Image: Visual guide to the one-step folding mechanism and compact storage capability of the PAISEEC R1 Rollator Walker.

only **22lbs**

Effortless mobility, all on your own.



Image: The lightweight design of the PAISEEC R1 allows for easy portability, even when navigating stairs.

6. MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your rollator walker.

- Periodically inspect all fasteners (screws, knobs) and tighten if necessary.
- Clean the frame with a damp cloth and mild detergent. Avoid abrasive cleaners.
- Check the wheels for wear and tear, and ensure they rotate freely.
- Inspect brake cables for fraying or damage and ensure brakes engage properly.
- Store the rollator in a dry place away from direct sunlight when not in use.

7. TROUBLESHOOTING

If you encounter any issues with your PAISEEC R1 Rollator Walker, refer to the following common troubleshooting tips:

- **Wheels not rolling smoothly:** Check for debris caught in the wheels or axles. Ensure all wheel fasteners are secure.
- **Brakes not engaging:** Inspect brake cables for damage. Ensure the brake levers are properly adjusted and can fully engage.
- **Unit feels unstable:** Verify that all assembly points are securely fastened. Ensure the frame is fully unfolded and locked.

For persistent issues or concerns, please contact PAISEEC customer support.

8. SPECIFICATIONS

Feature	Detail
Model Number	R1
Weight	22 lbs (10 kg)
Weight Capacity	300 lbs (136 kg)
User Height Range	5'1" to 6'7" (155 - 200 cm)
Frame Material	Aluminum alloy
Seat Dimensions	17.7" (W) x 9.1" (L)
Wheel Size	10" front / 8" rear
Storage Bag Capacity	20L (holds up to 11 lbs)
Color	Brown

Dimensions and Parameter Details



Image: Detailed dimensions and parameter specifications for the PAISEEC R1 Rollator Walker.

9. WARRANTY AND SUPPORT

The PAISEEC R1 Rollator Walker comes with a 1-year warranty for worry-free mobility. PAISEEC also offers lifetime after-sales support.

- **Warranty:** 1-year limited warranty on the product.
- **After-Sales Support:** Lifetime support is available for any questions or issues.
- **Contact Information:**
 - Phone Support: 1(833) 612-2678 (Mon-Sun 9:00AM-6:00PM EST)
 - Email Support: support@paiseec.com (Replies within 24 hours)



Image: PAISEEC customer support information, including phone number and email address.