

PH 01

URIT Hemoglobin Test Strips User Manual

Model: PH 01

INTRODUCTION

This manual provides essential instructions for the proper use and storage of URIT Hemoglobin Test Strips. These strips are designed for use with a compatible URIT Hemoglobin Test Meter (meter not included) to measure hemoglobin levels. Each package contains 50 test strips and 50 lancets, intended for home use.

Portable Hemoglobin Test Strip

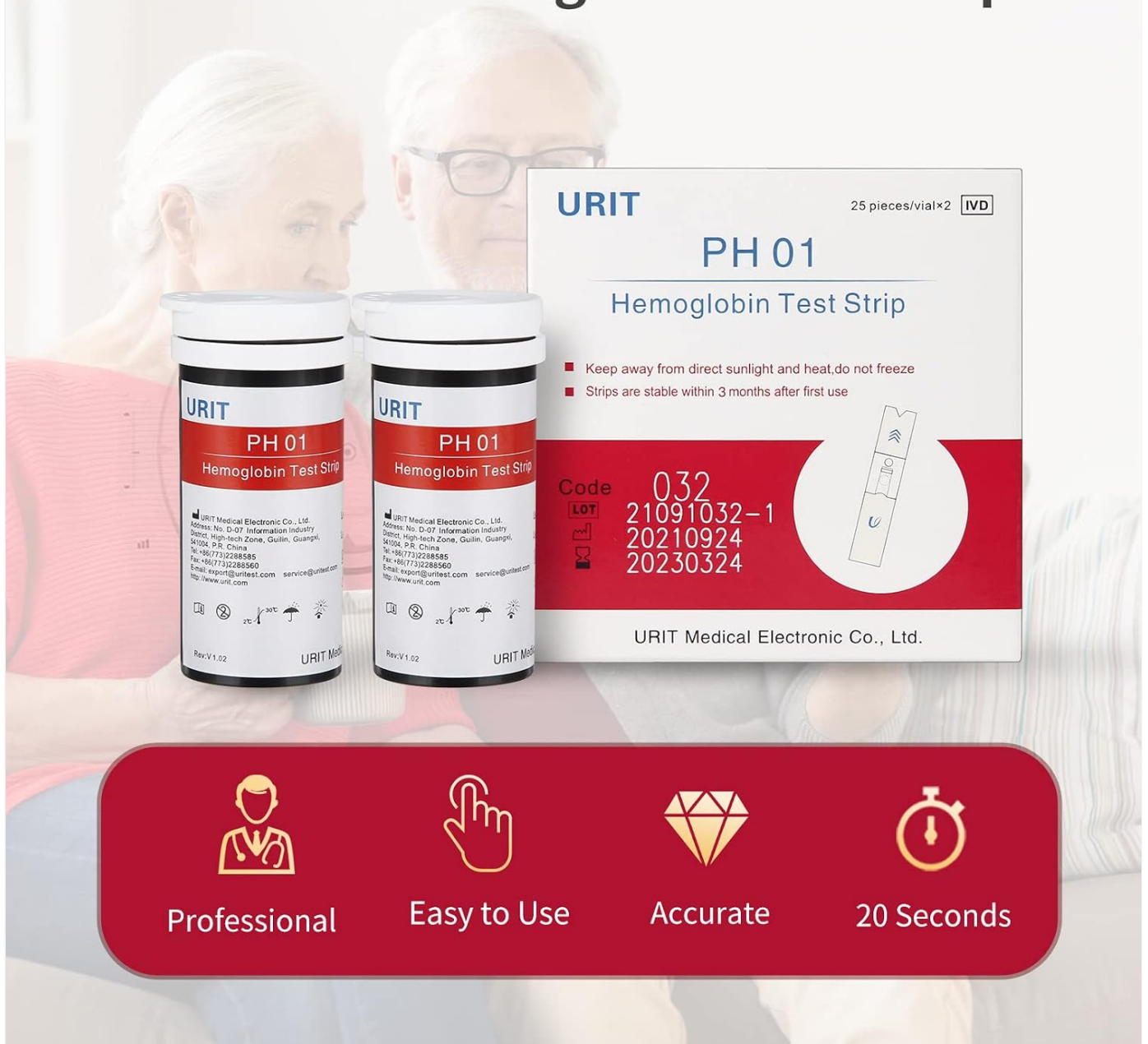


Image: URIT Hemoglobin Test Strips packaging, showing two vials of strips and a box with storage information. This image illustrates the product as it appears in its packaging.

PACKAGE CONTENTS

Verify that your package contains the following items:

- 50 Hemoglobin Test Strips
- 50 Lancets

PACKAGE INCLUDING



Test Strip × 50



Lancet × 50

Image: Contents of the package, showing one vial of 50 test strips and 50 individual lancets. This image confirms the items included with the product.

SETUP AND PREPARATION

1. **Ensure Compatibility:** These test strips are specifically designed for use with a URIT Hemoglobin Test Meter. Ensure your meter is compatible.
2. **Check Expiration Date:** Before use, always check the expiration date printed on the test strip vial. Do not use expired strips.
3. **Gather Supplies:** Have your URIT Hemoglobin Test Meter, a fresh test strip, a lancet, and an alcohol swab ready.
4. **Hand Hygiene:** Wash your hands thoroughly with soap and warm water, then dry them completely.

OPERATING INSTRUCTIONS

Follow these steps for accurate hemoglobin testing. Refer to your URIT Hemoglobin Test Meter's manual for specific instructions on meter operation.

1. **Prepare the Meter:** Turn on your URIT Hemoglobin Test Meter according to its instructions.
2. **Insert Test Strip:** Carefully remove one test strip from the vial. Immediately close the vial cap tightly to protect the remaining strips from moisture. Insert the test strip into the meter as directed by the meter's manual.
3. **Obtain Blood Sample:**
 - Use a fresh lancet to prick the side of your fingertip.
 - Gently squeeze your finger to form a small drop of blood.
4. **Apply Blood Sample:** Touch the tip of the test strip to the blood drop. The strip will draw the blood into the reaction area. Do not smear the blood.
5. **Read Results:** The meter will display your hemoglobin result after a short period (refer to your meter's manual for exact timing).
6. **Dispose Safely:** Dispose of the used test strip and lancet in a suitable sharps container.

MAINTENANCE AND STORAGE

Proper storage and handling are crucial for the accuracy and longevity of your test strips.

- **Storage Temperature:** Store strips in a dry place at a temperature between 39°F (4°C) and 86°F (30°C).
- **Avoid Extremes:** Do not expose test strips to heat, direct sunlight, or freezing temperatures.
- **Keep Vial Closed:** Always keep the test strip vial tightly capped immediately after removing a strip to prevent moisture exposure.
- **Usage After Opening:** Use test strips within 3 months after first opening the vial. Mark the opening date on the vial for reference.
- **Original Container:** Store test strips only in their original vial. Do not transfer them to other containers.

TROUBLESHOOTING

If you encounter issues during testing, consider the following:

- **Error Messages:** If your meter displays an error code, consult your URIT Hemoglobin Test Meter's user manual for specific troubleshooting steps related to the meter and strip interaction.
- **Inaccurate Results:**
 - Ensure strips are not expired.
 - Verify proper storage conditions have been maintained.
 - Confirm the test strip vial was tightly closed after each use.
 - Check if enough blood was applied to the strip.
 - Ensure the meter is clean and calibrated according to its manual.
- **No Blood Drawn:** Ensure the blood drop is sufficient and the strip is correctly touching the blood. Do not force blood onto the strip.

If problems persist, contact the manufacturer or a healthcare professional.

SPECIFICATIONS

Product Name	Hemoglobin Test Strips
---------------------	------------------------

Model	PH 01
Brand	Generic (Compatible with URIT Hemoglobin Test Meter)
Manufacturer	Bestar Lab
Unit Count	50 Test Strips, 50 Lancets
Test Type	Home Test
Storage Temperature	39°F - 86°F (4°C - 30°C)
Usage After Opening	Within 3 months

WARRANTY INFORMATION

Specific warranty details for these test strips are not provided in the product information. Please refer to the manufacturer's official website or contact their customer service for any warranty-related inquiries.

CUSTOMER SUPPORT

For technical assistance, questions about product usage, or to report any issues, please contact the manufacturer, Bestar Lab, or the seller.

Note: Always consult with a healthcare professional for medical advice and interpretation of test results.

