

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Fitiger](#) /

› [Fitiger EasyPumpVac Pro Portable Airway Suction Device Instruction Manual](#)

Fitiger EasyPumpVac Pro Home Kit 1

Fitiger EasyPumpVac Pro Portable Airway Suction Device

Instruction Manual - Model: EasyPumpVac Pro Home Kit 1

Brand: Fitiger

IMPORTANT SAFETY INFORMATION

This device is intended for emergency use only to clear airway obstructions. It is crucial to understand its proper operation before an emergency occurs. Always follow local first aid guidelines and seek professional medical attention immediately after using the device.

- Do not use if the device or any component appears damaged.
- Ensure a tight seal around the mouth and nose for effective suction.
- Practice with the included practice mask to familiarize yourself with the device.
- Keep out of reach of children when not in use.
- This device is not a substitute for professional medical care or CPR training.

PRODUCT COMPONENTS



Figure 1: The complete Fitiger EasyPumpVac Pro kit, including the suction device, three different sized masks (adult, child, baby), and a durable orange travel bag.

The Fitiger EasyPumpVac Pro Home Kit includes:

- **1x Suction Device:** The main unit that generates suction.
- **1x Adult Mask:** For individuals aged 12+ years.
- **1x Child Mask:** For children aged 3-12 years.
- **1x Baby Mask:** For infants aged 1-3 years.
- **1x Instruction Guide:** This manual.
- **1x Cloth Bag:** For convenient storage and portability.

CHOKING RESCUE DEVICE HOME KIT



Figure 2: Detailed view of all components included in the Fitiger Choking Rescue Device Home Kit, highlighting the different mask sizes for various age groups.

SETUP

1. **Unpack Components:** Remove all items from the orange travel bag.
2. **Select Appropriate Mask:** Choose the mask size (Adult, Child, or Baby) that best fits the individual's face. Ensure it covers both the nose and mouth completely.
3. **Attach Mask to Device:** Firmly press the selected mask onto the suction device's opening until it clicks securely into place.

UNIQUE DESIGN PREVENTS SECONDARY INJURY



Figure 3: The device features air release holes to safely release pressure during use, preventing the blockage from being pushed back into the throat and minimizing secondary injury.

OPERATING INSTRUCTIONS

In a choking emergency, rapid action is critical. The EasyPumpVac Pro is designed for quick and intuitive use.

1. **Position the Individual:** Lay the choking individual on their back.
2. **Apply the Mask:** Place the assembled device with the appropriate mask firmly over the individual's mouth and nose, ensuring a complete seal. The contoured design of the mask helps achieve this.
3. **Perform Suction:** With one hand holding the mask firmly against the face, use your other hand to quickly and forcefully pull the handle of the suction device upwards. This action creates a powerful vacuum to dislodge the obstruction.
4. **Check Airway:** After pulling, immediately remove the device and check if the airway is clear. If the obstruction is not removed, reposition the mask and repeat the suction process.

HOW DOES IT WORKING?

-----Design for Kids/Adults/Erderly Self-use without help



When you pull up the device and creat a vacuum environment, then the force will flow to the airway, not take much force to achieve apidly sucking out clogs.

Figure 4: The device is designed for self-use or assisted use. When pulled upwards, it creates a vacuum to clear airway blockages with minimal effort.


HIGH-EFFICIENT CHOKING RESCUE DEVICEN TO SAVE CRITICAL SECONDS

Fitiger Two-step operation quickly clears airway obstruction, saving precious seconds in emergencies.

Place



Pull

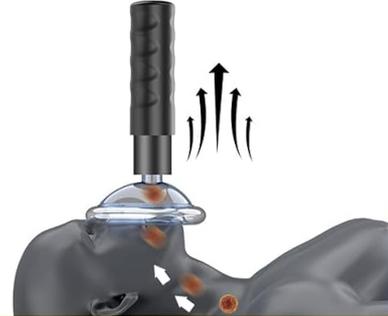


Other Brand

Three steps to remove airway obstruction will delay in emergency rescue.

Place



Press

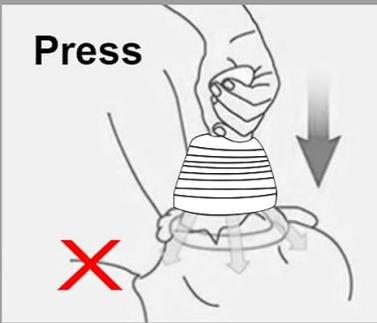

Pull

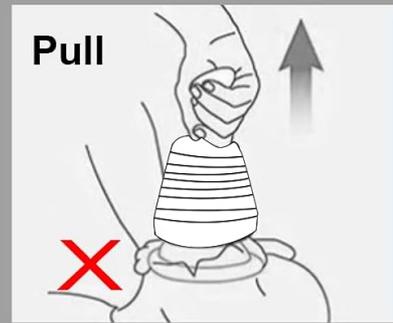


Figure 5: The Fitiger device offers a high-efficient, two-step operation (Place, Pull) to quickly clear airway obstructions, saving critical seconds in emergencies compared to multi-step devices.

MAINTENANCE

Proper maintenance ensures the device remains effective and hygienic.

- **Cleaning:** After each use, disassemble the mask from the suction device. Wash both components with warm, soapy water. Rinse thoroughly and allow to air dry completely before reassembling or storing.
- **Storage:** Store the device and masks in the provided cloth bag in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Mask Replacement:** It is recommended to inspect masks periodically for any signs of wear or damage. Replace masks every two to three years, or immediately if they show signs of degradation, to ensure optimal performance and seal.
- **Device Longevity:** The suction device itself does not require replacement unless it has been used in a choking rescue or shows signs of damage.

TROUBLESHOOTING

If the device does not perform as expected, consider the following:

- **No Suction:**

- Ensure the mask is firmly attached to the device.
- Verify a complete seal around the individual's mouth and nose. Any gaps will prevent effective suction.
- Check for any visible damage to the mask or device that might compromise the seal.

- **Difficulty Pulling Handle:**

- Ensure the device is clean and free of debris that might impede movement.
- The handle should move smoothly. If it feels stiff, ensure it's correctly assembled.

If issues persist, contact Fitiger customer support for assistance.

SPECIFICATIONS

Feature	Detail
Brand	Fitiger
Model	EasyPumpVac Pro Home Kit 1
Package Dimensions	12.28 x 8.03 x 3.7 inches
Item Weight	1.06 Pounds
Special Features	Compact, Portable, Reusable, Waterproof
Recommended Uses	Home, Travel, Emergency Preparedness

MEDICAL GRADE




Powerful Suction



Reliable Safety



Ergonomic Handle

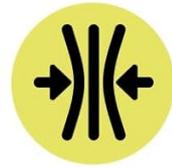

Prevent Bend/
Prevent Leak



Reusable Device



Rapid Rescue
Response



Figure 6: Key features of the medical-grade Fitiger EasyPumpVac Pro, emphasizing its reliability and ease of use.

WARRANTY AND SUPPORT

For warranty information, product support, or any inquiries, please contact Fitiger customer service through the retailer where the product was purchased or visit the official Fitiger website.

Manufacturer: Shenzhen Huayu Wantuo Technology Co., Ltd

Date First Available: April 5, 2025