

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

- › [Signia](#) /
- › [Signia Fun P Hearing Aid Instruction Manual](#)

## Signia Fun P

# Signia Fun P Hearing Aid Instruction Manual

Model: Fun P

## INTRODUCTION

The Signia Fun P is a Behind-The-Ear (BTE) hearing aid designed to provide clear sound and speech understanding for individuals with severe to profound hearing loss. It incorporates features for speech and noise management, making it suitable for senior citizens and others seeking enhanced auditory experiences in various environments.

This manual provides essential information on how to set up, operate, maintain, and troubleshoot your Signia Fun P hearing aid to ensure optimal performance and comfort.

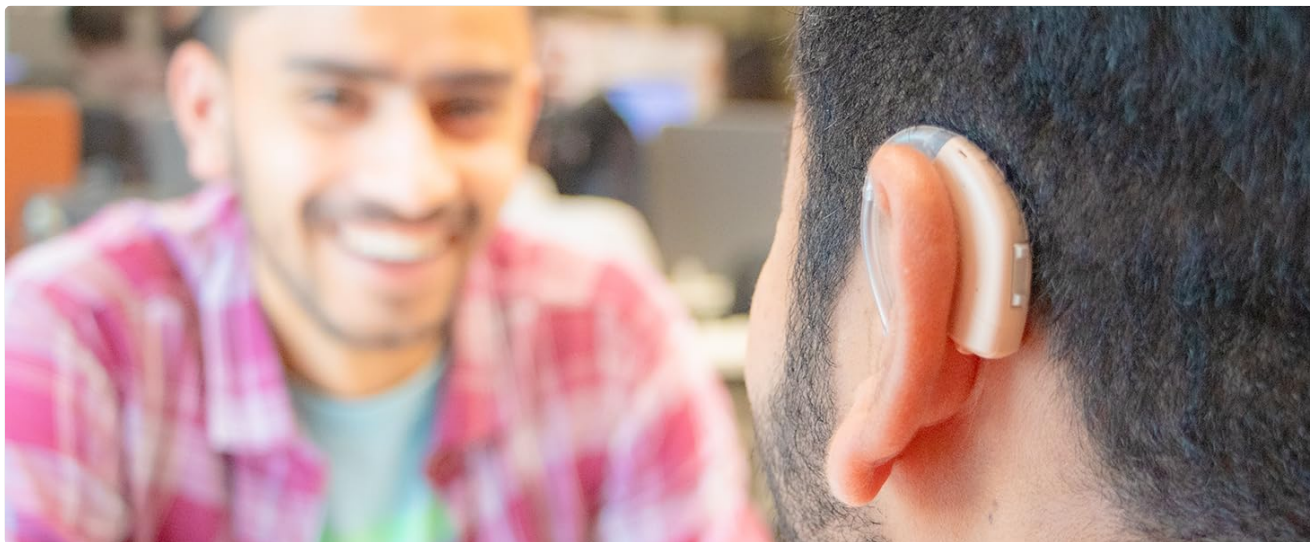


Image: A woman smiling, showcasing the discreet and comfortable design of the Signia Fun P hearing aid, designed to enhance daily interactions.

## PRODUCT COMPONENTS OVERVIEW

Familiarize yourself with the key components of your Signia Fun P hearing aid:

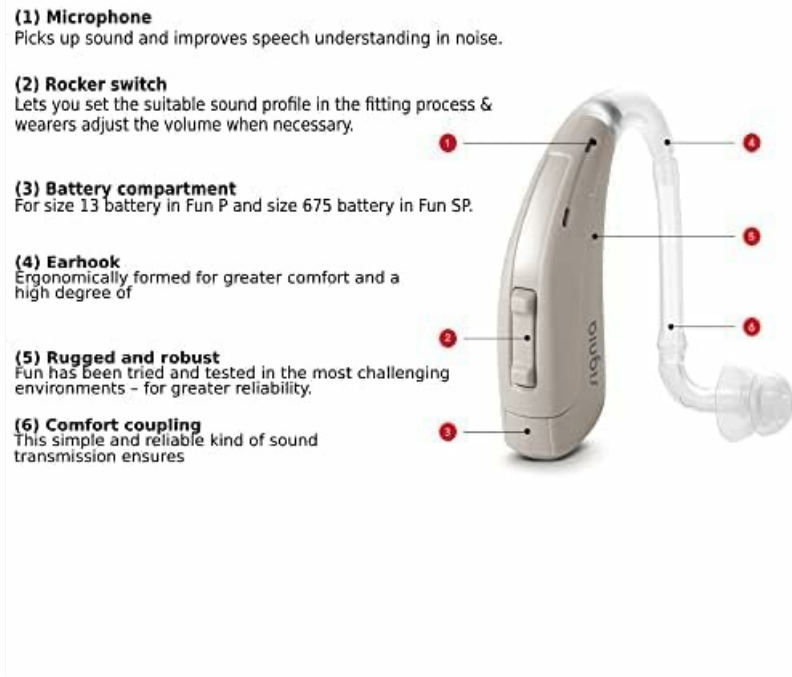


Image: Detailed diagram illustrating the main parts of the Signia Fun P hearing aid, including the microphone, rocker switch, battery compartment, earhook, rugged housing, and comfort coupling.

1. **Microphone:** Picks up sound and improves speech understanding, especially in noisy environments.
2. **Rocker Switch:** Used for easy volume control and adjusting sound profiles.
3. **Battery Compartment:** Houses the battery (Size 13 for Fun P).
4. **Earhook:** Ergonomically formed for comfort and high amplification.
5. **Rugged Housing:** Designed to be resistant to moisture, sweat, dust, and dirt.
6. **Comfort Coupling:** Ensures simple and reliable sound transmission.

## SETUP AND FIRST USE

---

### Unpacking and Inspection

Carefully remove your Signia Fun P hearing aid from its packaging. Inspect the device for any visible damage. Ensure all components, including the earhook and ear tip, are present.



Image: A user holding the Signia Fun P hearing aid and its accompanying case, demonstrating the compact design.

## Battery Insertion

The Signia Fun P uses a Size 13 battery. To insert the battery:

1. Gently open the battery compartment door (Component 3 in the diagram).
2. Remove the protective tab from the new battery. Wait for approximately one minute before inserting to allow the battery to activate fully.
3. Insert the battery with the positive (+) side facing the correct orientation, as indicated inside the compartment.
4. Close the battery compartment door firmly until it clicks into place.

## Wearing the Hearing Aid

The Fun P is designed for comfortable placement behind the ear.

1. Hold the hearing aid with the earhook (Component 4) pointing upwards.
2. Place the main body of the hearing aid behind your ear.
3. Gently guide the earhook over the top of your ear.
4. Insert the ear tip (connected to the comfort coupling) into your ear canal until it feels secure and comfortable. The



customizable tube length and super-soft tips ensure a precise fit.



Image: A Signia Fun P hearing aid worn behind the ear, demonstrating proper placement and fit.

## OPERATING YOUR HEARING AID

---

### Powering On/Off

To turn the hearing aid on, fully close the battery compartment door. To turn it off, open the battery compartment door slightly until it clicks into the 'off' position.

### Volume Control

Adjust the volume using the rocker switch (Component 2) located on the hearing aid. Press the upper part of the switch to increase volume and the lower part to decrease it.

### Speech and Noise Management

The Signia Fun P is equipped with 4 channels for advanced performance, allowing it to focus on conversations and manage background noise effectively. The microphone system (Component 1) is designed to capture and amplify speech while

reducing unwanted ambient sounds.

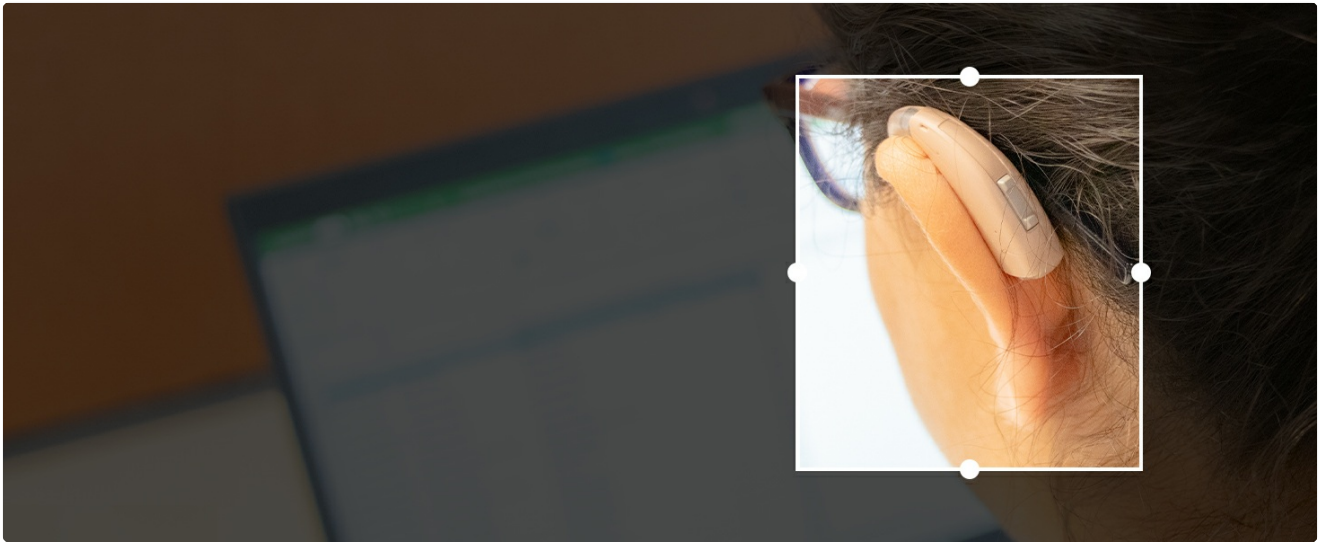


Image: A man wearing a Signia Fun P hearing aid, actively participating in a conversation, illustrating the device's ability to enhance communication.

## MAINTENANCE AND CARE

---

### Daily Cleaning

Regular cleaning is essential for optimal performance and hygiene. Use a soft, dry cloth to wipe the hearing aid daily. Ensure no moisture enters the device. The rugged housing (Component 5) is resistant to moisture, sweat, dust, and dirt, but proper cleaning extends its lifespan.

### Ear Tip and Earhook Care

Check the ear tip and earhook regularly for wax buildup or debris. Clean them gently with a specialized cleaning tool or a soft brush. Replace ear tips as recommended by your hearing care professional or when they show signs of wear.

### Battery Replacement

Replace the battery when you notice a decrease in sound quality or when the low battery indicator (if available) activates. Always use fresh Size 13 batteries.

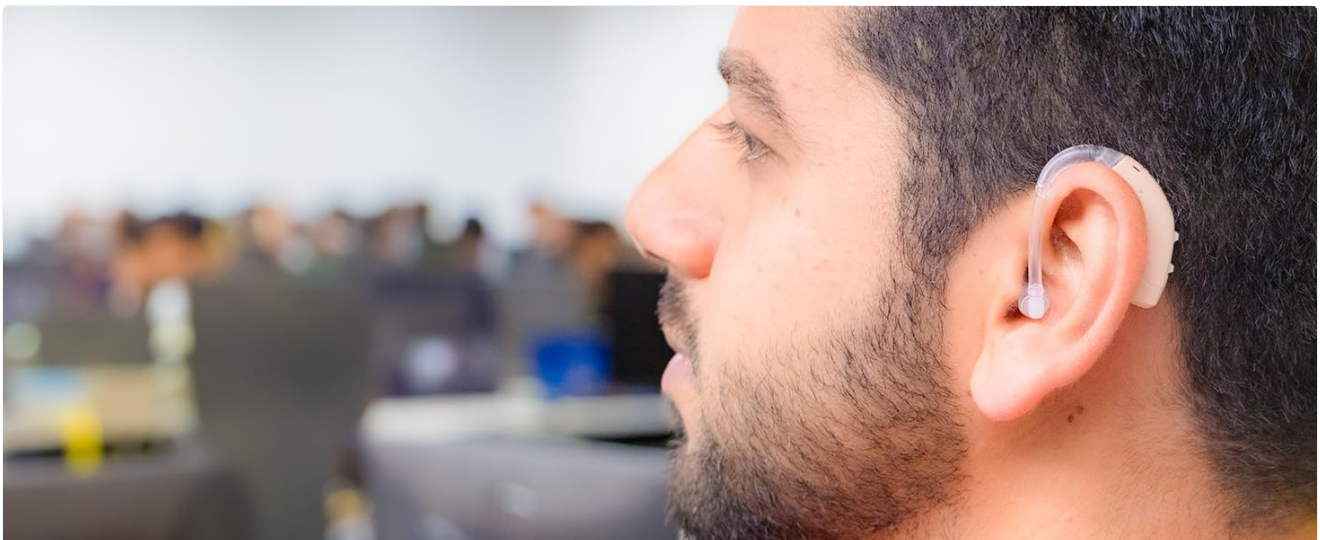


Image: Hands demonstrating the process of replacing the battery in a Signia hearing aid, showing the open battery compartment.



Image: A pack of Signia Size 13 hearing aid batteries, essential for powering the device.

## Storage

When not in use, store your hearing aid in its case in a cool, dry place. Open the battery compartment to conserve battery life and allow moisture to evaporate.

## TROUBLESHOOTING

If you experience issues with your Signia Fun P hearing aid, refer to the following common troubleshooting steps:

Problem	Possible Cause	Solution
No sound or weak sound	Dead or weak battery Blocked ear tip/earhook Volume too low	Replace battery Clean ear tip/earhook Increase volume using rocker switch



Problem	Possible Cause	Solution
Whistling or feedback	Improper fit Earwax buildup in ear canal Volume too high	Re-insert hearing aid for a snug fit Consult a hearing care professional for earwax removal Decrease volume
Sound is distorted or unclear	Low battery Moisture in microphone Damaged component	Replace battery Allow device to dry in a dehumidifier (if available) Contact support for inspection
Hearing aid not turning on	Battery compartment not fully closed Dead battery	Ensure battery compartment is fully closed Replace battery

If the problem persists after trying these steps, please contact customer support.

## SPECIFICATIONS

---

Key technical details of the Signia Fun P hearing aid:

- **Model:** Fun P
- **Type:** Behind-The-Ear (BTE)
- **Channels:** 4 Channels for enhanced performance and customization
- **Hearing Loss Level:** Suitable for Severe to Profound Hearing Loss
- **Special Feature:** Noise Cancellation, Speech and Noise Management
- **Battery Type:** Size 13
- **Dimensions (Product):** Approximately 10 x 10 x 10 cm
- **Item Weight:** Approximately 500 g
- **Manufacturer:** Signia GmbH
- **Country of Origin:** Germany

## WARRANTY AND SUPPORT

---


### Warranty Information

The Signia Fun P comes with a comprehensive **2-year warranty** provided by Soundbright. This warranty offers hassle-free doorstep service for repairs or replacements, ensuring reliable support for your device.

### Customer Support

For any concerns, technical assistance, or warranty claims, please contact our support team:

- **Phone:** +91 99108 87323
- **Email for Programming Support:** [store@hearstore.in](mailto:store@hearstore.in)
- **WhatsApp:** +91 99108 87323

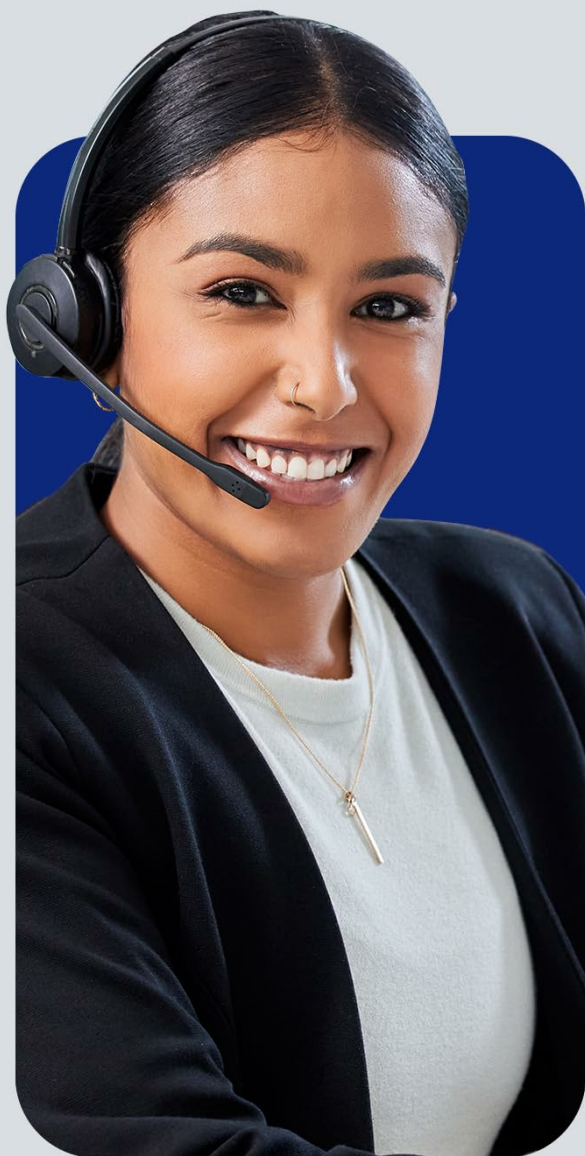

# Need Help?

Hearing Aid Experts Are Here For You!

For any concerns, call us at:

**+91 99108 87323**

Image: A customer service representative ready to assist, symbolizing dedicated support for your Signia hearing aid.



# Free Programming Support

Email your audiometry report to:  
**store@hearstore.in**

WhatsApp us at:  
**+91 9910887323**



Image: Information on free programming support, including contact details for email and WhatsApp assistance.