



Manuals.plus /

- › Starlite Smile /
- › Red Light Earbuds User Manual

Starlite Smile B0DXVCLX12

Red Light Earbuds User Manual

Model: B0DXVCLX12 | Brand: Starlite Smile

1. INTRODUCTION

Thank you for choosing the Starlite Smile Red Light Earbuds. This device is designed to direct 650nm red light to the ear canal, offering a convenient and portable solution for ear support. This manual provides essential information on how to set up, operate, maintain, and troubleshoot your new red light therapy device. Please read it thoroughly before use.



Image: The Red Light Earbuds with their USB connector, showing the red light emission. An inset image displays a user wearing the earbuds.

2. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before use and retain for future reference. Your safety is paramount.

- **Consult a Physician:** Always consult with your private physician or a qualified healthcare professional before using this device, especially if you have any pre-existing ear conditions, medical implants, or are undergoing other treatments.
- **FDA Disclaimer:** Statements regarding dietary supplements and health devices have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.
- **Children:** Keep this device out of reach of children.
- **Disassembly:** Do not attempt to disassemble, modify, or repair the device yourself. This may void the warranty and pose safety risks.
- **Water Exposure:** Do not immerse the earbuds or USB cable in water or other liquids. Avoid using in wet

environments.

- **Damage:** If the device appears damaged, discontinue use immediately and contact customer support.

3. PRODUCT COMPONENTS

The Starlite Smile Red Light Earbuds package includes the following items:

- One (1) pair of Red Light Earbuds with integrated USB cable
- Two (2) pairs of removable silicone ear canal covers (one pair pre-installed, one extra pair)



4. SETUP

The Red Light Earbuds are designed for simple plug-and-play operation. No charging is required.

4.1 Connecting the Device

To use the earbuds, simply plug the USB end into any compatible USB port. Compatible power sources include:

- Laptop or PC USB port
- Power strip with a USB port
- Power bank with a USB port



Ways to Use:

Plug USB end into any power bank.
(Power bank is NOT included, sold separately).



Plug USB end into any power strip with USB port.
(Power strip is NOT included, sold separately).



Plug USB end into any computer or laptop with USB port.
(Computer or laptop are NOT included, sold separately).

Image: Visual guide showing the USB connector being plugged into a power bank, a power strip, and a laptop.
Ensure the silicone ear canal covers are securely attached to the earbuds before insertion into your ears.

5. OPERATING INSTRUCTIONS

5.1 Usage Guidelines

Once plugged into a power source, the red light will activate automatically. Gently insert the earbuds into your ear canals, ensuring a comfortable fit.

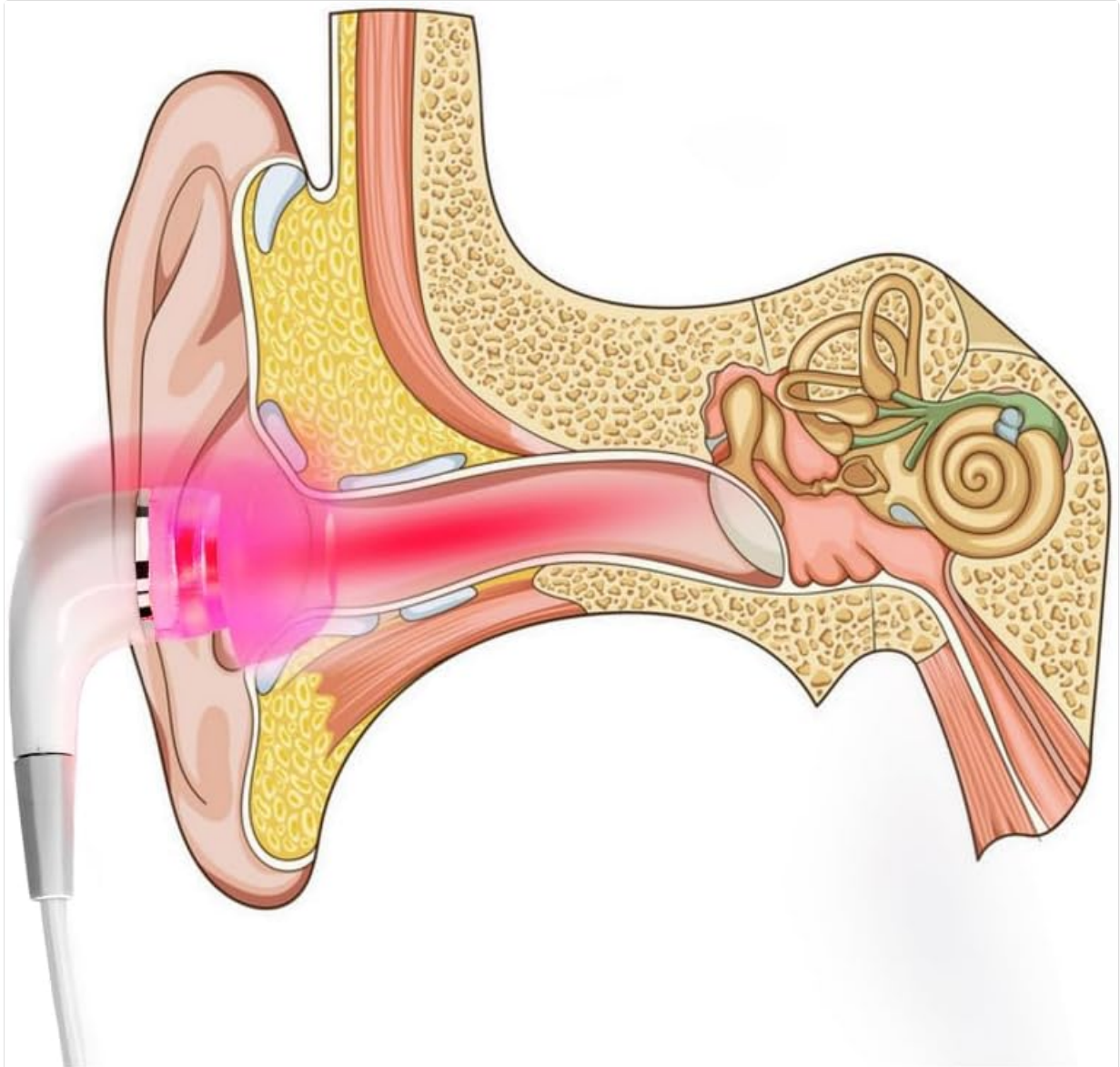


Image: An anatomical cross-section of the ear showing the red light from the earbud reaching into the ear canal.

For optimal use, it is generally recommended to use the red light earbuds 2 times a day for approximately 30 minutes each session. Always follow the specific advice and recommendations of your healthcare professional regarding usage duration and frequency.

5.2 Red Light Wavelength

These earbuds emit 650nm red light. This specific wavelength is chosen for its ability to penetrate through the thin skin of the ear canal.

6. MAINTENANCE AND CARE

6.1 Cleaning

For best ear hygiene, the silicone ear canal covers are removable and washable. Regular cleaning is recommended.

- Gently remove the silicone covers from the earbuds.
- Wash the silicone covers with mild soap and water.
- Rinse thoroughly and allow them to air dry completely before reattaching them to the earbuds.
- Do not immerse the main earbud unit or USB cable in water or any cleaning solutions.
- Wipe the main earbud unit with a soft, dry cloth if needed.

7. TROUBLESHOOTING

If you encounter any issues with your Red Light Earbuds, please refer to the table below for common problems and their solutions.

Problem	Possible Cause	Solution
Red light does not turn on.	Not properly plugged in or power source issue.	Ensure the USB connector is fully inserted into a working USB port. Try a different USB port or power source (e.g., another laptop, power bank, or power strip).
Earbuds feel uncomfortable or do not fit well.	Incorrect silicone cover size or improper insertion.	Ensure you are using the correct size silicone covers for your ears. Gently adjust the earbuds in your ear canal until a comfortable and secure fit is achieved.
Device appears physically damaged.	Accidental damage or manufacturing defect.	Discontinue use immediately. Refer to the warranty information in Section 9 for replacement options. Do not attempt to repair.

8. PRODUCT SPECIFICATIONS

- **Model:** B0DXVCLX12
- **Brand:** Starlite Smile
- **Light Wavelength:** 650nm Red Light
- **Power Source:** USB (No charging required)
- **Package Dimensions:** 3 x 3 x 1.5 inches

- **Weight:** 3 ounces
- **Country of Origin:** China
- **Date First Available:** February 20, 2025

9. WARRANTY AND SUPPORT

9.1 90-Day Replacement Window

This product comes with a 90-day replacement window. If you encounter any problems or defects with your device within 90 days of purchase, please contact Starlite Smile for assistance and a potential replacement.

9.2 Contact Information

For support, inquiries, or to claim your replacement, please visit the Starlite Smile store on Amazon or contact customer service through your purchase platform.

Starlite Smile Store: [Visit Store](#)