

## AlveyTech CTL1804039

# AlveyTech 10-Button 5-Pin Replacement Remote Hand Control (CTL1804039) - User Manual

For Pride VivaLift Lift Chairs (Standard & Early Version VivaLift 2 Models)

## 1. PRODUCT OVERVIEW

This manual provides instructions for the AlveyTech 10-Button 5-Pin Replacement Remote Hand Control, designed for use with specific Pride VivaLift Lift Chair models and other compatible power recliners. This remote serves as a direct replacement for the original hand control, offering full functionality for chair adjustments.

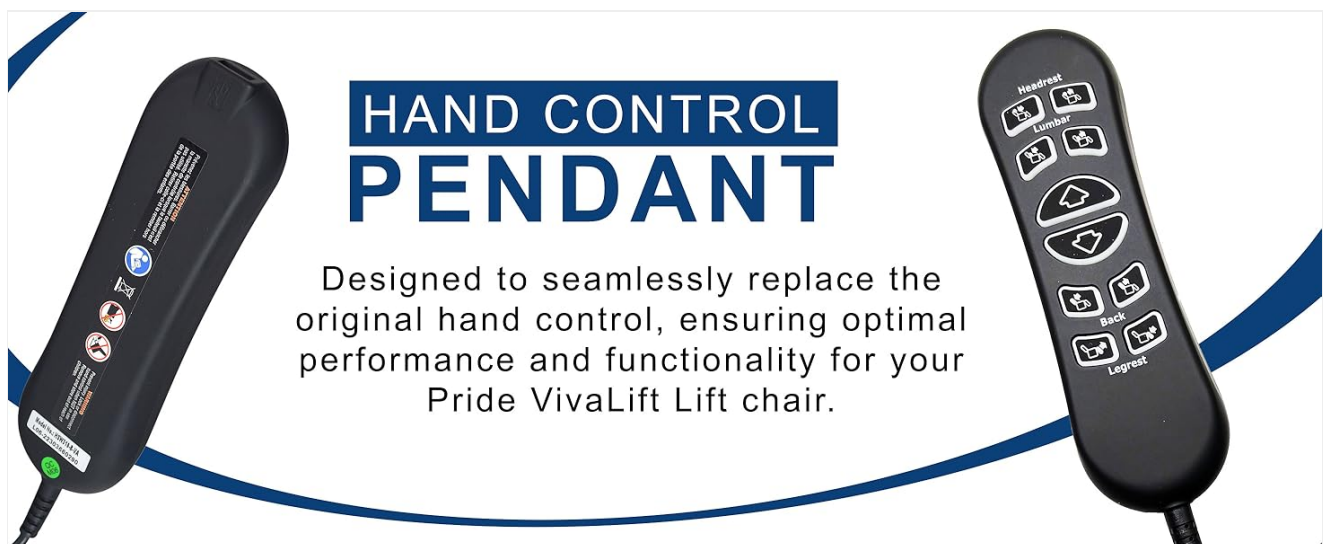


Image: The AlveyTech 10-Button 5-Pin Replacement Remote Hand Control.

## 2. COMPATIBILITY

This replacement hand control (Part Number CTL1804039) is manufactured by Pride Mobility and is compatible with the following models:

- Standard Pride VivaLift Lift Chairs
- Early Version VivaLift 2 Models (HDW310-KBA and HSW310-B-KBA)
- Mega Motion power recliners using the MM3702 dual motor assembly
- Jackson Catnapper Warner lift chairs
- Flexsteel Jenkins lift chair models

**Important Note:**

- This hand control is **NOT compatible** with Ultra, Radiant, or later VivaLift 2 models.
- Always verify the button configuration of your existing hand control prior to purchase to ensure correct functionality.

Compatible with  
**Standard VivaLift & Early VivaLift 2**  
models only



**\*Not Compatible with Radiant,  
Ultra, or Later VivaLift 2 models**

*Image: Illustration of compatible lift chair models. Note the explicit warning about incompatibility with Radiant, Ultra, or later VivaLift 2 models.*



Image: Visual guide for verifying button configuration. The correct remote has 10 buttons for Headrest, Lumbar, Back, and Legrest, plus up/down lift controls. The incorrect remote has additional 'M1' and 'M2' buttons.

### 3. KEY FEATURES

- **10-Button Control:** Dedicated buttons for Headrest, Lumbar, Back, Legrest, and Up/Down lift functions.
- **5-Pin Connector:** Ensures quick and easy plug-n-play installation.
- **USB Port:** Conveniently located at the top for charging mobile or other devices.
- **Ergonomic Design:** Shaped for comfortable handling.
- **Retractable Cord:** Provides flexibility and neat storage.
- **Upgraded Cable Joint:** Designed to prevent breakage in the latest models.


# Key Features



5-Pin connector



USB port

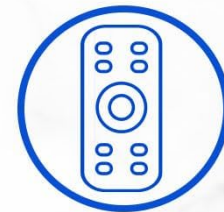

10 buttons

Image: The remote control with visual callouts for its 5-pin connector, USB port, and 10-button layout.



Ergonomic design perfectly fits in the hand



Made with high quality materials for durability



Retractable cord for better cable management

Image: The remote control demonstrating its ergonomic shape and retractable cord for ease of use and storage.

## 4. SETUP INSTRUCTIONS

---

Installation of the AlveyTech 10-Button 5-Pin Replacement Remote Hand Control is a simple plug-n-play process.

1. **Locate the Existing Connector:** On your lift chair, find the existing hand control cable and its connection point. This is typically located in a side pocket or underneath the chair.
2. **Disconnect Old Remote (if applicable):** If replacing an old remote, carefully disconnect its 5-pin connector from the chair's control system.
3. **Connect New Remote:** Align the 5-pin connector of the new AlveyTech remote with the corresponding port on your chair. Gently push until it clicks securely into place. Ensure the connection is firm.
4. **Test Functionality:** Once connected, test all buttons on the remote to ensure proper operation of the headrest, lumbar, back, legrest, and lift functions.



This hand control will only work with a **5-pin configuration**


*Image: A detailed view of the 5-pin connector, which ensures a secure and easy connection to your lift chair.*

*Video: A short demonstration of the 10-button hand control for Mega Motion & VivaLift Chairs, showing its features and ease of use.*

## 5. OPERATING INSTRUCTIONS

---

The AlveyTech 10-Button remote provides comprehensive control over your lift chair's various functions. Each button is clearly labeled with an icon indicating its purpose.

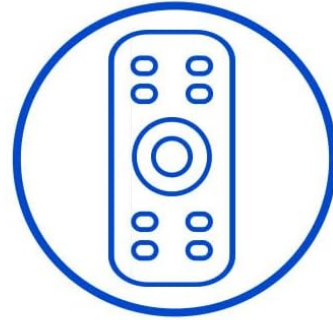

# 10 Buttons

- Headrest Controls
- Lumbar Controls
- Lift Controls
- Back Position
- Legrest Position

*Image: The remote control displaying its 10 buttons and their respective functions.*

## Button Functions

Button Label	Function
Headrest (Up/Down)	Adjusts the headrest position up or down.
Lumbar (Up/Down)	Adjusts the lumbar support position up or down.
Lift (Up/Down Arrows)	Moves the entire chair up for assistance in standing or down to a seated position.
Back (Recline/Upright)	Adjusts the backrest angle, reclining it backward or bringing it upright.

Button Label	Function
Legrest (Up/Down)	Raises or lowers the legrest.

**USB Charging:** To charge your mobile device, connect your device's USB cable to the USB port located at the top of the remote control.



USB  
port



## Charge your phone or other electronic devices while you relax

*Image: A close-up view of the USB port, allowing for convenient device charging.*

### 6. MAINTENANCE

---

- **Cleaning:** Wipe the remote control with a soft, dry cloth. Do not use liquid cleaners or abrasive materials, as these may damage the surface or internal components.
- **Storage:** When not in use, store the remote in a safe place where it will not be subjected to excessive force or extreme temperatures. The retractable cord helps keep it tidy.
- **Cable Inspection:** Periodically check the cable for any signs of wear, fraying, or damage. If damage is observed, discontinue use and seek a replacement.

### 7. TROUBLESHOOTING

---

If your lift chair is not responding correctly after installing the new remote, consider the following:

- **Check Connection:** Ensure the 5-pin connector is fully and securely plugged into the chair's control system. A loose connection can prevent the remote from functioning.
- **Power Supply:** Verify that the lift chair is properly plugged into a working power outlet and that the power supply is active.
- **Other Components:** If the hand control matches the one shown in the compatibility section and is correctly installed, but the issue persists, it is likely that another component of the lift chair, such as the control box or motor, may require repair or replacement.



*Image: The remote control shown with a control box and motor, illustrating other components that might need attention if the remote is not the source of the problem.*

## 8. PRODUCT SPECIFICATIONS

Feature	Detail
Brand	AlveyTech

Feature	Detail
Model Name	Hand Control
Part Number	CTL1804039
Color	Black
Material	Plastic
Connector Type	5-Pin
Buttons	10 (Headrest, Lumbar, Lift, Back, Legrest)
Additional Features	USB Charging Port, Retractable Cord
Item Dimensions	5.91 x 7.87 x 1.97 inches

## 9. IMPORTANT SAFETY INFORMATION

---

### WARNING:

- Always ensure you are replacing the correct part for your specific model. Incorrect replacement may cause the chair to malfunction, posing a safety risk.
- Consult your chair's manual or contact the manufacturer if you are unsure about compatibility.
- Keep the remote control out of reach of children to prevent accidental operation or damage.





*Image: The reverse side of the remote control, displaying important safety warnings and product identification.*

## 10. WARRANTY & SUPPORT

---

For information regarding warranty coverage, technical support, or replacement parts, please refer to the documentation provided with your original lift chair or contact AlveyTech customer service directly. Ensure you have your product model number (CTL1804039) and chair model information available when contacting support.