

manuals.plus /

- › [Invacare](#) /
- › [Invacare Essential SX3 Manual Wheelchair User Manual for 16-Inch Seat](#)

Invacare Essential SX3 Manual Wheelchair (ISX366DASA)

Invacare Essential SX3 Manual Wheelchair User Manual

Model: ISX366DASA

INTRODUCTION

This manual provides essential instructions for the safe and effective use, assembly, maintenance, and troubleshooting of your Invacare Essential SX3 Manual Wheelchair. Please read this manual thoroughly before using the wheelchair to ensure proper operation and to prevent injury.

PRODUCT OVERVIEW

The Invacare Essential SX3 Manual Wheelchair is designed for adults requiring mobility assistance. It features a lightweight, foldable frame, comfortable padded armrests, and adjustable footrests for enhanced user experience.



Figure 1: Invacare Essential SX3 Manual Wheelchair, front-side view.

Key Features:

- **Flip-Back Desk Length Armrests:** Padded armrests designed for easy access to desks and tables, and for simplified transfers.
- **Lightweight and Foldable Design:** Weighs approximately 38 lbs and folds to 12.3 inches wide for transport and storage.
- **Adjustable Seat Depth:** Standard 16-inch depth, extendable by 2 inches with an included kit.
- **Durable Vinyl Upholstery:** Easy to clean and resistant to wear.
- **Swing-Away and Adjustable Footrests:** Removable and height-adjustable for user comfort and easier transfers.
- **Smooth-Rolling Wheels:** 8-inch PVC front wheels and 24-inch rear mag PU tires with built-in wheel locks.

COMPONENTS

Familiarize yourself with the main components of your wheelchair:

1. **Frame:** Powder-coated steel frame.
2. **Seat and Backrest:** Vinyl upholstery.
3. **Armrests:** Padded, flip-back, desk-length.
4. **Footrests/Legrests:** Swing-away, height-adjustable.
5. **Front Casters:** 8-inch PVC wheels.
6. **Rear Wheels:** 24-inch mag PU tires with push-rims.
7. **Wheel Locks:** For securing the wheelchair.
8. **Push Handles:** For attendant propulsion.

Supportive and Accessible



Figure 2: Front view of the wheelchair, showing footrests and front casters.

SETUP AND ASSEMBLY

The Invacare Essential SX3 Wheelchair arrives fully assembled. Follow these steps for initial preparation:

1. **Unfolding the Wheelchair:**

- Place the folded wheelchair on a flat, stable surface.
- Grasp the push handles and gently pull them apart until the wheelchair frame begins to open.
- Push down firmly on both sides of the seat rails until the seat is fully flat and locked into position. Ensure there is an audible click.

2. Attaching/Adjusting Footrests:

- Align the footrest hangers with the receiver brackets on the front of the wheelchair frame.
- Swing the footrests inward until they lock into place.
- To adjust the height, press the release button on the legrest tube and slide the footplate up or down to the desired position. Release the button to lock.
- To remove, swing the footrests outward and lift them off the receiver brackets.

3. Adjusting Seat Depth (Optional):

- The wheelchair includes a kit to extend the seat depth by 2 inches. Refer to the separate instructions provided with the extension kit for detailed installation.

4. Inspecting Wheelchair: Before first use, ensure all components are securely attached, wheels roll freely, and wheel locks engage properly.

Multiple Sizes To Suit Your Needs

SEAT AVAILABLE IN THREE STANDARD SIZES:



Figure 3: Flip-back armrests allow seamless access to desks and tables.

OPERATING INSTRUCTIONS

Self-Propulsion:

- Grasp the push-rims on the rear wheels.
- Push forward on the rims to move forward, and pull back to slow down or reverse.
- To turn, apply more force to one push-rim than the other.


Attendant Propulsion:

- An attendant can push the wheelchair using the ergonomic non-slip hand grips on the back of the chair.
- Use the integrated steel step tube to assist in lifting the front casters over small obstacles like curbs or bumps.

On-the-go Convenience



Cross-frame design to ensure folding is simple and effortless

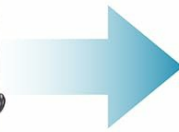

12.2"
Makes storage and transportation easy



Arrives fully assembled



Figure 4: Attendant pushing the wheelchair, demonstrating the use of the integrated step tube.

Engaging Wheel Locks:

- Always engage both wheel locks when transferring in or out of the wheelchair, or when stationary.
- Push the lever on each wheel lock forward until it clicks into place, securing the rear wheels.
- To disengage, pull the lever back.

Sturdy & Durable Design



Figure 5: Close-up of the secure locking wheel mechanism.

Folding the Wheelchair:

- Ensure footrests are swung away or removed.
- Grasp the center of the seat upholstery and lift upwards. The frame will begin to collapse inward.
- Continue lifting until the wheelchair is fully folded.

Effortless Control



Figure 6: The wheelchair folds compactly for storage and transport.

MAINTENANCE

Cleaning:

- Wipe down the vinyl upholstery and frame with a damp cloth and mild soap.
- Avoid abrasive cleaners or solvents that could damage the finish.
- Dry thoroughly after cleaning.

Inspection:

- Regularly check all nuts, bolts, and fasteners for tightness.
- Inspect tires for wear and proper inflation (if applicable for rear tires).
- Ensure wheel locks engage and disengage smoothly and securely.
- Check armrests and footrests for secure attachment and proper function.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Wheelchair does not roll smoothly.	Loose axle, debris in wheels, worn tires.	Check axle tightness, clean wheels, inspect and replace worn tires.
Wheel locks do not engage.	Misaligned brake pads, loose mechanism.	Adjust brake pads or tighten mechanism. Consult a qualified technician if unable to resolve.
Footrests are loose or difficult to attach.	Bent hanger, debris in receiver.	Inspect for damage, clean receiver brackets.

SPECIFICATIONS

- **Model Number:** ISX366DASA
- **Product Dimensions (Open):** 31.88 x 24.25 x 34 inches
- **Weight:** 38 Pounds
- **Seat Width:** 16 inches
- **Seat Depth:** 16 inches (extendable to 18 inches)
- **Front Wheels:** 8-inch PVC
- **Rear Wheels:** 24-inch Mag PU tires
- **Manufacturer:** Invacare

Transfer With Ease



Swing back or
detachable
armrests for simple
entry and exits



Figure 7: Key dimensions and maximum weight capacity of the wheelchair.

SAFETY INFORMATION

- Always engage wheel locks before transferring in or out of the wheelchair.
- Do not attempt to operate the wheelchair on steep inclines or declines without assistance.
- Ensure all parts are securely fastened before each use.
- Do not exceed the maximum weight capacity of 265 lbs.
- Avoid sudden movements or changes in direction that could cause instability.
- Regularly inspect the wheelchair for any signs of wear or damage. Do not use if damaged.

WARRANTY AND SUPPORT

For warranty information and customer support, please contact Invacare directly. Refer to the product packaging or the official Invacare website for specific warranty terms and contact details.





Manufacturer: Invacare



Return Policy: 30-day return policy (as per Amazon listing at time of purchase).



© 2024 Invacare. All rights reserved.

Related Documents - Essential SX3 Manual Wheelchair (ISX366DASA)

	<p>Invacare Kuschall Champion Folding Wheelchair - Configuration and Specifications</p> <p>Detailed product information, configuration options, and specifications for the Invacare Kuschall Champion folding wheelchair. Explore frame types, seating adjustments, backrest options, wheel choices, and accessories to customize your mobility solution.</p>
	<p>Invacare Kuschall Champion 2.0 Manual Active Wheelchair - Configuration and Price List</p> <p>Comprehensive configuration guide and price list for the Invacare Kuschall Champion 2.0 manual active wheelchair. Explore a wide range of customizable options for frame, seating, backrests, wheels, accessories, and more.</p>
	<p>Invacare Storm 4 X-plore Modulite Powered Wheelchair - Configuration and Price List</p> <p>Detailed configuration options, features, and pricing for the Invacare Storm 4 X-plore Modulite powered wheelchair. Explore seating, drive systems, accessories, and more.</p>
	<p>Invacare Action NG Wheelchair - Features, Specifications, and Options</p> <p>Detailed information on the Invacare Action NG wheelchair, including features, product characteristics, available options, technical specifications, and pricing. This guide covers various configurations, seating options, armrests, footrests, and accessories to customize the wheelchair.</p>

 <p>Invacare Kuschall Champion SK 2.0 Prescription Form and Specifications</p>	<p>Invacare Kuschall Champion SK 2.0 Prescription Form and Specifications</p> <p>Detailed prescription form and specifications for the Invacare Kuschall Champion SK 2.0 manual wheelchair, including features, dimensions, options, and pricing.</p>
 <p>Invacare Essentials EX1 Wheelchair User Manual</p>	<p>Invacare Essentials EX1 Wheelchair User Manual</p> <p>Comprehensive user manual for the Invacare Essentials EX1 Wheelchair. Covers safety guidelines, operating instructions, maintenance procedures, specifications, and warranty information for model numbers IEX166DASA, IEX186DASA, IEX106DASA, IEX166DAEL, IEX186DAEL, and IEX106DAEL.</p>