



Manuals.plus /

› GYENNO /

› GYENNO Parkinson Spoon and Fork Set for Elderly, Parkinson Utensils for Hand Tremor, Adaptive Utensils with Real-Time Tremor Assist User Manual

GYENNO B0F131QL1Z

GYENNO Parkinson Spoon and Fork Set User Manual

Model: B0F131QL1Z

Adaptive Utensils with Real-Time Tremor Assist for Elderly

1. INTRODUCTION

The GYENNO Parkinson Spoon and Fork Set is designed to assist individuals experiencing hand tremors, enabling them to dine independently and with confidence. This innovative utensil set utilizes real-time stabilization technology to adapt to unique tremor patterns, providing smoother eating experiences. The set includes interchangeable spoon and fork attachments to cater to various dining needs.



Figure 1: Overview of the GYENNO Parkinson Spoon and Fork Set.

2. SETUP AND INITIAL USE

2.1 Initial Charging

Before first use, the device must be charged for at least 3 hours. A charging cable is included in the package. Locate the Type-C charging port at the bottom of the device. Plug the Type-C cable into the device and connect the adapter to a power outlet.

- When the indicator light is **constantly yellow**, the device is charging.
- When the indicator light is **constantly green**, the device is fully charged.

Your browser does not support the video tag.

Video 1: Demonstrates initial charging and strap attachment. Ensure the device is fully charged before first use.

2.2 Attaching the Wrist Strap

Thread the wrist strap through the small hole located at the bottom of the GYENNO BRAVO Classic handle. This


strap helps prevent accidental drops during use.

3. OPERATING INSTRUCTIONS

3.1 Powering On/Off

Before turning on the device for the first time, ensure the charging cable is unplugged. To power on the device, press and hold the On/Off button for 2 seconds. The device will activate its tremor-assist function. To power off after use, press and hold the On/Off button for 1 second.

Extremely Easy to Use

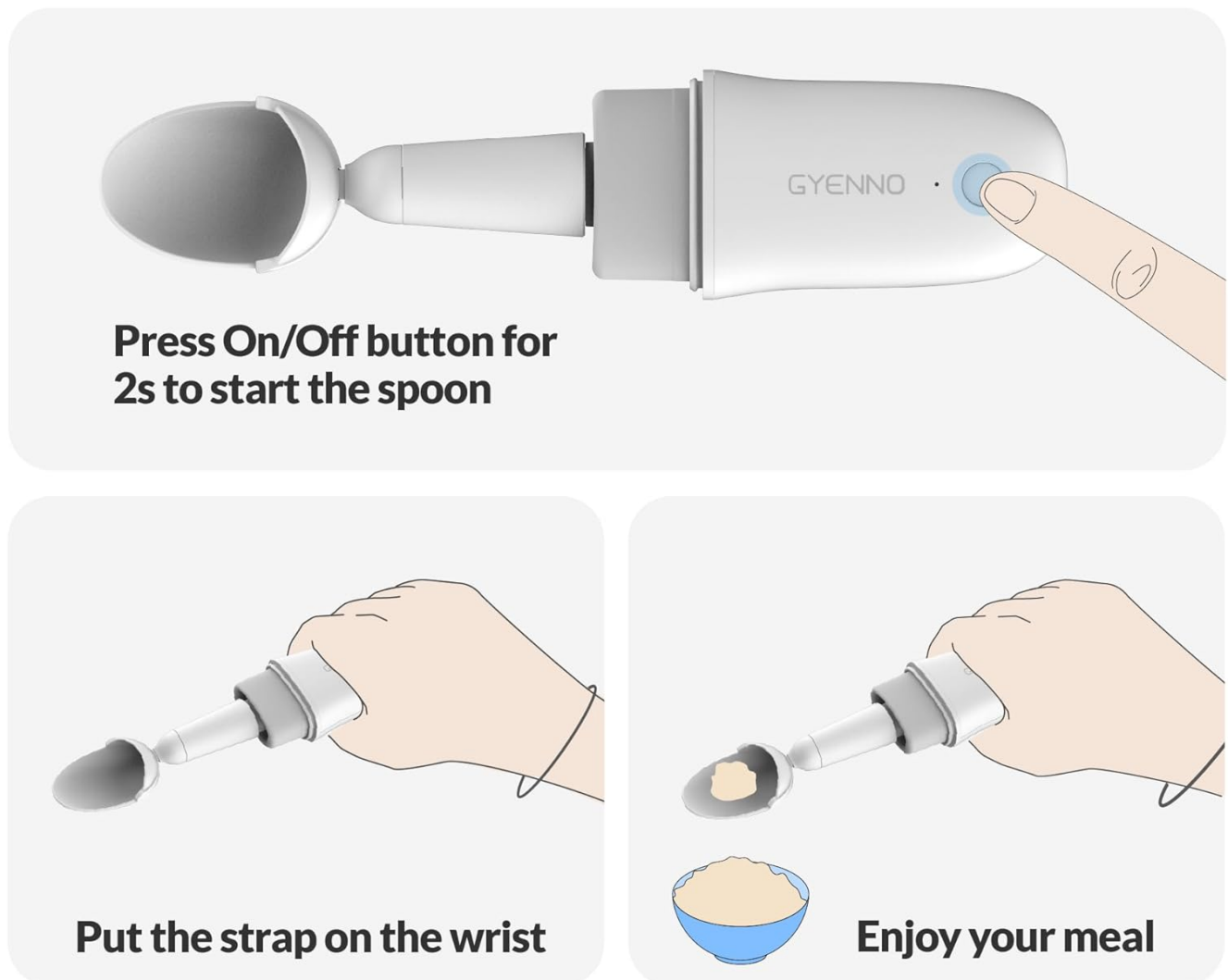


For additional instructions, please refer to the user manual and quick start guide found in the online document located at the center of the webpage.

Figure 2: Pressing the On/Off button to activate the device.

3.2 Using the Utensil

Once powered on, remove the device cover and place the wrist strap securely on your wrist. The real-time stabilization technology will begin to adjust to your tremor patterns. Simply use the spoon or fork as you normally would for eating. The device is designed to provide smoother and more controlled movements.

Hold It Any Way You Want

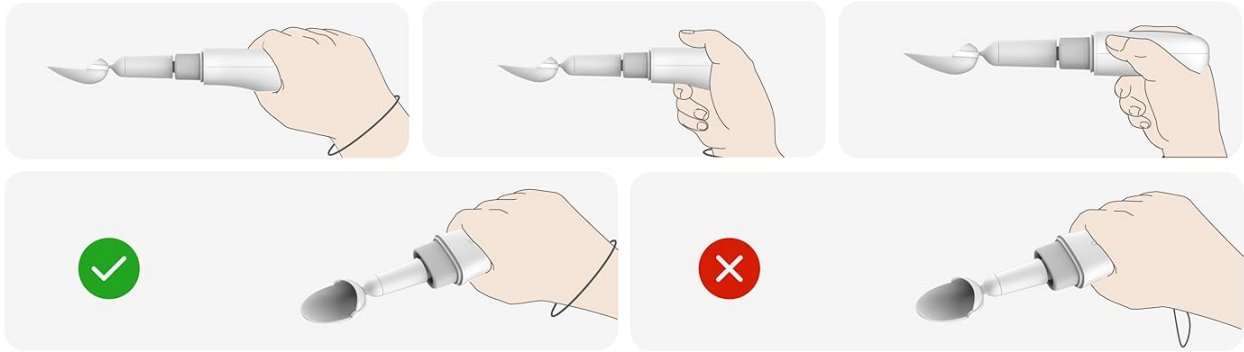

CAUTION: To prevent damage by accidental drops, make sure the strap is attached to your wrist when you use it.

Figure 3: The ergonomic design allows for various comfortable grips. Always use the wrist strap to prevent accidental drops.

4. REPLACING HEAD ATTACHMENTS

The GYENNO BRAVO Classic comes with two replaceable head attachments: a spoon and a fork. These can be interchanged based on your dining needs.

4.1 How to Replace

1. Before replacing an attachment, ensure the device is powered off by holding the On/Off button for 1 second.
2. The attachments connect magnetically to the handle. To remove an attachment, gently pull it away from the handle.
3. To attach a new head, align it with the magnetic connection point on the handle and push it into place until it clicks securely.
4. Once the replacement is complete, press the power button for 2 seconds to restart the device.

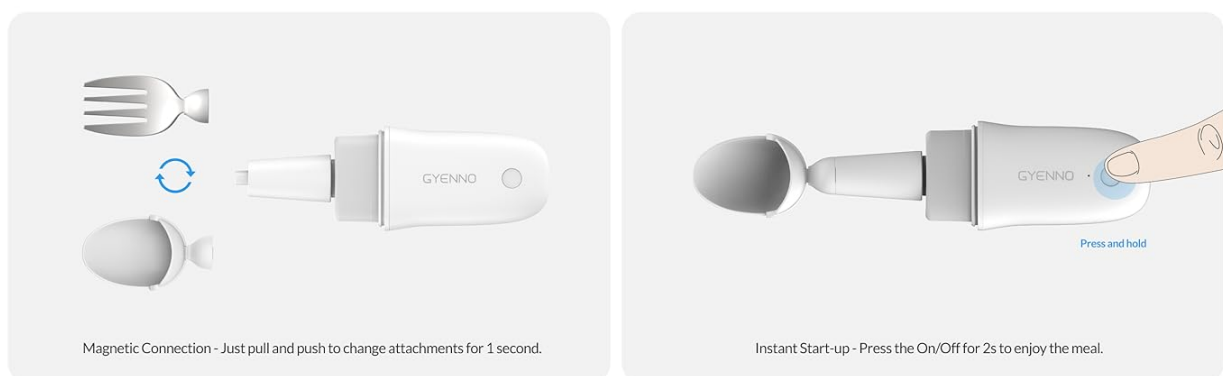


Figure 4: Magnetic connection allows for quick and easy swapping of spoon and fork attachments.

5. CLEANING AND STORAGE

5.1 Cleaning

After each use, ensure the device is powered off. The head attachments (spoon and fork) are dishwasher safe and can be cleaned with water. To clean stains on the handle part, use a moist towelette or alcohol pad to wipe, then dry

with a clean cloth.

- **IMPORTANT:** Make sure the head attachment has been removed from the handle before cleaning.
- **DO NOT** splash or immerse the handle part in water.
- **NEVER** clean the handle part with a dishwasher.
- **STRICTLY PROHIBIT** high-temperature sterilization or putting the handle part in the microwave.

Your browser does not support the video tag.

Video 2: Instructions for cleaning and proper storage of the utensil.

5.2 Storage

For daily storage, ensure both the head attachments and the handle are cleaned and completely dry before storing. Always place the protective cover on the device when not in use.

- Do not store the device in high temperature or high humidity environments.
- For prolonged periods of non-use, ensure the device is fully charged before storage, and charge the device each month to maintain battery health.

6. USABILITY TEST (OPTIONAL)


To evaluate the usability of the GYENNO BRAVO series more accurately for your specific needs, you can participate in a free usability testing program via the 'GYENNO SPOON' mobile application.

6.1 Using the 'GYENNO SPOON' App


1. Download and install the 'GYENNO SPOON' APP from Google Play or the App Store.
2. Register a new account and log in.
3. Enter the usability test section (either by clicking 'Test' on the main screen or navigating through 'Mine' -> 'Usability Test').
4. Read the test instructions carefully before proceeding. It is recommended to test during your daily mealtime for accurate results.
5. Select your dominant hand for dining and learn the test movement as instructed by the app.
6. Follow the on-screen prompts to perform the test. The test typically lasts for 30 seconds.
7. After the test, the result will be displayed, indicating if the product could improve your independent dining effectively.

How to Know If It Fits You

Step 1. Download 'GYENNO SPOON' APP

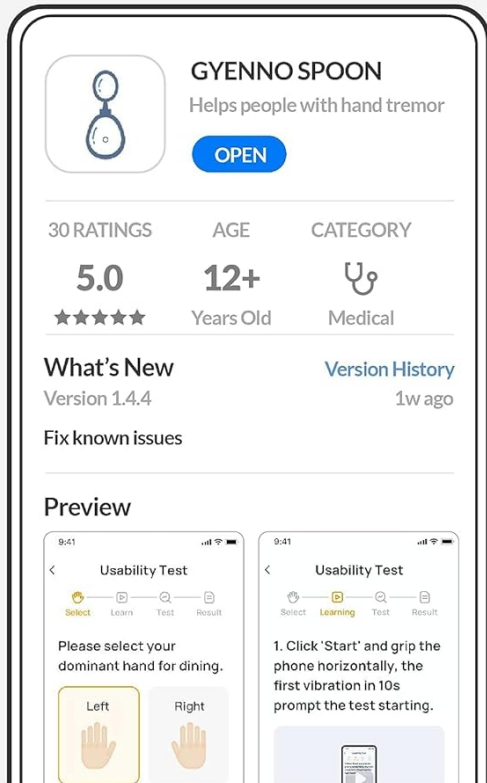

Step 2. Register & Login

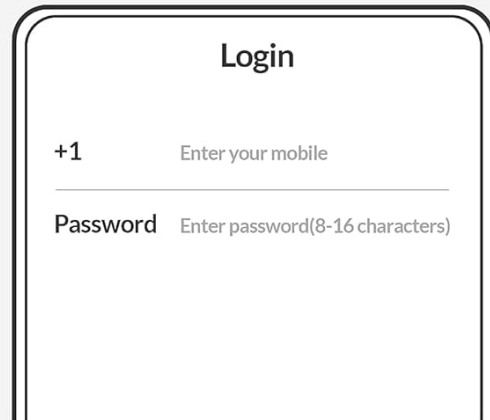

Step 3. Test & Get Results

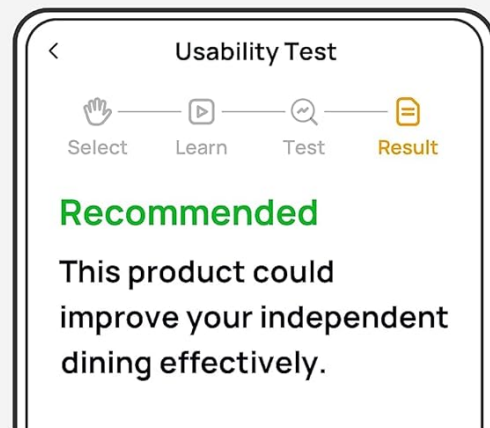


Figure 5: Steps to download and use the 'GYENNO SPOON' app for usability testing.

Your browser does not support the video tag.

Video 3: Guide on how to perform the usability test using the mobile application.

7. PRODUCT SPECIFICATIONS

Feature	Detail
Package Dimensions	9.84 x 3.43 x 2.68 inches
Item Weight	15.84 ounces
Batteries	5 Lithium Ion batteries (included)
Manufacturer	GYENNO
ASIN	B0F131QL1Z
Country of Origin	China

8. TROUBLESHOOTING

If you encounter any issues with your GYENNO Parkinson Spoon and Fork Set, please refer to the following common solutions:

- **Device not turning on:** Ensure the device is fully charged. Press and hold the power button for the specified duration (2 seconds for On).
- **Attachment not staying in place:** Verify that the attachment is properly aligned with the magnetic connection point and pushed firmly until it clicks. Ensure no debris is obstructing the connection.
- **Reduced tremor assistance:** Ensure the device is fully charged. If the issue persists, try restarting the device. Consider performing the usability test via the app to check device performance.
- **Cleaning concerns:** Always remove the head attachment before cleaning. Remember that the handle part is NOT waterproof and should only be wiped with a moist cloth. Do not use a dishwasher or high-temperature sterilization on the handle.

For further assistance, please contact customer support.

9. WARRANTY AND SUPPORT

GYENNO SCIENCE is committed to providing excellent support. While specific warranty details are not provided in this manual, typically products come with a manufacturer's warranty. Please refer to your purchase documentation or the official GYENNO website for detailed warranty information.

If you have any questions, require technical assistance, or need to report an issue, please contact GYENNO customer support. Contact information can usually be found on the product packaging, the official GYENNO website, or through the 'GYENNO SPOON' mobile application.



We Hear and We Care

GYENNO SCIENCE has an efficient and professional customer care team to provide you with support services.

If you have any questions, please feel free to contact us.

Figure 6: GYENNO's commitment to customer support.

10. IMPORTANT INFORMATION

Legal Disclaimer:

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

