

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [LINNER](#) /

> [LINNER Saturn New 2025 OTC Hearing Aids User Manual](#)

LINNER LINNER-SATURN-F703A

LINNER Saturn New 2025 OTC Hearing Aids User Manual

Model: LINNER-SATURN-F703A

INTRODUCTION

The LINNER Saturn New 2025 OTC Hearing Aids are designed to provide enhanced sound clarity and comfort for individuals with mild to moderate hearing loss. This updated model features an American ONSEMI processor, 16-channel digital processing, and improved battery life, offering a discreet and user-friendly experience. This manual provides essential information for the proper setup, operation, and maintenance of your hearing aids.



New 2025 Saturn

VS



Saturn

	New Saturn	Saturn
ONSEMI AMERICAN CHIP	✓	✗
16 WDRC Channels	✓	✗
Feedback shield Technology	Enhanced	Normal
Battery Life	24 hrs	20 hrs

Image: Comparison of the New 2025 Saturn model with the previous Saturn model, illustrating key technological advancements including the ONSEMI American chip, 16 WDRC (Wide Dynamic Range Compression) channels, enhanced feedback shield technology, and an extended battery life of 24 hours compared to 20 hours.

PACKAGE CONTENTS

Upon opening your LINNER Saturn New 2025 Hearing Aids package, you should find the following items:

- Charging Case
- Two (2) Hearing Aids (Left and Right)
- Quick Start Guide
- User Manual (this document)
- Type-C Charging Cable
- Replacement Wax Guards
- Various Sizes of Eartips

Package Contents



Image: All components included in the product package, neatly arranged. This includes the charging case, the two hearing aids, a quick start guide, the full user manual, a USB Type-C charging cable, replacement wax guards, and a selection of eartips in different sizes to ensure a comfortable fit.

SETUP

1. Charging the Hearing Aids

Before first use, ensure your hearing aids are fully charged. The compact charging case provides power for daily use and travel.

1. Place the hearing aids into their respective slots in the charging case. Ensure the charging points on the hearing aids align with the contacts in the case. The hearing aids are magnetically docked for secure placement.
2. Connect the provided USB Type-C charging cable to the charging port on the back of the case.
3. Plug the other end of the cable into a compatible USB power adapter (not included).
4. A green indicator light on the case will illuminate, indicating that charging is in progress. A full charge takes approximately 2 hours and provides up to 24 hours of continuous use.

3 Steps to Proper Charging



Image: Visual guide for charging, showing the correct placement of hearing aids in the case, the magnetic alignment of charging ports, and the meaning of the charging indicator light (green for charging, off for fully charged).

2. Selecting and Inserting Ear Tips


Choose the ear tips that provide the most comfortable and secure fit for your ears. A proper fit is crucial for optimal sound quality and to prevent feedback.

- Select an ear tip size from the provided assortment.
- Gently attach the chosen ear tip to the hearing aid.
- Identify the left (L) and right (R) hearing aids using their color-coded indicators (typically red for right, blue for left).
- Carefully insert each hearing aid into the corresponding ear canal until it feels snug and comfortable.

Effortless to Use

L/R Indicator

Effortless left and right identification

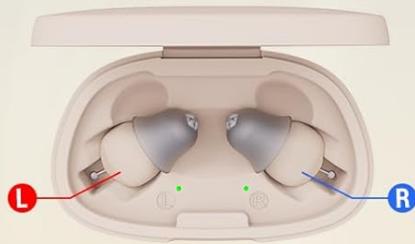

Magnetic Docking

Snap into place for easy charging

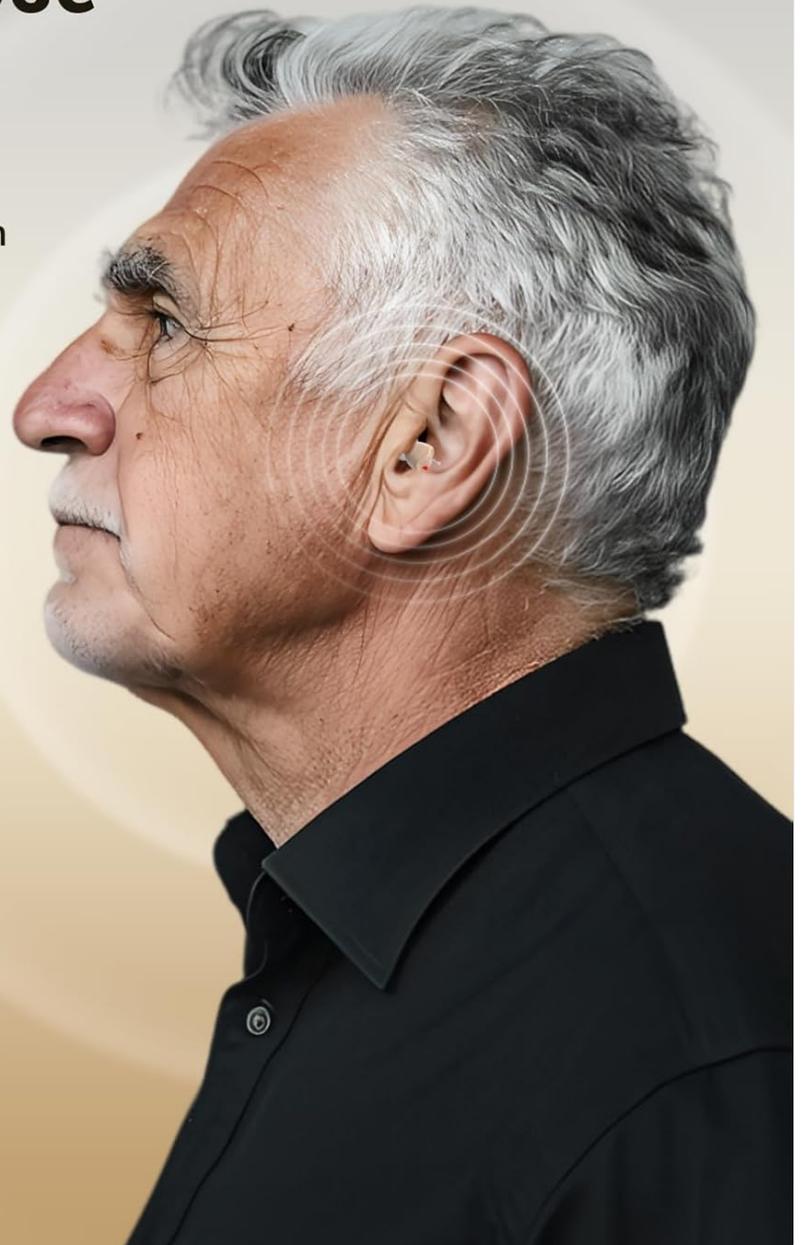


Image: The hearing aids feature clear L/R indicators for correct placement and magnetic docking for effortless charging. The discreet design ensures they are nearly invisible when worn.

3. Wearing the Hearing Aids

Once inserted, the hearing aids should sit discreetly and comfortably within your ear. Their lightweight and ergonomic design ensures a secure fit even during extended wear.

LINNER

Redefining Sound, Restoring Connection

Founded in 2016, LINNER emerged as a pioneering force dedicated to transforming the way the world experiences sound. We are driven by a singular mission: to seamlessly blend state-of-the-art auditory technology with elegant, ergonomic design, empowering you to rediscover the richness of every moment.

LINNER's commitment to excellence is underscored by FDA Approval, unwavering endorsements from over 160 Professional Audiologists, and the glowing testimonials of countless satisfied users. This trust fuels our continuous innovation as we lead the charge in the hearing enhancement landscape.

Image: A close-up view of the hearing aid worn in the ear, showcasing its compact and discreet design.

OPERATING INSTRUCTIONS

1. Power On/Off

The hearing aids automatically power on when removed from the charging case and power off when placed back into the case and the lid is closed. There is an 8-second power-on delay, so ensure you insert them into your ears immediately after removing them from the case.

2. Volume Adjustment

To adjust the volume or clarity, perform a short press on the control button located on the hearing aid. Experiment to find your preferred listening level.

3. Switching Programs/Modes

To switch between different listening programs or modes (e.g., for different environments), press and hold the control button on the hearing aid. Refer to the Quick Start Guide for specific program details.

4. LINNER HA V2.5 App Features

Enhance your experience with the fully upgraded LINNER HA V2.5 app. The app allows for personalized quick hearing screening and easy adjustments tailored to your needs. Please note that music streaming and app control are not supported.

2025 Sound Upgrade: Hear Life Again

Enhanced with ONSEMI chip and 16-channel tech for clearer, cleaner sound.



- 16-Channel DSP System
- ONSEMI AMERICAN CHIP



Image: An illustration detailing the advanced sound features of the 2025 model, such as enhanced noise cancellation, an improved feedback shield, a 16-channel Digital Signal Processing (DSP) system, and the integration of the ONSEMI American chip for clearer sound.

Product Demonstration Video

Your browser does not support the video tag.

Video: A comprehensive demonstration of the LINNEN New Saturn Hearing Aids, including unboxing, setup, charging, and a visual guide on how to wear and operate the devices. This video provides practical insights into the product's features and ease of use.

MAINTENANCE

1. Cleaning the Hearing Aids

Regular cleaning is essential to maintain optimal performance and hygiene. Use the provided wax guard tool to clean any wax buildup.

- Gently wipe the surface of the hearing aids with a soft, dry cloth.
- Use the wax guard tool to carefully remove any earwax or debris from the ear tips and sound outlets.
- Do not use water or cleaning solutions directly on the hearing aids.

2. Storage

When not in use, store the hearing aids in their charging case to protect them from dust, moisture, and physical damage. The case also keeps them charged and ready for use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No Sound / Low Volume	Low battery; blocked ear tip/sound outlet; incorrect volume setting.	Charge the hearing aids; clean ear tip/sound outlet with wax guard tool; adjust volume.

Problem	Possible Cause	Solution
Feedback (Whistling Sound)	Improper fit; earwax buildup; hearing aid not fully inserted.	Re-insert the hearing aid for a snug fit; try a different size ear tip; clean ear canal/hearing aid.
Hearing Aid Not Charging	Improper placement in case; cable not connected; faulty power source.	Ensure proper magnetic alignment; check cable connection; try a different USB power adapter.
App Connectivity Issues	Bluetooth not enabled; app outdated; device not paired.	Ensure Bluetooth is on; update LINNER HA app; re-pair device if necessary.

SPECIFICATIONS

- **Model Name:** LINNER-SATURN-F703A
- **Product Dimensions:** 0.8 x 0.5 x 0.1 inches; 0.24 ounces
- **Battery Type:** 3 Lithium Ion batteries (included)
- **Battery Life:** Up to 24 hours per charge
- **Charging Time:** Approximately 2 hours for full charge
- **Material:** Acrylonitrile Butadiene Styrene (ABS)
- **Color:** Beige
- **Style:** Discreet, In-Ear
- **Processor:** ONSEMI American processor
- **Channels:** 16-channel digital processing
- **Noise Cancellation:** Advanced noise cancellation
- **Feedback Suppression:** Enhanced feedback shield technology
- **App Compatibility:** LINNER HA V2.5 app (Note: Music Streaming / App Control Not Supported)
- **Manufacturer:** AUSTAR HEARING SCIENCE & TECHNOLOGY (XIAMEN) CO., LTD

WARRANTY AND SUPPORT

For any questions, concerns, or assistance with your LINNER Saturn New 2025 Hearing Aids, please utilize the available customer support channels.

- **Online Guidance:** Access detailed guides and FAQs through the LINNER HA V2.5 app.
- **Online Customer Service:** The LINNER customer support team is available 24/7 to assist you. Refer to the Quick Start Guide or the LINNER HA app for contact information.
- **Warranty:** Please refer to the product packaging or the official LINNER website for detailed warranty information.