

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

- › [BEBIRD](#) /
- › [BEBIRD Earsight Plus Ear Wax Removal Tool Camera User Manual](#)

## BEBIRD BP033\_BE\_ZE

# BEBIRD Earsight Plus Ear Wax Removal Tool Camera User Manual

Brand: BEBIRD | Model: BP033\_BE\_ZE

## INTRODUCTION TO YOUR EARSIGHT PLUS DEVICE

The BEBIRD Earsight Plus is an advanced ear care device designed for safe and effective earwax removal and general ear, nose, and throat (ENT) examinations. It features a high-definition visual otoscope with smart WiFi connectivity, allowing you to view clear, real-time visuals of the ear canal directly on your smartphone or tablet. The device comes with a variety of safe, flexible silicone scoops and tweezers, ensuring a comfortable and precise cleaning experience for both adults and children.

Its innovative design includes an adjustable-angle earpick handle for enhanced flexibility, a smart tweezer functionality for precise wax removal, and replaceable earpicks and ENT accessories for versatile use. Safety features such as a spiral cleaning head that remains securely fixed and a built-in safety scale with a protective stopper prevent over-insertion, safeguarding the eardrum during use.



Image: The BEBIRD Earsight Plus device, its charging case, various ear tips, and a visual representation of its app interface.

## SETUP GUIDE

### 1. Unboxing and Initial Charge

Upon receiving your BEBIRD Earsight Plus, carefully unbox all components. The package includes the main device, a charging case, various silicone ear tips, ENT accessories, and a USB-C charging cable. Before first use, ensure the device is fully charged. Connect the USB-C cable to the charging port on the device or its case and plug it into a compatible USB power source.

## What's in the Box



Image: A detailed view of all items included in the BEBIRD Earsight Plus package, including the device, case, and various tips.

## 2. App Download and Connection

To utilize the visual features of your Earsight Plus, download the official BEBIRD app from your device's app store (available for iOS and Android). Once installed:

1. Turn on the BEBIRD Earsight Plus device by pressing the power button. A light will illuminate on the device.
2. Open the BEBIRD app on your smartphone or tablet.
3. Follow the in-app instructions to connect the device via WiFi. The app will guide you through selecting the BEBIRD device's WiFi network from your phone's settings.
4. Once connected, you will see a real-time view from the device's camera on your screen.



The App 'Bebird' can work with iPhone, iPad and Android smart phones/tablets.

Image: Screenshots of the BEBIRD app on a smartphone and tablet, demonstrating the connection process.

## 3. Attaching and Changing Tips

The Earsight Plus comes with various interchangeable tips for different cleaning and inspection needs. To attach a tip, gently twist it onto the end of the device until it is secure. To remove, twist in the opposite direction. Ensure the tip is firmly attached before use.



### Four Classic Cleaning Tips

Dry ear, oil ear, ear skin, and ear debris can be easily removed




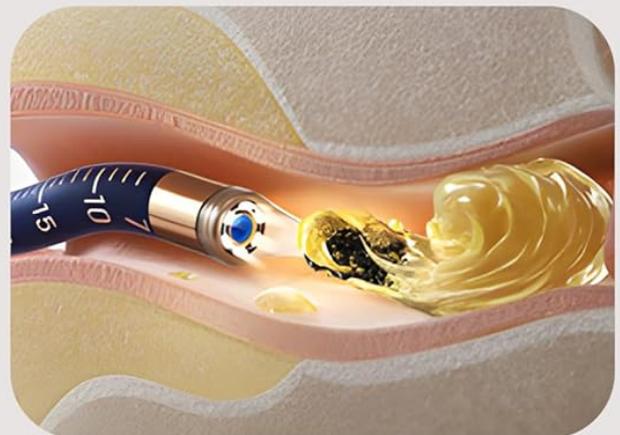


Image: Illustration of four classic cleaning tips designed for removing dry earwax, oily earwax, ear skin, and ear debris.

## OPERATING INSTRUCTIONS

### 1. Earwax Removal

For earwax removal, select the appropriate silicone scoop or tweezer tip. Gently insert the device into the ear canal while observing the live feed on your connected device. The adjustable-angle earpick handle allows for precise maneuvering to match the natural curves of your ear canal. Use the smart tweezer functionality to safely grip and remove stubborn earwax.

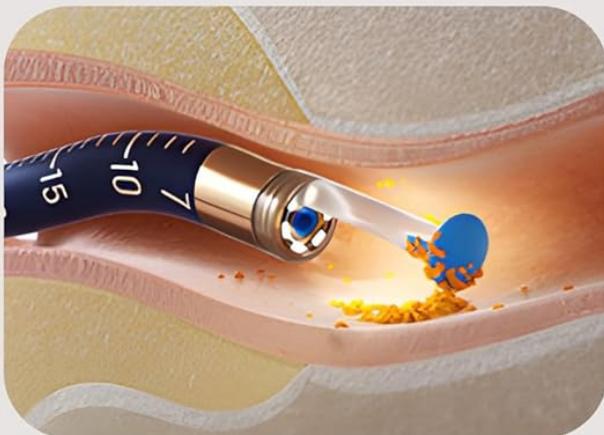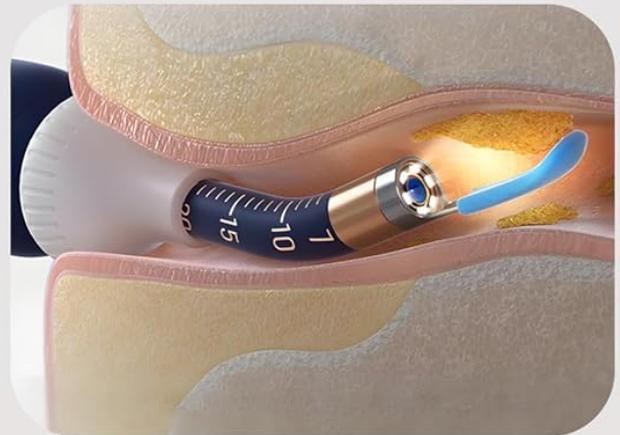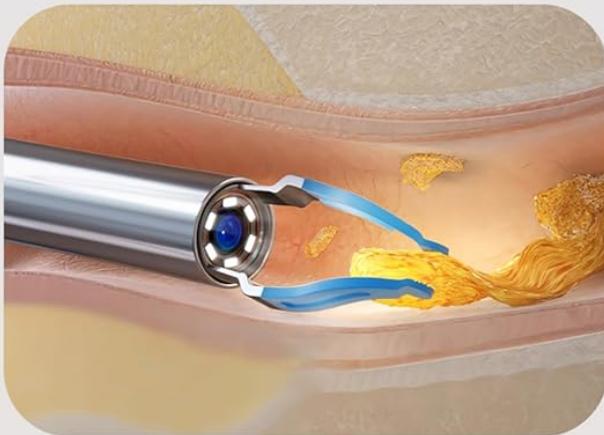

### **Tweezers**

Remove stubborn earwax



### **Scooper**

Removes regular earwax

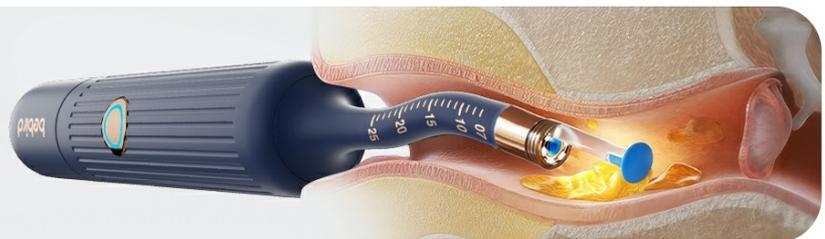


Image: Visual guide demonstrating the use of the tweezer tip for stubborn earwax and the scooper tip for regular earwax removal.



### **Scraper**

Extracts debris



### **Sponge**

Removes wet earwax



Image: Visual guide demonstrating the use of the scraper tip for debris extraction and the sponge tip for wet earwax removal.

## **2. Safety Features**

The device incorporates several safety features to prevent injury:

- **Safety Stopper:** An adjustable stopper can be placed on the earpick to prevent accidental over-insertion into the ear canal, protecting the eardrum.
- **Secure Tips:** The spiral cleaning head and other tips are designed to remain securely fixed during use, eliminating the risk of detachment.
- **Temperature Control:** A dual temperature control system ensures the device remains at a comfortable temperature during use.

### Safety Stopper Keeps It Safe

The adjustable topper can prevent the ear pick tool from accidentally going too deep into the ear canal, making it safe for use with your family.



Image: Illustration showing how the adjustable safety stopper prevents the ear pick tool from going too deep into the ear canal.

### 3. Versatile Home Check-Ups

Beyond earwax removal, the BEBIRD Earsight Plus can be used for general inspections of the ear, nose, and throat. Attach the appropriate ENT accessories for a clear visual of these areas, aiding in personal health monitoring.

# Home Check-Up Must-Have

Check ears, nose, tooth and throat for your family



Image: Demonstrations of the BEBIRD Earsight Plus being used for ear, nasal, and oral examinations, showing the corresponding app views.

## MAINTENANCE AND CARE

### 1. Cleaning the Device and Tips

After each use, it is crucial to clean the silicone tips and the camera lens. Gently remove the used tip and clean it with a soft cloth and mild disinfectant. For the camera lens, use a soft, dry cotton swab or lens cleaning cloth to wipe away any debris. Do not use abrasive cleaners or immerse the main device in water.

## Clean and Hygienic Storage

All-in-one, portable and convenient to use



Image: The BEBIRD Earsight Plus device and its accessories neatly stored in its portable and hygienic charging case.

## 2. Storage

Store the BEBIRD Earsight Plus and its accessories in the provided charging case when not in use. This keeps the components organized, protected, and hygienic. Store in a cool, dry place away from direct sunlight and extreme temperatures.

## TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. Connect it to the charger and wait for the charging indicator to show a full charge.
- **Cannot connect to WiFi:**
  - Make sure the device is powered on.
  - Check your phone's WiFi settings to ensure you are connected to the BEBIRD device's specific WiFi network.
  - Restart both the BEBIRD device and your smartphone/tablet.
  - Ensure the BEBIRD app has necessary permissions (e.g., location, local network access) on your phone.
- **Unclear image/Poor visibility:**
  - Clean the camera lens carefully with a soft, dry cloth.
  - Ensure there is sufficient lighting in your environment, although the device has a built-in light.
  - Check if the tip is properly attached and not obstructing the camera view.
- **Tips not staying on:** Ensure the tips are twisted on firmly until they click into place. If a tip is damaged, replace it with a new one.

## PRODUCT SPECIFICATIONS

<b>Model Number</b>	BP033_BE_ZE
<b>Package Dimensions</b>	7.91 x 5.28 x 1.46 inches; 10.05 ounces
<b>Batteries</b>	1 Lithium Ion batteries required (included)
<b>Manufacturer</b>	Heifeng Zhizao (Shenzhen) Technology Co., Ltd
<b>Country of Origin</b>	China
<b>Camera Resolution</b>	1080P HD
<b>Connectivity</b>	Smart WiFi (iOS & Android compatible)
<b>Lens Diameter</b>	4.1mm micro diameter
<b>Battery Life</b>	Up to 90 minutes of use, up to 30 days standby

## WARRANTY AND SUPPORT

The BEBIRD Earsight Plus typically comes with a 1-year warranty from the date of purchase, covering manufacturing defects. For specific warranty terms and conditions, please refer to the warranty card included in your product packaging or visit the official BEBIRD website.

For technical support, troubleshooting assistance, or inquiries regarding your device, please contact BEBIRD customer service. Contact information (email, phone, or support portal) can usually be found on the product packaging, the official BEBIRD website, or within the BEBIRD mobile application.

Always ensure you are using official BEBIRD accessories and following the instructions in this manual to maintain product validity and safety.

© 2025 BEBIRD. All rights reserved.

This manual is for informational purposes only. Product specifications are subject to change without notice.