

VOCIC AX15

VOCIC AX15 Electric Floor Lift Instruction Manual

Model: AX15 | Brand: VOCIC

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your VOCIC AX15 Electric Floor Lift. Please read all instructions carefully before assembly and operation. This device is designed to assist individuals with limited mobility in transitioning from the floor to a seated position, supporting up to 440 pounds.

2. WHAT'S INCLUDED

The VOCIC AX15 Electric Floor Lift package includes the following components:

- AX15 Main Unit x 1
- Long Base Foot x 2
- Short Base Feet x 1
- Charger x 1
- Hand Controller x 1
- Cushion x 1

What's Included?



AX15 x 1



Long Base Foot x 2



Charger x 1



Short Base Feet x 1



Hand Controller x 1

Image: All components included in the VOCIC AX15 package. This includes the main lift unit, various base feet, a charger, a hand controller, and a cushion.

3. SETUP AND ASSEMBLY

The VOCIC AX15 features an upgraded one-piece construction for simplified setup. The unit arrives largely pre-assembled. Follow these steps for final assembly:

1. **Unpack the Unit:** Carefully remove all components from the packaging.
2. **Unfold the Chair:** Flip the seat up and rotate the backrest into place until it locks securely.
3. **Attach Base Rails:** Attach and lock the base rails underneath the main unit. Ensure they are securely fastened.
4. **Connect Hand Controller:** Connect the hand controller to the designated port on the chair.
5. **Charge Battery:** Connect the charger to the hand controller and a power outlet. Ensure the battery is fully charged before first use.

Details



Image: The VOCIC AX15 in both its upright and folded configurations, highlighting its compact design and dimensions for storage.

The upgraded design allows for one-piece storage without disassembly of the feet, making it quick to set up and put away.



Image: Visual comparison illustrating the AX15's convenient one-piece storage versus a competitor's model requiring disassembly.

Setup Video

Your browser does not support the video tag.

Video: This video demonstrates the quick and easy setup process of the VOCIC AX15 Electric Floor Lift, showing how to unfold the unit and prepare it for use.

4. OPERATING INSTRUCTIONS

The VOCIC AX15 is designed for effortless operation using the hand controller. It allows for smooth and controlled lifting and lowering.

1. **Positioning:** Position the chair near the individual requiring assistance. The optimized wheel design allows for easier movement across various surfaces.
2. **Assisted Transfer to Seat:** Gently assist the individual to slide onto the seat of the lift chair while it is in its lowest position.
3. **Secure with Seatbelt:** Fasten the safety belt around the individual's waist to ensure secure transfers.
4. **Raise the Chair:** Use the hand controller to slowly raise the seat to the desired height. The ergonomic armrests provide additional support during the lift.
5. **Lower the Chair:** To lower, use the hand controller to gently descend the seat.

Operation Video

Your browser does not support the video tag.

Video: This video demonstrates the safe and controlled operation of the VOCIC AX15 Electric Floor Lift, showing how to assist an individual from the floor to a seated position.

The chair's movement is slow and steady, providing a sense of security and reducing the risk of sudden movements.

5. SAFETY GUIDELINES

- Always ensure the base feet are fully extended and locked before use to maintain stability.
- Use the safety belt for all transfers to prevent falls or sliding.
- Ensure the individual's feet are flat on the floor or supported during transfers.
- Do not exceed the maximum weight recommendation of 440 pounds.
- Operate the chair on a level, stable surface.
- Keep the hand controller dry and away from water.
- Regularly check all connections and moving parts for wear or damage.

Upgraded Stabilized Base



Image: The AX15 features an upgraded T-stabilized carbon steel base, designed for enhanced stability and reduced risk of tipping compared to less stable structures.

6. MAINTENANCE

- **Cleaning:** Wipe down the chair's surfaces with a damp cloth and mild detergent. Avoid abrasive cleaners.
- **Battery Care:** Recharge the battery regularly, even if the chair is not in frequent use, to prolong battery life.
- **Inspection:** Periodically inspect all bolts, connections, and moving parts for tightness and proper function.
- **Storage:** When not in use, fold the chair for compact storage in a dry, safe place.

7. TROUBLESHOOTING

- **Chair Not Responding:** Check if the hand controller is securely connected and if the battery is charged.
- **Unstable Movement:** Ensure all base feet are fully extended and locked. Verify the chair is on a flat, stable surface.
- **Unusual Noises:** Discontinue use immediately and contact customer support.

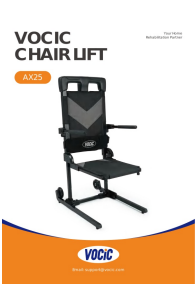
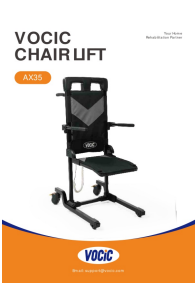
8. SPECIFICATIONS

Brand	VOCIC
Model Name	AX15
Color	Orange
Product Dimensions	16.5"D x 17.3"W x 22.4"H
Maximum Weight Recommendation	440 Pounds
Item Weight	45 pounds
Material	Carbon Steel
Fill Material	Memory Foam
Indoor/Outdoor Usage	Indoor
UPC	810190250044

9. WARRANTY AND SUPPORT

The VOCIC AX15 Electric Floor Lift comes with an industry-leading 1-year warranty. For any setup assistance, customization advice, or technical support, please contact our care team. Support is available 24/7.

Related Documents - AX15

	<p>VOCIC AX25 Chair Lift User Manual and Technical Specifications</p> <p>Comprehensive user manual and technical specifications for the VOCIC AX25 Chair Lift, designed for home rehabilitation and mobility assistance. Includes safety instructions, operating procedures, troubleshooting, and warranty information.</p>
	<p>VOCIC AX35 Chair Lift User Manual and Product Information</p> <p>Comprehensive guide for the VOCIC AX35 Chair Lift, covering safety instructions, product information, installation, operation, troubleshooting, and after-sales service. Designed for elderly and disabled individuals to assist with safe lifting and lowering.</p>
	<p>VOCIC Rollator Z51 Instruction Manual</p> <p>Comprehensive instruction manual for the VOCIC Rollator Z51, detailing safety precautions, product information, assembly, operation, and after-sales service.</p>



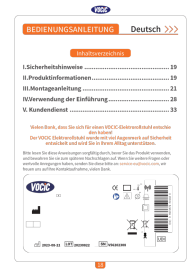
[VOCIC RLD-100W03B Electric Wheelchair Instruction Manual](#)

Comprehensive instruction manual for the VOCIC RLD-100W03B Electric Wheelchair, covering safety, product information, assembly, operation, charging, maintenance, and warranty.



[VOCIC Electric Wheelchair RLD-100W03B User Manual and Technical Specifications](#)

This document provides comprehensive instructions and technical specifications for the VOCIC Electric Wheelchair, model RLD-100W03B. It covers safety guidelines, product details, assembly, operation, charging, maintenance, and customer support.



[VOCIC Electric Wheelchair User Manual](#)

User manual for the VOCIC Electric Wheelchair Model RLD-100W01A, SKU V-EW-V56-E-SV. Includes safety instructions, product information, assembly, usage, charging, maintenance, and customer service details.