

TTQ FAID

TTQ FAID Hearing Aids User Manual

Model: FAID | Brand: TTQ

INTRODUCTION

The TTQ FAID Hearing Aids are advanced over-the-counter (OTC) devices designed to enhance auditory experience for individuals with moderate hearing loss. Equipped with intelligent DSP chip technology, these hearing aids not only amplify sound but also intelligently identify and enhance human voice while reducing environmental noise. They also function as wireless Bluetooth headphones, allowing for private audio entertainment and clear phone calls. The compact and ergonomic design ensures comfortable long-term wear, and the rechargeable charging box includes a dehumidification function for extended device longevity.

WHAT'S IN THE BOX

- TTQ FAID Hearing Aids (1 pair)
- Charging Box
- User Manual
- Charging Cable
- Ear Domes (6 pairs of different sizes: S, M, L, XL)

Help You Hear the World Better.

Suitable for different scenarios.



Conversation



Meeting



Natural



Sports

Figure 1: Package contents of the TTQ FAID Hearing Aids. This image displays the two hearing aids, their charging case, the instruction manual, a charging cable, and various sizes of ear domes (S, M, L, XL) to ensure a comfortable fit.

PRODUCT FEATURES

- **Clear and Natural Sound:** Built-in intelligent DSP chip enhances human voice and lowers environmental noise, reducing whistling and squeaking for a clearer hearing experience.
- **Wireless Bluetooth Headphones:** Doubles as Bluetooth earbuds for private entertainment and crystal-clear phone calls with intelligent noise cancellation.
- **Compact Design & Comfortable Wearing:** Ergonomically designed, each hearing aid weighs only 0.01 lb, reducing burden during long-term wear. Includes 6 ear domes for a better fit.
- **Longer Lasting & Dehumidification Charging Box:** Provides 18 hours of hearing aid use or 13 hours of Bluetooth headphone use with a 2-hour quick charge. The charger includes overcharge protection and a dehumidification function to extend device life. An LED display shows remaining power.
- **Easy to Use & Customized Settings in APP:** Simple one-button control for volume and noise-cancellation modes.

A dedicated app allows for a hearing test and automatic adjustment of equalizer settings for personalized sound.

Hear What Matters. Hear Everything Better.

The advertisement features a central image of an elderly man and woman smiling and talking. The man is wearing a black behind-the-ear hearing aid. Overlaid on this image is a diagram illustrating the 'Smart Noise Reduction' feature. The diagram is divided into three sections: 'Before', 'Smart Noise Reduction', and 'After'. In the 'Before' section, a large, complex sound wave is shown. In the 'Smart Noise Reduction' section, a hearing aid is shown with a smaller, more focused sound wave passing through it. In the 'After' section, the sound wave is significantly reduced in amplitude, indicating noise reduction. Below the diagram are four icons representing key features: an advanced chip, three different modes, less distortion, and less whistling.

Before **Smart Noise Reduction** **After**

Advanced Chip 3 Different Modes Less Distortion Less Whistling

Amplify human voice(max 35db) and reduce environment noise.

Figure 2: The TTQ FAID hearing aids provide clear and natural sound, enhancing the user's ability to hear conversations and environmental sounds.

Long Duration & Rechargeable Charging Box



Figure 3: The rechargeable charging box offers extended battery life and includes a dehumidification port to protect the hearing aids.

Easy to Use. Comfortable to Wear.



Figure 4: Detailed view of the hearing aid's components, emphasizing its ergonomic design and user-friendly features.

SETUP INSTRUCTIONS

- 1. Charge the Device:** Before first use, place both hearing aids into the charging box. Connect the charging cable to the box and a power source. The LED display on the charging box will indicate the battery level. Ensure a full charge for optimal performance.
- 2. Select Ear Domes:** Choose the ear dome size (S, M, L, XL) that provides the most comfortable and secure fit in your ear canal. A proper fit is crucial for sound quality and to prevent feedback (whistling).
- 3. Insert Hearing Aids:** Gently insert the hearing aid into your ear, ensuring the earhook rests comfortably over the back of your ear and the ear dome is securely seated in your ear canal.
- 4. Power On:** The hearing aids typically power on automatically when removed from the charging case. If not, refer to the "Operating Instructions" for manual power-on.
- 5. Initial Volume Adjustment:** Start with a low volume setting and gradually increase it until you achieve a comfortable listening level.



Figure 5: The ergonomic earhooks ensure a compact and lightweight fit, designed for comfort and stability during various activities.

OPERATING INSTRUCTIONS

Basic Controls (Function Button)

- **Volume Adjustment:** Short press the function button to cycle through volume levels.
- **Noise-Cancellation Modes:** Long press the function button to switch between different noise-cancellation modes (e.g., quiet, noisy, outdoor).
- **Power On/Off:** The devices power on automatically when removed from the case. To power off, place them back in the charging case.

Bluetooth Connectivity

1. **Pairing Mode:** When removed from the charging case, the hearing aids automatically enter pairing mode.
2. **Connect to Device:** On your smartphone or Bluetooth-enabled device, go to Bluetooth settings and select "F-Aid" from the list of available devices to connect.

3. **Phone Calls:** Once connected, you can answer or reject phone calls directly through the hearing aids using the function button (single tap to answer/end, long press to reject).
4. **Audio Playback:** Stream music, podcasts, or other audio directly to your hearing aids.

Help You Hear the World Better.

Suitable for different scenarios.

The central image shows two black TTQ FAID hearing aids with white accents, positioned in a white circle. Surrounding this central image are four circular inset images, each connected to the center by a thin orange line. The top-left inset shows a young man and an older man looking at a smartphone together, labeled 'Conversation'. The top-right inset shows a woman and a man looking at a tablet together, labeled 'Meeting'. The bottom-left inset shows a man meditating in a lush green field, labeled 'Natural'. The bottom-right inset shows an older man and a woman lifting dumbbells in a gym, labeled 'Sports'.

Conversation

Meeting

Natural

Sports

Figure 6: The TTQ FAID hearing aids offer seamless Bluetooth connectivity for audio streaming and phone calls.

Using the TTQ APP (Optional)

For customized settings and a personalized hearing experience, download the official TTQ app. The app allows you to:

- Conduct a hearing test to determine your specific hearing profile.
- Automatically adjust equalizer settings based on your hearing test results.
- Fine-tune volume and noise-cancellation modes with greater precision.

MAINTENANCE

- **Cleaning:** Regularly clean the ear domes and the body of the hearing aids with a soft, dry cloth. Do not use water or cleaning solutions directly on the device.
- **Charging Case Dehumidification:** The charging box features a dehumidification function. Always store the hearing aids in the charging box when not in use to help vaporize sweat and moisture, extending their lifespan.
- **Storage:** Store the hearing aids and charging case in a cool, dry place away from direct sunlight and extreme temperatures.
- **Battery Care:** Avoid fully depleting the battery frequently. Charge the devices regularly to maintain battery health. The charging box has overcharge protection.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No Sound/Low Volume	Low battery; incorrect volume setting; clogged ear dome.	Charge the hearing aids; adjust volume using the function button or app; clean or replace ear domes.
Whistling/Feedback	Improper fit; ear dome too small; earwax buildup.	Ensure proper insertion; try a larger ear dome size; consult a professional for earwax removal if necessary.
Bluetooth Disconnection	Out of range; interference; pairing issue.	Move closer to the connected device; restart both devices and re-pair.
Charging Issue	Loose connection; dirty charging contacts; faulty cable/adaptor.	Ensure secure connection; clean charging contacts on aids and case; try a different charging cable or adaptor.

SPECIFICATIONS

Feature	Detail
Brand	TTQ
Model Name	FAID
Color	Black
Material	Acrylonitrile Butadiene Styrene (ABS)
Style	Both Ear
Battery Life (Hearing Aid)	Up to 18 hours
Battery Life (Bluetooth Headphones)	Up to 13 hours
Quick Charge	2 hours for full charge
Noise Cancellation	Intelligent DSP chip, 3-level adjustable noise reduction
Bluetooth Version	Bluetooth 5.3 (implied from A+ content image)
Package Dimensions	7.48 x 4.76 x 2.2 inches
Item Weight	0.16 ounces (each hearing aid)

OFFICIAL PRODUCT VIDEO

Your browser does not support the video tag.

Video 1: An honest review of the TTQ OTC Hearing Aids for Seniors, demonstrating their features and performance. This video provides a user's perspective on the product's functionality and benefits.


WARRANTY AND SUPPORT

For warranty information, product support, or any inquiries, please refer to the contact details provided in the included User Manual or visit the official TTQ brand store on Amazon: [TTQ Brand Store](#).

Optional protection plans are available for purchase to extend coverage beyond the standard manufacturer's warranty. Please check the product listing for details on 2-Year and 3-Year protection plans.

© 2025 TTQ. All rights reserved.

Documents - TTQ – FAID

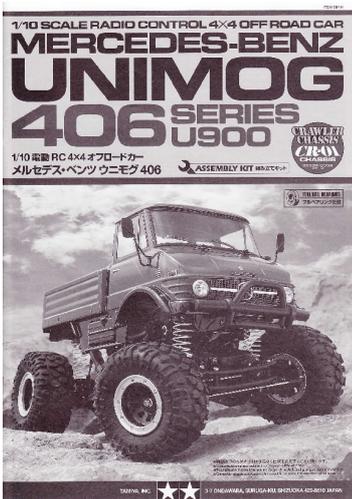

[\[pdf\]](#) User Manual

Page 1 ITEM 58414 10 SCALE RADIO CONTROL 4X4 OFF ROAD CAR MERCEDES BENZ UNIMOG SERIES PU900 EMAKULTET CRAWLER CHASSIS CR 01 DOSYON 110 RC 70 Ft XJLEFZ AW DOES 406 ASSEMBLY KIT REZITEYA 74955 EN Manual lindinger at media 74 07 87 1629195624 1/10s ,e,,ALE BADIO CONTFIOL 4X4 OFF FICAtrt CAFI MEFICEDESi.EIENZ 4oiCi= Eb== 1tr/1lr0ti5tzR dC74 ... 7 r o t. L r j ;:3 :XEP ;;Z: di - t-E4Fatr1.:E* L a7iri hL m r dMor nrmbds Tamry. P.rnr coroE. **Faid** b.dy rro omr aDi*s rrd6 using z eiicnhrsonforgplboltycd.lobonraamte,Ya a c6 ,rsne nd que aErenc... lang:en score:12 filesize: 10.64 M page_count: 28 document date: 2011-05-31

