

## **littneo SF20250706B0F7XZDX36**

# **Littneo Bedside Commode with Wheels, Shower Wheelchair SF20250706B0F7XZDX36 User Manual**

Model: SF20250706B0F7XZDX36

## **1. INTRODUCTION**

---

This manual provides essential instructions for the safe and effective use, assembly, and maintenance of your Littneo Bedside Commode with Wheels, Shower Wheelchair. Please read this manual thoroughly before initial use and retain it for future reference. This product is designed to assist individuals with limited mobility, offering versatility as a bedside commode, shower chair, and transport chair.



*Image 1.1: The Littneo Bedside Commode with Wheels, Shower Wheelchair, shown with its commode bucket and seat cover.*

## 2. SAFETY INFORMATION

Always prioritize safety when using this product. Failure to follow these guidelines may result in injury or damage to the product.

- **Weight Capacity:** Do not exceed the maximum weight capacity of 300 pounds (136 kg).
- **Wheel Locks:** Always engage all four wheel locks when the chair is stationary, especially during transfers, sitting, or standing.
- **Transfers:** Ensure the user is stable and the chair is locked before attempting transfers to or from the commode.
- **Standing:** Do not stand on the commode chair. It is not designed as a standing aid.
- **Obstacles:** When navigating obstacles, pull the chair backward rather than pushing it forward to maintain stability and prevent tipping.
- **Assembly:** Ensure all components are securely assembled according to the instructions before use. Regularly check for loose parts.
- **Supervision:** Users with significant mobility limitations or cognitive impairments should always be

supervised.

## ***How to Unfold?***



***Press Down Easy to Unfold***

*Image 2.1: Visual safety guidelines illustrating correct and incorrect methods for getting up/sitting down and overcoming obstacles.*

### **3. PACKAGE CONTENTS**

Verify that all components are present and undamaged upon opening the package. If any parts are missing or damaged, contact customer support immediately.

- Main Frame (pre-assembled sections)
- Backrest
- Seat Cushion with U-shaped opening
- Footrest Bar
- Commode Bucket with Lid
- Wheel Assemblies (4 with locks)
- Fasteners (wing nuts, screws, rubber plugs)
- Instruction Manual

## 4. SETUP AND ASSEMBLY

Assembly typically takes approximately 10 minutes and requires no tools. Follow these steps carefully:

1. **Unpack Components:** Carefully remove all parts from the packaging.
2. **Attach Backrest:** Insert the backrest into the designated slots on the main frame. Secure it using the provided wing nuts.
3. **Install Seat Cushion:** Place the seat cushion onto the frame. Ensure it is properly aligned.
4. **Attach Footrest:** Slide the footrest bar into the lower front openings of the frame until it clicks into place.
5. **Insert Commode Bucket:** Slide the commode bucket into the rails beneath the seat. Ensure the handle is positioned correctly for stability.
6. **Secure Rubber Plugs:** If switching between caregiver-assisted and self-operated modes, use the small rubber plugs to cover any unused screw holes for safety and to prevent rust. Store unused plugs in a secure location.

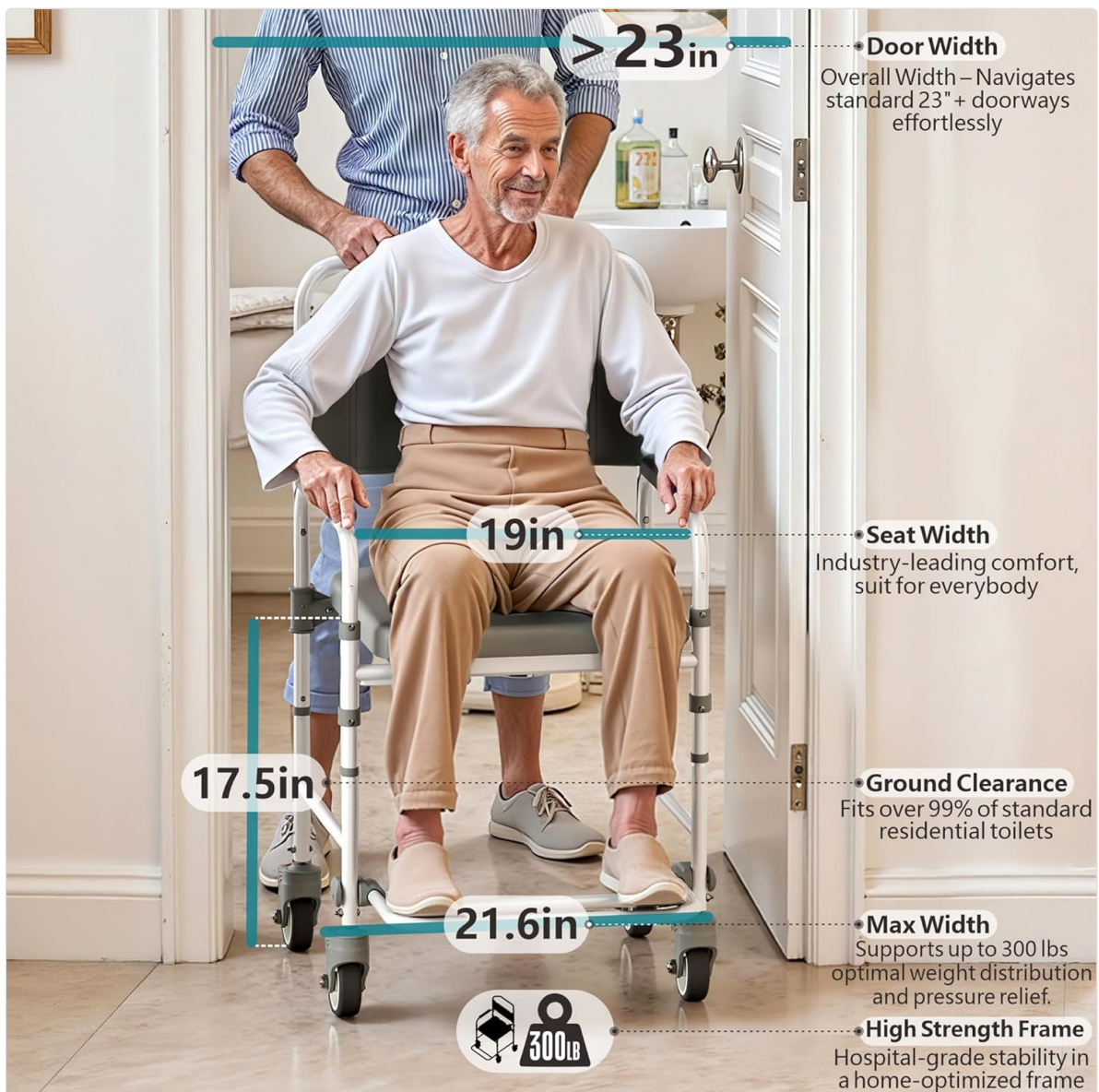


Image 4.1: Illustration of the commode bucket being inserted and the 360-degree swivel wheels with foot-lever locks.

## 5. OPERATING INSTRUCTIONS

## 5.1 Dual-Operation Modes

The commode chair features two operational modes, easily switchable by reorienting the backrest and seat cushion:

- **Caregiver-Assisted Mode:** The cushion opening faces rearward, providing easier access for caregivers to assist with hygiene. This mode is ideal for individuals requiring full-time caregiver support.
- **Self-Operated Mode:** The cushion opening faces forward, allowing for independent use and self-cleaning. This mode is suitable for individuals with limited mobility who can manage personal hygiene independently.

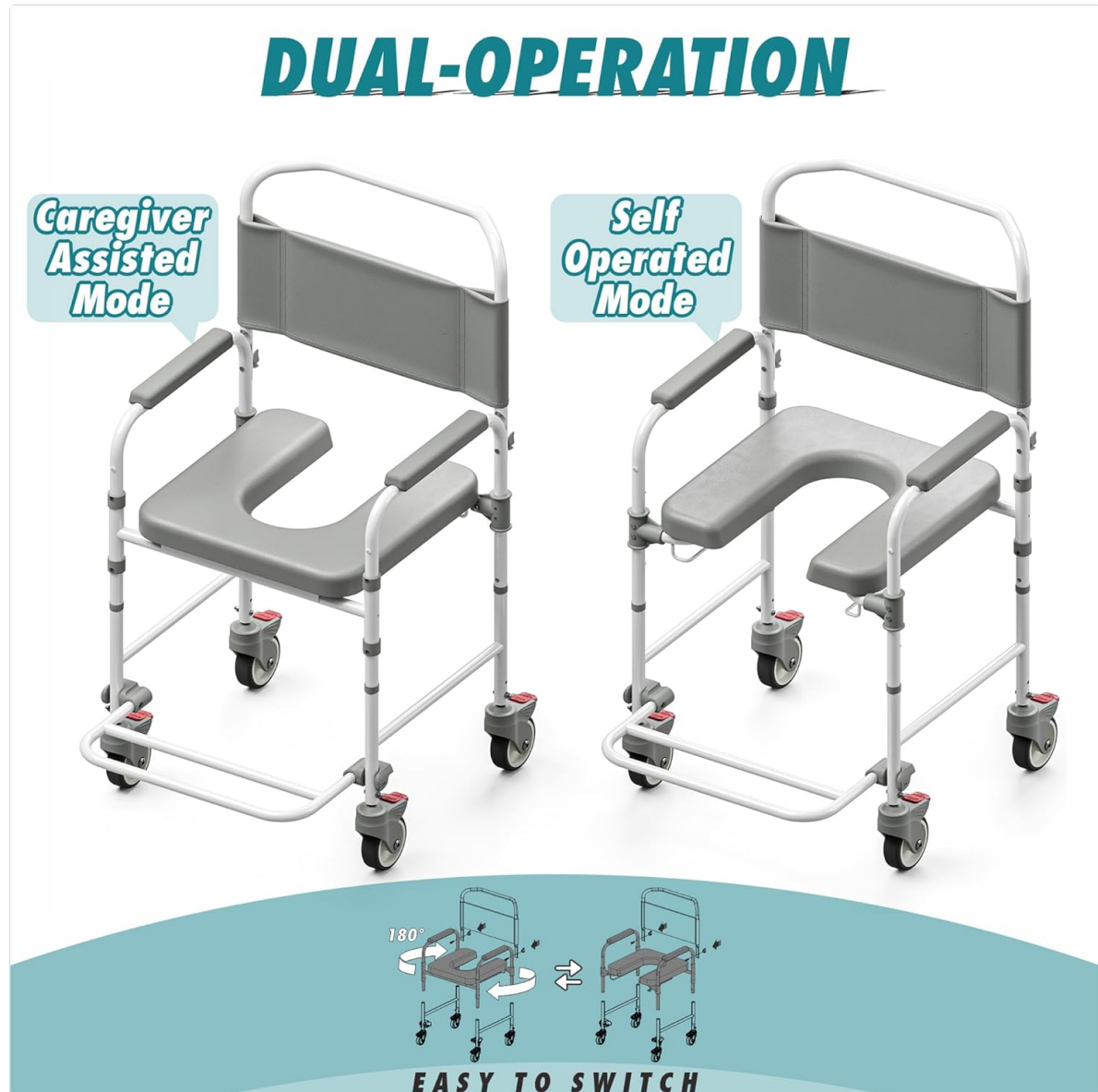


Image 5.1: Diagram illustrating the two operational modes and how to switch between them.





Image 5.2: The commode chair in Caregiver-Assisted Mode, showing a caregiver providing assistance.



Image 5.3: The commode chair in Self-Operated Mode, positioned over a standard toilet for independent use.

## 5.2 Using the Commode Bucket

The included 10-inch commode bucket can be used when a toilet is not accessible. It is designed for easy installation and removal.

1. Ensure the commode bucket is fully unfolded before use. To unfold, press down on the center of the collapsed basin.
2. Slide the bucket into the designated rails beneath the seat.
3. After use, remove the bucket, empty its contents, and clean thoroughly.



Image 5.4: Visual guide on how to unfold the collapsible commode bucket.

### 5.3 Mobility and Wheel Locks

The chair is equipped with four 3-inch swivel wheels, each with an individual foot-lever lock.

- **Moving the Chair:** Release all four wheel locks by pressing the lever upwards. The 360-degree swivel wheels allow for easy maneuverability.
- **Securing the Chair:** To prevent movement, press the foot-lever downwards on all four wheels until they are locked. Always lock the wheels before a user sits down, stands up, or during any stationary use.

## 6. MAINTENANCE AND CARE

Regular cleaning and maintenance will prolong the life of your commode chair.

- **Cleaning the Frame:** Wipe the carbon steel frame with a damp cloth and mild detergent. Dry thoroughly to prevent rust.
- **Cleaning the Seat and Backrest:** Use a mild disinfectant cleaner on the high-density foam cushion and non-slip armrests. Ensure all surfaces are dry before next use.
- **Cleaning the Commode Bucket:** Empty and rinse the commode bucket after each use. Clean with a disinfectant and allow to air dry.
- **Inspections:** Periodically check all fasteners, especially the wing nuts and wheel locks, to ensure they are secure. Inspect the wheels for any damage or excessive wear.
- **Storage:** Store the commode chair in a dry, clean environment when not in use.

## 7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Chair is difficult to move.	Wheel locks are engaged. Wheels are obstructed.	Disengage all wheel locks. Check wheels for debris and clear any obstructions.
Chair feels unstable.	Wheel locks not engaged. Loose assembly. Exceeding weight capacity.	Engage all wheel locks. Check all fasteners and tighten if necessary. Ensure user weight is within 300 lbs capacity.



Problem	Possible Cause	Solution
Commode bucket wobbles.	Bucket not fully inserted or handle misaligned.	Ensure the bucket is fully slid into the rails and the handle is correctly positioned within the metal rack.

## 8. SPECIFICATIONS

---

- **Model:** SF20250706B0F7XZDX36
- **Brand:** littneo
- **Material:** Carbon Steel (Frame), High-Density Foam (Cushion), Polypropylene (PP) (Seat Material Type)
- **Color:** White & Grey
- **Maximum Weight Capacity:** 300 lbs (136 kg)
- **Item Weight:** Approximately 19.81 lbs (9 kg)
- **Package Dimensions:** 21.5 x 18 x 9.5 inches (54.6 x 45.7 x 24.1 cm)
- **Seat Width (between armrests):** 18.9 inches (48 cm)
- **Overall Width (outer edges of armrests):** 21.875 inches (55.5 cm)
- **Height (floor to top of seat):** 21.75 inches (55.2 cm)
- **Wheel Size:** 3-inch swivel wheels with locks



Image 8.1: Diagram showing key dimensions and the 300 lbs weight capacity of the commode chair.

## 9. WARRANTY AND SUPPORT

For warranty information, product support, or to inquire about replacement parts, please contact Litneo customer service. Refer to your purchase documentation for specific warranty terms and contact details.

**Manufacturer:** litneo

**ASIN:** B0F7XZDX36