

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Cearvol](#) /
- › [Cearvol-Nano OTC Hearing Aids Instruction Manual](#)

## Cearvol Nano

# Cearvol-Nano OTC Hearing Aids Instruction Manual

Model: Nano | Brand: Cearvol

## 1. INTRODUCTION

The Cearvol-Nano OTC Hearing Aids are designed to provide clear and comfortable listening for individuals with mild to moderately severe hearing loss. This manual provides essential information for the proper setup, operation, and maintenance of your device to ensure optimal performance and longevity.

## 2. WHAT'S IN THE BOX

Upon unboxing your Cearvol-Nano OTC Hearing Aids, please verify that all the following items are included:

- Charging Case
- Nano Hearing Aids (Left & Right)
- Ear Tips (Small, Medium, Large)
- Type-C Charging Cable
- Cleaning Brush
- User Manual
- Quick Start Guide



**150+**

150+ Certified Patents



**95%**

95% Satisfaction Rate



**20+**

20+ Years of Audio Expertise

Image: All components included in the Cearvol-Nano package.

### 3. PRODUCT FEATURES

#### 3.1 All-in-One Case Control

The charging case serves as a remote control, allowing direct adjustment of volume levels, hearing programs, and environment modes without needing a smartphone app. This design minimizes whistling often associated with on-device adjustments.

**CEARVÖL**

**Stylish Hearing, Vibrant Living**

With nearly 20 years of acoustic expertise, we've grown from crafting high-performance headphones to redefining hearing aids. Driven by our vision "Technology Reshaping Confidence," we blend advanced sound technology with elegant design to create solutions that are lightweight, reliable, and accessible to more people.

Embracing our philosophy "Stylish Hearing, Vibrant Living," we see hearing aids not just as devices to restore hearing, but as companions that help people regain confidence and live with style and vitality.

Image: The Cearvol-Nano charging case with visible control buttons.

#### 3.2 Whistle-Free Comfort

Equipped with Advanced Feedback Cancellation (AFC Tech), these hearing aids effectively eliminate annoying whistling sounds, providing a comfortable and natural listening experience.



Image: Illustration of Advanced Feedback Cancellation technology.

### 3.3 Enhanced Speech Clarity

Advanced sound processing reduces background noise and amplifies voices, leading to clearer and more natural conversations. This feature is beneficial for users with mild to moderately severe hearing loss.

### 3.4 Personalized Hearing Programs

Three preset hearing programs (Mild, Moderate, Moderately Severe) are available directly from the charging case, allowing for customized sound amplification to match specific hearing needs.



Image: Hearing Loss Coverage Chart illustrating preset programs.

### 3.5 Three Environment Modes

Optimize your listening experience in any setting by switching between Indoor, Outdoor, and Restaurant modes directly on the case. This adapts the hearing aids to noisy or quiet environments for improved hearing.

## More Than a Hearing Aid

A stylish sound companion — confident, connected, and effortless, as natural as glasses.



Image: Visual representation of the three environment modes.

### 3.6 Hygienic Charging Case

The charging case features a built-in hygiene light that cleans your rechargeable hearing aids while they charge, ensuring they remain fresh and clean for optimal performance and hygiene.

## Invisible comfort design

Powerful sound in a tiny size. So small and discreet, you'll forget it's there.



Image: The UV hygiene charging case in action.

### 3.7 Extended All-Day Power & Fast Charge

Enjoy up to 17 hours of continuous use on a single charge. The portable charging case provides a total of 54+ hours of power. A quick 15-minute charge in the case offers an additional 3 hours of use.



## Vision

Technology Reshaping Confidence

## Mission

Ensuring people with hearing loss enjoy equality in health, technology, and accessibility

## Core Values

Innovative, Full of Love

Image: Battery life and fast charging capabilities.

### 3.8 Durable & Discreet Fit

The hearing aids feature an IPX5 water and dust-resistant design, protecting them from sweat and rain. Their compact, low-profile design ensures they sit nearly invisibly in the ear for confident, everyday wear.

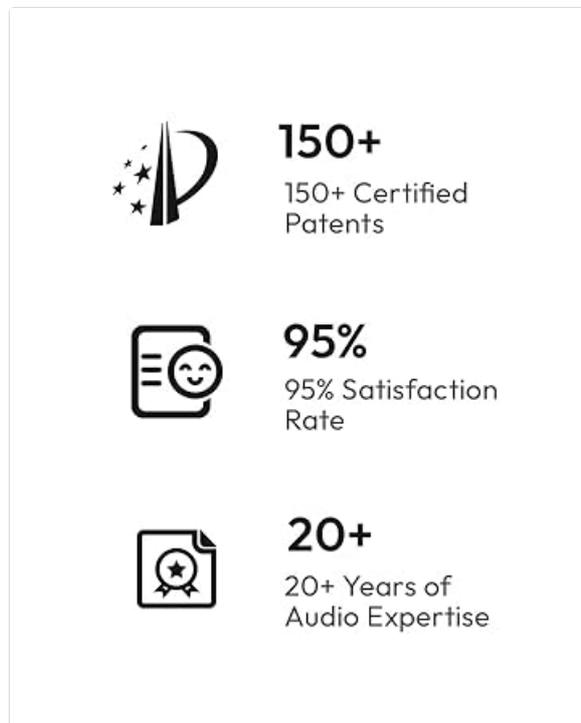


Image: Discreet design of the Cearvol-Nano hearing aids.

## 4. SETUP GUIDE

### 4.1 Initial Charging

Before first use, place both hearing aids into the charging case and connect the Type-C charging cable to a power source. Ensure the case is fully charged. The indicator lights on the case will show

the charging status.

## 4.2 Choosing Ear Tips

Select the ear tip size (S, M, L) that provides the most comfortable and secure fit in your ear canal. A proper fit is crucial for optimal sound quality and to prevent feedback.

## 4.3 Wearing and Removal

Gently pinch the body of the hearing aid with your index finger and thumb. Carefully insert it into your ear canal and twist it slightly until it feels comfortable and secure. To remove, gently pull the assistive pull line at the end of the hearing aid.

# 5. OPERATING INSTRUCTIONS

---

## 5.1 Power On/Off

The hearing aids automatically power on when removed from the charging case and power off when placed back into the case.

## 5.2 Volume Adjustment

Use the '+' and '-' buttons on the charging case to adjust the volume. There are 4 volume levels available. Adjust the volume for the left or right hearing aid independently.

## 5.3 Switching Hearing Programs

Press the 'M' button on the right side of the charging case to cycle through the 3 preset hearing programs: Mild, Moderate, and Moderately Severe. Each change is accompanied by an audible indicator.

## 5.4 Switching Environment Modes

Press the 'M' button on the left side of the charging case to cycle through the 3 environment modes: Indoor, Outdoor, and Restaurant. This allows the hearing aids to adapt to different acoustic settings.

Your browser does not support the video tag.

Video: Demonstrating the Cearvol-Nano's whistle-free technology and ease of use in various environments.

# 6. MAINTENANCE

---

## 6.1 Cleaning

Regularly clean your hearing aids using the provided cleaning brush to remove any earwax or debris. Pay special attention to the ear tips and sound outlets. The charging case's built-in UV light will sanitize the hearing aids during charging.

## 6.2 Storage

When not in use, always store the hearing aids in their charging case to protect them and ensure they are charged and sanitized for your next use.

# 7. TROUBLESHOOTING

---

If you experience any issues with your Cearvol-Nano OTC Hearing Aids, please refer to the following common solutions:

- **No Sound:** Ensure the hearing aids are charged and properly inserted. Check the volume level.
- **Whistling/Feedback:** Ensure ear tips are correctly sized and seated snugly in your ear canal. Adjust the volume or switch to a different hearing program.
- **Distorted Sound/Static:** Ensure the hearing aids are clean and free of earwax. Try switching environment modes.
- **Difficulty Hearing in Noisy Environments:** Switch to the 'Restaurant' or 'Outdoor' environment mode for enhanced speech clarity and noise reduction.
- **Short Battery Life:** Ensure the charging case is fully charged. Confirm hearing aids are properly seated in the case for charging.

If problems persist, please contact customer support for further assistance.

## 8. SPECIFICATIONS

---

- **Item Model Number:** Nano
- **Manufacturer:** Cearvol
- **Battery Type:** Li-ion Rechargeable (2x for hearing aids, 1x for charging case)
- **Continuous Use:** Up to 17 hours per charge
- **Total Power with Case:** 54+ hours
- **Fast Charge:** 15 minutes for 3 hours of use
- **Water Resistance:** IPX5 (sweat and rain resistant)
- **Gain:** Up to 30dB
- **Features:** Advanced Feedback Cancellation (AFC Tech), Digital Sound Processing (DSP), UV Hygiene Light
- **Package Dimensions:** 4.8 x 4.49 x 2.64 inches
- **Item Weight:** 1.62 ounces

## 9. WARRANTY AND SUPPORT

---

For manufacturer's return policy, please call Customer Care at 1(959)-206-0088 or email to [support@cearvol.com](mailto:support@cearvol.com). We work to ensure that both your return and credit refund are processed accurately.

For additional information and other labeling, including the user instructional brochure, these are available on the internet at [www.cearvol.com](http://www.cearvol.com).