



Manuals.plus /

- › Wuvee /
- › Wuvee OTC Hearing Aids Model TY CN430S User Manual

Wuvee TY CN430S

Wuvee OTC Hearing Aids Model TY CN430S User Manual

Your guide to clear hearing and optimal device performance.

1. INTRODUCTION

Thank you for choosing Wuvee OTC Hearing Aids Model TY CN430S. This manual provides essential information for the safe and effective use of your new hearing aids. Please read it thoroughly before first use and retain it for future reference.



OTC

Hearing Aids

Model: TY CN430S

Hear the Sound of Life



Easy operate



Comfort



Super-fit



Image: Wupee OTC Hearing Aids Model TY CN430S in charging case.

2. PACKAGE CONTENTS

Verify that all items are present in your package:

- Wupee Hearing Aids (Left and Right)
- Charging Box with LED Display
- Earplugs (Small, Medium, Large sizes)
- Wax Guards
- Cleaning Brush
- Charging Cable (USB Type-C)
- User Manual
- Screwdriver (for maintenance, if applicable)

What's in the box

- 1 Hearing Aid
- 2 Charging Box
- 3 Earplugs
- 4 Wax Guard
- 5 Cleaning Brush
- 6 Charging Cable
- 7 User Manual



Image: All components included in the Wuvée hearing aid package.

3. DEVICE COMPONENTS AND CONTROLS

Familiarize yourself with the different parts of your hearing aids and charging case:

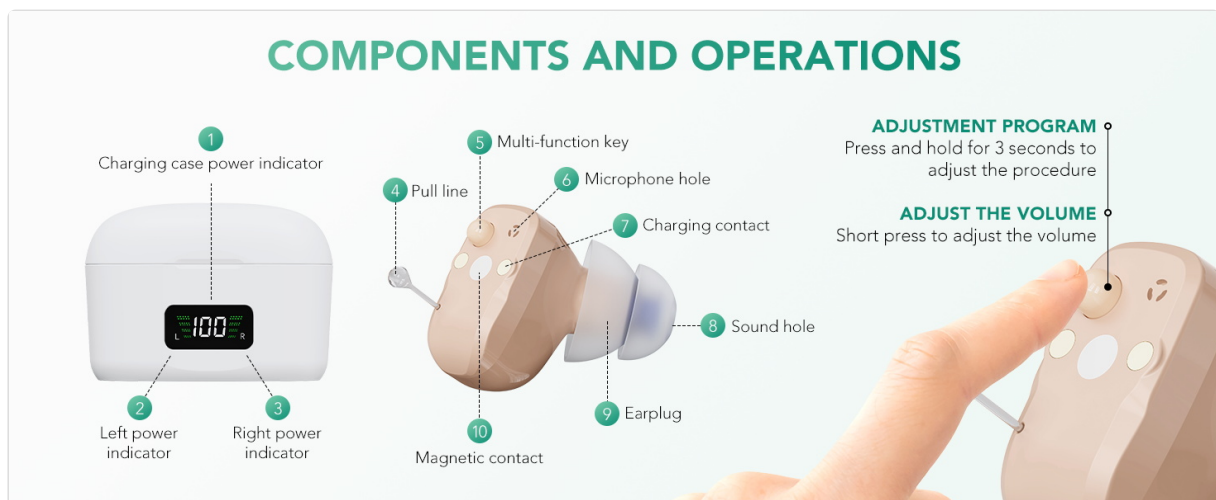


Image: Detailed diagram of hearing aid and charging case components.

1. **Charging Case Power Indicator:** Displays the charging status of the case.
2. **Left Power Indicator:** Shows the battery level of the left hearing aid.
3. **Right Power Indicator:** Shows the battery level of the right hearing aid.
4. **Pull Line:** For easy removal of the hearing aid from the ear.
5. **Multi-function Key:** Used for adjusting volume and switching programs.
6. **Microphone Hole:** Captures sound from the environment.
7. **Charging Contact:** Connects to the charging case for power.
8. **Sound Hole:** Delivers amplified sound into the ear canal.
9. **Earplug:** Soft silicone tip for comfortable fit and sound delivery.
10. **Magnetic Contact:** Secures the hearing aid in the charging case.

4. SETUP AND FIRST USE

4.1. Initial Charging

Before first use, fully charge the hearing aids and the charging case.

- Connect the charging cable to the charging case and a power source.
- The LED display on the charging case will indicate the charging status.
- A full charge for the case takes approximately 2 hours, providing 5-6 full charges for the hearing aids.
- The hearing aids will automatically charge when placed in the powered charging case.

FAST CHARGING & LONG-LASTING BATTERY LIFE

≥3 Times
Hearing Devices
Charging

≥80 hrs
Working Time

≤2 hrs
Charging Time

LED Power Display

Type-c Charging

Image: Charging the Wuvee hearing aids and case.

4.2. Selecting and Attaching Earplugs

Choose the earplug size (S, M, L) that provides a comfortable and secure fit in your ear canal. A proper fit is crucial for optimal sound quality and to prevent whistling.

4.3. Inserting the Hearing Aids

Follow these steps for proper insertion:

1. Select the correct size earplug for your ear.
2. Gently place the hearing aid inside your ear canal.

3. Adjust the hearing aid to a comfortable and secure position.



Image: Steps for selecting earplugs and inserting hearing aids.

5. OPERATING INSTRUCTIONS

5.1. Power On/Off

The hearing aids feature an automatic on/off function. They will power on when removed from the charging case and power off when placed back into the case.

5.2. Adjusting Volume and Programs

Use the multi-function key on the hearing aid to control settings:

- **Adjust Volume:** Short press the multi-function key to cycle through volume levels.
- **Change Program:** Press and hold the multi-function key for 3 seconds to switch between hearing programs.

The Wuvee hearing aids offer three distinct programs:

- **Quiet Mode (Program 1):** Ideal for daily conversations and watching television.
- **Noise Reduction Mode 1 (Program 2):** Suitable for environments with slight indoor background noise.

- **Noise Reduction Mode 2 (Program 3):** Designed for severe and complex noise environments.

3 PROGRAMS FOR ALL SITUATIONS



<p>Quiet Mode For Daily</p> <p>Conversations and watching TV</p>	<p>Noise Reduction Mode 1</p> <p>For indoor background noise, slight noise</p>	<p>Noise Reduction Mode 2</p> <p>For severe and complex noise environment</p>
---	---	--

Image: Overview of the three hearing programs.

5.3. LED Display Indicators

The charging case features an LED display to visualize the remaining power of the charging case and the individual hearing aids, helping you avoid unexpected power depletion.

6. MAINTENANCE

6.1. Cleaning the Hearing Aids

Regular cleaning is essential for hygiene and optimal performance. Use the provided cleaning brush to remove earwax and debris from the earplugs and sound holes.

6.2. Replacing the Earwax Guard

The earwax guard protects the hearing aid from earwax buildup. Replace it periodically or if you notice a decrease in sound quality.

1. Insert the "removal hook side" of the ear guard holder into the old earwax guard.
2. Pull out the earplug bracket along with the old earwax guard.
3. Rotate the ear cap holder to expose the new earwax guard.

4. Push the new earwax guard into the amplifier.
5. Pull out the ear guard holder straight out.



Image: Step-by-step guide for earwax guard replacement.

6.3. Charging Best Practices

- Always use the provided charging cable and adapter.
- Avoid overcharging by disconnecting once fully charged.
- Store the hearing aids in their charging case when not in use to keep them charged and protected.

7. TROUBLESHOOTING

If you encounter issues with your Wuvee hearing aids, refer to the following common problems and solutions:

- **Whistling or Feedback:**
 - Ensure the earplug is correctly sized and securely inserted into your ear canal. A loose fit can cause sound to leak and be picked up by the microphone, creating feedback.
 - Try a different size earplug for a better seal.
 - Check for earwax buildup in the ear canal or on the hearing aid.

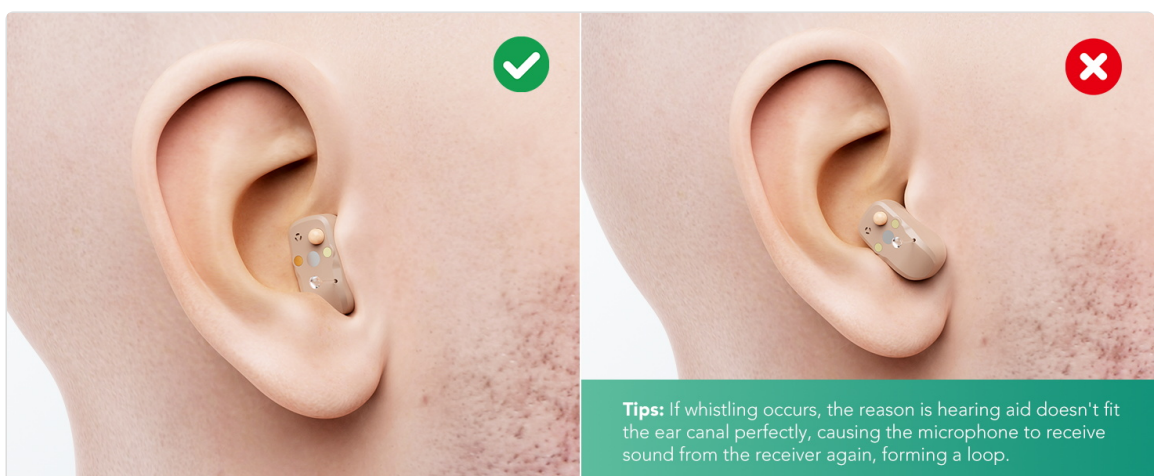


Image: Correct and incorrect hearing aid fit to prevent whistling.

- **No Sound or Low Sound:**
 - Check if the hearing aid is fully charged.

- Ensure the volume is adjusted to an audible level.
- Clean the earplug and sound hole for any blockages.
- Replace the earwax guard if it appears clogged.

- **Charging Issues:**

- Verify the charging cable is securely connected to both the case and the power source.
- Ensure the charging contacts on both the hearing aids and the case are clean and free of debris.
- Check if the power source is functional.

- **Uncomfortable Fit:**

- Experiment with different earplug sizes to find the most comfortable fit.
- Ensure the hearing aid is not inserted too deeply or at an awkward angle.

8. SPECIFICATIONS

Feature	Specification
Model Name	TY CN430S
Brand	Wuvec
Style	Noise Cancellation
Material	Resin
Color	White
Batteries	1 Lithium Ion battery required (included)
Charging Time	Approx. 2 hours for charging case
Standby Power	More than 80 hours (from fully charged case)
Package Dimensions	4.88 x 3.66 x 1.65 inches
Item Weight	5.29 ounces

9. SAFETY INFORMATION

This device is intended as a hearing aid. Consult a medical professional if you experience sudden hearing loss, ear pain, or discharge before using this product.

10. WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the documentation included with your purchase or contact Wuvec customer service directly. Contact details can typically be found on the product packaging or the official Wuvec website.