

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Loop](#) /

› [Loop Dream Silicone/Memory Foam Earplugs User Manual - Model en-dr-blk-03](#)

### Loop en-dr-blk-03

# Loop Dream Silicone/Memory Foam Earplugs Instruction Manual

Model: en-dr-blk-03

[Introduction](#)   [Setup](#)   [Operating](#)

[Instructions](#)

[Maintenance](#)

[Troubleshooting](#)

[Specifications](#)

[Warranty & Support](#)

[Safety Information](#)

## 1. INTRODUCTION

The Loop Dream earplugs are designed to provide comfortable and effective noise reduction for improved sleep quality. Featuring an innovative silicone body and a choice of ear tip styles, these earplugs are engineered to muffle disturbances such as snoring and nighttime noises, ensuring a more restful night's sleep.

### What's Included:

- 1 set of Loop Dream earplugs
- 4 sets of Loop Dream Foam-Silicone Ear Tips (sizes XS, S, M, L — M comes pre-installed)
- 4 sets of Loop Dream Double Tips (sizes XS, S, M, L)
- 1 Loop Dream Carry Case



Image: The Loop Dream earplugs packaging, displaying the earplugs, various tip sizes, and the compact carry case.

## 2. SETUP: FINDING YOUR PERFECT FIT

---

Achieving the correct fit is crucial for optimal comfort and noise reduction. Loop Dream earplugs come with two types of ear tips and multiple sizes to ensure a customized fit for every ear.

### 2.1 Choosing Your Ear Tips

Your Loop Dream earplugs include two styles of ear tips:

- **Foam-Silicone Tips:** These tips combine the softness of silicone with the conforming properties of memory foam for a snug and comfortable seal.
- **Double Tips:** These offer an alternative fit, providing extra options for a secure seal.

Experiment with both types and all provided sizes (XS, S, M, L) to determine which offers the best comfort and noise isolation for your ears. The Medium Foam-Silicone Tips come pre-installed.

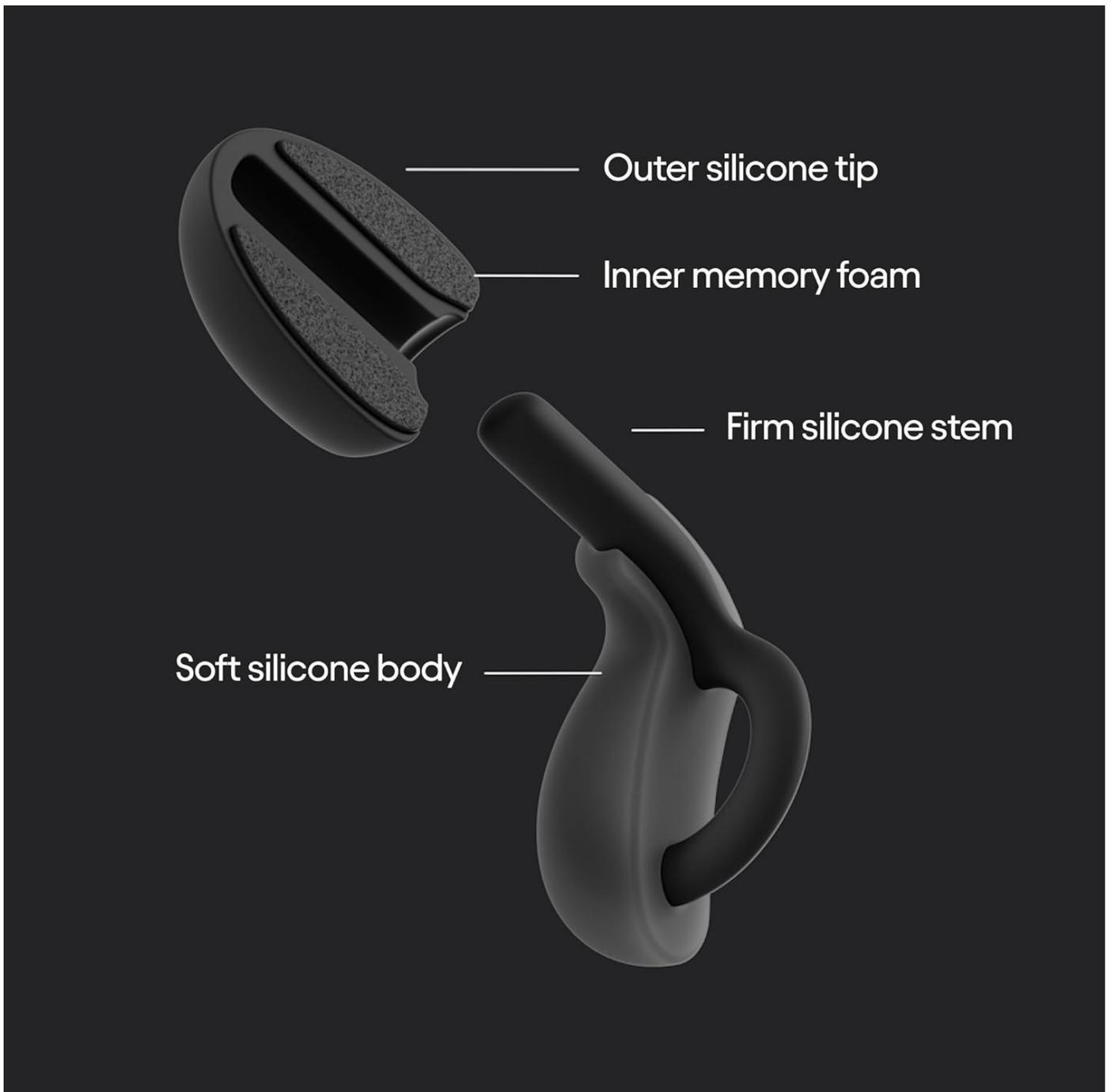


Image: An exploded view of the Loop Dream earplug, showing the outer silicone tip, inner memory foam, firm silicone stem, and soft silicone body.

## 2.2 Inserting the Earplugs

1. **Insert and Hold:** Gently pull your ear back to open up the ear canal. Insert the chosen ear tip into your ear canal.
2. **Secure the Seal:** Hold the earplug in place for 5 seconds to allow the tip to expand and create a secure seal.
3. **Check the Fit:** The earplugs should sit snugly in your ears without sticking out or being visible from the front. They should feel comfortable, even when lying on your side.

# Find your perfect fit

## 1. Insert and hold

Pull your ear back to open up the canal. Insert the ear tip and hold for 5 seconds to seal.

The earplugs should never stick out or be visible from the front.



## 2. Try all ear tip sizes

If they're sticking out or uncomfortable, try a smaller size.

If they don't feel secure, size up

For an alternative fit, try the Dream Double Tips

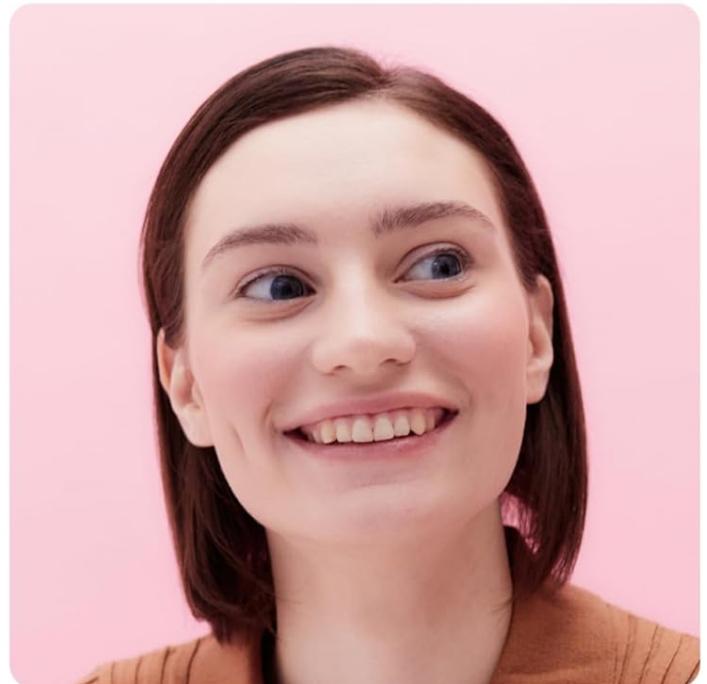


Image: A visual guide demonstrating the correct method for inserting Loop Dream earplugs, emphasizing pulling the ear back and holding for 5 seconds.

### 2.3 Adjusting the Size

- If the earplugs are sticking out or feel uncomfortable, try a smaller size.
- If they don't feel secure or fall out easily, try a larger size.
- For an alternative fit, experiment with the Dream Double Tips.


#### Video: How To Find Your Size

Your browser does not support the video tag.

Video Description: This video demonstrates the process of selecting the correct ear tip size for Loop Dream earplugs, showing how to test different sizes for optimal comfort and seal.

## 3. OPERATING INSTRUCTIONS

Once properly fitted, your Loop Dream earplugs will provide 27 dB (SNR) of noise reduction, effectively muffling ambient sounds to help you achieve a deeper and more restful sleep.

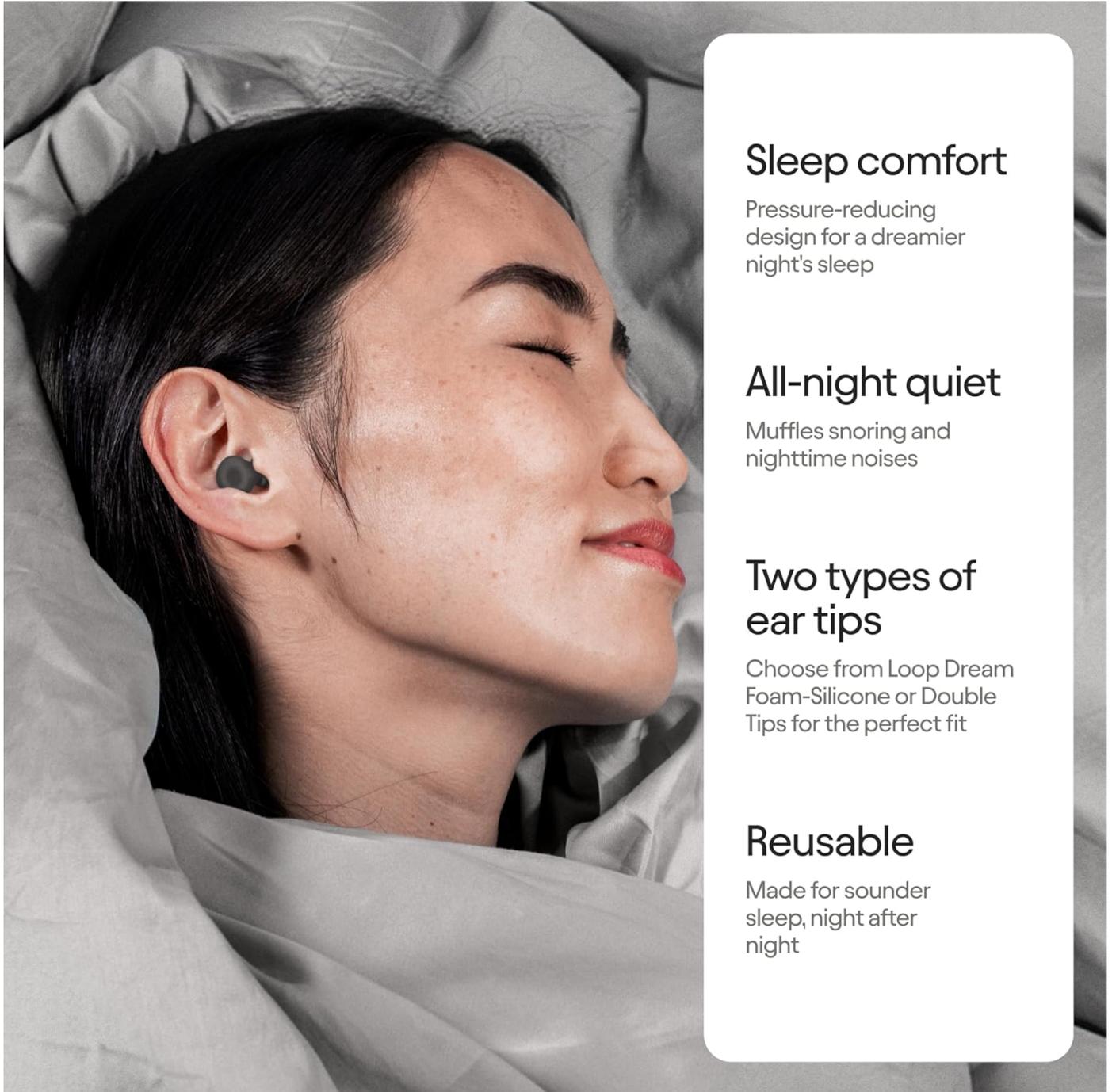

## Sleep comfort

Pressure-reducing design for a dreamier night's sleep

## All-night quiet

Muffles snoring and nighttime noises

## Two types of ear tips

Choose from Loop Dream Foam-Silicone or Double Tips for the perfect fit

## Reusable

Made for sounder sleep, night after night

Image: A woman sleeping peacefully on her side with Loop Dream earplugs inserted, illustrating comfort during sleep.



Image: A man sleeping peacefully on his side with Loop Dream earplugs inserted, illustrating comfort during sleep. The innovative silicone body is designed for comfort, especially for side sleepers, reducing pressure on the ear. The oval shape of the Foam-Silicone Tips follows your ear's natural contours for a secure fit throughout the night.

**Video: Loop Dream Overview**

Your browser does not support the video tag.

Video Description: This video provides a brief overview of the Loop Dream earplugs, highlighting their design for sleep comfort and noise reduction capabilities.

## 4. MAINTENANCE

---

To ensure the longevity and hygiene of your Loop Dream earplugs, regular cleaning is recommended.

### **Cleaning Instructions:**

- Gently remove the ear tips from the earplug body.

- Wash the silicone earplugs and tips with mild soap and warm water.
- Rinse thoroughly to remove all soap residue.
- Allow them to air dry completely before reassembling or storing them in the carry case.
- Avoid using harsh chemicals or abrasive materials, as these can damage the silicone.

## Storage:

Store your Loop Dream earplugs in the provided carry case when not in use. This protects them from dust, dirt, and damage, ensuring they remain clean and ready for your next use. The carry case features a no-slip grip bottom and a sizeable compartment for easy access at your bedside.

## 5. TROUBLESHOOTING

Problem	Possible Cause	Solution
Insufficient Noise Reduction	Incorrect ear tip size or improper insertion.	Ensure earplugs are inserted correctly and create a tight seal. Try different ear tip sizes (XS, S, M, L) or the Double Tips for a better fit.
Discomfort or Pain	Ear tip size is too large or earplugs are inserted too deeply.	Try a smaller ear tip size. Ensure the earplugs are not pushed too far into the ear canal. They should sit comfortably without causing pressure.
Earplugs Fall Out Easily	Ear tip size is too small or improper insertion.	Try a larger ear tip size to achieve a more secure fit. Ensure the earplugs are inserted correctly and held for 5 seconds to create a proper seal.
Reduced Hygiene	Lack of regular cleaning.	Clean earplugs and tips regularly with mild soap and warm water. Store them in the provided carry case when not in use.

## 6. SPECIFICATIONS

- **Model Name:** en-dr-blk-03
- **Material Type:** Silicone
- **Noise Attenuation (SNR):** 27 dB
- **Reusability:** Reusable
- **Product Dimensions:** 3.54 x 3.54 x 0.73 inches
- **Item Weight:** 30.5 Grams (1.08 ounces)
- **Color:** Black (other colors available separately)
- **Included Ear Tips:** 4 sets of Foam-Silicone Tips (XS, S, M, L), 4 sets of Double Tips (XS, S, M, L)

## 7. WARRANTY & SUPPORT

For any questions, concerns, or support needs regarding your Loop Dream earplugs, please refer to the user guide included in the box or contact Loop customer support.

- **Return Policy:** 30-day refund/replacement duration.
- **Manufacturer:** Loop
- **Contact:** Refer to the official Loop website or the contact information provided in your product packaging for

support.

## 8. SAFETY INFORMATION

---

**Choking Hazard Warning: This product contains small parts and may present a choking hazard for children under 3 years of age. Keep out of reach of young children.**

For detailed instructions on how to use the earplugs, please scan the QR code included in the packaging. The code will provide step-by-step guidance to ensure proper fit and optimal protection.



© 2026 Loop. All rights reserved.