



Manuals.plus /

› Momcozy /

› Momcozy BreezyClear Pro Electric Nasal Aspirator User Manual (Model BN008)

Momcozy BN008

Momcozy BreezyClear Pro Electric Nasal Aspirator User Manual

Model: BN008

Brand: Momcozy

INTRODUCTION

The Momcozy BreezyClear Pro is an innovative 2-in-1 electric nasal aspirator designed to provide gentle mucus-loosening and hospital-grade suction for babies and toddlers. This manual provides essential information for the safe and effective operation and maintenance of your device.



Image: Momcozy BreezyClear Pro Electric Nasal Aspirator and its retail packaging.

PRODUCT FEATURES

- **2-in-1 Design:** Combines gentle mucus-loosening with hospital-grade suction.
- **Effective & Safe Loosening:** Soft silicone mask for targeted moistening, automatically shuts off after 6 seconds.
- **Hospital-Grade, Adjustable Suction:** Up to 75 kPa suction power with 9 adjustable levels. Includes two different-sized soft silicone tips.
- **Easy to Clean:** Disassembles easily for cleaning by hand or in the dishwasher. A cleaning kit is included.
- **Low Noise & Soothing Sounds:** Operates at 59 dB and plays calming music to keep the baby relaxed.
- **Convenient Storage:** Built-in compartment for small parts and includes a travel pouch.

Clears All Congestion Levels With Ease

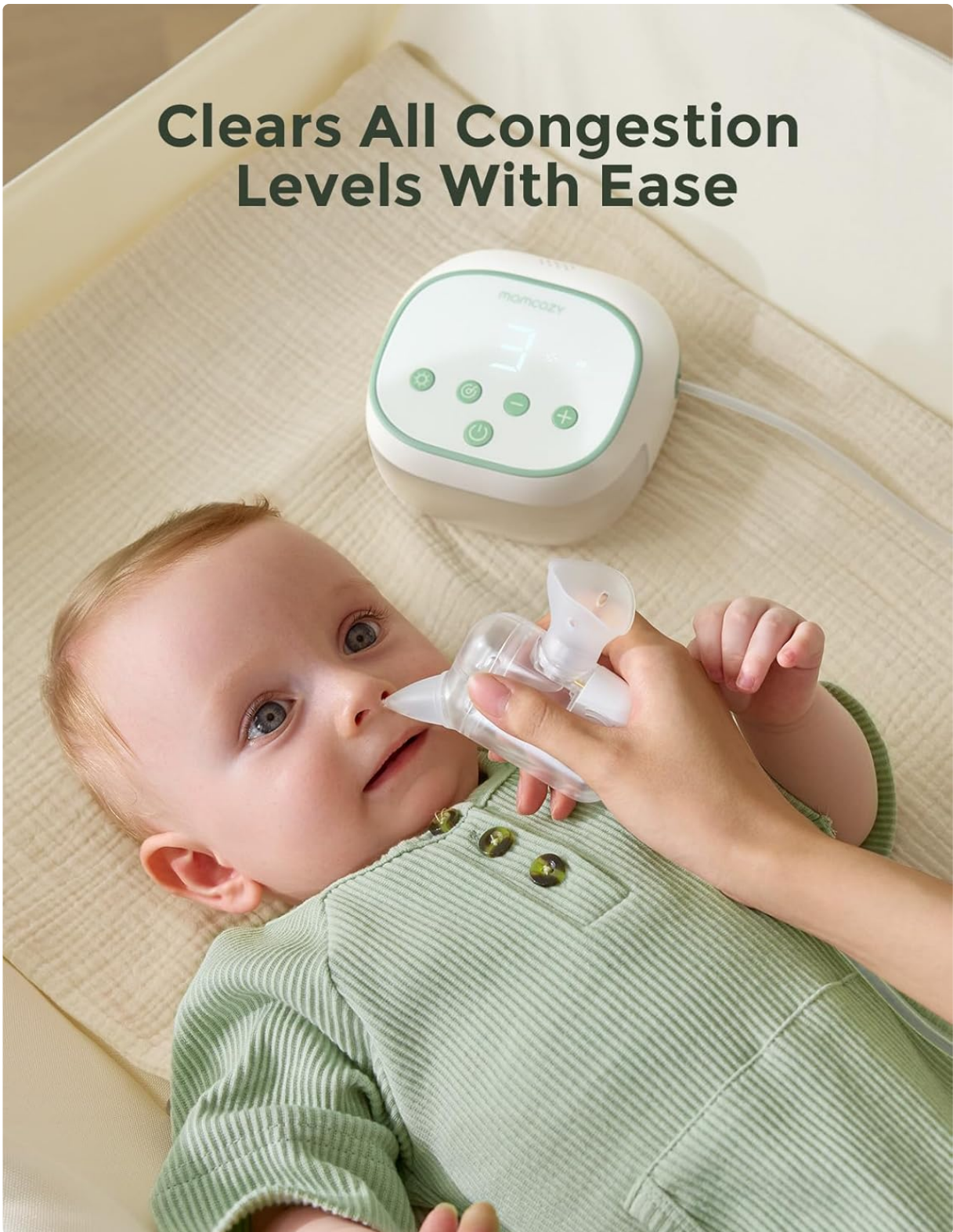


Image: The Momcozy BreezyClear Pro with its built-in storage compartment for accessories.

SETUP

1. **Unpack Components:** Carefully remove all items from the packaging. Verify that all parts listed in the "What's in the Box" section are present.
2. **Initial Cleaning:** Before first use, wash all mucus-contact parts (silicone tips, collection cup, silicone mask, silicone tube) with warm soapy water or place them in the dishwasher. Rinse thoroughly and allow to air dry completely.
3. **Assemble the Device:**
 - Attach the desired silicone tip to the collection cup.

- Connect the silicone tube to the collection cup and the main motor unit.
 - If using the moistening function, attach the soft silicone mask to the designated port on the collection cup.
4. **Charge the Device:** Connect the USB-C cable to the charging port on the main unit and plug it into a suitable USB power adapter. Ensure the device is fully charged before first use.

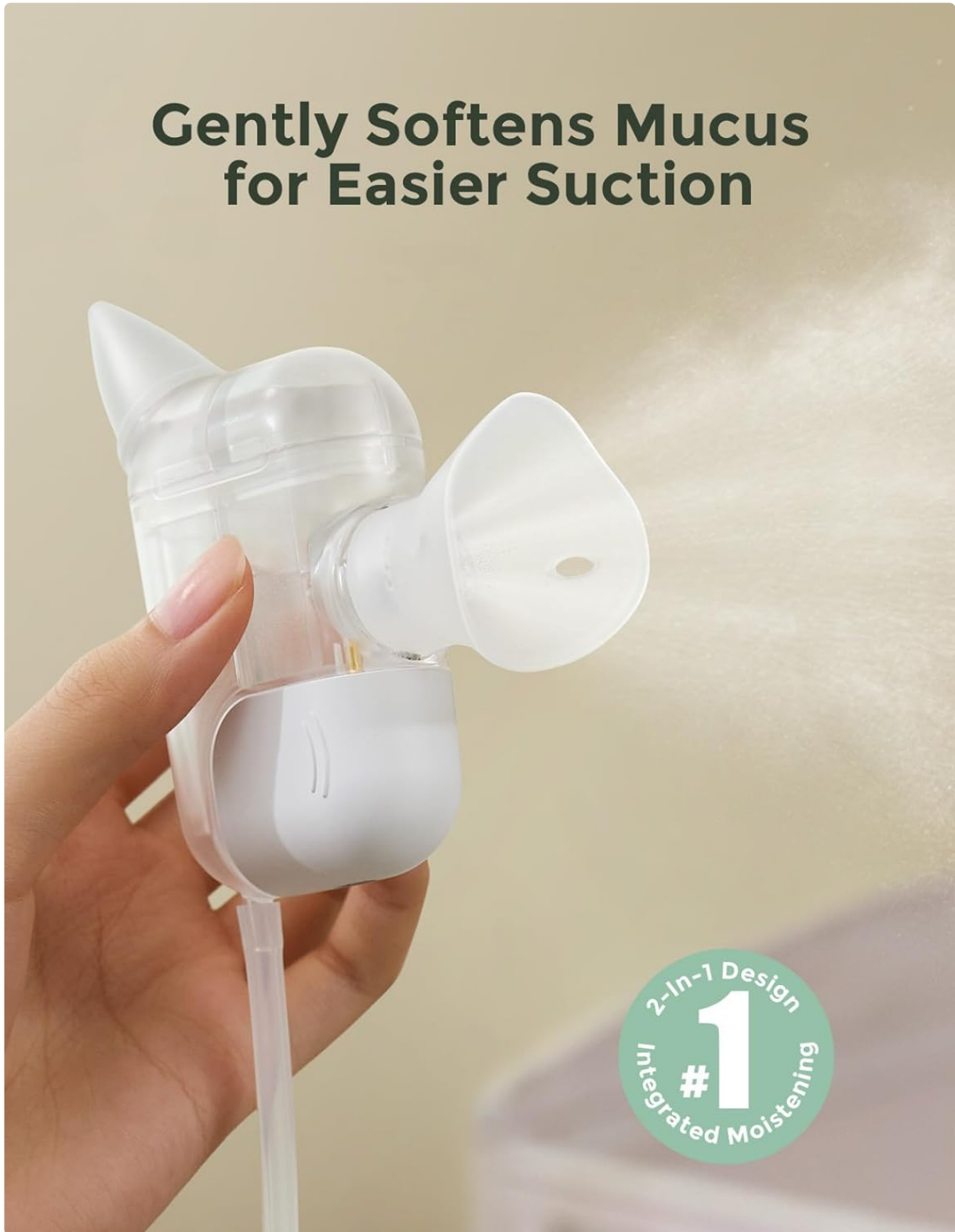


Image: All components of the Momcozy BreezyClear Pro laid out for assembly.

OPERATING INSTRUCTIONS

Power On/Off

- Press and hold the power button for 2 seconds to turn the device on (green light will illuminate).
- Press and hold the power button again for 2 seconds to turn the device off.

Using the Moistening Function

1. Pour purified or saline water into the water tank located on the collection cup.
2. Ensure the silicone mask is properly attached.
3. With the device powered on, click the moistening button to start moisturizing the nasal cavity. The feature will automatically stop after 6 seconds. Click again if additional moistening is needed.



Image: The moistening function being used to gently prepare a baby's nasal passages.

Using the Suction Function

1. Select the appropriate silicone tip for your child's age and congestion level.
2. Turn on the device.
3. Adjust the suction level using the '+' and '-' buttons. There are 9 adjustable levels.
4. Gently place the silicone tip into your child's nostril, ensuring a snug fit to create a seal.
5. Activate suction. Move the tip as needed to clear mucus.
6. Repeat for the other nostril as necessary.

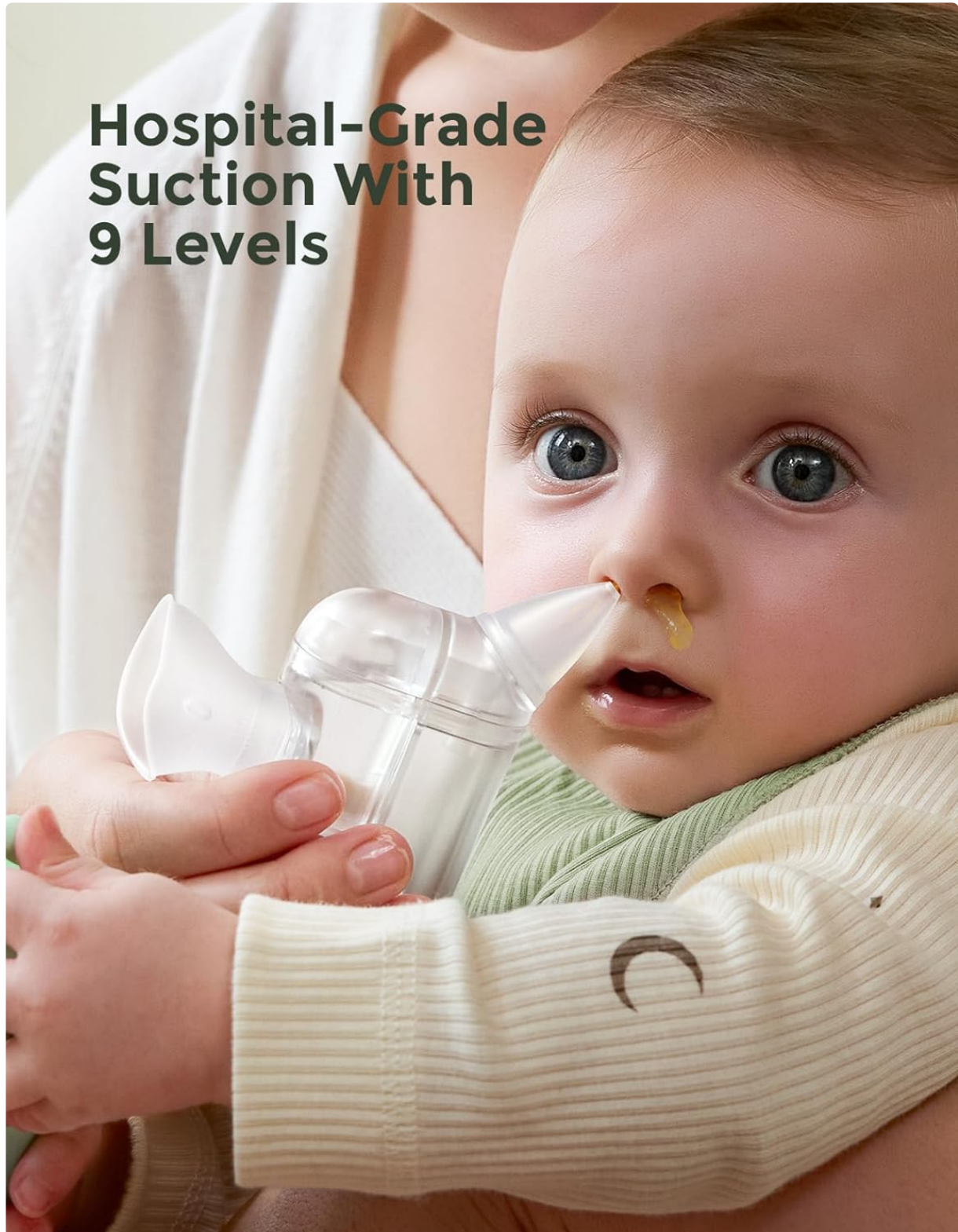


Image: The suction function of the Momcozy BreezyClear Pro being used to clear a baby's nasal congestion.

Music and Light Functions

- Press the music note button to cycle through 6 calming lullabies.
- Press the light bulb button to activate the night light feature, which has multiple brightness settings.

Your browser does not support the video tag.

Video: Demonstration of the Momcozy BreezyClear Pro's features, including suction levels, music, and light functions.

MAINTENANCE

Cleaning

- All mucus-contact parts are fully detachable for thorough cleaning.
- Wash parts by hand with warm soapy water and the included cleaning brush, or place them in the dishwasher.
- Rinse all parts thoroughly to remove any soap residue.
- Allow all parts to air dry completely before reassembling or storing to prevent mold or bacterial growth.
- Do not immerse the main motor unit in water. Wipe it clean with a damp cloth.



Image: Detachable parts of the Momcozy BreezyClear Pro being cleaned in a dishwasher basket.

Storage

- Store the clean and dry device in its built-in compartment or the provided travel pouch.
- Keep the device in a cool, dry place away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No suction or weak suction	<ul style="list-style-type: none"> • Improper assembly • Clogged tip/tube • Low battery 	<ul style="list-style-type: none"> • Ensure all parts are securely connected. • Clean the silicone tip and tube thoroughly. • Charge the device fully.
Moistening function not working	<ul style="list-style-type: none"> • Water tank empty • Silicone mask not attached correctly 	<ul style="list-style-type: none"> • Fill the water tank with purified or saline water. • Ensure the silicone mask is securely in place.
Device not turning on	<ul style="list-style-type: none"> • Low battery • Power button not held long enough 	<ul style="list-style-type: none"> • Charge the device. • Press and hold the power button for 2 seconds.

SPECIFICATIONS

- **Model Number:** BN008
- **Suction Power:** Up to 75 kPa
- **Suction Levels:** 9 adjustable levels
- **Noise Level:** 59 dB
- **Battery:** 2500mAh rechargeable (USB-C charging)
- **Usage Time per Charge:** Up to 180 minutes
- **Product Dimensions:** 8.7 x 5.5 x 2.4 inches
- **Weight:** 0.63 ounces
- **Included Tips:** 2 soft silicone tips (different sizes)

Honored by Awards and Media Recognition

Serving over 3,600,000* Moms Worldwide



*Data as of December 2024 from Amazon official platform.



Image: Visual summary of Momcozy BreezyClear Pro's key features and specifications.

WARRANTY AND SUPPORT

Momcozy products are designed with quality and customer satisfaction in mind. For warranty information, technical support, or any questions regarding your BreezyClear Pro Nasal Aspirator, please refer to the official Momcozy website or contact customer service directly.

Please visit the Momcozy Official Store for more information and product registration.