

Kortyio EG9AWE10W

Kortyio OrocLEAR Ear Cleaner Instruction Manual

Model: EG9AWE10W

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use, maintenance, and troubleshooting of your Kortyio OrocLEAR Ear Cleaner, Model EG9AWE10W. This device is designed to assist in the removal of earwax through gentle vibration and suction.

2. SAFETY INFORMATION

- **Do not insert the device too deeply into the ear canal.** The flexible silicone tip is designed to prevent over-insertion, but always use with caution.
- **Consult a medical professional** before use if you have ear infections, ear drum perforations, ear surgery, or any other ear conditions.
- **Keep out of reach of children.** Adult supervision is required if used by children.
- **Do not use excessive force.** The device operates with gentle vibrations and suction; forcing it can cause injury.
- **Ensure the silicone tip is securely attached** before each use.
- **Clean the device and tips thoroughly** after each use to prevent bacterial growth.
- **This device is for personal use only.** Do not share tips between users without proper sterilization.

3. PRODUCT COMPONENTS

The Kortyio OrocLEAR Ear Cleaner package includes the following items:

- OrocLEAR Ear Cleaner Device (main unit)
- Replacement Silicone Tips (multiple sizes/types)
- Cleaning Brush

- Storage Case
- AA Battery (may be included or sold separately)



Figure 3.1: Kortyio Oroclear Ear Cleaner and included accessories.



Figure 3.2: Detailed view of the device components.

4. SETUP

1. Battery Installation:

- Locate the battery compartment on the main unit.
- Open the compartment cover.
- Insert one AA battery, ensuring correct polarity (+/-).
- Close the battery compartment cover securely.


2. Attaching the Silicone Tip:

- Select a clean, appropriate-sized silicone tip.
- Align the tip with the device's nozzle.
- Gently push and twist the tip until it is firmly seated and secure.

5. OPERATING INSTRUCTIONS

Follow these steps for effective earwax removal:

1. **Preparation:** Ensure the device is clean and a fresh battery is installed. Select a comfortable silicone tip.
2. **Power On:** Press the power button located on the device. You will feel a gentle vibration and hear the suction motor activate.
3. **Insertion:** Gently insert the silicone tip into the ear canal. Do not force the tip. The flexible material is designed to fit various ear sizes comfortably.
4. **Movement:** Slowly and gently rotate the device within the outer part of the ear canal. Allow the gentle vibrations and suction to work. Avoid pushing the device deep into the ear.
5. **Duration:** Use for a short period, typically 15-30 seconds per ear. Discontinue use if you experience any discomfort or pain.
6. **Removal:** Gently withdraw the device from the ear canal.
7. **Power Off:** Press the power button again to turn off the device.

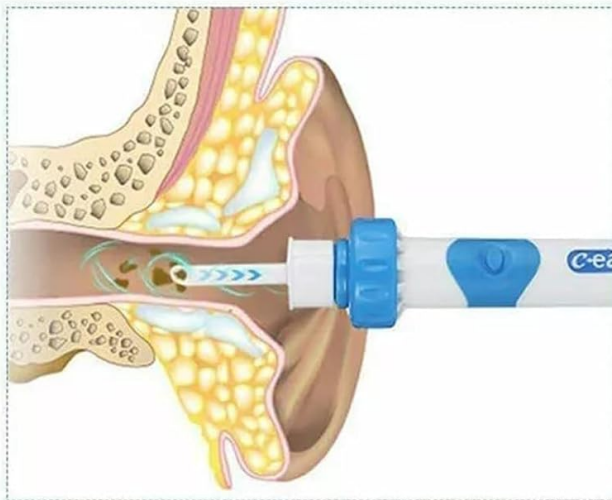

SAFE VACUUM EAR CLEANER

—
This removal cleans the earwax with light vibration, easily sucking away the powder and blocky earwax.

DANGEROUS

Cotton increases the chances of a ruptured eardrum and further hearing loss by pushing wax down the ear canal.


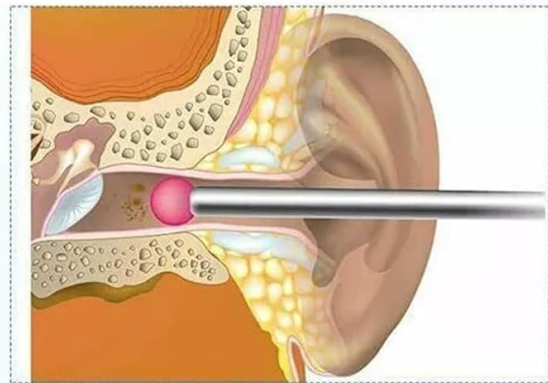


Figure 5.1: Illustration of earwax removal via suction compared to cotton swab use.

SUITABLE FOR ALL EAR HOLE SIZES
SOFT MATERIAL

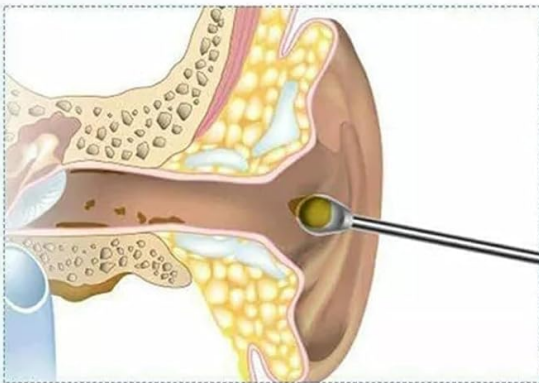

Made of 2.5cm soft and bendable material
This tool protects your ear canal and eardrum



Figure 5.2: Proper insertion of the soft silicone tip into the ear.

6. MAINTENANCE AND CLEANING

Regular cleaning ensures hygiene and extends the life of your device.

1. **Remove the Silicone Tip:** Gently twist and pull the silicone tip off the device.
2. **Clean the Tip:** Rinse the silicone tip under warm running water. Use mild soap if necessary. Ensure all earwax residue is removed.
3. **Clean the Device:** Use the provided cleaning brush to remove any earwax from the device's nozzle and internal components. Do not submerge the main unit in water. Wipe the exterior with a damp cloth.
4. **Dry Components:** Allow all components to air dry completely before reassembling or storing.
5. **Storage:** Store the clean, dry device and accessories in the provided storage case in a cool, dry place.



Figure 6.1: Removing accumulated earwax from the suction head.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Battery is dead or incorrectly installed.	Replace the battery or ensure it is inserted with correct polarity.
No suction or weak suction.	Battery low; tip is clogged with earwax; tip is not securely attached.	Replace battery; clean the silicone tip and device nozzle; reattach tip securely.
Device makes unusual noise.	Internal obstruction or malfunction.	Turn off the device and check for any visible obstructions. If the issue persists, contact customer support.
Discomfort during use.	Improper insertion; sensitive ear canal.	Ensure gentle insertion and movement. Discontinue use if discomfort persists and consult a medical professional.

8. SPECIFICATIONS

- **Model:** EG9AWE10W
- **Brand:** Kortyio
- **Material:** Acrylonitrile Butadiene Styrene (ABS)
- **Power Source:** Battery Powered (1x AA battery)
- **Tip Material:** Soft Silicone

9. CUSTOMER SUPPORT

For further assistance, questions, or concerns regarding your Kortyio OrocLEAR Ear Cleaner, please contact Kortyio customer support. Refer to your purchase documentation for specific contact details or visit the official Kortyio website.

Manufacturer: Kortyio