

[Manuals.plus](#) /

> [XIAOSHUNER](#) /

> XIAOSHUNER 4-in-1 Rolling Shower Chair User Manual

XIAOSHUNER SC8088A

XIAOSHUNER 4-in-1 Rolling Shower Chair User Manual

Model: SC8088A

1. INTRODUCTION

This XIAOSHUNER 4-in-1 Rolling Shower Chair is designed to provide versatile mobility and support for individuals with limited mobility, including the elderly and disabled. It functions as a shower chair, toilet commode, toilet riser, and transport wheelchair, offering enhanced safety and convenience in various settings.

2. SAFETY INFORMATION

Please read and understand all safety instructions before using this product. Failure to follow these instructions may result in injury.

- **Maximum Weight Capacity:** This chair supports a maximum user weight of 350 lbs (158 kg). Do not exceed this limit.
 - **Locking Wheels:** Always engage the wheel locks when the chair is stationary, especially during transfers, showering, or using the commode function, to prevent unintended movement.
 - **Non-Tip Design:** The chair features a non-tip design for stability. Ensure the chair is on a level surface before use.
 - **Secure Transfers:** Utilize the flip-up armrests for safer side transfers. Ensure the user is properly seated and balanced before moving the chair.
 - **Regular Inspection:** Periodically check all components for wear, damage, or loose fasteners. Do not use the chair if any part is damaged or missing.
 - **Supervision:** Users with significant mobility impairments should always be supervised during use.
-

3. PACKAGE CONTENTS

Verify that all components are present before assembly:

- Main Frame
- Seat Plate with Cushion
- Backrest
- Armrests (2)
- Footrests (2)
- Locking Casters (4)
- Commode Bucket with Lid
- Assembly Hardware (pins, knobs)

4. SETUP AND ASSEMBLY

The XIAOSHUNER 4-in-1 Rolling Shower Chair is designed for tool-free assembly. Follow these steps:

1. **Attach the Seat Plate:** Snap the seat plate into the main frame. Ensure it is level and firmly attached before proceeding.
2. **Install the Backrest:** Insert the backrest into the designated slots behind the seat. Push down firmly until it clicks into place.
3. **Secure Armrests:** Attach the armrests to the side of the frame. Lock them using the provided pins or knobs from the hardware pack.
4. **Install Footrests:** Slide the footrests into their respective slots on the lower front frame. Ensure they are securely attached.
5. **Insert Commode Bucket:** Slide the commode bucket into the opening beneath the seat from the rear.

For a visual guide, please refer to the installation video:

Your browser does not support the video tag.

Video: Xiaoshuner Installation Steps of The Toilet Chair. This video demonstrates the step-by-step assembly process for the shower chair, including attaching the seat, backrest, and other components.

Secure & Stable: Locking 5" casters & non-tip design

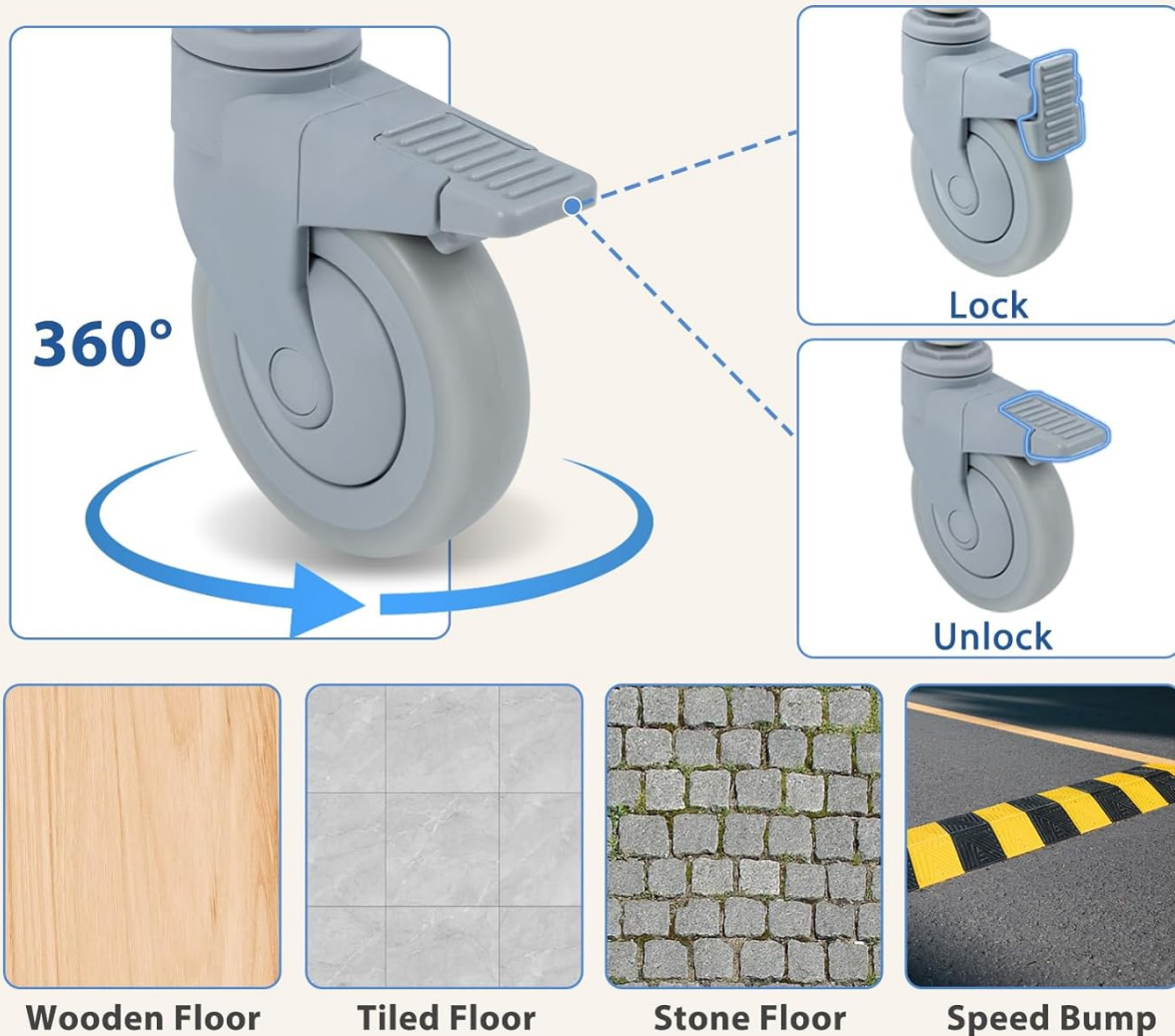


Image: Thoughtful Design. This image highlights key design features such as the ability to lift armrests for easier transfers, the hidden commode bucket, and foldable foot pedals for compact storage or user convenience.

5. OPERATING INSTRUCTIONS

This versatile chair offers four primary functions:

5.1. Shower Chair Function

- Ensure all four casters are locked before transferring the user into the chair.
- Carefully transfer the user onto the chair. The flip-up armrests can be raised to facilitate side transfers.
- Roll the chair into the shower area. Unlock the casters for movement and re-lock them once positioned.
- After showering, unlock the casters to remove the chair from the shower.

5.2. Toilet Commode Function

- Position the chair next to the bed or in a convenient location.

- Ensure the commode bucket is correctly placed under the seat.
- Lock all casters before the user sits on the chair.
- After use, remove the commode bucket from the rear for emptying and cleaning.

5.3. Toilet Riser Function

- Roll the chair over a standard toilet. Ensure the chair is centered and stable.
- Lock all casters to secure the chair in place over the toilet.
- The chair provides an elevated seating surface, making it easier for users to sit down and stand up from the toilet.

5.4. Transport Wheelchair Function

- Ensure the user is comfortably seated with footrests in the down position.
- Unlock all casters for smooth movement.
- The chair can be pushed by a caregiver to transport the user between rooms.
- Always lock the casters when stopping or leaving the user unattended.

5.5. Using Locking Casters

- To **Lock** a wheel: Press down firmly on the lever located on each caster until it clicks into the locked position.
- To **Unlock** a wheel: Lift the lever upwards until it clicks into the unlocked position.

4-in-1 Multifunctional Usage

Bathroom Wheelchair



Toilet Riser



Bedside Commode



Transport Wheelchair



Image: Secure & Stable Casters. This image illustrates the 360-degree swivel wheels with their locking and unlocking mechanisms, demonstrating their suitability for various floor types including wooden, tiled, and stone floors, as well as navigating over speed bumps.

NAVIGATES STANDARD 24' DOORWAYS EFFORTLESSLY



Image: 4-in-1 Multifunctional Usage. This image visually explains the four ways the chair can be used: as a bathroom wheelchair, a toilet riser, a bedside commode, and a transport wheelchair, showcasing its versatility.

Comfort Focused:

Cushioned seat/Rust-proof material



Image: Navigates Standard Doorways. This image shows the shower chair being easily maneuvered through a standard 24-inch doorway, highlighting its compact design for tight spaces.

6. MAINTENANCE AND CLEANING

Proper maintenance ensures the longevity and hygiene of your shower chair.

- **Cleaning:** The chair features a waterproof, wipe-clean surface. Use a mild soap and water solution to clean all surfaces. Avoid abrasive cleaners.
- **Commode Bucket:** The detachable commode bucket should be emptied and cleaned after each use. It can be washed with soap and water or a disinfectant solution.
- **Cushion:** The soft padded cushion is waterproof and can be wiped clean.
- **Drying:** After cleaning, ensure all parts are thoroughly dried to prevent water spots or residue buildup.
- **Storage:** Store the chair in a dry, clean environment when not in use.



Image: Comfort Focused Design. This image emphasizes the comfortable, anti-slip seat cushion and the durable, rust-proof plastic waterproof material used in the chair's construction, ensuring both user comfort and product longevity.

7. TROUBLESHOOTING

If you encounter any issues with your XIAOSHUNER 4-in-1 Rolling Shower Chair, refer to the following common solutions:

- **Chair feels unstable:** Ensure all assembly components are securely fastened. Verify that all four casters are locked when the chair is stationary. Check for uneven floor surfaces.
- **Wheels not rolling smoothly:** Check for debris or hair caught in the caster wheels. Clean the wheels if necessary. Ensure casters are in the unlocked position for movement.
- **Commode bucket difficult to remove/insert:** Ensure the bucket is aligned correctly with the opening. Check for any obstructions.
- **Armrests not flipping up/down easily:** Inspect the hinge mechanisms for any obstructions or damage. Ensure the locking pins (if applicable) are disengaged.

If the problem persists, please contact customer support.

8. SPECIFICATIONS

Feature	Detail
Brand	XIAOSHUNER
Model Number	SC8088A
Material	PP Plastic
Color	Lightblue
Item Weight	32 Pounds (14.5 kg)
Product Dimensions (D x W x H)	35" x 23.3" x 38.3" (88.9 cm x 59.2 cm x 97.3 cm)
Seat Depth	18 inches (45.7 cm)
Seat Width	18 inches (45.7 cm)
Seat Interior Width	17.4 inches (44.2 cm)
Maximum Weight Recommendation	350 Pounds (158 kg)

THOUGHTFUL DESIGN SUITS YOUR COMFORT

Can be lifted up the armrests



Hidden toilet



Foldable pedal

Image: Product Dimensions. This image provides a detailed diagram of the shower chair with all key dimensions labeled, including overall height, width, depth, and specific seat measurements, along with the 350 lbs weight capacity.

9. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the product packaging or contact XIAOSHUNER directly through their official channels. You can visit the [XIAOSHUNER Store](#) for more information.

