



[Manuals.plus](#) /

› [Mujocooker](#) /

› Upgraded Transport Wheelchair Lightweight Foldable User Manual

Mujocooker M-Z80

Mujocooker M-Z80 Transport Wheelchair User Manual

Lightweight Foldable Design

INTRODUCTION

This manual provides essential information for the safe and effective use, maintenance, and troubleshooting of your Mujocooker M-Z80 Lightweight Foldable Transport Wheelchair. Please read this manual thoroughly before operating the wheelchair to ensure proper function and user safety. Retain this manual for future reference.


SAFETY INFORMATION

Always prioritize safety when using the transport wheelchair. Adhere to the following guidelines:

- Ensure all locking mechanisms are engaged before use.
- Do not exceed the maximum weight capacity of 220 lbs (approximately 100 kg).
- Always use the anti-fall safety belt provided.
- Avoid using the wheelchair on steep inclines or uneven terrain without assistance.
- Regularly inspect the brakes, wheels, and frame for any signs of wear or damage.
- Keep hands and fingers clear of moving parts during folding and unfolding.



1 SAFTY BUCKLE
Safty Buckle Keep The Wheelchair
Fix And Not Easy To Deformation

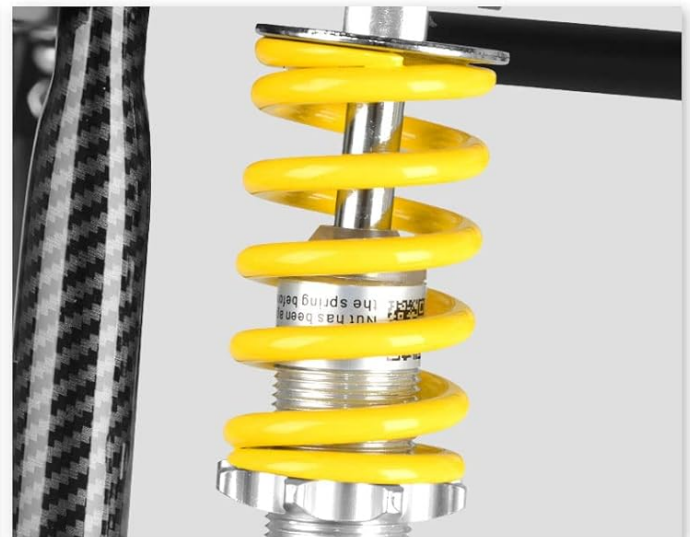

**2 SPRING SHOCK
ABSORPTION DESIGN**
Smooth and Comfortable Driving



3 STABLE STRUCTURE
Dispersed Stress Structure Is Stable



4 ANTI-FALL SAFETY BELT
Ensure The Safety of Users During use

Figure 1: Key safety features of the Mujocooker M-Z80 wheelchair, including the safety buckle, spring shock absorption design, stable stress structure, and anti-fall safety belt, all designed to enhance user security and comfort.

PRODUCT OVERVIEW AND COMPONENTS

The Mujocooker M-Z80 is a lightweight, foldable transport wheelchair designed for ease of use and portability. Key components include:

- **Frame:** Constructed from carbon fiber for durability and lightweight design.
- **Seat and Backrest:** Made from waterproof material for easy cleaning.
- **Wheels:** 10-inch rear honeycomb solid wheels and 6-inch 360-degree swivel rubber front wheels.
- **Suspension System:** Oversized spring dampening system on rear wheels and spring shock absorbers on front wheels for comfort on various terrains.
- **Handbrakes:** Dual handbrake system for flexible operation and quick response.
- **Storage Pockets:** Includes a travel bag, a luxury back pocket, and an under-seat storage pocket.

- **Telescopic Handle:** Located under the seat for easy transport when folded.



Figure 2: The Mujocooker M-Z80 Transport Wheelchair, illustrating its unfolded configuration ready for use and its compact folded state for storage or transport.



Figure 3: The Mujocooker M-Z80 wheelchair highlighting its integrated storage solutions, including a travel bag, a luxury back pocket, and a convenient under-seat storage pocket.

SETUP AND UNFOLDING

1. **Unpack:** Carefully remove the wheelchair from its packaging.
2. **Unfold Frame:** Gently pull the armrests outwards to begin unfolding the frame. The seat will start to drop into place.
3. **Secure Seat:** Press down firmly on both sides of the seat rails until the seat is fully flat and locked into position. Ensure there is no wobble.
4. **Adjust Footrests:** Lower the footrests to their operational position. They should click into place.
5. **Inspect:** Before first use, visually inspect all components to ensure they are securely attached and free from damage.




Figure 4: Visual representation of the Mujocooker M-Z80 transitioning from its folded, compact form to its fully unfolded, ready-to-use state.

OPERATING INSTRUCTIONS

Folding the Wheelchair

1. **Raise Footrests:** Lift the footrests to their upright position.
2. **Lift Seat:** Locate the folding mechanism (often a strap or handle) under the seat and pull upwards. The seat will begin to fold.
3. **Collapse Frame:** Push the armrests inwards to fully collapse the frame until it is compact.
4. **Engage Telescopic Handle:** Extend the telescopic handle located under the seat for easy pulling, similar to a suitcase.

Telescopic Handle

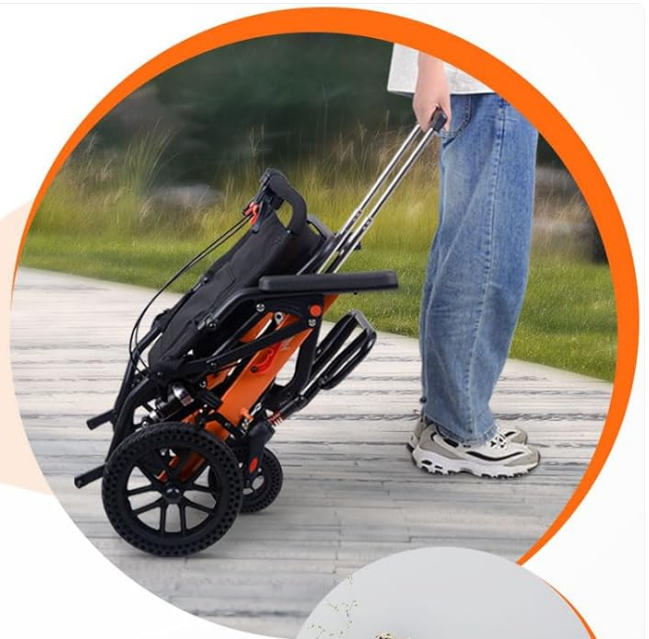

Lightweight
and Compact



Small and
flexible



Easy for storage
and transport



ENJOY YOUR TRAVEL!

Figure 5: The Mujocooker M-Z80 in its folded state, being transported effortlessly using its integrated telescopic handle, emphasizing its enhanced portability.

Braking System

The M-Z80 features an upgraded dual handbrake system for reliable stopping and parking.

- **To Engage Brakes (Stop/Park):** Hold the brake handle up and fix the braking by pushing the lever forward until it locks into the parking position. This will prevent the wheels from moving.
- **To Release Brakes:** Press the trigger release button on the brake handle to disengage the brake and allow the wheels to move freely.

Double Handbrake

Flexible operation and quick response!



Figure 6: Close-up view of the Mujocooker M-Z80's double handbrake mechanism, illustrating how to engage the brake by holding the handle up and fixing it, and how to release it by pressing the trigger.

DOUBLE HANDBRAKE

Flexible Operation and Quick Response!



① Hold The Brake Handle Up And Fix The Braking

② Press The Trigger Release The Braking

Figure 7: Another perspective of the handbrake, emphasizing the brake handle (1) and trigger (2) for intuitive control of the wheelchair's movement and parking.

Maneuvering

The 360-degree swivel front wheels allow for nimble cornering and easy navigation in various environments. Push the handles to move forward, and gently steer to change direction.

Using Storage Pockets

- **Travel Bag:** This bag is designed to store the folded wheelchair for convenient transport.
- **Luxury Back Pocket:** Ideal for carrying larger items such as computers, books, or documents.
- **Under-Seat Storage Pocket:** Suitable for smaller essentials like chargers, cups, or tissues.

MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your wheelchair.

- **Cleaning:** The seat and backrest are made of waterproof material and can be easily cleaned with a damp cloth and mild soap. Avoid abrasive cleaners.
- **Wheel Inspection:** Periodically check the wheels for debris, wear, or damage. Ensure they spin freely.
- **Brake Check:** Test the handbrakes regularly to ensure they engage and release smoothly and effectively.
- **Frame Inspection:** Check the carbon fiber frame for any cracks, bends, or loose fasteners. Tighten any loose screws as needed.
- **Storage:** Store the wheelchair in a dry, cool place away from direct sunlight when not in use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Wheelchair does not unfold completely.	Seat not fully pushed down or frame not fully extended.	Ensure the seat is pressed down firmly until it locks. Check for obstructions.
Brakes are not holding or are too loose.	Brake mechanism needs adjustment or is worn.	Inspect the brake cables and levers. If issues persist, contact customer support.
Wheels are not rolling smoothly.	Debris in wheels or worn bearings.	Remove any debris. If noise or resistance continues, professional inspection may be required.
Wheelchair feels unstable.	Loose fasteners or damaged frame component.	Check all bolts and screws for tightness. Discontinue use if frame damage is suspected and contact support.

SPECIFICATIONS

Feature	Detail
Model	M-Z80
Weight Capacity	220 lbs (approx. 100 kg)
Frame Material	Carbon Fiber
Rear Wheels	10-inch honeycomb solid wheels
Front Wheels	6-inch 360-degree swivel rubber wheels
Unfolded Dimensions (L x W x H)	Approx. 33.5" x 21" x 37.4"
Folded Dimensions (L x W x H)	Approx. 11.4" x 28.3" (height varies slightly with handle)
Seat/Backrest Material	Waterproof fabric

PRODUCT DIMENSION

Weight Capacity **220lbs**



Figure 8: Detailed product dimensions for the Mujocooker M-Z80, illustrating its size when unfolded and its compact dimensions when folded, along with its 220 lbs weight capacity.

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the product packaging or contact Mujocooker customer service directly. Keep your purchase receipt as proof of purchase.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a medical device and should be used as directed.