

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [iBstone](#) /

› [iBstone Nova Hearing Aids User Manual](#)

iBstone Nova

iBstone Nova Hearing Aids User Manual

Model: Nova

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your iBstone Nova Hearing Aids. Please read this manual thoroughly before using the device to ensure proper operation and to maximize your hearing experience. The iBstone Nova Hearing Aids are designed to provide clear sound with smart noise cancelling technology, offering comfort and ease of use for adults and seniors.

PACKAGE CONTENTS

Verify that all items listed below are included in your package:

What Included

Full Package • Real-human Support • Operation Video



Charging Case



Hearing Aids



Anti-lost Strap



Cleaning Brush



Wax Guard Filters




USB Cable



User Manual



Silicone Ear Domes

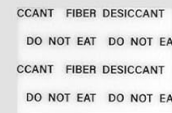

Dry Tablets

Image: The package contents include the charging case, two hearing aids, an anti-lost strap, a cleaning brush, wax guard filters, a USB cable, a user manual, silicone ear domes, and dry tablets.

- iBstone Nova Hearing Aids (Left and Right)
- Portable Charging Case
- USB Charging Cable
- Silicone Ear Domes (various sizes)
- Wax Guard Filters
- Cleaning Brush
- Anti-Lost Lanyard
- Dry Tablets (for dehumidifying capsule)
- User Manual (this document)

SETUP

1. Initial Charging

Before first use, fully charge the hearing aids and the charging case.

1. Connect the USB charging cable to the charging case and a power source.

2. Place the hearing aids into their respective slots in the charging case (L for left, R for right).
3. The indicator lights on the case will illuminate to show charging status. A full charge typically takes about 1.5 hours for the hearing aids. The charging case provides approximately 3 full recharges for the hearing aids.

Fast Charging & Long Battery Life

1.5H Fast Charge = **80H** Use



Charging time ≤ 1.5 hrs
Using time ≥ 20 hrs



Charging Case Backup Power ≥ 80 hrs
Charging Case Provides 3 Times recharging for a pair



Light up=Charging
No light up = Fully charge
1 light = 25% remaining power in the case

Image: The hearing aids are shown charging in their case, connected via a USB cable. Indicator lights on the case display the battery level. One light indicates 25% remaining power in the case.

2. Selecting Ear Domes

Choose the ear dome size that provides the most comfortable and secure fit, and minimizes sound leakage or whistling.

- Experiment with the different sizes and styles of silicone ear domes provided.
- A proper fit is crucial for optimal sound quality and to prevent feedback.



Image: A four-step guide illustrates how to fit the hearing aids: 1. Choose the right size of ear domes. 2. Place the device into the ear canal correctly, noting L/R markings. 3. Adjust for comfort and security. 4. Press the button to adjust volume/mode.

3. Inserting the Hearing Aids

Identify the left (L) and right (R) hearing aids. Gently insert each aid into the corresponding ear canal until it feels secure and comfortable.

OPERATING INSTRUCTIONS

1. Power On/Off

- **Automatic On/Off:** The hearing aids automatically turn on when removed from the charging case and turn off when placed back into the case.

2. Adjusting Volume and Modes

The iBstone Nova features a single multi-function button for easy control.

One Button, All Functions

Auto On/Off · 9 Volume Levels · 3 Listening Modes · Memory Function



Image: A close-up of the iBstone Nova hearing aid highlights its single function button. Text describes its functions: quick press for volume, 2-second press for modes, auto on/off, and auto memory function.

- **Volume Adjustment:**

- A **quick press** of the function button will cycle through 9 volume levels.
- Adjust to a comfortable listening level.

- **Mode Selection:**

- A **2-second press** of the function button will cycle through 3 listening modes:
 - **Quiet Mode:** For calm environments.
 - **Noisy Mode:** Enhances speech in environments with background noise.
 - **Outdoor Mode:** Optimized for open-air settings.
- The device has an auto-memory function, retaining your last setting when turned off and on.

MAINTENANCE

1. Daily Cleaning


- Use the provided cleaning brush to gently remove any earwax or debris from the ear domes and the hearing aid body.
- Wipe the hearing aids with a soft, dry cloth. Do not use water or cleaning solutions directly on the devices.

2. Dehumidifying


The charging case includes a built-in dehumidifying capsule to protect your hearing aids from moisture.

- Place the dry tablets into the designated compartment within the charging case.
- Store your hearing aids in the charging case with the dry tablets when not in use to extend their lifespan.


Extra Design Features




Anti-Lost Strap
Never worry about losing



Built-in Dehumidifying Capsule
Extent the lifetime of your devices



Portable Charging Case
Small case, big power



Abundant Ear domes
get your best fitting

Image: This image displays several features, including the anti-lost strap, the built-in dehumidifying capsule in the charging case, the portable charging case, and various ear domes for optimal fitting.

3. Replacing Wax Guard Filters

Regularly check the wax guard filters and replace them if they appear clogged with earwax. Refer to the detailed instructions in the full user manual for replacement procedures.

TROUBLESHOOTING

If you experience any issues with your iBstone Nova Hearing Aids, please refer to the following common solutions:

Problem	Possible Cause	Solution
No Sound or Weak Sound	Low battery, clogged ear dome/wax guard, incorrect volume setting.	Charge the hearing aids. Clean or replace the ear dome/wax guard. Increase volume.
Whistling or Feedback	Improper fit, ear dome too small, earwax blockage.	Re-insert the hearing aid. Try a larger ear dome. Clean ear canal or hearing aid.
Intermittent Sound	Low battery, loose connection in charging case, moisture.	Ensure full charge. Check placement in charging case. Use dehumidifying capsule.
Hearing Aid Not Charging	Improper placement in case, faulty cable/adaptor, dirty charging contacts.	Ensure proper alignment in case. Try a different cable/adaptor. Clean charging contacts.

If the problem persists, please contact customer support.

SPECIFICATIONS

- **Model:** Nova
- **Manufacturer:** iBstone
- **Battery Type:** Rechargeable (1 A battery included)
- **Charging Time:** Approximately 1.5 hours for hearing aids
- **Usage Time:** Up to 80 hours per charge (with charging case backup)
- **Volume Levels:** 9
- **Listening Modes:** 3 (Quiet, Noisy, Outdoor)
- **Dimensions (Package):** 6.73 x 4.21 x 2.13 inches
- **Weight (Package):** 11.36 ounces


WARRANTY AND SUPPORT

iBstone stands behind its products. For any questions or assistance, please contact our customer support team.


- **Customer Support:** U.S.-based phone support is available to assist with any after-sales problems.
- **Warranty Information:** Please refer to the warranty card included in your package or contact customer support for detailed warranty terms.

Lifetime Support from iBstone Team


Any question about our hearing aids, contact us by text, email, message or phone call



Real-human Instant Service



US Based Replacement & Repairment



Comprehensive Manual & Video Instruction


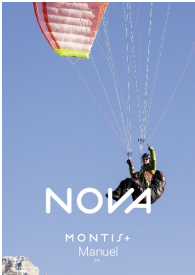







Image: A group of professionals in lab coats represents the iBstone support team, emphasizing real-human instant service, US-based replacement & repairment, and comprehensive manual & video instruction.

For further assistance, please visit the official iBstone website or contact the support team directly.

© 2025 iBstone. All rights reserved.

Related Documents - Nova

	<p>NOVA MONTIS+: Manuel pour Harnais de Parapente Hike & Fly</p> <p>Découvrez le harnais de parapente NOVA MONTIS+, idéal pour le hike & fly. Ce manuel détaillé couvre les caractéristiques, l'utilisation, la sécurité et l'entretien de votre équipement NOVA.</p>
	<p>Project Guide: How to Turn a Champagne Flute with NOVA Tools</p> <p>Learn how to craft an elegant champagne flute with this step-by-step woodworking project guide. Ideal for honing spindle turning techniques, this project emphasizes precision and balance, resulting in a functional and beautifully crafted piece perfect for celebrations. Features NOVA tools and expert instructions by Jeff Hornung.</p>
	<p>Manuale d'Uso e Manutenzione: Nova 15-22 Spezzatrice Manuale</p> <p>Scopri il manuale d'uso e manutenzione per la Nova 15-22, una spezzatrice manuale professionale progettata per panetterie, pasticcerie e pizzerie. Guida completa per l'installazione, l'uso sicuro e la manutenzione.</p>
	<p>NOVA Eco-NOVA Bidet Seat User Manual and Installation Guide</p> <p>This document provides comprehensive instructions for the NOVA Eco-NOVA Bidet Seat, including installation requirements, operation details, maintenance tips, troubleshooting guide, and warranty information.</p>

	<p>Nova Folding Cup Holder Assembly and Warranty Information</p> <p>Detailed assembly instructions and warranty information for the Nova Folding Cup Holder (Item # CH-1101), designed for Nova Medical rollators. Learn how to attach the cup holder to Express and Star series rollators.</p>
	<p>NOVA Traveler Rolling Walker User Guide</p> <p>Comprehensive user guide for the NOVA Traveler Rolling Walker (Model 4900). Includes essential safety warnings, step-by-step assembly instructions, guidance on attaching accessories, operating the walker, adjusting hand brakes, and details on the product warranty.</p>

Documents - iBStone – Nova

no relevant documents