

Eumspo JH-A490

Eumspo JH-A490 In-Ear Rechargeable Hearing Aids User Manual

Model: JH-A490

1. INTRODUCTION

Thank you for choosing Eumspo JH-A490 In-Ear Rechargeable Hearing Aids. These devices are designed to enhance auditory clarity for individuals with mild to moderate hearing loss. Featuring an ergonomic in-ear design, automatic power on/off, intelligent noise reduction, and anti-howling technology, they aim to provide a comfortable and effective listening experience. Please read this manual thoroughly before using your hearing aids to ensure proper operation and maintenance.

2. PACKAGE CONTENTS

Verify that all items are present in your package:

- Hearing Aids x 2 (Left and Right)
- Charging Case x 1
- Eartips x 6 (various sizes)
- USB Cable x 1
- Cleaning Tool x 1
- Instruction Manual x 1 (this document)
- Quick Start Guide x 1

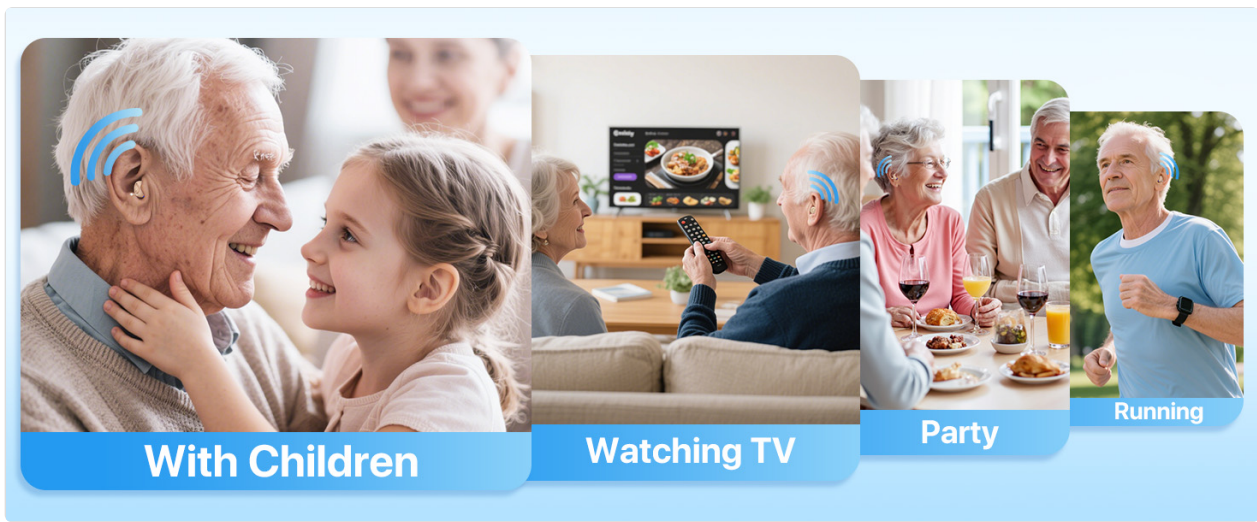


Figure 2.1: Package Contents Overview

3. PRODUCT OVERVIEW

Familiarize yourself with the components of your Eumspo JH-A490 hearing aids and charging case.



Figure 3.1: Hearing Aid and Charging Case Components

- **Eartip:** Soft silicone tips for comfortable fit and sound delivery. Available in multiple sizes.
- **Charging Contacts:** Connect the hearing aids to the charging case.
- **Sound Inlet:** Where sound enters the hearing aid.
- **Multi-function Button:** Used for volume adjustment and mode switching.
- **Charging Slots:** Designated spaces in the case for left (L) and right (R) hearing aids.
- **Charging Indicators:** LEDs on the case show charging status for the case and individual hearing aids.
- **Charging Port:** USB-C port on the case for external power.

4. CHARGING

Before first use, fully charge the hearing aids and the charging case.

1. **Charge the Charging Case:** Connect the provided USB cable to the charging port on the case and a suitable power adapter (not included). The charging indicators on the case will illuminate to show charging progress.
2. **Charge the Hearing Aids:** Place the left (blue) and right (red) hearing aids into their respective charging slots in the case. Ensure the charging contacts align. The hearing aids will begin charging automatically. A full charge takes approximately 2 hours.

Once fully charged, the hearing aids provide up to 16 hours of continuous use. The portable charging case allows for multiple recharges on the go.



Figure 4.1: Charging and Battery Duration

5. FITTING AND WEARING

Proper fit is crucial for comfort and optimal sound quality, and to prevent whistling.

1. **Select Eartips:** Choose the eartip size (S, M, L) that provides a snug yet comfortable seal in your ear canal. Experiment with different sizes to find the best fit.
2. **Identify Left and Right:** The left hearing aid is typically marked with blue, and the right with red.
3. **Insert Hearing Aid:** Hold the hearing aid with two fingers and gently insert it into your ear canal. Rotate slightly until it feels secure and comfortable. Ensure the eartip forms a good seal.



Comfortable In The Ear

no pain after wearing for a long time smooth communication and clear listening

-  Auto on /off
-  Made of medical -Grade materials
-  5-Level volume One-key to adjust
-  Red for right Blue for left
-  Ideal for seniors

Figure 5.1: How to Use and Fit the Hearing Aids

Note: Hearing aids may produce a whistling sound if not worn correctly. If this occurs, reinsert them into the ear canal, ensuring a proper seal.

6. OPERATION

The Eumspo JH-A490 hearing aids are designed for simple, one-button control.

1. **Power On/Off:** The hearing aids automatically power on within 10-15 seconds when removed from the charging case. They power off automatically when placed back into the charging case.
2. **Adjust Volume:** Short press the multi-function button to cycle through 5 volume levels. Each press increases the volume. The device will remember your last setting.
3. **Switch Modes:** Quickly press the multi-function button twice to switch between listening modes. The default mode at startup is Standard Mode. The second mode is Restaurant Mode.
 - **Standard Mode:** Enhances daily conversations and general listening.
 - **Restaurant Mode:** Reduces background noise to improve speech clarity in noisy environments.

Easy to Operate

one button control



Power On:

when the device is removed from the charging case, the hearing aid will automatically turn on

Adjust Mode:

short press twice, the device will turn on the default standard mode

Adjust Volume:

short press

Figure 6.1: Easy One-Button Control



HOW TO USE?

1 Choose right size of earplugs for your ears.

2 When the device is removed from the charging case, the hearing aids will automatically turn on within 10 s.

3 Hold the device with two fingers and insert it gently into your ear.

4 Press the button to cycle volume levels. x1

5 Quickly press the button twice to switch scene modes, standard mode and restaurant mode. (Default standard mode at startup) x2

Figure 6.2: Volume Adjustment Levels

7. MAINTENANCE

Regular cleaning and proper care will ensure the longevity and optimal performance of your hearing aids.

- **Daily Cleaning:** After each use, gently wipe the hearing aids with a soft, dry cloth. Use the provided cleaning tool to remove any earwax or debris from the eartips and sound inlet.
- **Eartip Care:** Periodically remove the eartips and clean them with a mild soap solution. Rinse thoroughly and ensure they are completely dry before reattaching.
- **Storage:** When not in use, store the hearing aids in their charging case to protect them from dust and damage.
- **Avoid Moisture:** Keep the hearing aids and charging case away from water and excessive humidity. Do not submerge them in liquid.

8. TROUBLESHOOTING

If you encounter any issues, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
No sound or low sound	Low battery; Volume too low; Eartip blocked by earwax; Device not powered on.	Charge the hearing aids; Increase volume; Clean eartip and sound inlet; Ensure device is removed from charging case.
Whistling/Feedback sound	Improper fit; Eartip too small or loose; Earwax buildup in ear canal.	Reinsert the hearing aid for a snug fit; Try a larger eartip size; Consult a healthcare professional for earwax removal if necessary.
Hearing aid not charging	Poor contact with charging pins; Charging case battery low; Faulty USB cable/adaptor.	Ensure hearing aids are correctly seated in the case; Charge the charging case; Try a different USB cable or power adapter.
Difficulty hearing in noisy environments	Incorrect listening mode.	Switch to Restaurant Mode by double-pressing the multi-function button.

If the problem persists after trying these solutions, please contact customer support.

9. SPECIFICATIONS

Feature	Detail
Model Number	JH-A490
Design Type	In-Ear (ITE)
Target Hearing Loss	Mild to Moderate
Charging Time	Approximately 2 hours
Working Time (Single Charge)	Up to 16 hours
Noise Reduction	Smart Noise Reduction Chip, Dual Modes (Standard, Restaurant)
Feedback Control	Intelligent Anti-Howling Technology
Volume Levels	5 adjustable levels

Feature	Detail
Automatic Features	Auto On/Off
Manufacturer	LIKELINKS TRADING LIMITED
Country of Origin	China

10. IMPORTANT SAFETY INFORMATION

- Consult a healthcare professional before using hearing aids if you experience sudden hearing loss, ear pain, discharge, or dizziness.
- Do not attempt to repair or modify the device yourself. Contact customer support for assistance.
- Keep out of reach of children and pets. Small parts may pose a choking hazard.
- Avoid exposing the devices to extreme temperatures, direct sunlight, or strong magnetic fields.
- Clean your ears regularly to prevent earwax buildup, which can affect hearing aid performance.




11. CUSTOMER SUPPORT




For any questions, technical assistance, or warranty inquiries regarding your Eumspo JH-A490 Hearing Aids, please contact our dedicated customer service team. Refer to the contact information provided in your Quick Start Guide or on the product packaging.

We are committed to providing detailed guidance and solutions to ensure your complete satisfaction.

© 2025 Eumspo. All rights reserved.

Related Documents - JH-A490

	<p>Eumspo G3 Wireless Headset User Manual</p> <p>User manual for the Eumspo G3 Wireless Headset, covering product lists, detailed description, pairing instructions, multi-function button operations, charging methods, how to wear the earpiece, troubleshooting, and safety precautions.</p>
	<p>Jinghao Hearing Aids JH-A51 User Manual and Instructions</p> <p>Comprehensive user manual for the Jinghao Hearing Aids JH-A51 model, covering setup, operation, safety, and maintenance for optimal hearing assistance.</p>
	<p>JINGHAO JH-A39 Hearing Aid: Instruction Manual and User Guide</p> <p>Comprehensive instruction manual for the JINGHAO JH-A39 hearing aid. Learn about safety, usage, maintenance, troubleshooting, and technical specifications for optimal hearing assistance.</p>

	<p>newgen medicals ZX-7702-675 Medizinische IdO-Hörgeräte Bedienungsanleitung</p> <p>Diese Bedienungsanleitung bietet detaillierte Informationen zur sicheren und effektiven Nutzung des newgen medicals ZX-7702-675 medizinischen IdO-Hörgeräts, einschließlich Einrichtung, Pflege und Fehlerbehebung für optimale Hörunterstützung.</p>
	<p>OTC Hearing Aid JH-DW2ET Instruction Manual</p> <p>Comprehensive user manual for the OTC Hearing Aid Model JH-DW2ET, covering setup, usage, care, troubleshooting, technical specifications, and regulatory information.</p>
	<p>JH-A17 Hearing Amplifier Instruction Manual</p> <p>Comprehensive instruction manual for the JH-A17 Hearing Amplifier by Huizhou Jinghao Medical Technology Co., Ltd., covering safety instructions, usage guidelines, maintenance, troubleshooting, and technical specifications.</p>