

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Coniler](#) /

› [Coniler A-26 Mini Digital Hearing Aids Instruction Manual](#)

Coniler A-26

Coniler A-26 Mini Digital Hearing Aids Instruction Manual

Model: A-26

1. INTRODUCTION

Thank you for choosing the Coniler A-26 Mini Digital Hearing Aids. These advanced hearing devices are designed to provide clear sound quality and comfortable wear, enhancing your listening experience in various environments. Featuring a 16-channel digital chip, noise cancellation, and multiple modes, the A-26 offers personalized sound amplification. Please read this manual thoroughly before use to ensure proper operation and maintenance.



Image 1: Coniler A-26 Mini Digital Hearing Aids and Charging Case.

2. PRODUCT FEATURES

- **16-Channel Intelligent Digital Chip:** Provides superior performance by actively identifying, eliminating, and suppressing intermittent noise for stable and clear sound.
- **Exclusive AI Dynamic Noise Reduction Algorithm:** Focuses on the sound source, reduces wind noise, filters out background noise, and accurately locates sound sources for scene recognition.
- **Super Clear Sound Quality:** Ensures smooth, uninterrupted, and crystal-clear sound for conversations, TV, and environmental sounds.
- **Lightweight, Portable, and Comfortable:** Ergonomically designed silicone earplugs offer long-lasting, pain-free wear and fit snugly in most ear shapes.
- **One-Touch Adjustment of Four Modes:** Easily switch between Normal, Outdoor, Noise Cancelling, and Indoor modes with a single press.
- **Adjustable Volume:** Five volume levels can be adjusted with push-button controls.
- **Long-Lasting Battery Life:** Up to 12 hours of use on a single charge, with the portable charging case providing up to 6 additional charges for a total of 72 hours of listening time.
- **Type-C Charging Port:** Ensures fast and convenient charging.

Top ten performances Zero-Noise Clarity

-  New digital chip
-  16 channels
-  AI dynamic noise
-  Third gear selection
-  WDRD Dynamics
-  Long-lasting battery life
-  Simple operation
-  Charging compartment
-  Multiple scenarios



Image 2: Key features of the Coniler A-26 hearing aids.

Powerful Features Because the Strong "Core" is Dirty

Digital chip

- ✓ Millisecond reduction
- ✓ Whistling inhibition
- ✓ Speech tracking
- ✓ Algorithm

Analog chip

- Bad noise reduction x
- Howling current x
- Recognition ambiguity x
- Can't hear clearly x

WDRD
Wide Dynamic
Compression

Stereo
Human Ear Experience

Noise Reduction
Filter Noise

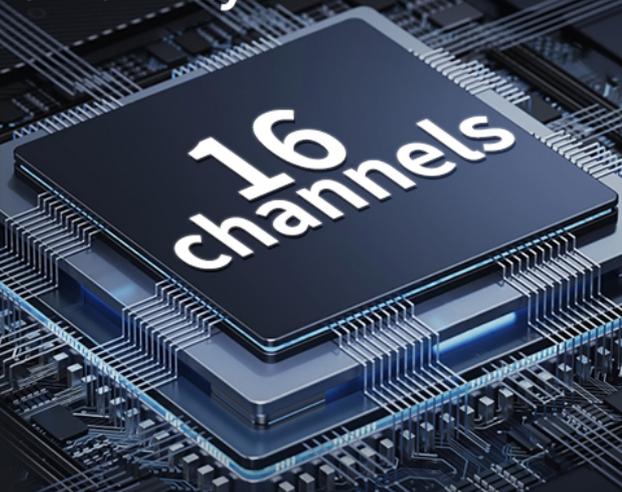


Image 3: Digital chip technology overview.

Fine intelligent noise reduction Multi-channel digital chip



Image 4: Intelligent noise reduction capabilities.

3. PACKAGE CONTENTS

The Coniler A-26 package includes:

- 1 Pair of Coniler A-26 Mini Digital Hearing Aids
- 1 Portable Charging Case
- Multiple sizes of silicone earplugs
- 1 Type-C Charging Cable
- 1 Cleaning Brush
- User Manual

4. SETUP

4.1 Charging the Device

1. Before first use, fully charge the hearing aids and the charging case.
2. Connect the provided Type-C charging cable to the charging case and a suitable power source (e.g., USB wall adapter, computer USB port).
3. The LED indicators on the charging case will show the charging status. Once fully charged, the indicators will typically show a solid light or turn off, depending on the model.
4. Place the hearing aids into their respective slots in the charging case. They will begin charging automatically.

4.2 Selecting and Installing Earplugs

1. Choose the silicone earplug size that provides the most comfortable and secure fit in your ear canal. A proper fit is crucial for optimal sound quality and to prevent whistling.
2. Gently attach the selected earplug to the nozzle of the hearing aid.

4.3 Wearing the Hearing Aids

1. Carefully insert the hearing aid into your ear canal until it feels snug and comfortable.
2. Ensure the hearing aid is securely seated to prevent it from falling out and to maximize sound clarity.

5. OPERATING INSTRUCTIONS

5.1 Power On/Off

- The hearing aids typically power on automatically when removed from the charging case and power off when placed back into the case.
- Refer to the specific instructions in the included user manual for manual power control if available.

5.2 Adjusting Volume

The Coniler A-26 features push-button controls for volume adjustment.

- Press the button once to cycle through the 5 available volume levels.
- Select the volume level that provides optimal listening comfort and clarity for your environment.



Image 5: Simple volume operation and controls.

5.3 Switching Modes

The hearing aids offer four distinct listening modes to adapt to different environments:

- **Normal Mode:** For everyday listening in quiet to moderately noisy environments.
- **Outdoor Mode:** Optimized for outdoor use, reducing wind noise and enhancing distant sounds.
- **Noise Cancelling Mode:** Actively reduces background noise for clearer conversations in loud settings.
- **Indoor Mode:** Tailored for indoor environments, such as watching TV or small group discussions.

To switch between modes, press and hold the control button on the hearing aid until the desired mode is

activated. The device may provide an audible cue or indicator light change to confirm the mode selection.

Suitable for a Variety of Scenarios Multiple Mode Options

Adaptive Environment
Automatic Frequency Modulation




Bus station



Normal indoor
communication



The sounds of
nature

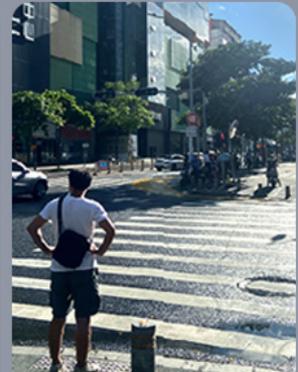

Highway
passing

Image 6: Multiple mode options for various scenarios.

6. MAINTENANCE

6.1 Cleaning

- Regularly clean your hearing aids to ensure optimal performance and hygiene.
- Use the provided cleaning brush to remove any earwax or debris from the earplugs and sound output.
- Wipe the body of the hearing aids with a soft, dry cloth. Do not use liquids or harsh chemicals.

6.2 Storage

- When not in use, store the hearing aids in their charging case to protect them and keep them charged.
- Store in a cool, dry place away from direct sunlight and extreme temperatures.

7. TROUBLESHOOTING

If you encounter any issues with your Coniler A-26 hearing aids, please refer to the following common solutions:

- **No Sound or Weak Sound:**
 - Ensure the hearing aids are fully charged.
 - Check if the volume is set to an audible level.
 - Clean the earplugs and sound output to remove any blockages.
 - Ensure the earplugs are correctly seated in your ear canal.
- **Whistling or Feedback:**
 - This often indicates an improper fit. Try a different size earplug or re-insert the hearing aid to

ensure a snug seal.

- Ensure there are no obstructions near the microphone.

- **Inconsistent Sound:**

- Ensure the hearing aids are fully charged.
- Clean the devices thoroughly.
- Try resetting the devices by placing them back in the charging case for a few seconds, then removing them.

- **Charging Issues:**

- Verify the Type-C cable is securely connected to both the case and the power source.
- Ensure the charging contacts on both the hearing aids and the case are clean and free of debris.

If the issue persists, please contact customer support for further assistance.

8. SPECIFICATIONS

Feature	Specification
Model Number	A-26
Digital Chip	16 Channel
Material	Acrylonitrile Butadiene Styrene (ABS)
Style	Invisible CIC (Completely-in-Canal)
Power Source	Battery Powered
Single Charge Usage	Up to 12 hours
Total Usage with Charging Case	Up to 72 hours
Charging Port	Type-C
Dimensions (Package)	4.96 x 3.78 x 2.67 inches
Weight (Package)	4.8 ounces

9. WARRANTY AND SUPPORT

For warranty information, technical support, or any product-related inquiries, please refer to the documentation included in your product packaging or contact the manufacturer directly. Keep your proof of purchase for warranty claims.