



[Manuals.plus](#) /

> [VITASY](#) /

> VITASY Ingrown Toenail Treatment Kit with Upgraded 6-Angle Shape Instruction Manual

## VITASY B0FPJM97QV

# VITASY Ingrown Toenail Treatment Kit with Upgraded 6-Angle Shape Instruction Manual

Model: B0FPJM97QV

## INTRODUCTION

---

The VITASY Ingrown Toenail Treatment Kit provides an effective solution for ingrown toenails. This kit features an upgraded 6-angle titanium alloy design, engineered to increase correction and adhesion strength, significantly shortening the treatment period. It is designed to offer painless correction and relief from discomfort associated with ingrown toenails, including unilateral or bilateral paronychia and curved nails.

## SAFETY INFORMATION

---

This product is for external use only. Keep out of reach of infants. If you are diabetic or have circulation disorders, please consult a professional doctor before use. Observe the treated area for 2 hours after application; if you feel uncomfortable, adjust or remove the corrector. Do not exercise vigorously within 24 hours of application. If the glue is accidentally smeared on your hands or other body parts, do not be nervous. Wipe with paper towels, apply the thinner to the place where the glue is, wait about 2 minutes, and gently scrub the glue with your hands.

## PRODUCT COMPONENTS

---

The VITASY Ingrown Toenail Treatment Kit includes the following items:

- Ingrown Toenail Correction Patches (4 multi-sizes: 0.51in, 0.59in, 0.67in, 0.79in - 2 pieces per size)
- Diluent (8ml)
- Glue Adhesive (8ml)
- Clean Pads (8 pieces)
- Fingertip Gloves (6 pieces)
- Auxiliary Patch (10 pieces)
- Polishing Strip (2 pieces)
- Anti-friction Stickers (1 sheet)

- Tweezers (1 piece)




Figure 1: Complete contents of the VITASY Ingrown Toenail Treatment Kit.

## SETUP AND PREPARATION


1. **Clean the Toenail:** Use an alcohol pad (Clean Pad) to thoroughly clean the surface of the affected toenail and the surrounding area. Ensure the nail is dry and free of oils or debris.
2. **Roughen the Nail Surface:** Gently polish the entire area where the correction patch will be applied using the provided polishing strip. This helps improve adhesion.
3. **Select Correct Patch Size:** Choose an ingrown toenail correction patch that best fits the width of your toenail. The kit includes four optional lengths (0.51in, 0.59in, 0.67in, 0.79in).

# 4 OPTIONAL LENGTHS

Adaptable to different nail widths

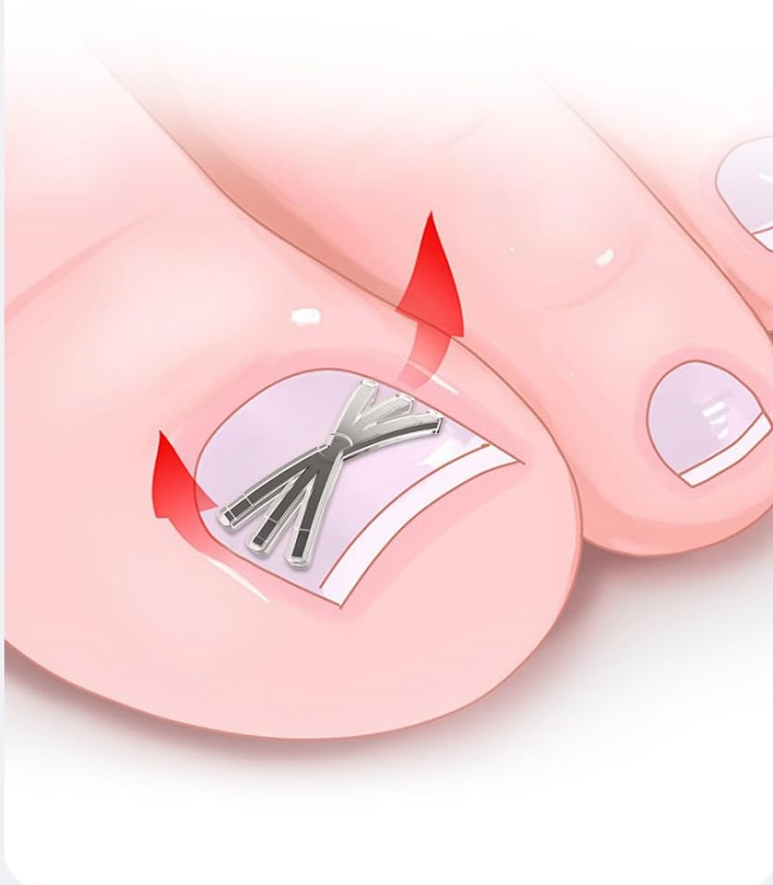


2 PCS

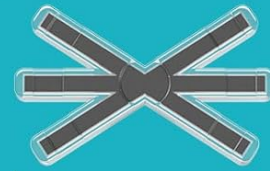


0.51in

2 PCS

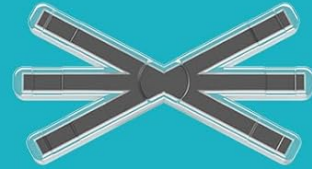


0.59in

2 PCS

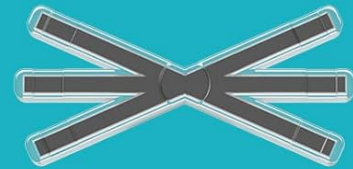

0.67in

2 PCS

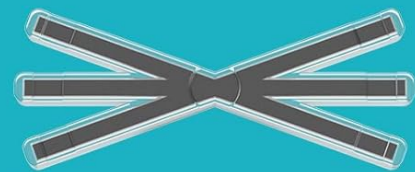

0.79in

Figure 2: Available patch sizes and proper placement on the toenail.

## OPERATING INSTRUCTIONS (APPLICATION)

Follow these steps carefully for optimal results:

1. **Prepare the Patch:** Remove one correction patch from its adhesive strip.
2. **Apply Glue to Patch:** Apply a thin, even layer of the Glue Adhesive to the center of the correction patch.
3. **Position the Patch:** Carefully place the patch onto the center of the affected toenail, ensuring it covers the curved area.
4. **Secure the Patch:** Press down firmly on both ends of the correction patch for 15-20 seconds, then slowly release. This helps the patch adhere securely and begin its corrective action.
5. **Reinforce Adhesion:** Apply additional Glue Adhesive 3-5 times over the patch to further reinforce its bond with the toenail.
6. **Apply Anti-friction Stickers:** For improved comfort and to reduce friction with socks, place the anti-friction stickers over the applied correction patch.

Your browser does not support the video tag.

Video 1: Step-by-step guide on applying the VITASY Ingrown Toenail Correction Patch. This video demonstrates cleaning the nail, applying the patch with adhesive, and pressing to secure it. It highlights that both inlaid toenails operate the same way.



Figure 3: Detailed usage steps for the ingrown toenail correction patch.

## MAINTENANCE

---

The correction patches are designed to be durable and waterproof, allowing for normal daily activities such as walking, sports, bathing, and showering. It is recommended to replace the correction patch once a month for continuous treatment. The kit's high-elasticity TPE silicone layer and memory function titanium alloy ensure the corrector maintains its shape and effectiveness over time.

# Significant Effects

Premium memory titanium alloy patches relieve ingrown toenails pain within 1~3 days and completely correct nail shape within 3~5 weeks.

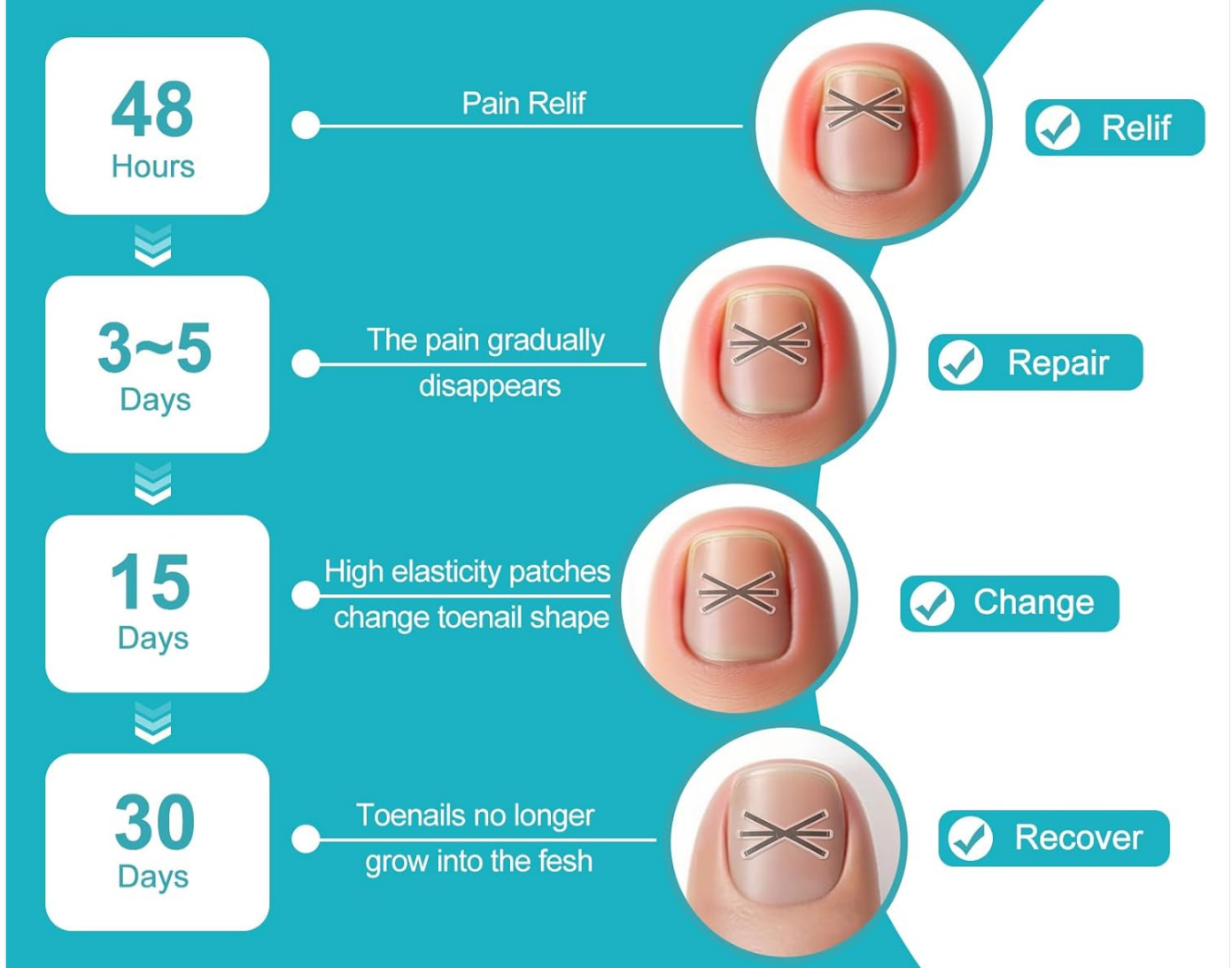


Figure 4: Expected timeline for pain relief and toenail correction.

## TROUBLESHOOTING

- **Patch not sticking:** Ensure the toenail surface is thoroughly cleaned with an alcohol pad and roughened with the polishing strip before applying the glue. Apply sufficient glue and press firmly for 15-20 seconds.
- **Discomfort after application:** If you experience discomfort, observe the area for 2 hours. If discomfort persists, adjust the patch's position or remove it. Avoid vigorous exercise for 24 hours after application.
- **Glue on skin:** If glue gets on your skin, do not panic. Wipe with paper towels, apply the Diluent to the affected area, wait 2 minutes, and gently scrub to remove.
- **No immediate results:** Toenail correction is a gradual process. Significant relief is typically observed within 15 days, with complete correction in about 30 days. Consistent application is key.

## SPECIFICATIONS

Feature	Detail
Brand	VITASY
Model	B0FPJM97QV
Correction Patch Material	Memory Function Titanium Alloy, High-elasticity TPE Silicone Layer
Correction Patch Sizes	0.51in, 0.59in, 0.67in, 0.79in (2 pieces per size)
Package Dimensions	5.35 x 3.82 x 1.26 inches
Item Weight	2.89 ounces
Unit Count	1.0 Count (Kit)

## WARRANTY AND SUPPORT

---

For any questions or concerns regarding the VITASY Ingrown Toenail Treatment Kit, please contact VITASY customer service. While specific warranty details are not provided, the manufacturer is committed to solving foot problems and offers support for product usage.