



Manuals.plus /

› Marsyu /

› Marsyu Smartwatch User Manual: ECG, Blood Pressure, Uric Acid, Lipids, BMI, 1.97" AMOLED Display

Marsyu B0FRLVPJP8

Marsyu Smartwatch User Manual

Model: B0FRLVPJP8

1. INTRODUCTION

This manual provides detailed instructions for the setup, operation, and maintenance of your Marsyu Smartwatch. Please read this manual thoroughly to ensure proper use and to maximize the features of your device.



Image 1.1: Marsyu Smartwatch with both strap options.

2. PACKAGE CONTENTS

Verify that all items are present in your package:

- Smartwatch x 1
- Charging Cable x 1
- User Manual x 1
- Pink Silicone Strap x 1
- Gold Milanese Steel Strap x 1

Your browser does not support the video tag.

Video 2.1: Unboxing the Smartwatch and its contents.

3. INITIAL SETUP

3.1 Charging the Device

Before first use, fully charge your smartwatch. Connect the provided charging cable to the charging port on the back of the watch and to a USB power source. A full charge typically takes approximately 2.5 hours. The watch has a 350 mAh battery capacity.

3.2 App Installation and Pairing

To access all features and synchronize data, download the **H Band** application on your smartphone. The app is compatible with iOS 10.0+ and Android 5.0+ devices.

1. Scan the QR code in the manual or search for "H Band" in your smartphone's app store.
2. Install and open the H Band app.
3. Ensure Bluetooth is enabled on your smartphone.
4. In the app, follow the instructions to search for and pair with your Marsyu Smartwatch. The device uses Bluetooth 5.0 for stable connection.

4. BASIC OPERATION

4.1 Touchscreen and Buttons

The smartwatch features a 1.97-inch AMOLED HD touchscreen with a resolution of 390x450 pixels. Navigate by swiping across the screen and tapping icons. The side button can be rotated for navigation and pressed for selection or to return to the home screen.

4.2 Bluetooth Calls and Notifications

Once paired with your smartphone via Bluetooth, the smartwatch supports stable calls and notifications:

- **Making/Receiving Calls:** You can answer, make calls, add contacts, and view call history directly from the watch.
- **Message Notifications:** Receive notifications from various applications directly on your watch.
- **SOS Function:** Press and hold the side button for three seconds to dial a pre-set emergency contact in the H Band app.

4.3 Smart Features

The smartwatch includes several convenient features:

- **Voice Assistant:** Interact with your phone's voice assistant directly from the watch.
- **Music Control:** Control music playback on your smartphone.
- **Camera Control:** Use the watch as a remote shutter for your phone's camera.
- **Find Phone:** Locate your misplaced smartphone by triggering an alert from your watch.
- **Weather Forecast:** View current weather conditions and forecasts.
- **Alarm Clock:** Set alarms directly on the watch.
- **Calculator:** Basic calculator functionality.
- **Sedentary Reminder:** Get alerts to move after periods of inactivity.
- **Water Reminder:** Receive reminders to drink water throughout the day.

5. HEALTH MONITORING

The Marsyu Smartwatch offers comprehensive health monitoring features. **Note: The data provided by this device is for informational purposes only and should not be used for medical diagnosis or treatment. Consult a medical professional for any health concerns.**

5.1 ECG Monitoring

Equipped with an ECG button, the watch can generate an ECG report in 60 seconds with a single touch. The data can be analyzed and viewed in the H Band app, providing insights into heart health and detecting abnormal heart rate alerts.

Your browser does not support the video tag.

Video 5.1: Demonstration of ECG and other health monitoring features (functionality similar to Marsyu Smartwatch).



Image 5.1: ECG monitoring and app integration.

5.2 Non-Invasive Uric Acid and Lipids Monitoring

The smartwatch supports non-invasive monitoring of uric acid and blood lipid levels, providing real-time tracking for comprehensive blood health insights. This feature is for informational purposes only and not for medical use.



Image 5.2: Uric acid and lipid monitoring.

5.3 BMI, Emotion, and Fatigue Measurement

The watch can measure body composition, including BMI, by pressing the ECG button for 30 seconds. It also features emotion and fatigue measurement functions, allowing you to monitor your emotional state, stress levels, and fatigue degree effectively within 30 seconds.

Monitorización de ECG

Recopile señales de ECG a través de tecnología de onda de pulso inteligente, genere un informe de ECG en 60 segundos y vea el ECG y el estado cardíaco en la aplicación móvil.



Image 5.3: BMI and body composition measurement.

5.4 24/7 Heart Rate Monitoring

The smartwatch continuously monitors your heart rate throughout the day, providing timely reminders if any abnormalities are detected.



Image 5.4: 24-hour heart rate monitoring.

5.5 Body Temperature Monitoring

Automatic 24-hour body temperature monitoring generates a report of temperature changes, helping you track your body's thermal trends.

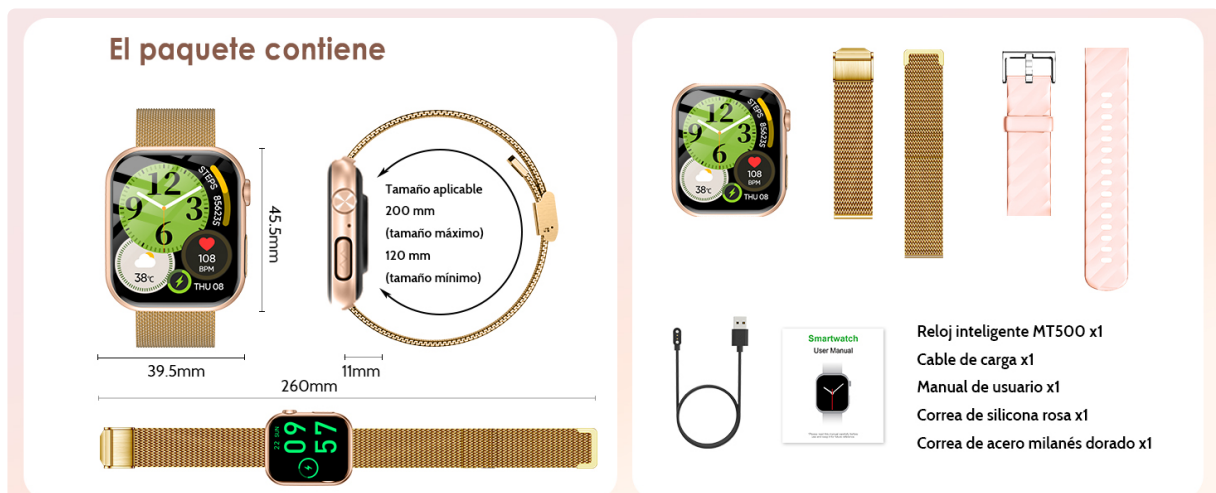


Image 5.5: Body temperature monitoring.

5.6 Sleep Monitoring

The watch registers deep sleep, light sleep, and nocturnal blood oxygen levels, quantifying your sleep quality to help you understand and improve your rest patterns.



Image 5.6: Sleep monitoring features.

5.7 Female Health Tracking

The smartwatch provides detailed data for female health, including menstrual period tracking, prediction, and ovulation tracking, helping users stay informed about their body's changes.



Image 5.7: Female health tracking.

6. FITNESS MODES AND WATER RESISTANCE

6.1 Sports Modes

The smartwatch includes over 50 sports modes. It accurately records daily steps, mileage, calorie consumption, and heart rate during exercise. Examples of supported activities include cycling and various other exercises.



Image 6.1: Multiple sports modes and activity tracking.

6.2 Water Resistance (IP67)

The device is IP67 water-resistant, suitable for daily use and workouts. **It is not recommended for use during bathing, prolonged underwater activities, or in high-temperature environments like saunas.**

7. CUSTOMIZATION

7.1 Watch Faces

Through the H Band app, you can download over 100 watch faces or customize exclusive dials using your own photos, allowing you to change styles at any time.



Image 7.1: Customizable watch faces.

7.2 Interchangeable Straps

The smartwatch comes with two interchangeable straps: an elegant gold stainless steel strap and a casual pink silicone strap, allowing you to switch between classic and sporty styles.



Image 7.2: Pink silicone strap and gold Milanese strap options.

8. SPECIFICATIONS

Brand	Marsyu
Model Number	B0FRLVJP8
Display Size	1.97 Inches
Display Type	AMOLED
Resolution	390 x 450 pixels
Battery Capacity	350 mAh
Battery Composition	Lithium Polymer
Connectivity Technology	Bluetooth 5.0
Memory Capacity	128 MB
Water Resistance	IP67
Compatible Devices	Smartphone (iOS 10.0+, Android 5.0+)
Product Weight	190 g

9. MAINTENANCE

To ensure the longevity and optimal performance of your Marsyu Smartwatch, follow these maintenance guidelines:

- **Cleaning:** Regularly wipe the watch screen and straps with a soft, dry, lint-free cloth. For stubborn dirt, slightly dampen the cloth with water. Avoid using harsh chemicals or abrasive materials.
- **Charging:** Use only the provided charging cable. Ensure the charging contacts on the watch and cable are clean and dry before charging.
- **Storage:** Store the watch in a cool, dry place away from direct sunlight and extreme temperatures when not in use.
- **Avoid Impacts:** Protect the watch from strong impacts or drops, which can damage the screen or internal components.

- **Water Exposure:** While IP67 rated, avoid prolonged submersion, hot water, or chemical exposure (e.g., soaps, detergents) which can compromise the water resistance seals.

10. TROUBLESHOOTING

If you encounter issues with your smartwatch, refer to the following common solutions:

- **Watch not turning on:** Ensure the watch is fully charged. Connect it to the charger and wait a few minutes before attempting to turn it on again.
- **Unable to pair with phone:**
 - a. Ensure Bluetooth is enabled on your phone and the watch.
 - b. Make sure the watch is within the Bluetooth range of your phone.
 - c. Restart both your phone and the smartwatch.
 - d. Check if the H Band app is updated to the latest version.
 - e. If previously paired, try unpairing and re-pairing the device through the app.
- **Inaccurate health data:**
 - a. Ensure the watch is worn snugly on your wrist, not too tight or too loose.
 - b. Clean the sensors on the back of the watch.
 - c. Remember that data is for informational purposes only and not medical-grade.
- **Notifications not appearing:**
 - a. Verify that notifications are enabled in both your phone's settings and the H Band app.
 - b. Ensure the watch is connected via Bluetooth.
 - c. Check if your phone is in 'Do Not Disturb' mode.

11. WARRANTY AND SUPPORT

For warranty information, please refer to the documentation included with your purchase or contact the retailer. If you require further assistance or have questions not covered in this manual, please contact Marsyu customer support through the contact information provided on the product packaging or the official Marsyu website.