

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [OMAGISM](#) /

› [OMAGISM ZF-M909 Rechargeable Behind-the-Ear Digital Hearing Aids User Manual](#)

OMAGISM ZF-M909

OMAGISM ZF-M909 Rechargeable Behind-the-Ear Digital Hearing Aids User Manual

Model: ZF-M909 | Brand: OMAGISM

1. INTRODUCTION

Thank you for choosing the OMAGISM ZF-M909 Rechargeable Behind-the-Ear Digital Hearing Aids. This manual provides essential information for the proper setup, operation, and maintenance of your new hearing aids. Please read this manual thoroughly before using the device to ensure optimal performance and longevity.

These hearing aids are designed to provide clear and natural sound amplification with advanced noise-cancelling technology, offering a comfortable and discreet solution for improved hearing.




Image: The OMAGISM ZF-M909 hearing aids with their portable charging case.

2. PACKAGE CONTENTS

Please verify that all items listed below are included in your package:

- 2 x OMAGISM ZF-M909 Hearing Aids
- 1 x Portable Charging Case
- 1 x USB Charging Cable
- Multiple Pairs of Silicone Earplugs (various sizes)
- 2 x Cleaning Brushes
- 1 x User Manual (this document)

What's In The Package

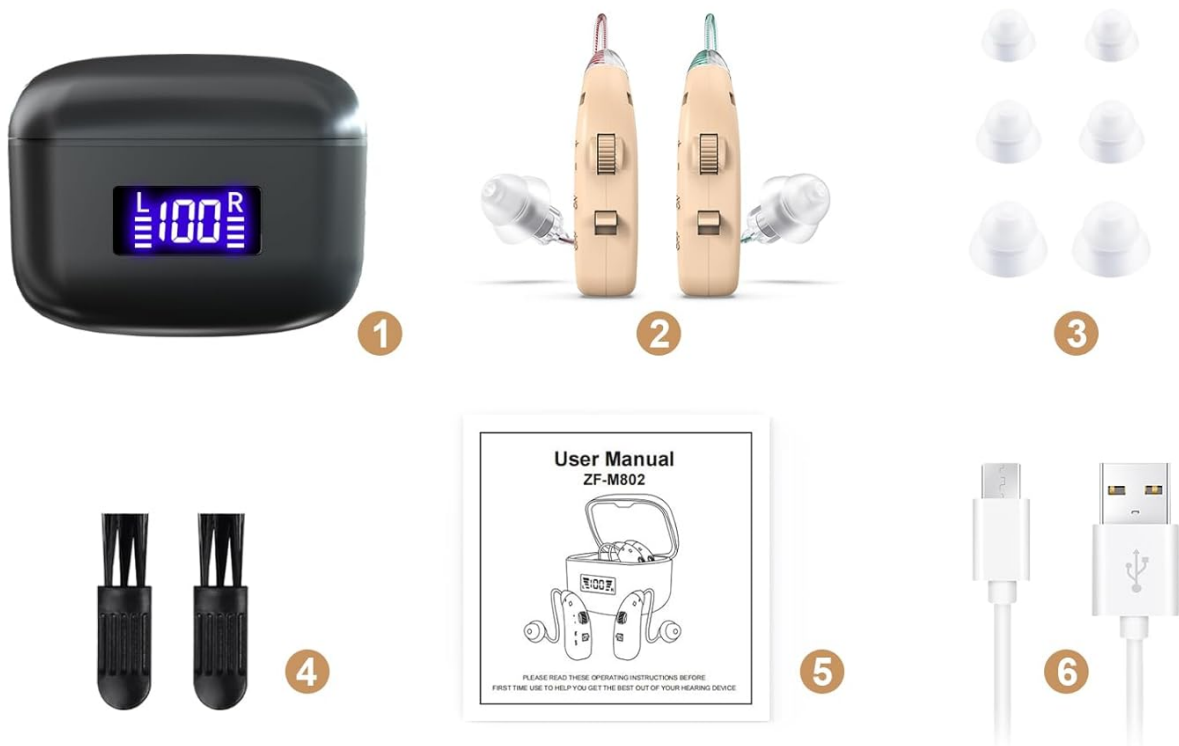

1. Charging Case

2. Hearing Aids

3. Earplugs

4. Cleaning brush

5. Instructions

6. USB Cable

Image: A visual representation of all components included in the product package.

3. PRODUCT OVERVIEW

Familiarize yourself with the different parts of your OMAGISM ZF-M909 hearing aids and charging case:

Rechargeable Behind-The-Ear Hearing aids



Image: Detailed diagram illustrating the sound tube, earplug, microphone, volume control, ON/OFF switch on the hearing aid, and the charging indicators (Left, Right, and overall) on the charging case.

- **Sound Tube:** Connects the hearing aid to the earplug.
- **Earplug:** Sits inside the ear canal to deliver amplified sound.
- **Microphone:** Captures ambient sounds.
- **Volume Control:** A scroll wheel for adjusting sound amplification levels.
- **ON/OFF Switch:** Powers the hearing aid on or off.
- **Charging Point:** Connects the hearing aid to the charging case.
- **Charging Case:** Stores and recharges the hearing aids, with indicators for battery status.

4. SETUP

Follow these steps for initial setup and wearing your hearing aids:

1. **Initial Charge:** Before first use, fully charge the hearing aids and the charging case. Refer to the 'Charging' section for details.

2. **Select Ear Tips:** Choose the silicone earplug size that provides the most comfortable and secure fit in your ear canal. A good seal is crucial for optimal sound quality and to prevent whistling.
3. **Place in Ear:** Gently place the hearing aid behind your ear. Insert the earplug firmly into your ear canal. Ensure the sound tube is positioned correctly.
4. **Power On:** Locate the ON/OFF switch on the hearing aid and slide it to the 'ON' position.
5. **Adjust Volume:** Use the scroll wheel volume control to adjust the sound to a comfortable level. Start with a lower volume and gradually increase it.



Simple And Intuitive

Designed for Ease



Remove from Case:
Ready for use.



Select Ear Tips:
Find your perfect fit.



Place in Ear:
Comfortable and secure fit.



Turn on/Adjust Volume:
Use the dials for perfect sound.

Image: A visual guide demonstrating the simple and intuitive steps for using the hearing aids, from removing them from the case to adjusting the volume.

5. CHARGING

Your OMAGISM ZF-M909 hearing aids come with a portable charging case for convenient recharging.

Charging the Hearing Aids:

- Place the hearing aids into their respective slots in the charging case. Ensure they are correctly seated.
- The charging case display will show the charging status for each hearing aid (L for left, R for right) and the case's battery level.
- A full charge takes approximately 3 hours. Once fully charged, the hearing aids can operate continuously for up to 48 hours.

Charging the Charging Case:

- Connect the provided USB charging cable to the Type-C port on the charging case and to a standard USB power adapter (not included) or a computer USB port.
- The charging case can provide up to 140 hours of additional battery life for the hearing aids, allowing for multiple recharges on the go.
- The digital display on the case indicates its current battery percentage.

Fast Charging & Long Lasting

Stay Connected All-Day Long

 24Hrs Working Time	 96Hrs Powered By Charging Case	 3Hrs Charging Time
--	--	---

Connector
Type-C



Image: The charging case connected to a power source, highlighting its fast charging capabilities and extended battery life.

6. OPERATION

Power On/Off:

- To turn on, slide the ON/OFF switch to the 'ON' position.
- To turn off, slide the ON/OFF switch to the 'OFF' position. Turning off the hearing aids when not in use conserves battery life.

Volume Adjustment:

- Use the small scroll wheel on the hearing aid to adjust the volume.
- Roll the wheel upwards to increase volume and downwards to decrease volume.
- Adjust the volume to a comfortable level that allows you to hear clearly without discomfort.

Sound Modes:

The hearing aids offer three different modes to adapt to various environments. While specific details on mode selection are not provided, typically these modes are cycled through by a short press of a button (often the ON/OFF switch or a dedicated mode button if present). Experiment to find the best mode for your current listening situation:

- **Normal Mode:** General amplification for everyday environments.
- **Noise Reduction Mode:** Enhanced filtering for noisy environments like restaurants or crowded places.
- **Outdoor Mode:** Optimized for open spaces, reducing wind noise and enhancing distant sounds.

Easy To Use

Clear and intuitive on off switch, and volume control via scroll wheel, offering a clear and easy way to check device status



Greater capacity



One-Button Adjustment



Wear Comfortable



Operation Simple

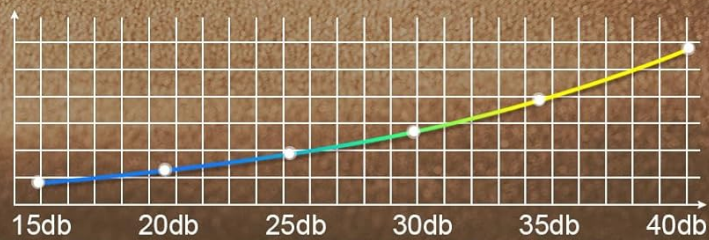


Image: The hearing aid's user-friendly design, emphasizing its capacity, ease of adjustment, comfort, and simple operation, alongside a graph showing sound amplification levels.

7. MAINTENANCE

Regular cleaning and proper care will extend the life and performance of your hearing aids.

Daily Cleaning:

- After each use, gently wipe the hearing aids with a soft, dry cloth to remove any earwax or debris.
- Use the provided cleaning brush to carefully clear any wax or dirt from the earplug and microphone openings.
- Ensure the sound tube and earplug are free from blockages.

Storage:

- Always store your hearing aids in the charging case when not in use. This protects them from damage and keeps them charged.
- Store in a cool, dry place, away from direct sunlight and extreme temperatures.

General Care:

- Avoid exposing the hearing aids to water, moisture, or excessive humidity.
- Do not use alcohol, solvents, or abrasive cleaners.
- Keep out of reach of children and pets.

8. TROUBLESHOOTING

If you encounter issues with your OMAGISM ZF-M909 hearing aids, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
No sound or very faint sound	<ul style="list-style-type: none">• Hearing aid is off• Low battery• Volume too low• Earplug blocked by earwax• Sound tube blocked	<ul style="list-style-type: none">• Turn on the hearing aid• Charge the hearing aid• Increase volume• Clean the earplug• Clean the sound tube
Whistling or feedback sound	<ul style="list-style-type: none">• Improper fit of earplug• Earplug too small• Earwax buildup in ear canal• Hand or object too close to hearing aid	<ul style="list-style-type: none">• Re-insert the earplug firmly• Try a larger earplug size• Clean your ear canal• Ensure nothing is obstructing the microphone
Sound is distorted or unclear	<ul style="list-style-type: none">• Low battery• Microphone blocked• Moisture in the hearing aid	<ul style="list-style-type: none">• Charge the hearing aid• Clean the microphone opening• Place in a dry environment or use a drying kit
Charging case not charging	<ul style="list-style-type: none">• USB cable not connected properly• Faulty USB cable or power adapter	<ul style="list-style-type: none">• Ensure secure connection• Try a different USB cable or power adapter

If the problem persists after trying these solutions, please contact customer support for further assistance.

9. SPECIFICATIONS

Feature	Detail
---------	--------

Feature	Detail
Model Number	ZF-M909 (Model Name: ZF-M802)
Dimensions (Hearing Aid)	Approx. 1.38 x 0.61 x 0.39 inches (3.5 x 1.55 x 1 cm)
Weight (Hearing Aid)	Approx. 4.5 g
Battery Type	Lithium Ion (included)
Charging Time	Approx. 3 hours for full charge
Continuous Working Time	Up to 48 hours (per charge)
Charging Case Battery Life	Provides up to 140 hours of additional hearing aid life
Peak Gain	Up to 38 dB
Noise Reduction	Automatic noise reduction and feedback suppression
Power Source	Battery Powered
Color	Skin

10. SAFETY INFORMATION

Please observe the following safety precautions:

- Consult a medical professional if you experience sudden hearing loss, ear pain, or discharge before using hearing aids.
- Do not attempt to modify or repair the device yourself. Contact customer support for service.
- Avoid using the hearing aids in extremely loud environments without proper ear protection, as this may cause further hearing damage.
- Keep the charging case and hearing aids away from water and liquids.
- Dispose of batteries and electronic components according to local regulations.
- If you experience any discomfort or adverse reactions while using the hearing aids, discontinue use and consult a medical professional.

11. WARRANTY AND SUPPORT

OMAGISM is committed to providing high-quality products and customer satisfaction. If you encounter any problems with your hearing aid, please do not hesitate to contact us. Our customer service team is ready to assist you with troubleshooting, warranty claims, or any questions you may have.

Please refer to your purchase documentation or the seller's contact information for details on warranty coverage and how to reach customer support.

For assistance, please visit our official website or contact the seller directly through the platform where you purchased the product.



© 2025 OMAGISM. All rights reserved.