

LAMYOO 915

LAMYOO Digital Rechargeable Invisible Hearing Aids Model 915 User Manual

Model: 915

INTRODUCTION

This manual provides detailed instructions for the proper use, care, and maintenance of your LAMYOO Digital Rechargeable Invisible Hearing Aids, Model 915. Please read this manual thoroughly before using the device to ensure optimal performance and longevity.

PACKAGE CONTENTS

Verify that all items are present in your package:

- LAMYOO Digital Rechargeable Invisible Hearing Aids (Left and Right)
- Portable Charging Case
- USB Charging Cable
- Various Sizes of Ear Domes
- Cleaning Brush
- User Manual (this document)



Figure 1: LAMYOO Digital Rechargeable Invisible Hearing Aids and Charging Case.

SETUP

1. Initial Charging

Before first use, fully charge the hearing aids and the charging case. The portable charging case provides power for the hearing aids and displays battery levels.

1. Connect the USB charging cable to the charging case and a power source.
2. Place the hearing aids into their respective slots in the charging case. Ensure proper alignment.
3. The LED display on the charging case will indicate the charging status. A full charge for the hearing aids takes approximately 100 minutes, providing up to 100 hours of continuous use. The charging case holds enough power for five full charges of the hearing aids.

Extra-Long Battery Life

Better Hearing Throughout the Day

No more fumbling with tiny button batteries. Hearing aids are ideal for those with dexterity issues or who are simply tired of constantly buying and replacing hearing aid batteries.

2-Hour

Charge Time

30 Hours

Of Use on a Single Charge

150 Hours

With Charging Case



Figure 2: Hearing aids in charging case, illustrating battery life and charging times.

2. Selecting Ear Domes

Choose the ear dome size that provides a comfortable and secure fit in your ear canal. A proper fit is crucial for optimal sound quality and to prevent feedback.

3. Inserting the Hearing Aids

Gently insert the hearing aid into your ear canal. Ensure it sits snugly and comfortably. The device is designed to be discreet and lightweight for all-day wear.

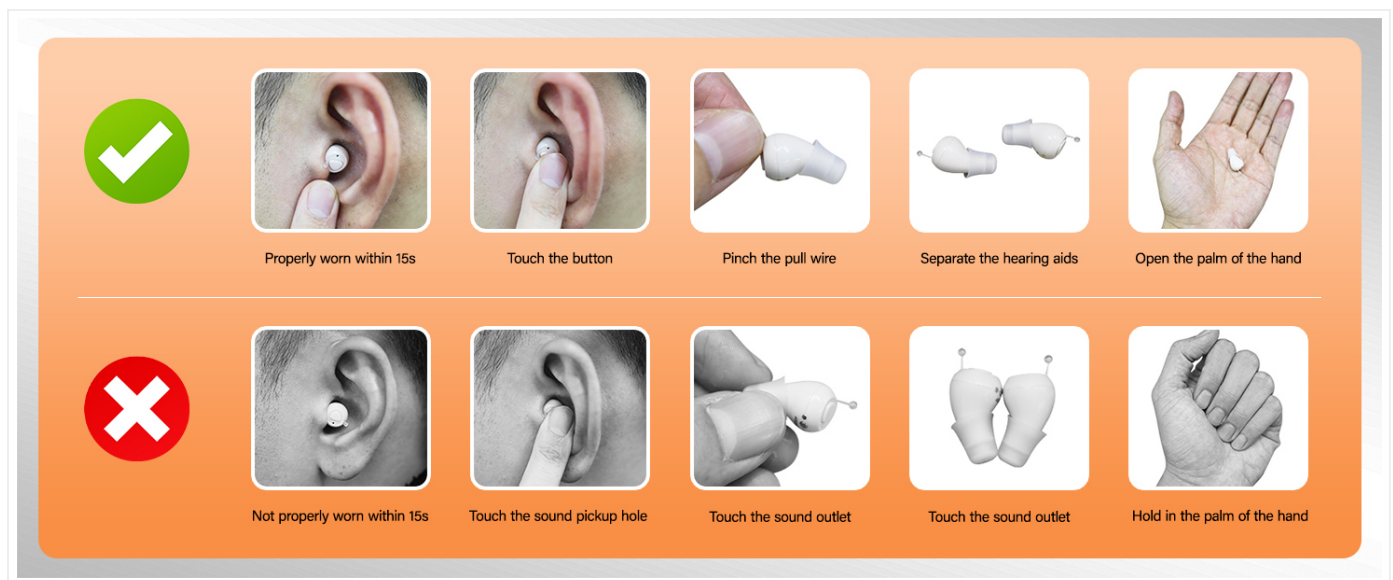


Figure 3: Visual guide for proper insertion and handling of the hearing aids.

OPERATING INSTRUCTIONS

1. Power On/Off

The hearing aids feature a simple one-button control. Press and hold the button on each hearing aid to power it on or off.

2. Volume Adjustment

Briefly press the control button to cycle through the four adjustable volume settings. The device includes a memory function that retains your last selected volume level.

3. Intelligent Noise Reduction

The advanced digital chip technology in your LAMYOO hearing aids automatically minimizes background noise, enhancing clarity for conversations and other audio experiences. This feature operates continuously to provide clear sound in various environments.

Intelligent Background Noise Reduction



Designed with advanced chip technology, it optimizes human speech and background noise, providing clear and distinguishable sound even in the most challenging acoustic environments.

Figure 4: The intelligent noise reduction feature helps users hear clearly in various settings.

MAINTENANCE

1. Daily Cleaning

After each use, gently wipe the hearing aids with a soft, dry cloth. Use the provided cleaning brush to remove any earwax or debris from the ear domes and sound outlets. Regular cleaning prevents blockages and maintains sound quality.

2. Storage

When not in use, store the hearing aids in their charging case. This protects them from damage and keeps them charged for your next use. Store in a cool, dry place away from direct sunlight and extreme temperatures.

3. Battery Care

To maximize battery life, avoid fully discharging the hearing aids and charging case frequently. Recharge them regularly, even if not completely depleted.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No Sound / Weak Sound	<ul style="list-style-type: none">• Low battery• Ear dome blocked by earwax• Volume too low• Improper fit	<ul style="list-style-type: none">• Recharge the hearing aids.• Clean the ear dome and sound outlet with the brush.• Increase the volume using the control button.• Adjust the hearing aid for a secure fit or try a different ear dome size.
Whistling / Feedback Sound	<ul style="list-style-type: none">• Improper fit in ear canal• Volume too high	<ul style="list-style-type: none">• Re-insert the hearing aid for a tighter seal.• Try a different size ear dome.• Lower the volume.
Hearing Aid Not Charging	<ul style="list-style-type: none">• Charging cable not properly connected• Charging contacts are dirty• Charging case battery is depleted	<ul style="list-style-type: none">• Ensure the USB cable is securely connected to both the case and power source.• Clean the charging contacts on both the hearing aids and the case.• Charge the portable charging case.

SPECIFICATIONS

Feature	Detail
Model Name	915
Brand	LAMYOO
Color	Carbon Black
Material	Acrylonitrile Butadiene Styrene (ABS)
Special Features	Adjustable, Lightweight, Noise Cancellation, Rechargeable, Volume Control
Battery Type	Lithium Polymer (included)
Battery Life (Hearing Aids)	Up to 100 hours per charge
Charging Time (Hearing Aids)	Approximately 100 minutes
Charging Case Capacity	Provides 5 full charges for hearing aids
Package Dimensions	4.92 x 3.66 x 1.73 inches
Item Weight	5.93 ounces (package)

WARRANTY AND SUPPORT

LAMYOO is committed to providing quality products and customer satisfaction. For any issues or inquiries regarding your

hearing aids, please contact our customer service team. We aim to respond and assist within 12 hours.
Please refer to your purchase documentation for specific warranty terms and conditions.

