

Dolonm 3 Motor Electric Chair

Dolonm 3 Motor Electric Chair Instruction Manual

Model: 3 Motor Electric Chair

1. SAFETY INFORMATION

Please read all safety instructions carefully before operating the Dolonm 3 Motor Electric Chair. Failure to follow these instructions may result in injury or damage to the product.

- Always assemble according to the instructions provided.
- Keep fingers and loose clothing away from moving parts during operation.
- Do not exceed the maximum weight limit of 159 pounds (72 kg).
- Use caution when adjusting the chair to prevent tip-over. Ensure the chair is stable before use.
- Keep children and pets away from the chair when in use or during assembly.
- Ensure the power cord is not pinched or damaged. Disconnect power before cleaning or maintenance.

2. PRODUCT OVERVIEW

The Dolonm 3 Motor Electric Chair is a versatile and adjustable treatment chair designed for various professional settings, including medical exams, facial treatments, spa services, eyelash extensions, tattooing, podiatry, and massage. It features a robust 3-motor system for independent control of the backrest, leg rest, and seat height, along with an upgraded hand remote for easy operation.



Figure 2.1: Dolom 3 Motor Electric Chair, showcasing its design and functionality.

Key Features:

- **Efficient 3-Motor Adjustable System:** Independently controls backrest angle, leg rest angle, and seat height for customized positioning.
- **Automatic One-Touch Control:** Upgraded hand remote allows effortless adjustments with a single press to start and stop movement.
- **Versatile Applications:** Suitable for facials, eyelash extensions, tattoos, massages, and podiatry care.
- **Sturdy Structure & Comfortable Surface:** Reinforced steel base for stability, soft PU leather, and high-density sponge for comfort. Water-resistant and easy to clean.
- **Removable Pillow:** The headrest features a removable pillow with a face hole, allowing for adjustable head distance and optimal client positioning.

3. SETUP INSTRUCTIONS

The Dolonm 3 Motor Electric Chair is shipped in two secure packages: one for the bed surface and one for the base. Follow these steps for quick and easy assembly.

3.1 Unpacking

1. Carefully open both packages. Ensure all components are present and undamaged.
2. Remove all packing materials. Keep them for future transport if needed.



Figure 3.1: The Dolonm Electric Chair is shipped in two separate packages for safe delivery.

3.2 Assembly

1. Place the base of the chair on a flat, stable surface.
2. Carefully align the bed surface with the connection points on the base.
3. Secure the bed surface to the base using the provided hardware. Ensure all connections are tight and secure.
4. Connect the hand remote control cable to the designated port on the chair.
5. Plug the power cord into a suitable electrical outlet.

3.3 Leveling the Chair

The chair is equipped with four high-load leveling feet to ensure stability on uneven surfaces.

- To lower the height of a corner, turn the leveling foot to the right (clockwise).
- To raise the height of a corner, turn the leveling foot to the left (counter-clockwise).
- Adjust each foot until the chair is perfectly level and stable. This is crucial for safety and proper operation.



Figure 3.2: Adjusting the 4 high-load leveling feet to ensure chair stability.

4. OPERATING INSTRUCTIONS

The Dolonm 3 Motor Electric Chair features an intuitive hand remote for precise control over its three independent motors.

4.1 Hand Remote Control

The upgraded hand remote allows for one-touch adjustments. Simply press a button once to initiate movement and press it again to stop at the desired angle or height.



The thickened PU leather is water- and oil-resistant, easy to clean, and remains soft and elegant after long-term use, ideal for daily spa and beauty care environments.

Figure 4.1: The hand remote provides one-touch control for automatic position changes.

4.2 Adjusting the Backrest

The backrest can be adjusted to an angle of up to 75 degrees, allowing for various client positions from upright to reclined.



Figure 4.2: Backrest adjustable function, demonstrating the range of motion.

4.3 Adjusting the Leg Rest

The leg rest can be adjusted up to 45 degrees, providing comfortable support for the client's legs.



Figure 4.3: Footrest adjustable function, illustrating its range of movement.

4.4 Adjusting the Height

The chair's height can be adjusted from 23 inches to 32.5 inches, allowing practitioners to maintain an ergonomic posture during treatments.



Figure 4.4: Height adjustable function, detailing the chair's dimensions and height range.

4.5 Removable Pillow and Headrest Adjustment

The chair features a removable pillow and an adjustable headrest to accommodate various client needs and treatment types.

- The headrest distance can be adjusted by sliding it in or out.
- The removable pillow reveals a face hole, ideal for massage or other prone treatments.



EASY MAINTENANCE

The removable armrests are designed for effortless cleaning and maintenance. Simply detach them for thorough cleaning after each use. The premium PU surface is smooth, waterproof, and easy to wipe down, helping maintain a hygienic and professional workspace.



Figure 4.5: The removable pillow and adjustable headrest provide versatile client positioning.

4.6 Usage Examples

The chair's adjustability makes it suitable for a wide range of professional applications:



Figure 4.6: Configuration for Facial Spa treatments, allowing clients to lie in a lower position for better access.

Smart One-Touch Auto Adjustment-Designed for professional efficiency the hand remote activates automatic position changes with a single press no need for long holding. Simply press again to stop at the desired angle. This feature is ideal for facial, medical, tattoo, and podiatry treatments where quick and smooth adjustments enhance workflow and client comfort.



Figure 4.7: Configuration for Tattoo application, with independent backrest and leg rest adjustments.

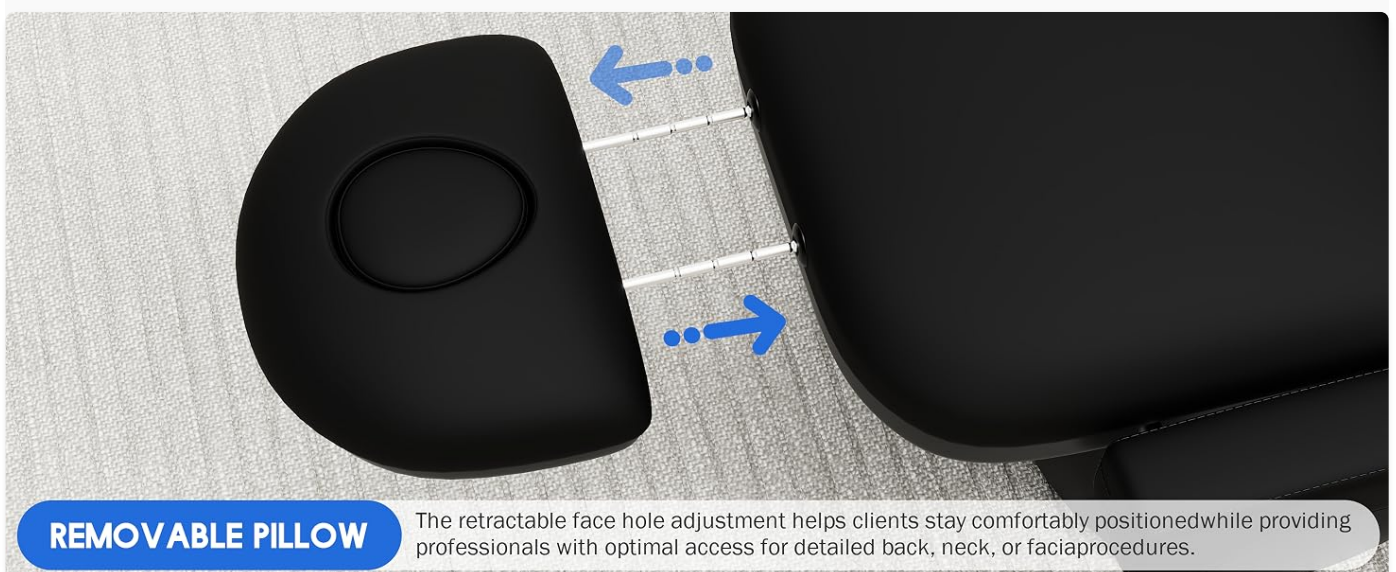


Figure 4.8: Configuration as a Treatment Exam Chair, supporting ergonomic posture for practitioners and comfort for patients.

4 HIGH LOAD LEVELING FEET

They provide a solid and stable base for the chair. Whether the floor is slightly uneven due to construction flaws, wear and tear, or other factors, the leveling feet can be adjusted to ensure that the chair stands firmly. This is especially important during treatments such as tattooing, eyelash extensions, or medical examinations where precision and a stable work surface are crucial. If the chair were to wobble, it could disrupt the procedure and potentially lead to mistakes or discomfort for the client.



Figure 4.9: Configuration for Podiatry Treatment, with a sturdy frame and smooth motorized lift system.

4.7 Operational Video (Reference Only)

This video demonstrates the general operation and adjustability of a Dolonm electric chair. **Please note: This video features a previous model. The remote control shown in the video has been upgraded to the new version, as displayed in the images above.**

Your browser does not support the video tag.

Video 4.1: Reference video demonstrating the adjustability of a Dolonm electric chair. (Duration: 0:46)

5. MAINTENANCE

Proper maintenance ensures the longevity and hygienic operation of your Dolonm 3 Motor Electric Chair.

5.1 Cleaning the Upholstery

- The chair's soft PU leather surface is water-resistant and oil-resistant.
- Clean regularly with a damp cloth and mild, non-abrasive cleaner.
- Avoid harsh chemicals or solvents that may damage the material.
- Wipe dry immediately after cleaning.

5.2 Removable Armrests

The armrests are designed to be removable for thorough cleaning and maintenance.

- Simply detach the armrests after each use for easy cleaning.
- Clean the armrests using the same method as the main upholstery.
- Ensure armrests are securely reattached before next use.



Figure 5.1: Removable armrests facilitate easy cleaning and maintenance.

5.3 General Care

- Regularly inspect all moving parts and connections for wear or looseness.
- Keep the chair away from direct sunlight and extreme temperatures to prevent material degradation.
- Do not attempt to lubricate internal motor components unless specified by the manufacturer.

6. TROUBLESHOOTING

If you encounter any issues with your Dolonm 3 Motor Electric Chair, refer to the following troubleshooting guide.

Problem	Possible Cause	Solution
Chair does not respond to remote control.	1. Power cord not connected. 2. Remote control cable loose or disconnected. 3. Power outlet issue.	1. Ensure power cord is securely plugged into the chair and a working outlet. 2. Check the remote control cable connection to the chair. 3. Test the power outlet with another device.
Motor makes unusual noise or movement is jerky.	1. Obstruction in moving parts. 2. Overload. 3. Loose connections.	1. Inspect for any obstructions and remove them. 2. Ensure the weight limit is not exceeded. 3. Check all assembly connections for tightness.
Chair is unstable or wobbles.	1. Uneven floor surface. 2. Leveling feet not adjusted correctly.	1. Move the chair to a flat, even surface. 2. Adjust the four high-load leveling feet until the chair is stable (refer to Section 3.3).

If the problem persists after attempting these solutions, please contact Dolonm customer support for further assistance.

7. SPECIFICATIONS

Detailed technical specifications for the Dolonm 3 Motor Electric Chair.

Feature	Specification
Brand	Dolonm
Model	3 Motor Electric Chair
Color	Black
Item Weight	159 Pounds (72 kg)
Material Type	Metal (Frame), PU Leather (Upholstery)
Item Dimensions (L x W x H)	76.8 x 24.5 x 32.5 inches (195 x 62.2 x 82.5 cm)
Adjustable Height Range	23 - 32.5 inches (58.4 - 82.5 cm)
Backrest Adjustment Angle	Up to 75 degrees
Leg Rest Adjustment Angle	Up to 45 degrees
Assembly Required	Yes

8. WARRANTY AND SUPPORT

Dolonm is committed to providing high-quality products and customer satisfaction.

8.1 Warranty Information

For specific warranty details regarding your Dolonm 3 Motor Electric Chair, please refer to the warranty card included with your purchase or visit the official Dolonm website. The warranty typically covers manufacturing defects under normal use conditions.

8.2 Customer Support

If you have any questions, require technical assistance, or need to report an issue, please contact Dolonm customer support:

- **Website:** Visit the official Dolonm Store on Amazon or their corporate website for FAQs and support resources.
- **Email:** Refer to your product packaging or warranty documentation for customer service email addresses.
- **Phone:** Contact numbers may be available on the official website or product documentation.

When contacting support, please have your product model (3 Motor Electric Chair) and purchase information readily available.



