

[Manuals.plus](#) /

› [Medcursor](#) /

› Medcursor TENS Unit Replacement Pads SM9187 Instruction Manual

Medcursor SM9187

Medcursor TENS Unit Replacement Pads

Model: SM9187

IMPORTANT SAFETY INFORMATION

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.


These replacement pads are designed for use with Medcursor TENS units only. Do not use with other brands.

PRODUCT OVERVIEW

The Medcursor TENS Unit Replacement Pads are designed to provide effective and comfortable electrical stimulation when used with compatible Medcursor TENS units. These pads feature an upgraded biocompatible hydrogel for reliable conductivity and strong adhesion, ensuring optimal performance during therapy sessions.

Upgraded Gel Formula - More reuses

MedPulse pads feature a biocompatible hydrogel that provides superior conductivity, strong adhesion, and a gentle, skin-friendly feel.

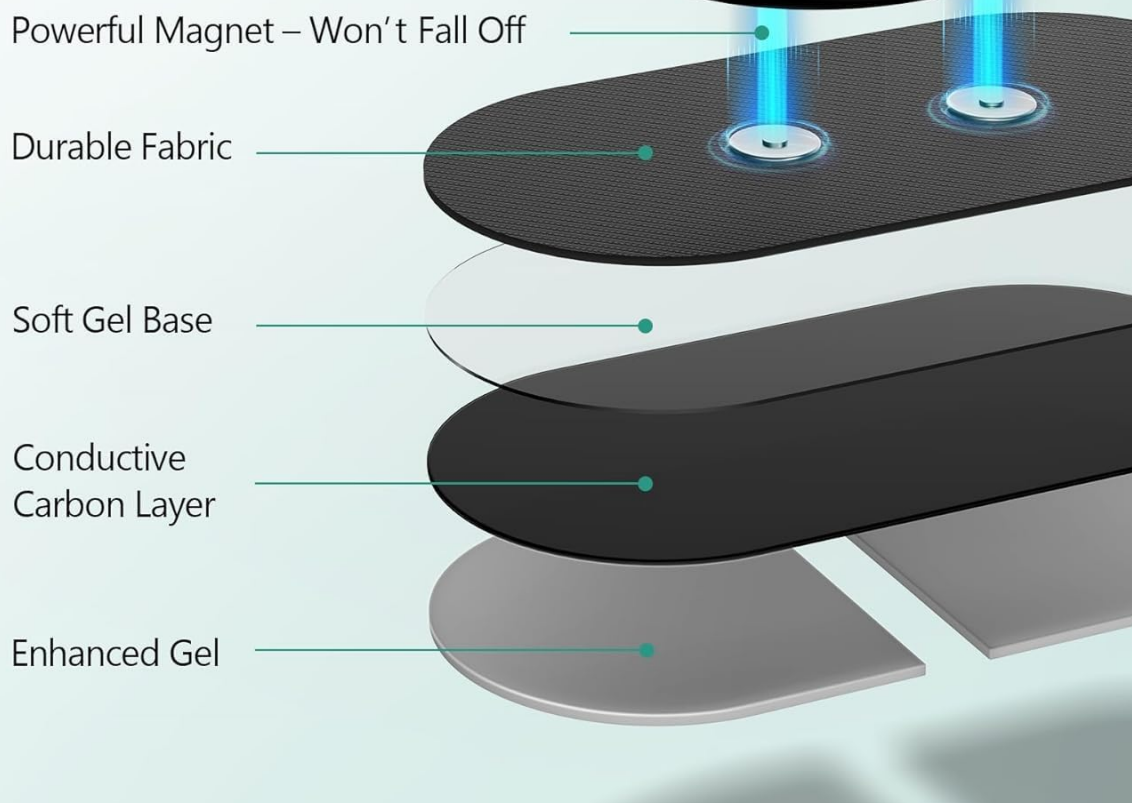

This diagram illustrates the advanced construction of Medcursor TENS unit electrode pads. It highlights the enhanced gel for superior conductivity, the conductive carbon layer, a soft gel base, and a durable fabric backing. The powerful magnetic snap ensures a secure connection to your TENS unit.

Key Features:

- **Upgraded Biocompatible Hydrogel:** Delivers reliable conductivity and strong, skin-friendly adhesion.
- **Powerful Magnetic Snap:** Ensures a secure connection to your TENS unit, maintaining stability even during movement.
- **Durable Construction:** Made from high-quality, latex-free, and skin-safe materials with a durable cloth backing.
- **Reusable:** With proper care, pads are reusable for over 30 sessions.

PACKAGE CONTENTS

Each package contains a value pack of 8 electrode pads, specifically:

- 16 Round Pads (2" x 2")
- 8 Rectangular Pads (2" x 3.5")

SETUP AND APPLICATION

Follow these steps for proper application of your Medcursor TENS Unit Replacement Pads:

1. **Prepare Skin:** Ensure the skin area where the pads will be applied is clean and dry. Moisture or sweat can affect adhesion.
2. **Attach Pads:** Firmly attach the electrode pads to flat skin. Ensure the magnetic snap connects securely to your Medcursor TENS unit leads.

Easy to Use & Longer-Lasting Pads

Proper use helps extend pad lifespan



1 Clean and dry your skin before use.



2 Attach pads firmly to flat skin.



3 After use, rinse pads with water and let them air dry.



4 Stick pads back on the plastic film and store in a cool, dry place.

This visual guide provides clear instructions for the proper use and care of your Medcursor TENS unit electrode pads. It demonstrates how to prepare your skin, apply the pads, clean them after use, and store them correctly to maximize their lifespan and performance.

OPERATING INSTRUCTIONS

Once the pads are securely attached to your skin and connected to your Medcursor TENS unit, operate your TENS unit according to its specific instruction manual. The powerful magnetic clasp on these pads is designed to maintain a secure connection throughout your therapy session, even during light movement.

MAINTENANCE AND STORAGE

Proper care extends the lifespan and performance of your electrode pads:

- **Restoring Stickiness:** If adhesion decreases, gently rinse the pads with water. Carefully remove any surface residue. Allow the pads to air-dry completely before re-applying or storing. Do not wipe with abrasive materials.
- **Storage:** After use and cleaning, stick the pads back onto their original plastic film. Store them in the resealable pouch in a cool, dry place.
- **Reusability:** With proper cleaning and storage, these pads are reusable for over 30 sessions.

TROUBLESHOOTING

Reduced Adhesion:

- Ensure skin is clean and dry before application.
- Rinse pads with water and air-dry to restore stickiness as described in the Maintenance section.
- Avoid applying pads to areas with excessive hair or lotions.

SPECIFICATIONS

Specification	Detail
Model Number	SM9187
Pad Material	Biocompatible Hydrogel, Durable Cloth Backing
Pad Sizes Included	16 Round (2" x 2"), 8 Rectangular (2" x 3.5")
Compatibility	Medcursor TENS Units Only
Reusability	Over 30 sessions with proper care
Package Dimensions	9.06 x 5.98 x 0.63 inches
Item Weight	7.37 ounces

SUPPORT AND WARRANTY

Medcursor stands by the quality of its TENS unit pads. For any questions or concerns regarding your replacement pads, please contact Medcursor customer service.

These replacement pads are HSA/FSA eligible.

Contact Information: Please refer to the Medcursor official website or your product packaging for the most current customer service contact details.