

VEVOR FBL650506

VEVOR FBL650506 Swivel Shower Chair Instruction Manual

Model: FBL650506 | Brand: VEVOR

1. INTRODUCTION

Thank you for choosing the VEVOR Swivel Shower Chair. This product is designed to provide a safe and comfortable bathing experience for individuals with limited mobility, including the elderly, pregnant women, those recovering from surgery, and people with reduced mobility. Please read this manual thoroughly before assembly and use to ensure proper installation and safe operation.

Safety Information

- Always ensure all components are securely fastened before each use.
- Do not exceed the maximum weight capacity of 100 kg (220 lbs).
- Ensure the chair is stable and properly seated on the bathtub edges before transferring onto it.
- Use caution when transferring to and from the chair. Seek assistance if necessary.
- Regularly inspect the chair for any signs of wear, damage, or loose parts. Discontinue use if any issues are found.
- Keep children away from the shower chair unless supervised by an adult.



Figure 1: VEVOR Swivel Shower Chair in use, demonstrating ease of entry and exit with a 220 lbs (100 kg) static load capacity.

2. PRODUCT OVERVIEW AND COMPONENTS

The VEVOR Swivel Shower Chair is constructed from reinforced aluminum alloy and durable plastic, designed for stability and comfort. It features a 360° swivel seat, armrests, and a backrest.

Package Contents:

- 1 x Swivel Shower Chair Frame with Seat and Backrest
- 4 x Rubber Sleeves (for feet)
- 4 x Slotted Knobs
- 4 x C-Shaped Clips
- 2 x Threaded Rods
- 1 x User Manual (this document)

Item Model Number: **FBL650506**

Net Weight: **9.3 lbs / 4.2 kg**

Product Dimensions:

31.1 × 21.3 × 18.5 in / 790 × 540 × 470 mm

Accessories:

4 x Rubber Sleeves

4 x Slotted Knobs

4 x C-Shaped Clips

2 x Threaded Rods

1 x User Manual

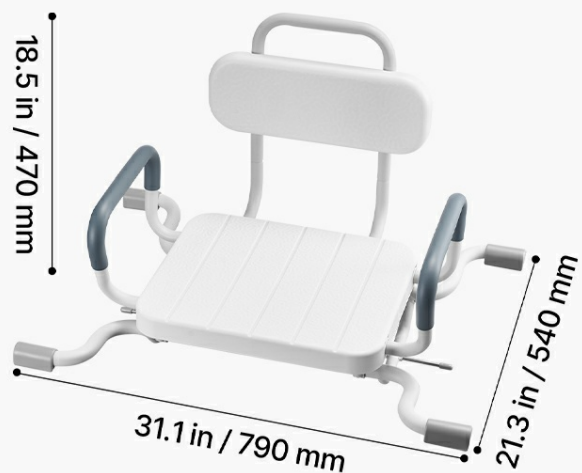


Figure 2: Product components and dimensions. Please verify all parts are present before assembly.

3. SETUP AND ASSEMBLY INSTRUCTIONS

Follow these steps to assemble and install your VEVOR Swivel Shower Chair securely.

1. **Unpack Components:** Carefully remove all parts from the packaging and verify against the package contents list (Section 2).
2. **Attach Rubber Sleeves:** Slide the four rubber sleeves onto the ends of the chair's legs to prevent slipping and protect your bathtub surface.
3. **Assemble Frame:** Connect the main frame components. Ensure all connections are firm.
4. **Install Seat and Backrest:** Secure the seat and backrest to the frame using the provided hardware (slotted knobs, C-shaped clips, threaded rods). Ensure they are tightly fastened.
5. **Adjust Bathtub Width:** The chair is designed to fit bathtubs with outer widths of 740 mm (29.1 inches) or 790 mm (31.1 inches). Adjust the chair's width by sliding the frame sections to match your bathtub's dimensions.
6. **Position in Bathtub:** Carefully place the assembled chair onto the bathtub edges. Ensure the four large non-slip feet are firmly seated on the tub's rim.
7. **Secure Locking Mechanism:** Locate the two clamps under the seat. These clamps press against the inner wall of

the bathtub. Tighten the locking mechanism to secure the chair firmly in place. This prevents movement during use.

8. **Final Check:** Before first use, gently push and pull on the chair to confirm it is stable and does not wobble. Retighten any loose connections if necessary.



Figure 3: Bathtub mounting and adjustment details. Ensure proper fit and secure locking.

4. OPERATING INSTRUCTIONS

Using your VEVOR Swivel Shower Chair is straightforward. Follow these guidelines for safe and comfortable use.

1. Transferring onto the Chair:

- Ensure the chair is securely locked in position (refer to Section 3).
- Use the armrests for support as you carefully sit down on the seat.
- Avoid sudden movements.

2. Using the Swivel Seat:

- The seat can rotate 360 degrees. To swivel, disengage the locking mechanism (if applicable, some models have a continuous swivel or a simple push-button release).
- Gently rotate the seat to the desired position. The swivel feature allows for easy transfer into and out of the tub without needing to shift your body.
- Once in the desired position, engage the locking mechanism to secure the seat.

3. Using Armrests and Backrest:

- The ergonomic armrests provide additional support when sitting down, standing up, or repositioning.
- The comfortable backrest offers lumbar support during bathing.

4. Exiting the Bathtub:

- Swivel the seat to face the exit direction.
- Use the armrests for support as you carefully stand up and transfer out of the tub.

Figure 4: The 360° swivel seat facilitates easy transfers.

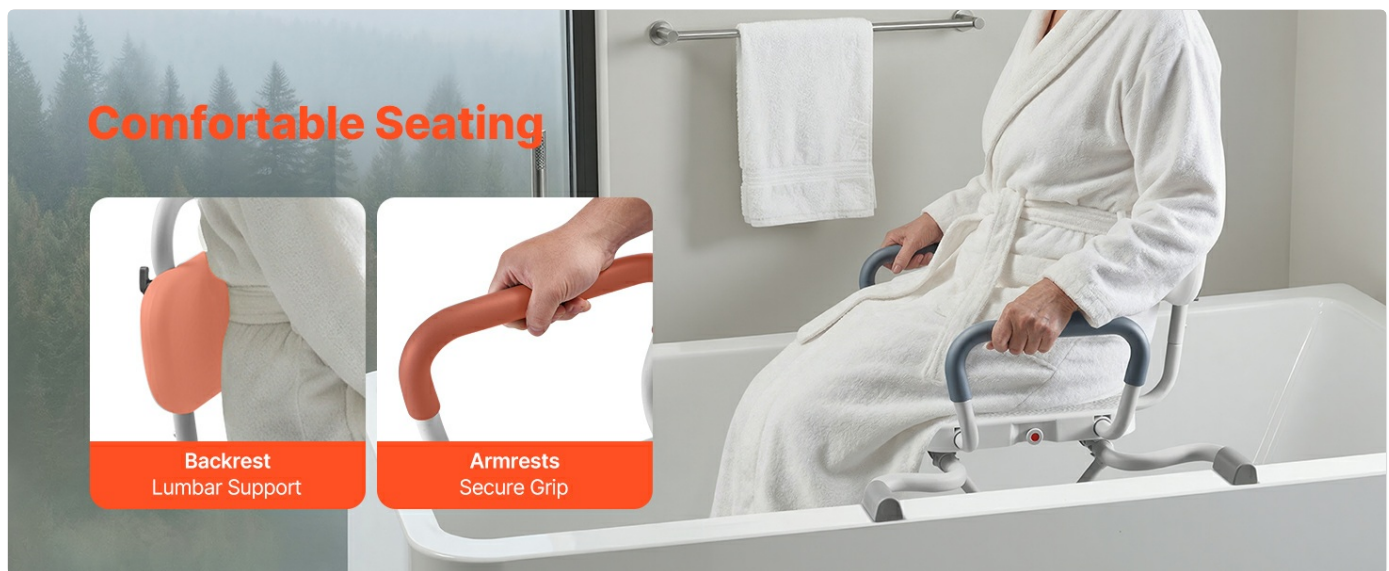


Figure 5: Ergonomic design for comfortable and secure seating.

5. MAINTENANCE

Proper maintenance ensures the longevity and safe operation of your VEVOR Swivel Shower Chair.

- **Cleaning:** Clean the chair regularly with mild soap and water. Avoid abrasive cleaners or harsh chemicals that could damage the materials. Wipe dry after cleaning to prevent water spots and mildew.
- **Inspection:** Periodically check all bolts, screws, and connections to ensure they are tight. Inspect the rubber sleeves for wear and tear. Replace any damaged parts immediately.
- **Storage:** If not in regular use, store the chair in a dry, cool place away from direct sunlight.

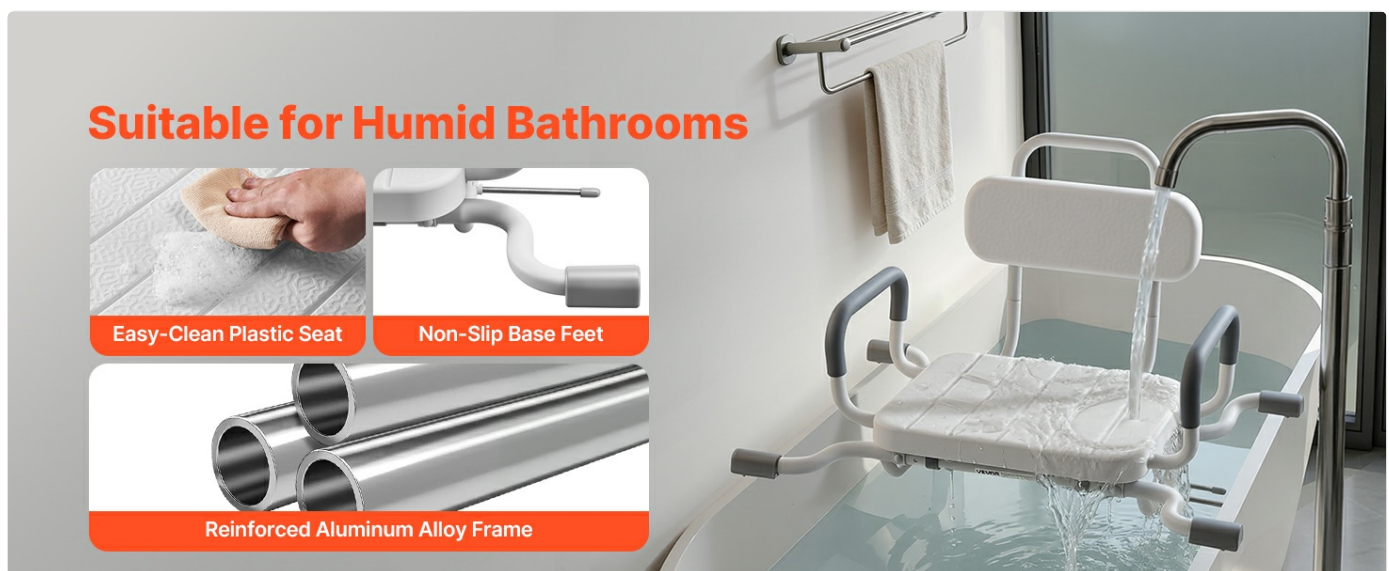


Figure 6: The chair's design is suitable for humid bathroom environments and is easy to clean.

6. TROUBLESHOOTING

If you encounter any issues with your VEVOR Swivel Shower Chair, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Chair feels unstable or wobbly.	Not properly seated on bathtub edges. Locking mechanism not engaged or loose. Bathtub width not correctly adjusted. Loose assembly bolts.	Ensure all four non-slip feet are firmly on the tub rim. Tighten the locking mechanism under the seat. Re-adjust the chair's width to match your bathtub (refer to Section 3). Check and tighten all assembly bolts and knobs.
Seat does not swivel smoothly.	Swivel mechanism is dirty or obstructed. Locking mechanism is engaged.	Clean the swivel mechanism area. Ensure the swivel locking mechanism is disengaged before attempting to rotate.
Chair slides on bathtub.	Non-slip feet are worn or missing. Bathtub surface is excessively slippery.	Replace worn rubber sleeves. Ensure the bathtub surface is clean and free of excessive soap residue.

If the problem persists or you encounter an issue not listed here, please contact VEVOR customer support for assistance.

7. SPECIFICATIONS

Feature	Detail
Model Number	FBL650506
Brand	VEVOR
Material	Reinforced Aluminum Alloy, Durable Plastic
Color	Grey + White
Net Weight	4.2 kg (9.26 lbs)
Static Load Capacity	Up to 100 kg (220 lbs)
Product Dimensions (L x W x H)	790 x 540 x 470 mm (31.1 x 21.3 x 18.5 inches)
Adjustable Bathtub Width	Fits outer widths of 740 mm (29.1 inches) or 790 mm (31.1 inches)
UPC	197988183655

Item Model Number: **FBL650506**

Net Weight: **9.3 lbs / 4.2 kg**

Product Dimensions:

31.1 × 21.3 × 18.5 in / 790 × 540 × 470 mm

Accessories:

4 x Rubber Sleeves

4 x Slotted Knobs

4 x C-Shaped Clips

2 x Threaded Rods

1 x User Manual

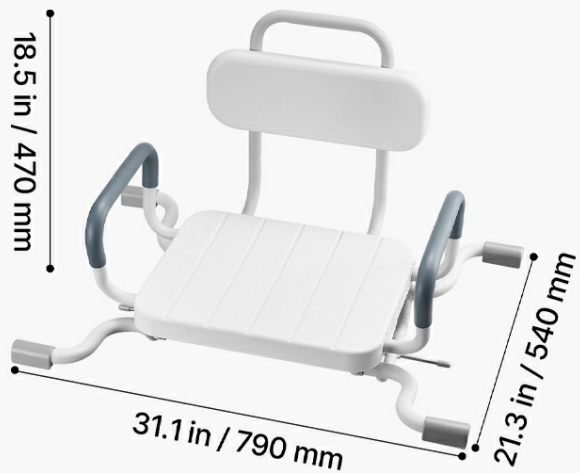


Figure 7: Detailed product dimensions.

8. WARRANTY AND SUPPORT

VEVOR is committed to providing high-quality products. For warranty information, product support, or any questions not covered in this manual, please contact VEVOR customer service through the retailer's platform or visit the official VEVOR website.

Please have your model number (FBL650506) and purchase date ready when contacting support.

© 2026 VEVOR. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question