

VOCIC Z35

VOCIC Z35 PackGo Rollator Instruction Manual

Model: Z35

1. INTRODUCTION

Thank you for choosing the VOCIC Z35 PackGo Rollator. This lightweight and portable rollator is designed to provide mobility assistance and support for seniors, offering stability and convenience for daily use and travel. Please read this manual thoroughly before assembly and use to ensure safe and proper operation.

2. SAFETY GUIDELINES

- **Weight Capacity:** Do not exceed the maximum user weight of 300 lbs (136 kg).
- **Brake System:** Always engage the foot brakes when sitting on the rollator or when stationary to prevent accidental movement.
- **Stability:** Ensure the rollator is fully unfolded and all locking mechanisms are secure before use.
- **Terrain:** Use the rollator on flat, stable surfaces. Avoid uneven terrain, stairs, or escalators.
- **Guidance:** The rollator rolls with ease. Guide it gently and maintain a steady, controlled pace.
- **Supervision:** Individuals with limited cognitive or physical abilities should use the rollator under supervision.
- **Maintenance:** Regularly inspect all parts for wear and tear. Replace damaged components immediately.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a mobility aid and not a medical treatment device.

3. PACKAGE CONTENTS

Verify that all components are present in the package:

- Z35 Main Frame * 1

- Front Wheel * 2
- Rear Wheel * 2
- VOCIC Backpack * 1 (Extra Accessory)

First Ever 4-Wheel Removal One-Click, Wheel Off



Figure 3.1: Included components of the VOCIC Z35 PackGo Rollator.

4. ASSEMBLY INSTRUCTIONS

4.1 Unfolding the Rollator

1. Unfold the front legs to a vertical position.
2. Position the rollator upright.
3. Push the lever right and pull the Z35 frame out to fully extend.
4. Press the seat down until the buckle fully engages, securing the frame in the open position.



Figure 4.1: Steps to unfold the rollator.

4.2 Wheel Installation

Attach the front and rear wheels to their respective positions on the rollator frame. Ensure they click securely into place.


Your browser does not support the video tag.

Video 4.1: This video demonstrates the complete assembly process for the Z35 PackGo Rollator, including wheel installation.

4.3 Handle Height Adjustment

The handles are adjustable to 7 different levels, accommodating users from 4'11" to 5'11" (150 cm to 180 cm). To adjust the handle height:

1. Locate the adjustment pins or buttons on the handle shafts.
2. Press the pins/buttons and slide the handles up or down to the desired height.
3. Ensure the pins/buttons fully engage in the adjustment holes to lock the handles securely in place.

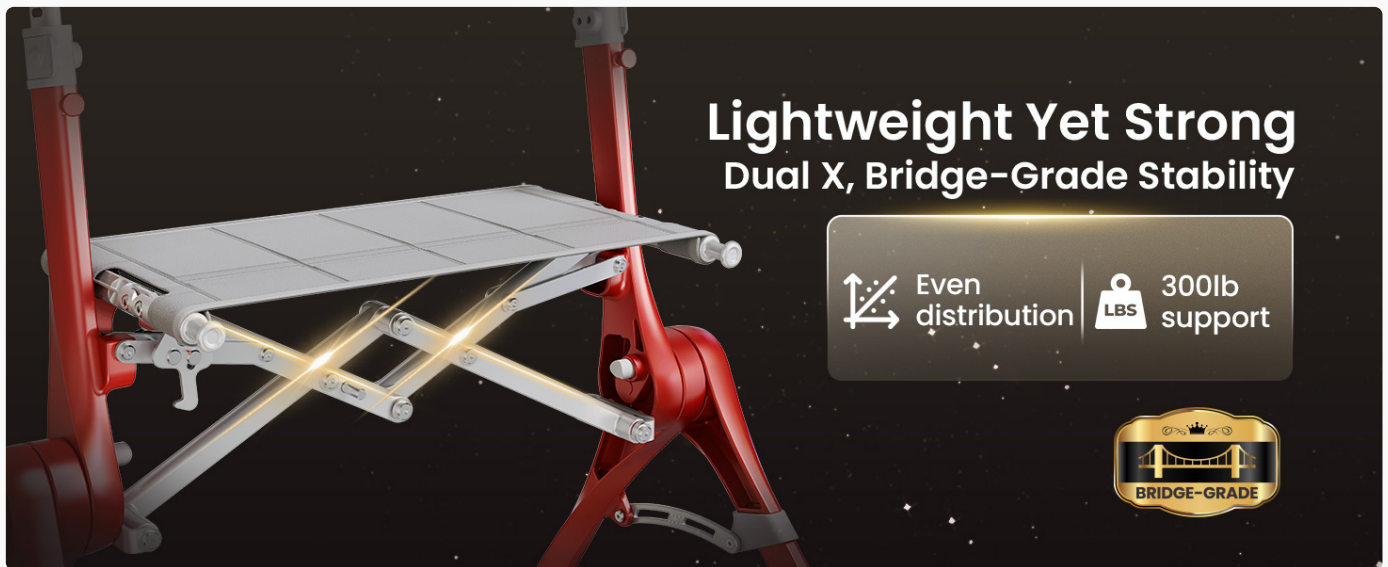

Frame Material Magnesium Alloy



Figure 8.1: Detailed technical specifications and dimensions.

9. WARRANTY AND SUPPORT

The VOCIC Z35 PackGo Rollator comes with 3 years of coverage, long-term support, and professional after-sales care. For any questions, concerns, or warranty claims, please contact VOCIC customer service through the retailer where you purchased the product or visit the official VOCIC website.



Figure 9.1: VOCIC offers comprehensive support and warranty for its products.

© 2026 VOCIC. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question